



**Secure
Care**
Products, LLC

**Please contact your
Distributor /
Installer for service ...**

Tel.: _____

Transmitter User Guide

Adult Wandering ID and Non-ID Transmitters

39 Chenell Drive
Concord, NH USA 03301-8501
Phone: (800) 451-7917/(603) 223-0745
Fax: (603) 227-0200
<http://www.securecare.com>

© 2021 Secure Care Products®, LLC.

CONTENT IS SUBJECT TO CHANGE WITHOUT NOTICE

TABLE OF CONTENTS

IMPORTANT NOTICES	5
1. Secure Care is Not Responsible for the Locks	5
2. Secure Care Is Not Responsible for The Computer Hardware.	5
3. Several Factors Outside the Secure Care System Can Affect its Performance	6
4. The Secure Care System Must be Properly Installed	6
5. Performance of the Secure Care System Software Depends on Proper Maintenance	6
6. Only a Qualified Service Technician Should Work on a Secure Care System	6
7. Only a Authorized Distributor or Installer can Install a Secure Care System	7
8. Any Work Must Comply with Electrical and Life Safety Codes	7
9. Immediately Have Replacements or Repairs Checked On-Site by a Qualified Service Technician	7
10. The Secure Care System is not a Substitute for Careful Identification and Monitoring by Professional Staff	7
11. Warranty	8
12. Acceptance of Order	8
POLICY AND PROCEDURES	9
Ordering New or Replacement Transmitters	9
Transmitter Warranty	9
Daily Transmitter Testing	9
Assembling the Strap and Transmitter	10
Attaching the Transmitter to the Resident	10
Testing Transmitter	10
Removing Transmitter from the Resident	11
Cleaning the Transmitter	11
Using a Transmitter with Secure Care Software	11

Radio Equipment Directive Declaration of Conformity (DoC)

Error! Bookmark not defined.

EU Declaration of Conformity

Error! Bookmark not defined.

1 Concord, New Hampshire 03301 USA
2 Telephone Number: 603.223.0745

Error! Bookmark not defined.
Error! Bookmark not defined.

GENERAL PRODUCT WARRANTY	16
1. Notices	19
2. Limited Warranty	21
3. Limitations of Liability	22
4. Governing Law and Arbitration	22
5. Severability	23
DAILY TRANSMITTER TESTING LOG	24

Figure 1 Transmitter Information.....	10
Figure 2 Strap Assembly	10
Figure 3 Attaching Strap	10
Figure 4 View of transmitter on wrist.....	10
Figure 5 Testing Transmitter	10
Figure 6 Testing 90 Degree Transmitter.....	10
Figure 7 Admit screen in the SCP Software.....	11

IMPORTANT NOTICES

PLEASE READ THIS MANUAL BEFORE BEGINNING THE INSTALLATION OF A SECURE CARE SYSTEM

This installation manual is provided for reference by purchasers and installers of Secure Care Products, LLC (“Secure Care’s”) systems. This manual is not intended as a catalog of warnings for the protection of anyone or as a substitute for obtaining professional training or assistance in the design of a facility’s security procedures and systems, or in the installation, set-up, testing, support, operation, maintenance, repair or use of Secure Care’s systems. Nothing in this manual modifies the terms of Secure Care’s General Product Warranty Statement or of any written agreement signed by Secure Care or creates further warranties or extends benefits of any sort to anyone beyond those already expressly established in Secure Care’s General Product Warranty Statement and in any written contract signed by Secure Care.

1. Secure Care is Not Responsible for the Locks

ALL LOCKS USED WITH SECURE CARE’S SYSTEM ARE DESIGNED, MANUFACTURED, LABELED AND DELIVERED SOLELY BY AN INDEPENDENT VENDOR OVER WHOM SECURE CARE HAS NO CONTROL AND FOR WHOSE ACTIONS OR FAILURES TO ACT SECURE CARE DISCLAIMS ALL RESPONSIBILITY. REGARDLESS OF WHETHER THE LOCKS CARRY SECURE CARE’S LOGO OR NAME OR ANY OTHER TRADEMARK, SERVICE MARK OR TRADE NAME USED OR CLAIMED BY SECURE CARE, SECURE CARE DISCLAIMS ALL WARRANTIES, EXPRESS OR IMPLIED, WITH RESPECT TO THE LOCKS AND/OR THEIR USE WITH OR OPERATION IN THE SECURE CARE SYSTEM, INCLUDING, WITHOUT LIMITATION, ALL IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, TITLE AND/OR NON-INFRINGEMENT. SECURE CARE ALSO DISCLAIMS ALL OBLIGATIONS WITH RESPECT TO THE LOCKS AND/OR THEIR USE WITH OR OPERATION IN THE SECURE CARE SYSTEM THAT MIGHT OTHERWISE ARISE OR BE IMPLIED FROM THE FACT THAT SUCH LOCKS CARRY SECURE CARE’S LOGO OR NAME OR ANY OTHER TRADEMARK, SERVICE MARK OR TRADE NAME USED OR CLAIMED BY SECURE CARE OR FROM THE DELIVERY OR INSTALLATION OF THE LOCKS WITH SECURE CARE SOFTWARE, PARTS AND/OR PRODUCTS OR FROM A COURSE OF DEALING OR USAGE IN TRADE. ALL RESPONSIBILITY FOR DESIGNING, MANUFACTURING, LABELING AND WARNING OF HIDDEN DEFECTS OR DANGERS IN THE LOCKS AND/OR THEIR USE WITH AND OPERATION IN THE SECURE CARE SYSTEM RESTS EXCLUSIVELY WITH THE INDEPENDENT VENDOR, AND ANY CLAIMS, COSTS, DAMAGES OR LIABILITIES ARISING FROM THE LOCKS AND/OR THEIR USE WITH OR OPERATION IN THE SECURE CARE SYSTEM SHALL BE MADE SOLELY AGAINST THE INDEPENDENT VENDOR.

2. Secure Care Is Not Responsible for The Computer Hardware.

IF YOU PURCHASE COMPUTER HARDWARE THROUGH SECURE CARE AND REQUEST THAT SECURE CARE SOFTWARE BE INSTALLED AND TESTED ON THAT HARDWARE AT THE FACTORY, SECURE CARE WARRANTS ONLY THAT THE HARDWARE AND THE SOFTWARE PACKAGES WERE INSTALLED, SET-UP AND TESTED PRIOR TO SHIPMENT IN ACCORDANCE WITH ALL SECURE CARE PRODUCT MANUALS AND THAT, AT THE TIME THE HARDWARE AND THE SOFTWARE PACKAGES WERE FINALLY INSPECTED AT THE FACTORY, THEY WERE PERFORMING (SUBJECT TO SECURE CARE’S SPECIFIED TOLERANCES) IN ACCORDANCE WITH SECURE CARE’S SPECIFICATIONS. SECURE CARE WILL NOT BE RESPONSIBLE FOR ANY DEFECTS IN OR PROBLEMS CAUSED BY THE HARDWARE, ALL CLAIMS FOR WHICH MUST BE MADE TO THE HARDWARE MANUFACTURER AND/OR VENDOR. SECURE CARE DISCLAIMS ALL WARRANTIES, EXPRESS OR IMPLIED, WITH RESPECT TO THE HARDWARE AND/OR ITS USE WITH OR OPERATION IN THE SECURE CARE SYSTEM, INCLUDING, WITHOUT LIMITATION, ALL IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, TITLE AND/OR NON-INFRINGEMENT. SECURE CARE ALSO DISCLAIMS ALL OBLIGATIONS WITH RESPECT TO THE HARDWARE AND/OR ITS USE WITH OR OPERATION IN THE SECURE CARE SYSTEM THAT MIGHT OTHERWISE ARISE OR BE IMPLIED FROM THE FACT THAT SUCH HARDWARE CARRIES SECURE CARE’S LOGO OR NAME OR ANY OTHER TRADEMARK, SERVICE MARK OR TRADE NAME USED OR CLAIMED BY SECURE CARE OR FROM

THE DELIVERY OR INSTALLATION OF THE HARDWARE WITH SECURE CARE SOFTWARE, PARTS AND/OR PRODUCTS OR FROM A COURSE OF DEALING OR USAGE IN TRADE. ALL RESPONSIBILITY FOR DESIGNING, MANUFACTURING, LABELING AND WARNING OF HIDDEN DEFECTS OR DANGERS IN THE HARDWARE AND/OR ITS USE WITH AND OPERATION IN THE SECURE CARE SYSTEM RESTS EXCLUSIVELY WITH THE HARDWARE MANUFACTURER AND/OR VENDOR, AND ANY CLAIMS, COSTS, DAMAGES OR LIABILITIES ARISING FROM THE HARDWARE AND/OR ITS USE WITH OR OPERATION IN THE SECURE CARE SYSTEM SHALL BE MADE SOLELY AGAINST THE HARDWARE MANUFACTURER AND/OR VENDOR.

3. Several Factors Outside the Secure Care System Can Affect its Performance

Secure Care's software, parts and products are designed for operation in a wireless system. However the range, performance, and predictability of any wireless system, including Secure Care's, is dependent on several factors, including, but not limited to, the following: building structure; environmental extremes (e.g., temperature, earth tremors, air pollution, etc.); the proximity of other wireless devices; the presence of variable speed products; sources of Radio Frequency Interference (RFI); physical orientation and positioning of the equipment; and sources of Electro Static Discharge (ESD). Secure Care is not responsible for the effect of these types of factors on operation of its software, parts and products and disclaims all responsibility for any claim relative thereto.

4. The Secure Care System Must be Properly Installed

Secure Care's system must be installed, set-up, tested, supported, operated, maintained, repaired and used only in accordance with all manuals and instructions (including the user, installation, technical and other manuals) issued by Secure Care (the "Product Manuals"). It is your responsibility to assure that any person who might be installing, setting-up, testing, supporting, maintaining or repairing the Secure Care system knows the contents of and has access to the Product Manuals and has successfully completed Secure Care technical training. It is also your responsibility to assure that any person who might be operating or using this Product knows the contents of and has access to the Product Manuals and has successfully completed Secure Care in-service training. Secure Care cannot be responsible for performance problems caused by a failure to follow prescribed and appropriate procedures for installation, set-up, testing, support, operation, maintenance, repair and use.

All adjustable features on new and repaired Secure Care software, parts and products are shipped with "factory default" settings. These "factory default" settings may not comply with building and life safety codes or other applicable laws and regulations in the location where they are installed or operated. Secure Care strongly recommends, therefore, that the settings on all Secure Care software, parts and products be checked and, if necessary, reset to comply with local building and life safety codes and other applicable laws and regulations at the time of any installation, set-up, testing, support, maintenance or repair.

5. Performance of the Secure Care System Software Depends on Proper Maintenance

Secure Care's system is driven by software. However, the performance and reliability of any software-driven system depends on adequately maintaining the recommended minimum configuration of computing platform, operating systems and applications programs and on regularly performing industry-standard and application-specific backup processes. If recommended minimum configurations of computing platform, operating systems, and applications programs are not adequately maintained, or if appropriate backups are not regularly performed, the software may not drive the system as intended. Secure Care is not responsible for operational problems caused by a failure to perform these maintenance and backup procedures and disclaims all responsibility for any claim relative thereto.

6. Only a Qualified Service Technician Should Work on a Secure Care System

Secure Care does not authorize, and strongly recommends against, any installation or field replacement of software, parts or products by untrained contractors or facility staff. Such work can be hazardous, can render the system ineffective and will void any Secure Care warranty or liability that might otherwise relate to the system.

Before any software, parts or products which have been designed and manufactured by Secure Care can be safely installed, set-up, tested, supported, maintained or repaired, technical training in accordance with standards established by Secure Care is required. Regardless of how Secure Care's software, parts or products are obtained, they should not be installed, set-up, tested, supported, maintained or repaired by any person who has not satisfactorily completed that technical training (a "qualified service technician".) When Secure Care's software, parts or products are sold separately from installation services, it is assumed that only a qualified service technician will conduct any installation, set-up, testing, support, maintenance or repair involving that software, part or products.

7. Only a Authorized Distributor or Installer can Install a Secure Care System

Secure Care Products LLC, requires all installations, upgrades or servicing of an existing installation of any and all SCP products or systems to be performed by factory certified Distributors and/or Installers with signed distributor or installer agreements. Customers that gain technical service certifications and maintain their certifications over time are allowed to install replacement hardware and service previous installations conducted by certified distributors for the facilities they own and operate only for the products they are certified to and for the sole purpose of maintenance and repairs. In doing so, they assume liability for those repairs and maintenance. Any individual working for and assisting a company who has not signed a distributor agreement and as a result is not a certified distributor of Secure Care Products LLC will assume all liability of the equipment/system in its entirety. Certified individuals no longer working for a certified Secure Care distributor is no longer considered certified to install/service Secure Care Product equipment/systems. Equipment/systems installed outside the above criteria will void any and all warranty given by Secure Care Products, LLC.

8. Any Work Must Comply with Electrical and Life Safety Codes

It is important that any installation, set-up, testing, support, operation, maintenance, repair or use involving the system comply with all local and national electrical and life safety codes. If you have any questions about compliance with those codes, please contact your local authorities.

9. Immediately Have Replacements or Repairs Checked On-Site by a Qualified Service Technician

Secure Care receives and responds to telephone and dial-in inquiries (the "Help Line") about its software, parts and products for the purpose of discussing users' experiences with Secure Care's system, helping users better understand how their systems work, and providing ideas about what may be causing difficulties. However, Secure Care cannot accurately diagnose the cause of any problems or give complete instructions on how to fix problems over the telephone or Internet. The only way to assure that software, parts or products are installed, set-up, tested, supported, maintained or repaired correctly or that a Secure Care system is functioning properly is to have it examined on site by a qualified service technician. In addition, Secure Care software, parts and products cannot be operated or used correctly by anyone who has not successfully completed Secure Care in-service training. Secure Care's Help Line is not a substitute for on-site diagnosis and servicing by a qualified service technician or for successful completion of Secure Care in-service training. Secure Care strongly recommends that any installation, set-up, testing, support, maintenance or repair of a system that is performed by a person who has not satisfactorily completed technical training in accordance with standards established by Secure Care be immediately checked on-site by a qualified service technician.

WARNING: EVEN SLIGHT MODIFICATIONS TO THE SYSTEM OR CHANGES IN THE OPERATING ENVIRONMENT MAY CAUSE SECURE CARE'S SYSTEM TO MALFUNCTION. THE ONLY WAY TO ASSURE THAT SECURE CARE'S SYSTEM HAS BEEN INSTALLED, SET-UP, TESTED, SUPPORTED, MAINTAINED AND REPAIRED CORRECTLY IS TO HAVE A QUALIFIED SERVICE TECHNICIAN DO THE WORK.

10. The Secure Care System is not a Substitute for Careful Identification and Monitoring by Professional Staff

Secure Care's software, parts and products have been designed to augment a facility's reasonable procedures for protecting residents, patients, and infants. However, no system or combination of procedures and equipment can eliminate all risk or assure complete security. Secure Care's system is not intended as a substitute for the careful identification and monitoring of residents, patients, and infants by a facility's professional staff.

11. Warranty

Departures, substitutes, or declining Secure Care Products supplied hardware will void the product/system warranty and the agency certification. You will assume all liability for any deviations.

12. Acceptance of Order

Acceptance of this order and its terms and conditions, also certifies that the equipment described above has been delivered to and received by the purchaser and is in good operating order and condition, and that the equipment is accepted by the purchaser for all purposes for which the equipment is intended by Secure Care Products. Revised 3/20/17

POLICY AND PROCEDURES

Ordering New or Replacement Transmitters

Transmitters and straps are purchased directly from Secure Care Products. Please call 1-800-451-7917 to order transmitters.

Transmitters and straps will typically be shipped within two weeks after receipt of your order. All transmitter sales are final because they are date sensitive products. Transmitter straps are sold by the dozen for a single use.

Transmitter Warranty

Just prior to shipping all transmitters are labeled with a warranty date indicating their expected end of life. The actual expiration date is the last day of the month labeled on the transmitter. If within this warranty period a transmitter is not performing to our specification, **CALL** our toll free number (1-800-451-7917), **for a TXRA (Transmitter Return Authorization) BEFORE returning any Transmitter.**

A replacement transmitter will be sent along with a postage paid label marked with a Return Authorization number. Simply replace the suspect transmitter with the new replacement and send the suspect transmitter back to Secure Care Products within 15 days with facility name and address, and the TXRA Number displayed, and a description of the transmitter problem. When the transmitter is returned to Secure Care Products, it is tested on a system similar to that in the facility to verify whether or not the transmitter is working to our specifications.

- If the transmitter fails testing, a credit may be given. The value of the credit is calculated by deducting the dollar value of the life remaining (in months) on the transmitter from the cost of the new transmitter. Example: if the transmitter is rate at 12 months and costs \$X, the transmitter cost per month is \$X/12. If the transmitter failed at the end 10 months a credit for 2 months or 2(\$X/12) may be applied.
- If the transmitter tests properly, and is in good working condition, the facility may not receive a credit for the transmitter and may be responsible for any freight charges incurred.
- If the transmitter is physically damaged, the facility will be charged the full replacement transmitter charge.
- If the suspect transmitter is not returned to us within the 15-day period, the facility will be charged the full transmitter amount.

Daily Transmitter Testing

Each day, the aide responsible for the care of the residents utilizing the Secure Care System must ensure that the transmitter is in place. This **must be done at each shift change**. Documentation of this check should be made by the nurse's aide. This should be recorded on the aide's daily checklist for the particular resident. Each transmitter should be tested **daily** to ensure working properly. The Hand Held Tester is available to purchase for fast, easy testing of all transmitters.

A documented test of each transmitter at the facility must be made **each day**. This testing should also include those transmitters not currently in use. The procedure involves using the transmitter Tester and/or the Exit Panel on the wall by the exit, and documenting the performance of the transmitter.

Towards the end of this manual is a daily transmitter testing log, which you may find useful.

Assembling the Strap and Transmitter

Verify the expiration date has not been reached. Route the strap through the slots in the transmitter.

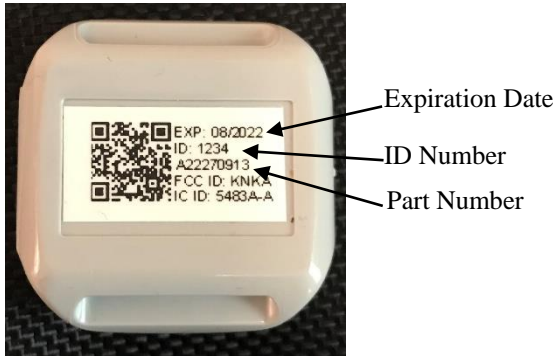


Figure 1 Transmitter Information



Figure 2 Strap Assembly

Attaching the Transmitter to the Resident

Attach the transmitter to the resident's wrist or ankle comfortably (not too tight – keep two finger tips between the band and the wrist or ankle). Insert the post of the strap's snap through the nearest hole in the strap and fold the snap's receptacle down onto the post. Press together until it snaps securely. Cut off excess strap.



Figure 3 Attaching Strap

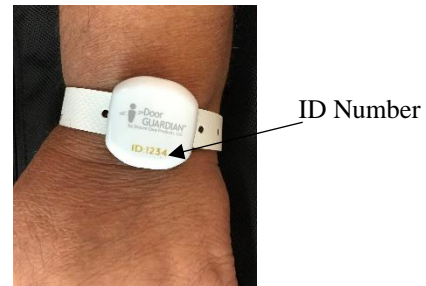


Figure 4 View of transmitter on wrist

Testing Transmitter

Test the Transmitter using the Tester to verify the ID code or if this is a non-ID transmitter the tester should read TX OK as shown. **Figure 5** shows the transmitter to tester orientation for a standard TX. **Figure 6** shows the transmitter to tester orientation for a 90 Degree TX



Figure 5 Testing Transmitter



Figure 6 Testing 90 Degree Transmitter

Removing Transmitter from the Resident

Using a pair of safe scissors, cut the strap and pull out of the transmitter's slot and discard. Clean the transmitter and store in a transmitter containment box for next use.

Cleaning the Transmitter

Secure Care recommends transmitters are cleaned after each use. To properly clean the transmitter use lint free IPA (Isopropyl Alcohol) wipes. With the transmitter removed from the resident, wipe the transmitter on all sides with the IPA wipe. When cleaning is finished store the transmitter in a transmitter containment box.

Using a Transmitter with Secure Care Software

If you have any questions please call your Distributor/Installer first, then if necessary call Secure Care Products Service Department at 1-800-451-7917 or email Order_39@securecare.com.

If you are using the Secure Care Software (SCS), please complete the following steps:


NOTE: New transmitters need to be added to the transmitter list before they will be displayed in the transmitter ID drop down menu.

- 1 Log in to Secure Care Software with your respective User Name and Password.
- 2 Click on "Admit" and fill out all of the required information following the hospital policy and procedure.
- 3 Click on the "Admit/Update" button.
- 4 Click the "Done" button.
- 5 Log out by clicking the "Log Out" button.

The screenshot shows the 'Admit' screen in the Secure Care Software. The interface is divided into a sidebar on the left and a main content area. The sidebar contains a logo at the top and a list of buttons: 'Logged In: admin', 'Log Out', 'Register' (highlighted), 'Release', 'Schedule Absence', 'Create Reports', 'Transmitter List', 'Clear Current Events', 'Staff Permissions', 'Review Current Events', and 'Device Status'. The main content area has a title bar 'Admit' and two tabs: 'Current' and 'View History'. Below the tabs, there are several form sections. The 'New person:' section includes 'Transmitter ID' (a dropdown menu showing 'NONE'), 'Pendant ID' (a dropdown menu showing 'NONE'), 'Last' (a text field), and 'First' (a text field). To the right of these is a 'Photo:' section with a camera icon. Below these are 'Floor' (a dropdown menu showing 'Office'), 'Room Number' (a text field), and 'Bed' (a dropdown menu). Further down are 'Admit Date' (a date/time picker showing '3/10/2021 3:58:03 PM'), 'Identifier (SSN / Phone)' (a text field), and 'Comment' (a text field). At the bottom of the main area is a table titled 'Currently Admitted (click name to select or deselect):'. The table has three columns: 'Name', 'ID', and 'Admit Date'. Below the table is a checkbox labeled 'Add Photo...' and a red asterisk followed by the word 'Required'. At the bottom right of the screen are three buttons: 'Reset', 'Add', and 'Done'. The Secure Care logo is visible in the bottom left corner of the main area.

Figure 7 Admit screen in the SCP Software

TRANSMITTER TROUBLE SHOOTING FLOWCHART

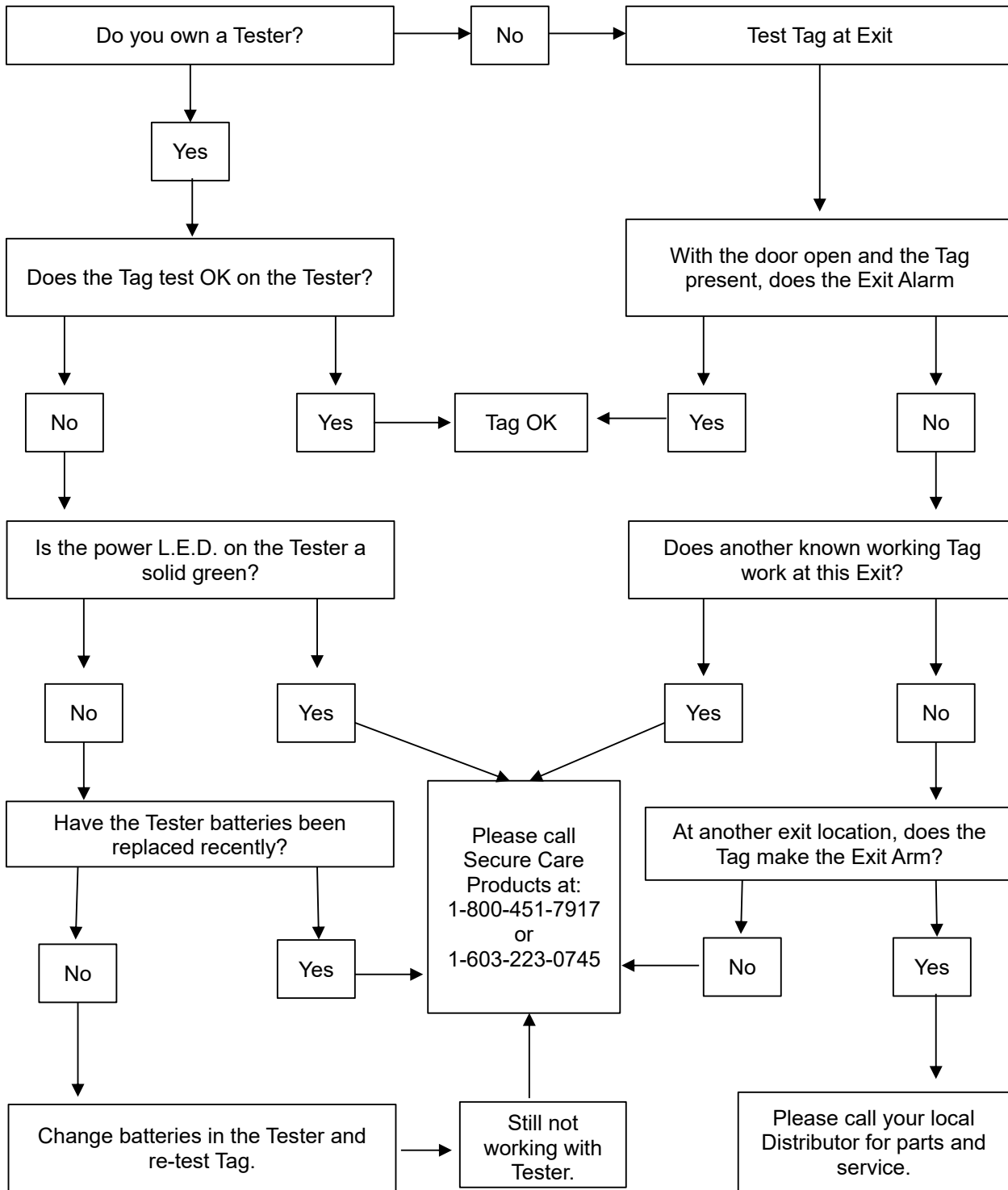

COMPLIANCE INFORMATION

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Radio Equipment Directive Declaration of Conformity (DoC)

Unique identification of the DoC: 13.56/433.92 MHz Transmitter

We, Secure Care Products LLC
39 Chenell Dr. Concord, NH 03301.

Declare under our sole responsibility that the product(s):

Product name(s): 13.56/433.92 Patient Transmitter
Trade name(s): RFID Transmitter
Assembly part#: A20180011

Type or model(s): A22274915, A20240900, A20240901, A20240903, A20270912, A21350910,
A22270912, A21350915, A20240905, A20240904, A20250903, A20270937, A20274935, A22260912,
A22270937, A22274935, A22274916, A22270938, A22270913, A22260913

Relevant supplementary information: Software version 7.154

To which this declaration relates is in conformity with the essential requirements and other relevant requirements of the RED Directive (2014/53/EU). The product(s) is/are in conformity with the following standards and/or other normative documents:

ETSI EN 300 220-1 V3.1.1: 2017
ETSI EN 300 220-2 V3.1.1: 2017
ETSI EN 330 330-1 V2.1.1 (2017-2)
ETSI EN 60950-1:2006 + A2:2013
ETSI EN 301 489-1
ETSI EN 301 489-17 (Article 3.1(b) of R&TTE Directive)

Limitations of validity (if any): N/A
Supplementary information:

Technical file held by:

Secure Care Products LLC
39 Chenell Dr.
Concord, NH 03301

Place and date of issue (of this DoC):

Secure Care Products LLC. 11/06/2017

Signed by or for the manufacturer:

Printed name: Michael Singer
Title: R&D Engineer
Signature: ON FILE

A00020756 Rev. G

EU DECLARATION OF CONFORMITY

THE EU DIRECTIVE COVERED BY THIS DECLARATION:

Radio Equipment Directive 2014/53/EU

PRODUCT COVERED BY THIS DECLARATION:

Name: **13.56 / 433.92 MHz ID STAT Transmitter**

Model: **135**

Part Number(s): **A20240905, A20274935, A21350915, A22274915,
A22274935**

THE BASIS ON WHICH CONFORMITY IS BEING DECLARED:

This device, Model 135 ID STAT Transmitter, complies with the essential requirements of the Radio Equipment Directive 2014/53/EU on the basis of conformance to the harmonized standards listed below:

EMC – EN 300 330 V2.1.1:2017
Radio – EN 300 220-1 V3.1.1:2017, EN 300 220-2 V3.1.1:2017
Safety – IEC/EN 60950-1

NOTIFIED BODY:

Retlif Testing Laboratories
101 New Boston Rd
Goffstown, NH 03045 USA

MANUFACTURER AND RESPONSIBLE PARTY:

Secure Care Products, Inc
39 Chenell Drive
Concord, New Hampshire 03301 USA

Printed Name: Mike Singer

Signed: Signature on File

Title: Electrical Engineer

Date: 3/23/2020

A00020536 Rev F

SECURE CARE PRODUCTS, LLC	
Name: Compliance Information Statement	F04-0018
Date: 04-06-15	Issue Version: 1
ECO #: 11563	Page 1 of 1
[Type here]	

Secure Care Products, LLC.	
Concord, New Hampshire USA	
<h1>Compliance Information Statement</h1>	
Conforming Product Name(s):	6.78MHz Transmitters
Conforming Model or Part Number(s):	A20180920, A20180921, A20180925, A22180900, A22180910, A21180900, A21180905, A21180910, A21180911, A22180901
Manufacturer and Responsible Party:	Secure Care Products, LLC. 39 Chenell Drive Concord, New Hampshire 03301 USA
Reference Document(s):	A00020674

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Name of Responsible Person	Company Title or Position
Mike Singer	R&D Engineer
Authorized Signature	Date
Signature on File	3/6/20

Secure Care Document PN A00020666 Rev M

Secure Care Products, LLC.	
Concord, New Hampshire USA	
<h1>Compliance Information Statement</h1>	
Conforming Product Name(s):	SC-40 Transmitters
Conforming Model or Part Number(s):	A20270909, A20270947, A20270948, A20270949, A21350935, A21350936, A21350940, A21350941, A21350942, A21350945, A21350946, A21350947, A21350948, A21350957, A21350958, A22270940, A22270941, A22270942, A22270944, A22350949,
Manufacturer and Responsible Party:	Secure Care Products, LLC 39 Chenell Drive Concord, New Hampshire 03301 USA Telephone Number: 603.223.0745
Reference Document(s):	A00020772, FCC ID: KNK-B

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Name of Responsible Person	Company Title or Position
Mike Singer	R&D Engineer
Authorized Signature	Date
Signature on File	05/05/19

SECURE CARE PRODUCTS, LLC	
Name: Compliance Information Statement	F04-0018
Date: 04-06-15	Issue Version: 1
ECO #: 11563	Page 1 of 1
[Type here]	

Secure Care Products, LLC.	
Concord, New Hampshire USA	
<h1 style="text-align: center;">Compliance Information Statement</h1>	
Conforming Product Name(s):	13.56/433.92MHz Transmitter(s)
Conforming Model or Part Number(s):	A20250903, A20240900, A20240901, A20240903, A20240904, A20240905, A22260913, A22270913, A20270937, A20274935, A21350910, A21350915, A22260912, A22270912, A22270937, A22274915, A22274916, A22274935
Manufacturer and Responsible Party:	Secure Care Products, LLC. 39 Chenell Drive Concord, New Hampshire 03301 USA
Reference Document(s):	A00020771, FCC ID: KNKA

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Name of Responsible Person	Company Title or Position
Mike Singer	R&D Engineer
Authorized Signature	Date
Signature on File	5/6/19

Secure Care Document PN A00020535 Rev H

GENERAL PRODUCT WARRANTY

SECURE CARE PRODUCTS, LLC

GENERAL PRODUCT WARRANTY STATEMENT

BY PERMITTING INSTALLATION OR BY MAKING USE OF ANY PRODUCT OR SERVICE DESIGNED OR MANUFACTURED BY SECURE CARE PRODUCTS, LLC (“SECURE CARE”) (INCLUDING SUPPORT SERVICES, MAINTAINED SOFTWARE AND MAJOR RELEASES, WHETHER OR NOT IT IS COVERED BY ANY SOFTWARE MAINTENANCE OR LICENSE AGREEMENT) (“THIS PRODUCT”), YOU ACKNOWLEDGE THAT YOU HAVE READ ALL THE TERMS AND CONDITIONS OF THIS GENERAL PRODUCT WARRANTY STATEMENT, THAT YOU UNDERSTAND THEM, AND THAT YOU AGREE TO BE BOUND BY THEM. YOU UNDERSTAND THAT, IF YOU PURCHASED THIS PRODUCT FROM ANY AUTHORIZED DISTRIBUTOR OF SECURE CARE, THAT DISTRIBUTOR IS NOT SECURE CARE’S AGENT AND IS NOT AUTHORIZED TO MAKE ANY REPRESENTATIONS OR WARRANTIES OR TO AGREE TO ANY TERMS OR CONDITIONS WHICH ARE DIFFERENT FROM ANYTHING EXPRESSLY SET FORTH IN THIS GENERAL PRODUCT WARRANTY STATEMENT.

If you do not agree to the terms and conditions of this General Product Warranty Statement, do not permit the installation or make use of this Product and promptly return this Product to the place where you obtained it for a full refund. If you have any difficulty obtaining a refund, please contact Secure Care at the telephone number provided in Section 2.B below.

1. Notices

A. ALL LOCKS USED WITH THE SECURE CARE SYSTEM ARE DESIGNED, MANUFACTURED, LABELED AND DELIVERED SOLELY BY AN INDEPENDENT VENDOR OVER WHOM SECURE CARE HAS NO CONTROL AND FOR WHOSE ACTIONS OR FAILURES TO ACT SECURE CARE DISCLAIMS ALL RESPONSIBILITY. REGARDLESS OF WHETHER THE LOCKS CARRY SECURE CARE’S LOGO OR NAME OR ANY OTHER TRADEMARK, SERVICE MARK OR TRADE NAME USED OR CLAIMED BY SECURE CARE, SECURE CARE DISCLAIMS ALL WARRANTIES, EXPRESS OR IMPLIED, WITH RESPECT TO THE LOCKS AND/OR THEIR USE WITH OR OPERATION IN THE SECURE CARE SYSTEM, INCLUDING, WITHOUT LIMITATION, ALL IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, TITLE AND/OR NON-INFRINGEMENT. SECURE CARE ALSO DISCLAIMS ALL OBLIGATIONS WITH RESPECT TO THE LOCKS AND/OR THEIR USE WITH OR OPERATION IN THE SECURE CARE SYSTEM THAT MIGHT OTHERWISE ARISE OR BE IMPLIED FROM THE FACT THAT SUCH LOCKS CARRY SECURE CARE’S LOGO OR NAME OR ANY OTHER TRADEMARK, SERVICE MARK OR TRADE NAME USED OR CLAIMED BY SECURE CARE OR FROM THE DELIVERY OR INSTALLATION OF THE LOCKS WITH SECURE CARE SOFTWARE, PARTS AND/OR PRODUCTS OR FROM A COURSE OF DEALING OR USAGE IN TRADE. ALL RESPONSIBILITY FOR DESIGNING, MANUFACTURING, LABELING AND WARNING OF HIDDEN DEFECTS OR DANGERS IN THE LOCKS AND/OR THEIR USE WITH AND OPERATION IN THE SECURE CARE SYSTEM RESTS EXCLUSIVELY WITH THE INDEPENDENT VENDOR, AND ANY CLAIMS, COSTS, DAMAGES OR LIABILITIES ARISING FROM THE LOCKS AND/OR THEIR USE WITH OR OPERATION IN THE SECURE CARE SYSTEM SHALL BE MADE SOLELY AGAINST THE INDEPENDENT VENDOR.

B. IF YOU PURCHASE COMPUTER HARDWARE THROUGH SECURE CARE AND REQUEST THAT SECURE CARE SOFTWARE BE INSTALLED AND TESTED ON THAT HARDWARE AT THE FACTORY, SECURE CARE WARRANTS ONLY THAT THE HARDWARE AND THE SOFTWARE PACKAGES WERE INSTALLED, SET-UP AND TESTED PRIOR TO SHIPMENT IN ACCORDANCE WITH ALL SECURE CARE PRODUCT MANUALS AND THAT, AT THE TIME THE HARDWARE AND THE SOFTWARE PACKAGES WERE FINALLY INSPECTED AT THE FACTORY, THEY WERE PERFORMING (SUBJECT TO SECURE CARE’S SPECIFIED TOLERANCES) IN ACCORDANCE WITH SECURE CARE’S SPECIFICATIONS. SECURE CARE WILL NOT BE RESPONSIBLE FOR ANY DEFECTS IN OR PROBLEMS CAUSED BY THE HARDWARE, ALL CLAIMS FOR WHICH MUST BE MADE

TO THE HARDWARE MANUFACTURER AND/OR VENDOR. SECURE CARE DISCLAIMS ALL WARRANTIES, EXPRESS OR IMPLIED, WITH RESPECT TO THE HARDWARE AND/OR ITS USE WITH OR OPERATION IN THE SECURE CARE SYSTEM, INCLUDING, WITHOUT LIMITATION, ALL IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, TITLE AND/OR NON-INFRINGEMENT. SECURE CARE ALSO DISCLAIMS ALL OBLIGATIONS WITH RESPECT TO THE HARDWARE AND/OR ITS USE WITH OR OPERATION IN THE SECURE CARE SYSTEM THAT MIGHT OTHERWISE ARISE OR BE IMPLIED FROM THE FACT THAT SUCH HARDWARE CARRIES SECURE CARE'S LOGO OR NAME OR ANY OTHER TRADEMARK, SERVICE MARK OR TRADE NAME USED OR CLAIMED BY SECURE CARE OR FROM THE DELIVERY OR INSTALLATION OF THE HARDWARE WITH SECURE CARE SOFTWARE, PARTS AND/OR PRODUCTS OR FROM A COURSE OF DEALING OR USAGE IN TRADE. ALL RESPONSIBILITY FOR DESIGNING, MANUFACTURING, LABELING AND WARNING OF HIDDEN DEFECTS OR DANGERS IN THE HARDWARE AND/OR ITS USE WITH AND OPERATION IN THE SECURE CARE SYSTEM RESTS EXCLUSIVELY WITH THE HARDWARE MANUFACTURER AND/OR VENDOR, AND ANY CLAIMS, COSTS, DAMAGES OR LIABILITIES ARISING FROM THE HARDWARE AND/OR ITS USE WITH OR OPERATION IN THE SECURE CARE SYSTEM SHALL BE MADE SOLELY AGAINST THE HARDWARE MANUFACTURER AND/OR VENDOR.

C. Secure Care's software, parts and products are designed for operation in a wireless system. However, the range, performance, and predictability of any wireless system, including Secure Care's, is dependent on several factors, including, but not limited to, the following: building structure; environmental extremes (e.g., temperature, earth tremors, air pollution, etc.); the proximity of other wireless devices; the presence of variable speed products; sources of Radio Frequency Interference (RFI); physical orientation and positioning of the equipment; and sources of Electro Static Discharge ("ESD"). Secure Care cannot be responsible for the effect of these types of factors on operation of its software, parts and products.

D. This Product must be installed, set-up, tested, supported, operated, maintained, repaired and used only in accordance with all manuals and instructions (including the user, installation, technical and other manuals) issued by Secure Care (the "Product Manuals"). It is your responsibility to assure that any person who might be installing, setting-up, testing, supporting, maintaining or repairing this Product knows the contents of and has access to the Product Manuals and has successfully completed Secure Care technical training. It is also your responsibility to assure that any person who might be operating or using this Product knows the contents of and has access to the Product Manuals and has successfully completed Secure Care in-service training. If you do not have the Product Manuals or if you have any questions regarding this Product and/or its installation, set-up, testing, support, operation, maintenance, repair or use, please call Secure Care at the telephone number provided in section 2.B below. Secure Care cannot be responsible for performance problems caused by a failure to follow published and appropriate procedures for installation, set-up, testing, support, operation, maintenance, repair and use.

All adjustable features on new and repaired Secure Care software, parts and products are shipped with "factory default" settings. These "factory default" settings may not comply with building and life safety codes or other applicable laws and regulations in the location where they are installed or operated. Secure Care strongly recommends, therefore, that the settings on all Secure Care software, parts and products be checked and, if necessary, reset to comply with local building and life safety codes and other applicable laws and regulations at the time of any installation, set-up, testing, support, maintenance or repair.

E. Secure Care's system is driven by software. However, the performance and reliability of any software-driven system depends on adequately maintaining the recommended minimum configuration of computing platform, operating systems and applications programs and on regularly performing industry-standard and application-specific backup processes. If recommended minimum configurations of computing platform, operating systems, and applications programs are not adequately maintained, or if appropriate backups are not regularly performed, the software may not drive the system as intended. Secure Care cannot be responsible for operational problems caused by a failure to perform these maintenance and backup procedures.

F. Secure Care does not authorize, and strongly recommends against, any installation or field replacement of software, parts or products by untrained contractors or facility staff. Such work can be hazardous, can render the system ineffective and will void any Secure Care warranty or liability that might otherwise relate to the system.

Before any software, parts or products, which have been designed and manufactured by Secure Care can be safely installed, set-up, tested, supported, maintained or repaired, technical training in accordance with standards established by Secure Care is required. Regardless of how Secure Care's software, parts or products are obtained, they should not be installed, set-up, tested, supported, maintained or repaired by any person who has not satisfactorily completed that technical training (a "qualified service technician".) When Secure Care's software, parts or products are sold separately from installation services, it is assumed that only a qualified service technician will conduct any installation, set-up, testing, support, maintenance or repair involving that software, part or products.

G. It is important that any installation, set-up, testing, support, operation, maintenance, repair or use involving the system comply with all local and national electrical and life safety codes. If you have any questions about compliance with those codes, please contact your local authorities.

H. Secure Care receives and responds to telephone and dial-in inquiries (the "Help Line") about its software, parts and products for the purpose of discussing users' experiences with Secure Care's system, helping users better understand how their systems work, and providing ideas about what may be causing difficulties. However, Secure Care cannot accurately diagnose the cause of any problems or give complete instructions on how to fix problems over the telephone or Internet. The only way to assure that software, parts or products are installed, set-up, tested, supported, maintained or repaired correctly or that a Secure Care system is functioning properly is to have it examined on site by a qualified service technician. In addition, Secure Care software, parts and products cannot be operated or used correctly by anyone who has not successfully completed Secure Care in-service training. Secure Care's Help Line is not a substitute for on-site diagnosis and servicing by a qualified service technician or for successful completion of Secure Care in-service training. Secure Care strongly recommends that any installation, set-up, testing, support, replacement, maintenance or repair of a system that is performed by a person who has not satisfactorily completed technical training in accordance with standards established by Secure Care be immediately checked on-site by a person who has completed that technical training.

WARNING: EVEN SLIGHT MODIFICATIONS TO THE SYSTEM OR CHANGES IN THE OPERATING ENVIRONMENT MAY CAUSE SECURE CARE'S SYSTEM TO MALFUNCTION. THE ONLY WAY TO ASSURE THAT SECURE CARE'S SYSTEM HAS BEEN INSTALLED, SET-UP, TESTED, SUPPORTED, MAINTAINED, AND REPAIRED CORRECTLY IS TO HAVE A QUALIFIED SERVICE TECHNICIAN DO THE WORK.

I. Secure Care's software, parts and products have been designed to augment a facility's reasonable procedures for protecting residents, patients, and infants. However, no system or combination of procedures and equipment can eliminate all risk or assure complete security. Secure Care's system is not intended as a substitute for the careful identification and monitoring of residents, patients, and infants by a facility's professional staff.

2. Limited Warranty

A. Subject to the limitations set forth in this general product warranty statement (as amended from time to time by Secure Care in its absolute discretion), and unless a different period is specified in writing by Secure Care for a particular product or service, Secure Care warrants that this product (subject to Secure Care's specified tolerances and excluding any expendable items), if sold by Secure Care to an authorized Secure Care distributor, shall conform to the specifications which accompany this product for a period of one (1) year from the date of delivery of this product by Secure Care to a common carrier, f.o.b. Secure Care's manufacturing facility in Concord, New Hampshire or, in the case of services, from the date of first provision of such services. This warranty does not extend to and is not for the benefit of any person other than an authorized Secure Care distributor who purchases this product from Secure Care, any sub-distributor thereof and/or the customer to whom this product is first provided for use, by Secure Care, an authorized Secure Care distributor or any sub-distributor thereof. In the event that this product does not comply with this warranty, Secure Care will, at its option, either repair or replace this product or refund the purchase price, provided that this product is returned as provided in section 2.b below. Replacement of this product under warranty will not extend the original warranty period.

Secure Care will also, at its option, either repair or replace this Product after the warranty has expired, for an additional charge, provided that this Product is returned as provided in Section 2.B below. If Secure Care repairs or replaces this Product after the warranty has expired, the terms of the warranty set forth in this Section 2.A for a new Product will apply to the repaired or replaced Product, with the exception that the term will run for ninety (90) days from the date of repair or replacement.

Repair may include the replacement of parts and products with functionally equivalent, reconditioned parts or products. Any part or product replaced by Secure Care will become the property of Secure Care upon replacement.

B. Warranty service is available by contacting Secure Care at 800-451-7917 and obtaining a Return Authorization Number. No Product may be returned to Secure Care without first obtaining a Return Authorization Number. When this Product is returned to Secure Care, please include the Return Authorization Number and a detailed written description of the problem. Issuance of a Return Authorization Number by Secure Care will not constitute an admission that there is a problem with the Product being returned, that any problem is covered by warranty or that Secure Care has any responsibility to repair, replace, make refunds for or pay claims, costs, damages or liabilities connected with the Product being returned.

NOTE: Transmitters are not repaired, nor is the warranty extended, beyond the expiration date.

If this Product is returned to Secure Care for any reason, you will retain title (unless and until a part or product is replaced, in which case you will obtain title to the replacement part or product at the time of replacement), the risk of loss, and the obligation to pay all costs of shipping, storage and other charges and obligations relating to this Product.

C. Except as stated in this section 2, Secure Care disclaims all warranties, express or implied, with respect to the whole or any part of this product, including, without limitation, all implied warranties of merchantability, fitness for a particular purpose, title and/or non-infringement. Secure Care also disclaims all obligations that might otherwise arise or be implied from a course of dealing or usage in trade.

3. Limitations of Liability

A. Regardless of the form of any claim or action, Secure Care's total liability to all persons, whether singly or together, for all occurrences combined, for claims, costs, damages or liabilities based on any cause whatsoever and arising from or in connection with this product, or the manufacture, distribution, promotion, sale, installation, set-up, testing, support, maintenance, operation, servicing, use or performance of this product, or from or in connection with any delay or failure in providing this product, shall not exceed the aggregate price (without interest) paid to Secure Care for this product.

B. In no event shall Secure Care be liable to anyone for any loss of data, loss of profits or loss of use of this product or any equipment, or for any special, incidental, consequential, exemplary, punitive, multiple, or other damages, arising from or in connection with the manufacture, distribution, promotion, sale, installation, set-up, testing, support, maintenance, operation, servicing, use or performance of this product or from or in connection with any delay or failure in providing or delivering this product.

C. In no event shall Secure Care be liable to anyone for any claims, costs, damages or liabilities caused by: (i) any distributor's failure to perform its obligations and responsibilities under a distributor agreement with Secure Care; (ii) improper or defective promotion, distribution, sale, installation, set-up, testing, support, maintenance or repair of this product, including work performed, without Secure Care's prior written consent in its absolute discretion, by a person who has not satisfactorily completed Secure Care technical training, or in a manner not consistent with Secure Care technical training; (iii) improper or defective operation or use of this product by a person who has not successfully completed Secure Care in-service training, or in a manner not consistent with Secure Care in-service training; (iv) supply of this product by a distributor for use in, or the use of this product in, any system or configuration not designed to Secure Care standards or in which a distributor or any third party has substituted materials and/or goods not specified by Secure Care; or (v) deterioration of this product during storage. This liability extends to voiding the product/system warranty and agency certifications. You will assume all liability for any deviations from the above.

D. You agree to indemnify and hold Secure Care harmless from all claims, costs, damages and liabilities asserted by anyone for any damages that are excluded and waived, or are intended to be excluded and waived, by this section 3, or which are imposed by law on behalf of anyone but which are not expressly stated in this general product warranty statement.

E. The exclusions, waivers and limitations on claims, costs, damages and liabilities and any rights of indemnification set forth in this section 3 shall be enforceable to the maximum extent allowed by law and shall not be expanded or negated in any respect by Secure Care's operation of a "help line" to receive and respond to telephone or dial-in inquiries about this product, by any communications through that "help line" or by any actions taken by anyone following communications with Secure Care over such "help line."

4. Governing Law and Arbitration

A. This General Product Warranty Statement, and all questions arising out of or relating to it, shall be governed by and construed in accordance with the laws of the State of New Hampshire, without giving effect to the conflict of laws provisions thereof, and excluding the United Nations Convention on contracts for the international sale of goods, the 1974 convention on the limitation period on the international sale of goods (the "1974 convention"), and the protocol amending the 1974 convention, done at Vienna April 11, 1980.

B. Any dispute, controversy or claim arising out of or relating to this general product warranty statement shall be resolved by arbitration. Regardless of the amount in dispute, the arbitration shall be conducted by a single arbitrator selected by the parties or, if they cannot agree, by a single arbitrator selected in accordance with the commercial arbitration rules of the American Arbitration Association without regard to the amount in dispute. The arbitration shall be conducted in English, in accordance with the commercial arbitration rules of the American Arbitration Association, in Concord, New Hampshire. The decision of the arbitrator shall be binding and enforceable by any state or federal court in New Hampshire, and you hereby consent to the personal jurisdiction of any state or federal court in New Hampshire for that purpose. The expense of the arbitration (excluding each side's own attorneys' fees, costs, and related expenses) shall initially be paid in equal shares

by each side, but the total of such expenses plus any award of attorneys' fees, cost and expenses shall finally be paid by the parties as the arbitrator determines. Nothing in this section 4.b shall preclude Secure Care from seeking provisional or equitable relief from any appropriate court to protect its rights prior to, pending or in the absence of such arbitration proceedings.

5. Severability

The invalidity or unenforceability of any provision of this General Product Warranty Statement shall not affect the validity or enforceability of any other provision hereof.

6. Waiver

No term or condition of this General Product Warranty Statement may be waived except in writing signed by Secure Care. A waiver on one or more occasions of any term or condition of this General Product Warranty Statement shall not constitute or be deemed to be a waiver of such term or condition on any other occasion. No delay or failure of Secure Care to exercise any right or remedy under this General Product Warranty Statement will operate as a waiver thereof; no failure to enforce or insist upon compliance with any provision of this General Product Warranty Statement on any one occasion shall be deemed to be a waiver of Secure Care's right to do so on another occasion; and no course of dealing will constitute a waiver, alteration, limitation or expansion of any of the parties' rights and obligations under this General Product Warranty Statement.

Revised 02/14/17



DAILY TRANSMITTER TESTING LOG

DAILY TRANSMITTER TESTING LOG DATE _____

Resident Last Name	Resident First Name	Transmitter ID #	Transmitter Test Pass/Fail	Transmitter Expiration Date	Comments

NOTE:
Test all transmitters in use daily.

NOTES: