
Model: R-A012 / R-A020

Solar-Powered Scale

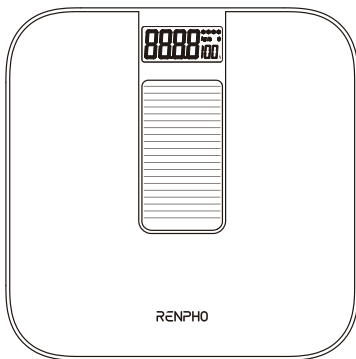RENPHO is dedicated to protecting the environment by reducing paper usage. Scan the QR Code or enter the URL below to view or download the full version of the electronic User Manual:

<https://renpho.com/pages/solar-powered-smart-scale>

CONTENTS

• ENGLISH	01-11
-----------	-------

• DEUTSCH	12-22
-----------	-------

• FRANÇAIS	23-33
------------	-------

• ITALIANO	34-44
------------	-------

• ESPAÑOL	45-55
-----------	-------

Important Safety Instructions

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

1. The scale is for household use only and not for medical or commercial use.
2. Individuals who are pregnant, have osteoporosis, a pacemaker, artificial joints, portable ECG or other metal implants should use the scale for weighing purposes only and make sure to wear shoes to prevent any electric current when stepping onto the scale.
3. All body composition data is for reference and should not be used as a substitute for advice from a licensed healthcare professional.
4. This device is a solar-powered scale. Charge it with direct sunlight or other light sources. It does not support USB charging.
5. Turn on the power switch before the first use. Switch the scale off if it is not used for a long time.
6. To activate the scale, step on it without shoes.
7. Do not use the scale in the environment with sudden temperature changes.
8. Use the scale on a dry, flat, solid and even surface.
9. Do not stand on the sides of the scale or jump on the scale to avoid any injury.
10. Do not overload the scale (Max. 180kg/400lb).
11. Place/use the scale in a well-lit area (At least 2 "  " appear on the display.) to ensure it receives sufficient light for prolonged use.
12. Place the scale in the area where 4 "  " appear on the display to obtain faster charging when it is in low battery.
13. Avoid splashing water on the solar panel on the scale.
14. Do not immerse the scale in water or use chemical cleaning agents. Clean the scale with a slightly damp cloth.
15. Do not attempt to disassemble the scale for any reason, otherwise, the warranty may be void.
16. Do not heat or dismantle the battery.
17. A lithium battery with a damaged container should not be exposed to water.

SAVE THESE INSTRUCTIONS

About the Product

1. Technical Specifications

Product Name: Smart Scale

Model: R-A012 / R-A020

Working Voltage: 3.7V

Working Current: $\leq 13.5\text{mA}$

Built-in Battery Capacity: 500mAh

Unit: kg/lb (Switch via the app)

Precision: $< 50\text{kg}$: $\pm 0.3\text{kg}$;

$50\text{-}100\text{kg}$: $\pm 0.4\text{kg}$;

$> 100\text{kg}$: $\pm 0.5\text{kg}$;

$< 110\text{lb}$: $\pm 0.66\text{lb}$

$110\text{lb}\text{-}220\text{lb}$: $\pm 0.88\text{lb}$

$> 220\text{lb}$: $\pm 1.1\text{lb}$


Division: 0.05kg / 0.2lb

Capacity: $0.2\text{-}180\text{kg}$ / $0.4\text{-}400\text{lb}$

Product Size: $280 \times 280 \times 23\text{mm}$ / $11 \times 11 \times 0.9\text{inch}$

Product Weight: 1.3kg / 2.9lb

2. Display

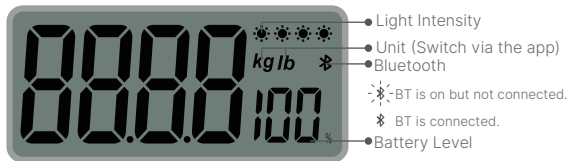- a. To ensure the scale receives sufficient light for prolonged use, place/use it in a well-lit area. (At least 2 "  " appear on the display.)
- b. To obtain faster power storage, place the scale in the area with direct sunlight (4 "  " appear on the display.) when it is in low battery.


3. Icons

 Overloaded (Max. 180kg/400lb)	 Low battery (Place the scale in the sun timely.)	 Measuring body fat	 Calibrating
---	--	---	--

Using the Product

1. Before Use

a. Place the scale on a dry, flat, solid and even surface. Avoid any wet, slippery, carpeting or uneven flooring.

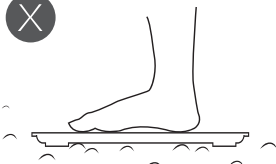b. Turn ON the power switch on the back of the scale.


c. Calibrate: Step on the scale platform gently and then step off. The scale's display will show "CAL" and then "0.00", indicating that calibrating is complete.

d. Start measuring.

2. Measuring

- Step on the scale with dry and bare feet and wait for 3-5 seconds.
- When the number shown on the scale's display flashes twice, measuring is complete.

* The measurements can be viewed via the Renpho Health App if the connection is successful.

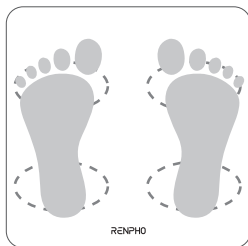Tips:

- To ensure accuracy, please try to measure at the same time and under the same conditions.
- The scale will be on the standby status if left idle for 10 seconds. Step on the scale gently to activate it.

1. Download the Renpho Health App

Scan the QR code below or search for “Renpho Health” from the App Store or Google Play to download the app.

System Requirements: iOS 12.0 or later, Android 7.0 or later.

Renpho Health

Note: Due to continuous updates and improvements, the "Renpho Health" App may appear slightly different.

2. Log in/Register

- Log in directly with an existing Renpho App account.
- If you don't have any Renpho App account, tap **"Sign up"**. Register a new account using your email address and complete your profile information, which is required for calculating the measurements only.



Note: You can also log in using Guest Mode if you do not wish to register an account. However, your measurements will not be saved if you log in on another device.

3. Add Device

- * Enable Bluetooth on your smartphone before adding a device.
(Location is required by Android System.)





Method 1:

- Step on the scale and a window will pop up.
- Tap "Pair" to pair your smartphone with the scale.

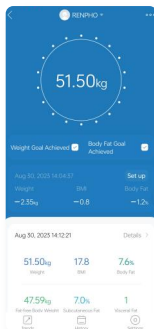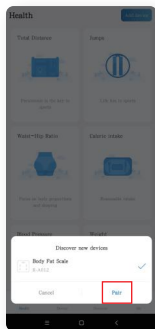
Method 2:

- Tap “Add device” on the top right corner of the “Health Page” to connect manually.
- Tap “Body Fat Scale”.
- Step on the scale following the in-app instructions until a window pops up on the screen.
- Tap “+” to connect. (The Bluetooth icon will turn solid on the scale’s display once the pairing is successful.)

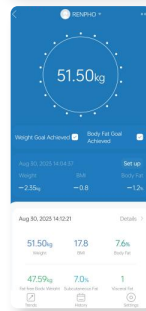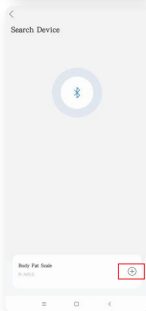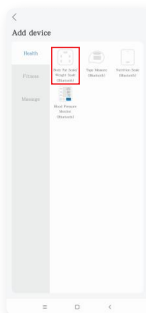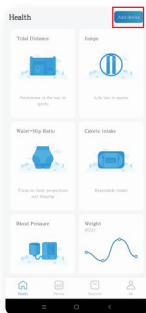4. Switch Units

To switch units, tap **“Settings”** and then tap **“Units”** to choose the unit according to your needs.

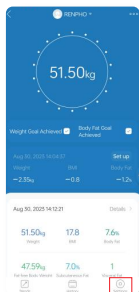1. How can I turn on the scale?

Turn ON the power switch on the back of the scale and step onto activate it without shoes.

2. Do I need to replace the battery?

No, you don't. The scale is powered by solar power and other light sources.

3. How can I charge the scale?

Please place the scale where light is noticeable. The stronger the light intensity, the faster the charge.

4. How can I know the current battery level?

Check the battery percentage on the scale's display. 100% means that the scale is at full charge.

5. How can I know the current light intensity?

Check the number of "  " on the scale's display. The more "  " there are, the stronger the light intensity is.

6. Why does the scale stop working?

The scale may have run out of power. Please place the scale in a well-lit place for charging.

Warranty Policy

Your RENPHO product purchase is covered by a one year limited manufacturer warranty from the date of delivery. For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not void the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **Tel:** +44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)

 **Email:** support-eu@renpho.com (UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

Notice:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

RF Exposure Compliance

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

ISED Regulatory compliance

This device contains licence-exempt transmitters that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s).

Operation is subject to the following two conditions:

(1) This device may not cause interference.

(2) This device must accept any interference, including interference that may cause undesired operation of the device.

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This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment.

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EU Importer: Renpho GmbH
Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.
UK REP Richers & Vision Ltd.
EU REP C&E Connection E-Commerce (DE) GmbH