arpha

EMS Titanium Ion Beauty Instrument

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Good skin starts with skincare.

Thank you for using arpha EMS Titanium Ion Beauty Instrument.

Before using this product, you must please read and follow all warnings, prohibitions, and safety reminders to ensure safety, and any failure to follow the instructions in this manual may cause the product to burn or damage your skin.

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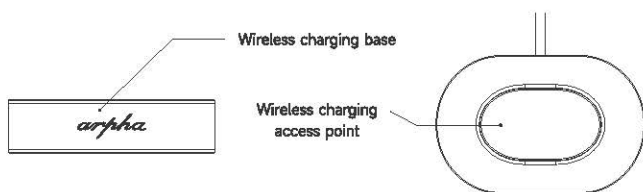
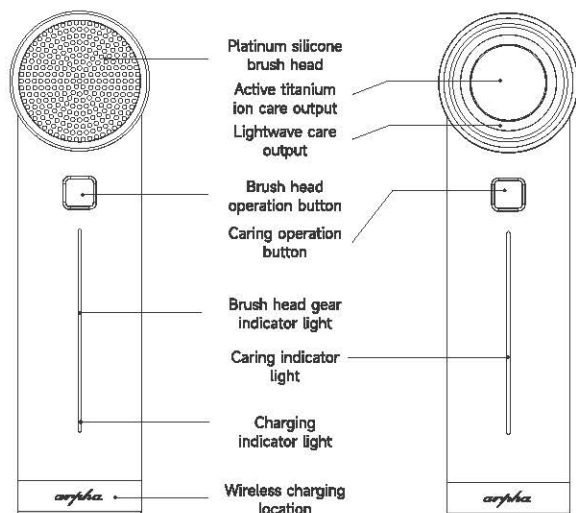
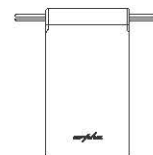
Portable flannel bag

Product Specification

Product size	144.5 x 40.1 x 50.5 (mm)
Battery capacity	≈ 1450mAh
Rated voltage	3.7V
Rated power	9.0W
Adapter specifications	USB-A / 5.0V-1.5A
Charging time	≈ 8 hours
Standby time	250 hours
Main material	food-grade ABS, platinum silicone, aviation aluminum, ion-activated titanium

Warnings

- The function of EMS is to inject micro-electric current into the muscles, so as to stimulate the muscles for firming effect. The actual use effect varies according to individual skin type, and there will be a warm tingling sensation for 5-30 minutes.
- It is recommended to use face cleaning / EMS mode every day, that do not use for more than 6 minutes continuously.
- After the cooling function of the skin beauty instrument (Brightening Mode), the treatment head will automatically dissipate heat, and burning hot is a normal phenomenon.
- Do not use this product on delicate and sensitive parts of the body skin such as eyelids, lips and ears.
- This product is intended for care of specified functions only. Do not use for any purpose other than as directed.
- This product can not be used as a toy.
- When charging, please use the original charger of this product for charging.
- This product is not suitable for infants and young children, please place the product out of the reach of infants young children.
- Keep hands completely dry while using this product.
- For gel, cream, lotion and other care products used with this product, please refer to their instructions.
- For people who are under 18 years old and those who have undergone physical surgery, please consult a medical professional before use.
- This product cannot be cleaned and sterilized with strong acidic and alkaline detergents or disinfectants.

Prohibitions

- Do not use the instrument on the thyroid or lymphatic area.
- Do not use the instrument on the throat, eyes, ears, upper eyelids, or temples.
- Do not place this product near fire or heat sources.
- Do not drop or strongly bump this product.
- Do not clean this product with water higher than 60°C.
- Do not clean this product with dishwashers or other mechanical devices.
- Do not immerse this product or the charger base in water or other liquids for long periods of time.
- Do not use this product when running or walking to avoid accidents.
- Do not store or use this product in places with strong magnetic fields, electric waves, or in dusty places.

Cautions

- Do not use the instrument on the thyroid or lymphatic area.
- Please be careful when using this product on thin skin areas such as the corners of the eyes, forehead, neck, shins, etc.
- A lower setting is recommended for more sensitive skin.
- Do not use skin care products containing alcohol, acetone, perfume, disinfectant and other flammable and volatile properties within one hour before using this product.
- If treatment is unpleasant or you feel uncomfortable during use, please stop using immediately and consult a medical professional.
- When using this product, it should be in close contact with your skin and you should choose the appropriate gear according to your skin's tolerance to temperature and vibration.
- The treatment area should be moistened with gel or lotion. Applying the product directly on dry skin may cause burns or other injuries.

- Storage of this product in extreme temperature conditions may damage the body material and reduce its service life.
- Please clean this product in time after use. When cleaning the instrument, please make sure to turn off the power of this product first.
- Please check this product for obvious damage before each use. If there is damage or the unit or charging base has been severely dropped, please do not use it.
- Before connecting the power supply, check that the local power supply voltage matches the voltage shown on the product.
- If you do not intend to use the product for a while, please unplug the wireless charging adapter and clean the product thoroughly, store it in a cool and dry place, and remove the charging power to avoid fire or other accidents.
- The instrument has the 30-second protection function. If the treatment side operates in one mode for more than 30 seconds, then switching to the next mode, the treatment side will automatically enter the protection mode and shut down. You can turn on the instrument again to continue using it.

The following people are prohibited to use the product

- Who suffer from heart disease, coronary heart disease, thyroid, nervous system, hypertension, diabetes, glaucoma, etc;
- Who currently have or have had cancer, especially skin cancer or precancerous lesions;
- Who equipped with electronic devices such as hearing aids, pacemakers or cardiac defibrillators;
- Who are pregnant, lactating or menstruating women;
- Who have tumors and abnormal wound healing;
- Who are undergoing hormone therapy;
- Who are sensitive to electric current;
- Who have impaired immune systems due to immunosuppressive diseases or who are taking immunosuppressive drugs;
- Whose care area has any heat irritant disease, such as recurrent herpes simplex virus;
- Whose care area has any skin diseases such as sores, phlegm, rashes, etc.;
- Whose care area has localized lesions such as severe lesions, scar growths and inflammation;
- Whose care area has severe sunburn;
- Whose care area has had Laser skin rejuvenation, chemical skin rejuvenation or wound healing ;
- Whose care area has permanent fillers, such as buried metal, plastic, silicone;
- Whose care area has temporary fillers, such as hyaluronic acid and Botox.

Product Features

- The platinum food-grade silicone brush head is soft and gentle, making it more suitable for sensitive and acne-prone skin, and the addition of sulfur in the brush head prevents dust mites from breeding.
- Microcurrent (EMS / Electric Muscle Stimulation) converts the nutritional essence of maintenance products into a water-soluble ionic state. Active titanium ion penetrates deep into epidermal cells and EMS micro-current penetrates sebum, eliminating facial flab from the root, penetrating directly into the skin's stratum corneum and stimulating collagen growth, making the face firmer and more three-dimensional in contour.
- The sonic vibrating brush head has hundreds of soft bristles with aligned contacts to ensure that all parts of the face are properly cleaned. The soft silicone brush head can cleanse the skin while gently dissolving dirt and excess sebum, cosmetic residue and removing the skin's aged and dulled cuticle.
- Thermostatic hot compresses accelerate metabolism, stimulate the internal circulation of the skin, repair damaged tissues, enlarge the pores during the hot compress, and accelerate the absorption of nutrients into the skin.
- The low-temperature ice compress shrinks pores. After the skin is exposed to moisturizing water, it quickly shrinks pores and makes the skin more hydrated, firm and bright.
- Red LED will stimulate blood flow, accelerate blood circulation, enhance cellular vitality, improve metabolism, promote collagen fiber synthesis, and increase skin elasticity and luster.
- Blue LED can inhibit oil secretion, diminish inflammation and sterilize, destroy acne bacillus, and treat acne and scars. Green LED has a balancing and stabilizing effect and can improve skin roughness, wrinkles, blackheads, acne and acne caused by fatigue and mental tension.

Charging the Beauty Instrument

⚡ Charging the battery

- Insert the USB-A port of the dedicated wireless charging base into the power charger (5.0V/1.5A is recommended), and then insert the charger into the power outlet. At this time, the charging base light will turn solid blue, indicating that the charging base is correctly connected to the socket.
- After inserting the beauty instrument into the wireless charging base, all the lights of the beauty instrument flash 5 times, representing the correct access to power; at this time, the charging indicator light at the bottom of the brush head of the beauty instrument will flash blue, indicating that the beauty machine is charging. Charging is completed in about 8 to 12 hours (charging time may vary depending on the amount of charging). When charging has completely, the charging indicator light will turn solid blue.
- In order to ensure safety and save energy, it is recommended to remove the power from the charger after charging is completed.

⚠ Attention

- While charging or charging just completed, it is a normal phenomenon that the bottom of the beauty instrument will be slightly hot.
- Charging time may change when first used or not used for a long time.
- If the beauty instrument is left unused for a long time, the battery of the beauty instrument will enter a protective state and the charging indicator light may not light up for a few minutes due to the battery wake-up process. Therefore, if it is continuously charged, the indicator will light up.
- When the beauty instrument is on and detects low battery, the charging indicator light will flash which shows low battery and the beauty instrument needs to be charged. The instrument can not be turned on and used during the charging process.
- It is recommended that the beauty instrument should be charged once a month, if it is not be used for a long time. It should be completely charged to ensure the battery life of the beauty instrument.

State of power indicator light		Explanation
	Gear indicator light flashes 3 times	Low battery needed charging
	Gear indicator light flashes 5 times	Correct power connection
	The charging indicator flashes as a breathing light	Charging
	Charging indicator light is always on	Charging completed, fully charged

How to Operate the Cleaning Function

- After cleaning your face with makeup remover, wash your face with water and moisten the skin of all areas of your face well.
- Apply an appropriate amount of cleanser evenly on the face and rub the face to produce foam.
- Place the brush head end of the beauty instrument close to the skin and press the brush head operation button to select the cleansing mode for the cleaning function depending on your own need.

⚠ Notice

- When using the instrument for skin cleansing, you need to keep the brush head moving and do not apply it on the same area of the facial skin for a long time, which would lead to part of the skin stratum corneum too thin, thus producing skin health problems.
- If the face has seborrheic dermatitis or local infection, the cleansing function should not be performed, which may lead to more serious skin problems
- If the face has subcutaneous implants or surgical stitches, the cleansing function should not be performed.
- After finishing use, you need to rinse the beauty instrument with water and wipe the water stains with the attached cleaning cloth or soft towel. Please wipe the instrument clean and store it in a dry and ventilated place. Make sure to keep the beauty instrument dry.

		Mode	State
During any state, press and hold the brush head button, and the beauty instrument is forced to shut down.	Off-state	Off-state Mode	None
	Press and hold the brush head operation button to turn on the machine	Mild Cleansing Mode	The brush head will be activated in Mild Cleansing Mode, and the instrument will run for two minutes and then shut down automatically.
	Press the brush head operation button once	Standard Cleansing Mode	The brush head will be activated in Standard Cleansing Mode, and the instrument will run for two minutes and then shut down automatically.
	Press the brush head operation button twice	Deep Cleansing Mode	The brush head will be activated in Deep Cleansing Mode, and the instrument will run for two minutes and then shut down automatically.
	Press the brush head operation button three times		

※ Press means pressing for 0.5 seconds and then release;
press and hold means pressing for more than 1 second and then release.

Steps	Button Operation	Model	Status of gear indicator	Function Description
1	Press and hold the brush head operation button	Mild cleansing		The brush head starts to vibrate in Mild Cleansing Mode with one "tick" and continues to run in "Mild Cleansing Mode" for 2 minutes before the instrument shuts down with another "tick".
2	Press the brush head operation button once	Standard cleansing		During the operation of Mild Cleansing Mode, press the brush head operation button once again, and then the brush head will vibrate to "Standard Cleansing Mode" with a "tick", and continue to run in "Standard Cleansing Mode" for 2 minutes before the instrument shuts down with another "tick".
3	Press the brush head operation button twice	Deep cleansing		During the operation of Deep Cleansing Mode, press the brush head operation button once again, and then the brush head will vibrate to Deep Cleaning Mode with a "tick", and continue to run in Deep Cleaning Mode for 2 minutes before the instrument shuts down with another "tick".
4	Press and hold the brush head operation button	Off-state		During any operation mode, press and hold the operation button of the brush head, with "tick" of the beauty instrument, the gear indicator light flashes 3 times, and the beauty instrument is forced to shut down.

※ Press means pressing for 0.5 seconds and then release;
press and hold means pressing for more than 1 second and then release.

Skin Caring Function Operation

1. Rinse your face with water to remove excess foam or other dirt, and then dry your face with a clean, soft towel.
2. Take out an appropriate amount of skin care essence or care gel and apply it evenly on the skin of the face.
3. Put the treatment head of the beauty instrument close to the skin and press the treatment operation button to turn on the "Positive Electrode Mode"; and it is necessary to slowly drag it to achieve best introduction effect.
4. After the "Positive Electrode Mode", wash off the residue from the face with water and dry the face with a clean and soft towel.
5. Spray an appropriate amount of moisturizing water evenly over the face, then activate the beauty instrument to "Brightening Mode" and hold the treatment head close to the skin, dragging it slowly and gently for the best results.

Cautions

1. During the operation of the treatment mode, it is necessary to keep the treatment end close to the skin throughout the whole process, and not to apply it on the same area for a long time.
2. This function is not available for pregnant women, those with heart conditions, those with pacemakers, those with subcutaneous implants or those with surgically closed wounds.
3. After the operation of the treatment mode is over, some dirt may remain on the care head of the beauty instrument, indicating that the instrument exports heavy metals and other residues from the skin, just use water to scrub the beauty instrument clean.
4. After the treatment mode, you need to rinse the beauty instrument with water and wipe the water stains with the attached cleaning cloth or soft towel. Please wipe the instrument clean and store it in a dry and ventilated place.

Steps	Button Operation	Mode	State of gear indicator light	Function Description
1	Press and hold the treatment button	Mild Positive Electrode Treatment		The beauty instrument starts with a short "tick" sound and enters the Mild Positive Electrode Treatment Mode, while the treatment indicator light flashes red.
2	Press the brush head treatment button once	Deep Positive Electrode Treatment		The beauty instrument starts with a short "tick" sound and enters the Deep Positive Electrode Treatment Mode and activates the micro-vibration massage function, while the treatment indicator light turns solid red.
3	Press the brush head treatment button twice	Brightening Mode		The beauty instrument starts with a short "tick" sound and enters the Skin Brightening and Rejuvenation Mode, while the treatment indicator light turns solid blue.
4	Press and hold the treatment button	Off-state		During any operation mode, press and hold the operation button of the treatment button with short "tick" sound of the beauty instrument, while the treatment indicator light flashes 3 times, and the beauty instrument is forced to shut down.

※ Press means pressing for 0.5 seconds and then release;
press and hold means pressing for more than 1 second and then release.

	Mode	State
During any operation mode, press the treatment button, and the beauty instrument is forced to shut down.	Off-state	None
Press and hold the treatment button to turn on the instrument	Mild Positive Electrode Treatment Mode	The brush head will be activated in Mild Cleansing Mode, and the instrument will run for two minutes and then shut down automatically.
Press the treatment button once	Deep Positive Electrode Treatment Mode	The brush head will be activated in Standard Cleansing Mode, and the instrument will run for two minutes and then shut down automatically.
Press the treatment button twice	Skin Brightening and Rejuvenation Mode	The brush head will be activated in Deep Cleansing Mode, and the instrument will run for two minutes and then shut down automatically.
Press the treatment button three times		

※ Press means pressing for 0.5 seconds and then release;
press and hold means pressing for more than 1 second and then release.

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This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator & your body.