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## Contents

Introduction .....	1
Product Specifications .....	1
Important Safety Instructions .....	1
User Manual Inserts For Part 15 .....	5
Product Diagram .....	6
Cleaning & Maintenance .....	7
Battery Disposal Instructions .....	8

## INTRODUCTION

Thank you for purchasing the Lyric. Please read this user guide carefully and keep it in a safe place for future reference.

## INCLUDED

1 x Lyric massager  
1 x Docking station  
4 x Massage attachments  
1 x Massage attachment tray  
1 x Extension handle  
1 x Charging cable  
1 x Quick Start Guide  
1 x User manual

## PRODUCT SPECIFICATIONS

Model	LYRC0087
Name	Massage Gun
Rated Input:	5V == 2A
Battery	Lithium-ion Battery Pack 2200mAh 7.4V
Charging time	Charging duration for device: 4 hours
Working time	Approx. 4 hours

## IMPORTANT SAFETY INSTRUCTIONS

**WARNING**-When using the massager, the following precautions should be followed to avoid damage to the device and to reduce the risk of burns, fire, electric shock or injury to persons:

1. Inspect the massager and read the user manual carefully before use.
2. Do not charge the massager in an environment where temperature is below 0 °C /32°F or above 40 °C /104°F. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
3. Do not use the massager in an environment where the temperature is higher than 50 °C /122°F.
4. Do not expose the massager to fire or excessive temperature. Exposure to fire or temperature

above 130°C/266°F may cause explosion.

5. Do not use the massager in an environment where it could be exposed to steam or moisture (such as a bathroom).
6. Under abusive conditions, liquid may be ejected from the battery; avoid contact. If contact accidentally occurs, flush with water. If liquid contacts eyes, seek medical help. Liquid ejected from the battery may cause irritation or burns.
7. Do not use the massager on wet skin.
8. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
9. If the massager is moved too quickly from a low-temperature environment to a warm environment, moisture may accumulate on internal metal parts, which might influence the normal operation of the product, or cause damage to it. To avoid this happening, avoid use for at least one hour after moving it between contrasting temperatures.
10. Do not use a massager that is damaged or modified. A damaged or modified device may exhibit unpredictable behavior resulting in fire, explosion or risk of injury.
11. If you encounter any of the following, please discontinue use of the massager and contact customer service (support@experiencecyril.com):
  - a. If it has a damaged charging cord, adapter or charging port.
  - b. If it is not working normally.
  - c. If it has been dropped or damaged.
  - d. If it has been submerged in water.
12. Do not modify or attempt to repair the massager except as indicated in the instructions for use and care.
13. Stop using immediately if you feel discomfort during use.
14. Avoid using on or near the head to avoid injury or health risks.
15. Do not use on or near bony prominences.
16. Only use on soft tissue.
17. Stop using if you experience pain or discomfort.
18. Do not use on areas that are hyper-sensitive or have diminished sensitivity due to nerve or circulatory issues.
19. If you fall into any of the following groups, please consult a physician before using the product:
  - a. If you are pregnant or experiencing a menstrual cycle.
  - b. If you have osteoporosis or osteopenia.
  - c. If you have cardiovascular disease.
  - d. If you have neuropathy.
  - e. If you have a pacemaker or other implanted medical devices.
  - f. If you have an acute disease.
  - g. If you have an active/ open skin lesion and/ or trauma to the skin.
  - h. If you are post-operative and still under the direction of a medical practitioner.
20. Avoid using if you are impaired.
21. Close supervision is necessary when this device is used by, on, or near children and physically challenged individuals.
22. Do not fall asleep while operating the massager.
23. Do not use this product while driving.

24. Use this massager only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
25. Do not carry this massager by its charging cord or use its cord as a handle.
26. Do not use outdoors.
27. Prevent unintentional powering of the device. Ensure the switch is in the OFF-position before connecting the device to a power source, picking it up or carrying it. Carrying the device while your finger or fingers are on the ON/OFF switch or connecting the device to a power source with the switch in the ON-position can lead to accidents.
28. **FCC RF Radiation Exposure Statement Caution:**  
To maintain compliance with the FCC's RF exposure guidelines, place the docking station at least 20cm from nearby persons.
29. **IC RF Radiation Exposure Statement Caution:**  
To maintain compliance with the RF exposure guidelines, place the docking station at least 10cm from nearby persons.

## **ANGER**

### **To reduce the risk of electric shock:**

1. Unplug from the electrical outlet before cleaning.
2. Do not reach for the massager if it falls into water. Unplug immediately.
3. Do not use while bathing or showering.
4. Do not place or store massager where it can fall or be pulled into a bathtub or sink. Do not place in or drop into water or other liquid.
5. Do not use massager in close proximity to loose clothing or jewelry.
6. Keep long hair away from massager while in use.
7. Do not tamper with or alter the unit in any way. There are no user-serviceable components or parts inside.
8. This unit contains a rechargeable Lithium-ion battery, which is neither replaceable nor user-serviceable.
9. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
10. Avoid inserting objects into any opening.
11. To disconnect, power down the Lyric (press trigger button until it powers down), then remove plug from outlet.
12. Recharge only with the charger specified by the manufacturer. A charger that is suitable for one type of battery pack may create a risk of fire when used with another battery pack.
13. Have servicing performed by a qualified repair person using only identical replacement parts. This will ensure that the safety of the product is maintained.
14. This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance (5V  $\overline{\text{---}}$  2A).

## SAVE THESE INSTRUCTIONS USER MANUAL FOR IC

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

## USER MANUAL INSERTS FOR PART 15

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the warranty and the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.


However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- ☐ Reorient or relocate the receiving antenna.
- ☐ Increase the separation between the equipment and receiver.
- ☐ Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- ☐ Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

## PRODUCT DIAGRAM
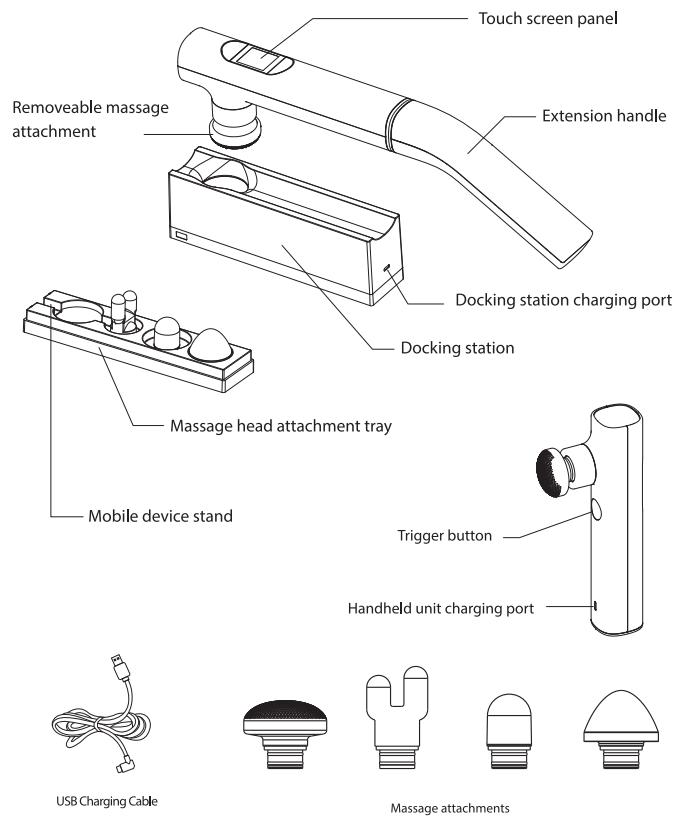
## CLEANING & MAINTENANCE

1. Battery level should be maintained above the medium level at all times. Regular charging and maintenance will help prolong battery life.
2. Keep the device clean. For daily maintenance, clean with neutral detergent and a soft cloth. Do not clean using organic solvents or other chemicals. Do not rinse in water.
3. Avoid high pressure washing and avoid hard objects scratching the product surface.

Certificate Charger List:  
Use only with Listed Class 2 Power Unit or Lyric  
LYRC0087 wireless charger

**BATTERY DISPOSAL INSTRUCTIONS**

The lithium battery in this device is a valuable resource that can be recycled. When the device is at the end of its service life, remove the lithium battery and dispose of it at a facility or establishment that offers battery recycling services. This device should not be discarded with household waste.

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