




beoka

INSTRUCTIONS

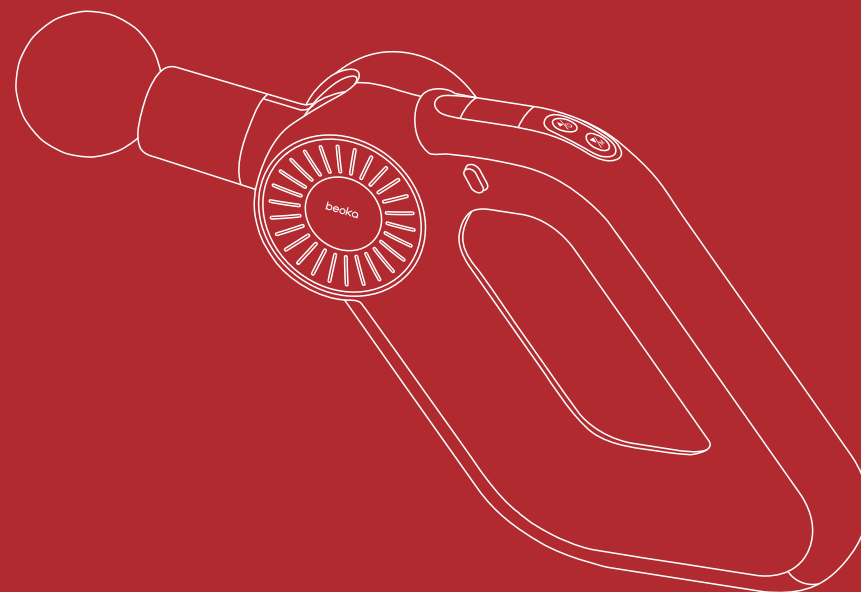

DEEP MUSCLE STIMULATOR

D6 Pro-A

- Thank you for buying Beoka deep muscle massager.
- Before you use this product, please read the instruction manual carefully and keep it properly!
- Execution standard: Q/91510108629517429X•3-2022
- In case of any deviation from the illustrations in the manual, the real object shall prevail



Sichuan Qianli-beoka Medical Technology Inc.

Address: Longtan Industrial Park 2nd Sec, east 3rd ring road, Chenghua district, Chengdu, Sichuan, China

Phone: 86-28-84215340 84215342

Fax: 86-28-84215341

Postcode: 610052

Hotline: 4000-8899-01 Website: www.qlbeoka.com

FCC Warning:

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.












This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Matters needing attention.....	01
Figure、 symbol, abbreviation.....	01
Safety notes	01
Outline.....	04
Instruction.....	04
Purpose.....	04
Main technical parameters.....	04
Main structure.....	05
Main structure of the device.....	05
Usage method.....	06
Preparation before use.....	06
Operation steps	06/07
Malfunction and warranty.....	08
Maintenance and repair.....	08
Notice when products are discarded.....	09
Installation, transportation and storage conditions of products.....	10
Product warranty certificate/warranty card.....	11

Matters needing attention

Figure、symbol, abbreviation

Read the manual carefully before use, and pay special attention to the warning points.

Figure, symbol, abbreviation			
	Direct current		Identification of reclamation and utilization of lithium batteries
	WEEE Recovery Identification		Do not roll
	Be careful! Access to Random Documents		Please read the instructions carefully.
	Fragile		Upward
	Prevent of rain		Paper recycling
	Stacking layer limit		

Safety notes

The following contents must be known before use, otherwise it may cause damage or danger!

About the User Population and Location

- For use only by healthy adults other than the elderly;
- Use the product only in the soft tissue of the body without pain or discomfort;
- Do not use on the head or any hard or bony part of the body;
- Silting may occur either under moderate pressure or under pressure. Check the massage area regularly and stop using it once there are signs of pain or discomfort;
- It can be used on the dry and clean skin surface of body parts, and gently press and move on the skin;
- Do not use near the skeleton with less muscle coverage, such as tibia, acrotarsium, the back of hand;
- Do not use in or around sensitive areas such as eyes, teeth, breast augmentation,

Matters needing attention

reproductive organs, prostheses, etc;

- Do not use the product at any possible fracture site;
- If the pain is persistent, stop using it immediately and consult your doctor;
- This product is not a medical device. It is not recommended for medical use and is not intended for medical use in hospitals;
- This product is not a medical device and is not recommended for medical use. It is not intended for medical use in hospitals;
- Please check the host and cable before use. If there is any damage, the manufacturer or its designated professionals must replace it in order to avoid danger;

About Instrument Operation

- Vibrating head which can produce the best effect should be used;
- Do not place any object in the ventilation hole of the product;
- Do not immerse in water or let water enter the ventilation holes of the product;
- Do not fall or improperly use;
- Before each use, we should carefully check the operation status of the main engine, batteries, chargers and so on to ensure that they can meet the normal use;
- Do not change or modify the product in any way;
- Never let the product run or charge in an unattended state;
- The special charger provided by the manufacturer must be used for charging. It is strictly forbidden to use other power supply equipment to charge the product;
- It is forbidden to start the equipment during charging;
- It is forbidden to plug the battery after the power switch is turned on;
- It is forbidden to put equipment into fire;
- Don't use it for more than one hour continuously. Let the equipment stop for 30 minutes before reusing.
- Note: When the device appears " ⚠ ", it means that the temperature of the device is too high, and it needs to stand for 10 minutes before restarting the device.


About Batteries

Note: The battery used in this product is designed to be non-removable by the user; please do not disassemble the battery without permission.

- Do not put batteries in water or wet them!
- It is forbidden to charge the battery under the condition of fire source or extreme heat. Do not use or store batteries near heat sources such as fires or heaters! If the battery leaks or smells, it should be removed immediately from the vicinity of the open fire;
- Do not put batteries into fire or heat them;
- It is forbidden to short circuit the positive and negative poles of batteries with wires or other metal objects. It is forbidden to transport or store batteries with necklaces, hairpins or other metal objects;

- Do not use nails or other sharp objects to pierce the battery shell. Do not hammer or pedal the battery!
- Do not hit or throw batteries!
- It is prohibited to decompose batteries in any way!
- Do not place batteries in microwave ovens or pressure vessels!
- If the battery emits odor, heat, distortion, discoloration or any other abnormal phenomena, it should not be used; if the battery is in use or charging, it should be immediately removed from the instrument or charger and stopped using!

About Charging

- The device must be charged using the power adapter provided with the product. Do not use the ordinary 5V power adapter to charge the device. If the icon “” appears on the display, it means the charging is abnormal, please replace the power adapter.
- Connect the output end of the charging cable to the TYPE-C charging port at the bottom of the massager, and then connect the input end of the charging cable to the power adapter provided with the product;
- Before the first use, charge for about three hours to fully charge;
- There is a charging animation on the display screen of the instrument, and the indicator light at the charging interface is solid green, indicating that it is charging;
- When the battery is fully charged on the display screen of the instrument, the indicator light at the charging interface is off, indicating that the battery is fully charged;
- The battery can be charged at any time and at any power level;
- When the power on the display of the instrument flashes, it means that the power is insufficient and needs to be charged;
- It is not recommended to completely drain the battery;
- The average operating time of the equipment after loading is about three hours, which depends on the vibration frequency and applied pressure during use.

Introduction


Deep Muscle Massage Instrument is a device for deep muscle and fascia relaxation. It is suitable for physical massage after exercise. It can reduce muscle stiffness and pain, alleviate muscle fatigue and accelerate blood circulation to achieve the purpose of improving the health of body soft tissue.

Purpose


- Physical massage after exercise

Main technical parameters

- Product Name: DEEP MUSCLE MASSAGER
- Product Model: D6 Pro-A
- Structure: It mainly consists of control part, massage head and lithium battery;
- The vertical vibration frequency ≤ 60 Hz;
- Amplitude: 16mm;
- Rated power: 25W;
- Instrument: rated voltage: 21.6V \approx battery capacity: 2500mAh
- Adapter: Input: 100-240V~50/60Hz 0.5A
Output: 5V \approx 3A/9V \approx 2.22A/12V \approx 1.67A
- Timing: 10 minutes intelligent timing
- The instrument has output control function switch;
- Massage head can be replaced as needed;
- Running mode: intermittent loading.

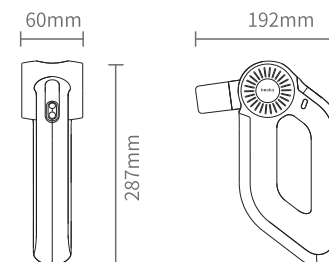

Main structure of the device

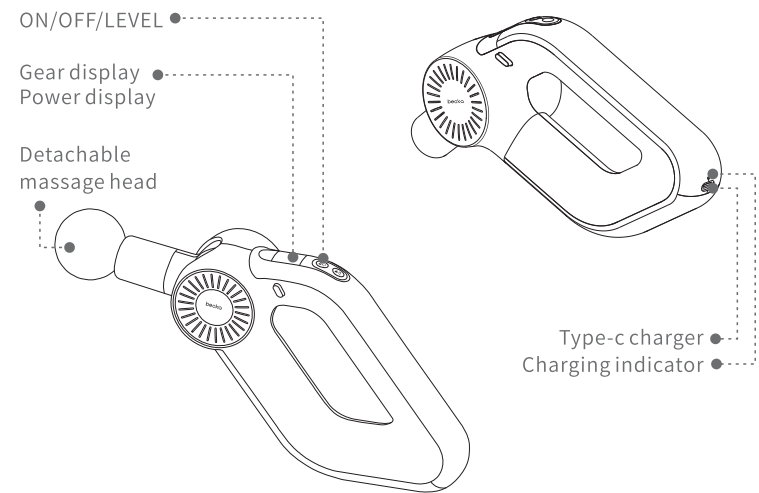


Preparation before use

- Open the packing box and check whether the main engine and accessories are in good condition. If there is any defect, please inform your supplier immediately;
- select the appropriate massage head to connect with the host.

Operation steps

1) Connect the device and the massage head

When the device is turn off, insert the massage head completely with pressure and pull it out when removing the vibrating head (note: the installation of the U-shaped head should be positioned. When the buckle is inserted into the slot, it means the installation is in place. If the massager head is not installed in place, there will be a risk of falling off while use.)

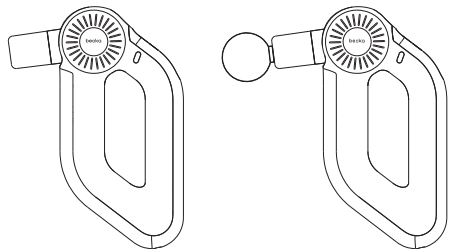


2) Button function description refer to the following table

Button	Operation	Function
	Long press	Turn on the device while it is off
		Turn off the device while it is on
	Press	In M1 mode or standby state, increase the gear
		In the standby state of M2/M3/M4 mode, the device works
	Long press	In standby mode, mode switching
		In the working state of M2/M3/M4 mode, enter the standby state
	Press	In the working state of M1 mode, reduce the gear
		In the working state of M2/M3/M4 mode, enter the standby state


Malfunction and warranty

Notice:

- 1 When the device is powered on, it defaults to M1 mode
 - 2 Switching modes can only be done in standby mode
 3. M1 is a common gear mode (0-6 gears can be adjusted), M2, M3, M4 are fixed modes and cannot be adjusted.
 4. The device will automatically enter the standby state after running for 10 minutes
 - 5 Automatic off after 60s without any operation in standby state
- 3) Display Description

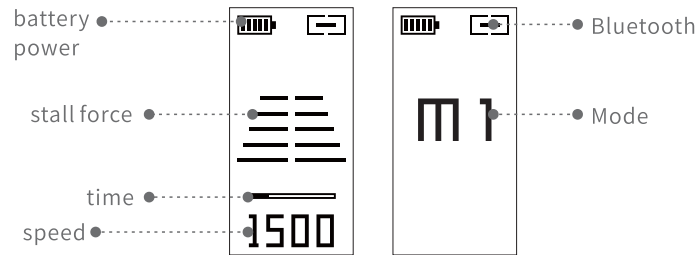

- 4) The direction of the elbow arm of this equipment can be adjusted, and the total adjustment range is 0-90°. Press and hold the elbow arm rotary switch on the left side of the device to rotate the elbow arm.

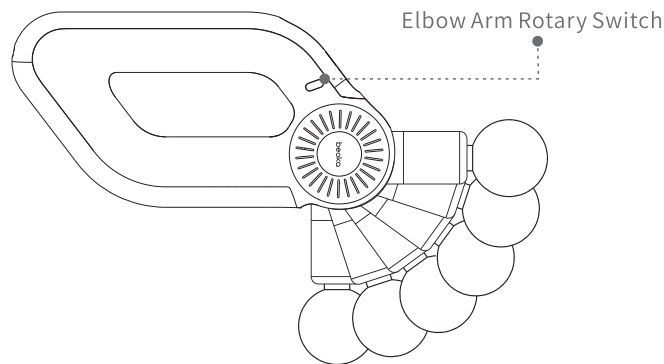

Malfunction and warranty

About Maintenance



- Clean the equipment with a slightly wet towel and dry it with a soft cloth;
- If you do not use the equipment for a long time, please separate the battery from the mainframe and store it in the box;
- Do not disassemble the screw or attempt to disassemble the equipment;
- In order to achieve the best use of new batteries in a short period of time, complete charging and discharging can be carried out three or more times.
- When the equipment is not used for a long time, the battery should be charged once a month.

About repair

- It is forbidden to disassemble the product without permission, and artificial ring breaking is not within the scope of warranty;

Notice when products are discarded

Please support the recycling of rechargeable batteries (lithium-ion batteries).

 Danger	<p>The battery used in this product is a special rechargeable battery, which shall not be used in other products, please do not take out Rechargeable battery for charging.</p> <p>The above conditions can cause heating, fire, rupture and other phenomena.</p> <ul style="list-style-type: none">• Do not throw it into the fire to overheat it;• Do not pierce with nails, impact with gravity or conduct decomposition and transformation;• Do not connect the positive and negative electrodes of the battery directly with metal;• Please do not store it with necklace, hairpin, etc.• Please do not charge, use or place in high temperature places such as near the fire or in the hot sun. <p>The above conditions can cause heating, fire, rupture and other phenomena.</p>
 Danger	<p>Please keep the removed rechargeable battery out of the reach of infants.</p> <ul style="list-style-type: none">• If swallowed by mistake, please see a doctor immediately.• When the battery liquid leaks, do not touch the liquid with your bare hands, and handle it as follows:• When the battery liquid enters into the eyes, there is a risk of blindness. Do not rub the eyes directly with your hands. Wash them with clear water immediately and seek medical advice. When the battery liquid touches the body or clothes, it may cause skin inflammation or injury. Please wash it with clear water. If there is still discomfort, seek medical advice immediately.• Only when the product is discarded, it can be disassembled. If it is not observed, it may cause fire, electric shock, injury and other accidents. If it needs to be repaired (except for battery replacement), please consult the original seller.

The lithium ion battery used in this product is a valuable resource that can be recycled and regenerated. Please remove the lithium ion battery and dispose of it to a professional organization that can perform rechargeable battery recycling and recycling.

Installation, transportation and storage conditions of products

Environmental conditions for product installation:

- Ambient temperature: +5°C~+40°C
- Relative humidity: ≤90%
- Atmospheric pressure: 70kPa~106kPa

Transportation

- Ambient temperature:-20°C~+50°C;
- Relative humidity: ≤90%;
- This product is suitable for automobiles, trains, ships, airplanes and other modes of transportation;
- During transportation, it should be fixed to prevent violent vibration, collision, roll and fall; and be upward, careful, light and rainproof;
- No more than 5 layers of massager stacked.

Storage

- Massage equipment should be stored in a well ventilated, non-corrosive gas warehouse, 250 mm off the ground;
- No more than 5 layers of massager stacked;
- When massage instrument is not used for a long time, it should be packed and stored after cleaning. It is better to take it out and recharge once a month to prevent dampness.

Malfunction and warranty

Product warranty certificate

1. From the date of purchase, beoka products can enjoy a free one-year warranty for the body of massage instrument with the purchase invoice.
2. We will not provide free warranty service for faults caused by the following user's personal reasons. Such as:
 - A) failure caused by unauthorized disassembly and refitting of the product;
 - B) failure caused by careless falling during use and handling;
 - C) failure caused by lack of reasonable maintenance;
 - D) failure caused by failure to operate in accordance with the correct instructions in the operation manual;
 - E) faults caused by improper repair of repair shop not authorized by beoka.
3. Repair services beyond the scope of warranty will be charged according to regulations.
4. Please call the customer service hotline when you need to provide warranty service.
5. During the warranty service, if necessary, we can provide the product circuit diagram and repairable component information to the qualified technical personnel identified by us.
6. Ensure to continue to provide repair parts within 3 years after the production stop.
7. Please refer to the host or packing box for the production date.

Warranty card

Product model: _____

Name: _____

Date of purchase: _____

Address: _____

Sales shop name: _____

Postal Code: _____

Address: _____

TEL: _____

remarks