

Visual blackhead suction instrument  
M/N:X15S



# USER MANUAL

# Table of Contents

App Operation.....	1
App Download.....	1
App Functions.....	2
Care Steps.....	7
Video for Use.....	14
Introduction.....	14
Functions.....	17
Attention.....	19
Charging.....	20
Cleaning and Safety.....	20
Specification.....	21
Precautions of Usage.....	22
Customer Service.....	22

# App Operation

The phone can easily connect to the WiFi of the device. You can more clearly watch the pore enlargement and blackheads sucking out on the App.

## App Download

- Scan the following QR code to download the App.
- Please search [ULOOK] on App Store/Google Play to download if you cannot scan the QR code.




- Available for iOS, Android 6.0, and above systems.
- Follow the latest guidelines of App for operation, maybe slightly different due to the upgrade and update of the App.

# App Functions

1.Device: Press the power button for two seconds to turn it on. Ensure the device is on, or the phone can't connect to the WiFi.

2.Phone: Turn on [Settings] → [WLAN] → Connect the WiFi [Vellcam-xxxxxx]. Note: [xxxxxx] has six different characters.

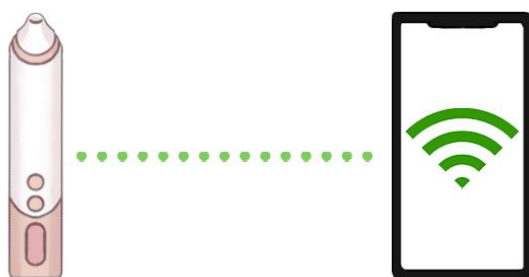

- A ✓ shows a successful connection.



- Please turn on the device again and reconnect if connect is unsuccessful or unstable.

# App Functions

3.The interface shows two buttons.


- Click [Start] to watch the image transmitted with the camera of the device.
- Without WiFi connected, it will turn to the App Access permission settings.



- Click the button to select the App system language.

4.The interface shows three available settings.



- Record: Click the button to save the video.

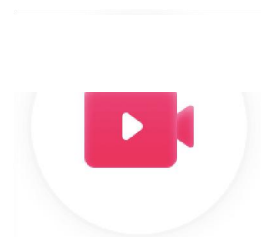

- Snapshot: Click the button to save the image transmitted by the camera .

# App Functions

- Zoom: Change the range watched and magnification to clean the pores as desired.

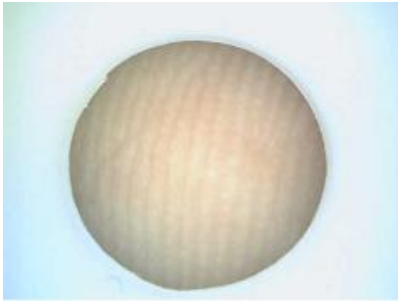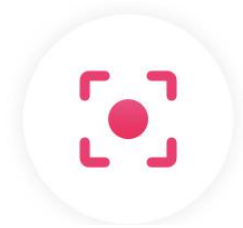

[General]



[Enlarge]

5.The App won't work if you reject the access permission.


- Read the [Terms of Use and Privacy Policy] carefully.
- Click the pink characters to know more details.

## Terms of Use and Privacy Policy


Please read carefully and fully understand the terms and conditions, including but not limited to: In order to provide you with basic services, we need to collect personal information such as your device information and operation logs. You can read [Terms of Use and Privacy Policy](#) for details. If you agree, please click "Agree" to start enjoying our services.

Disagree and exit **Agree**

# App Functions

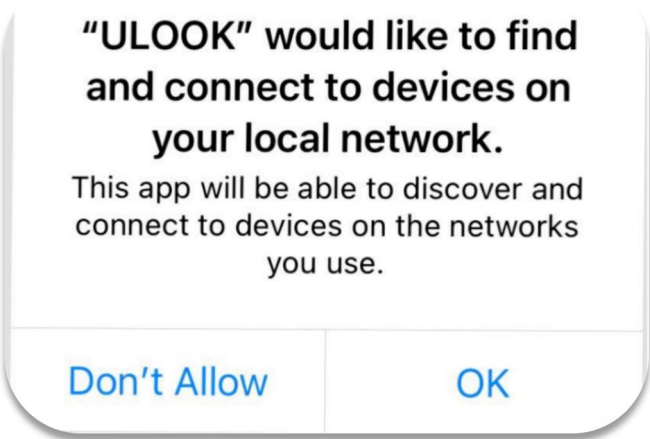

- The App will access the [Cellular] and [WLAN].

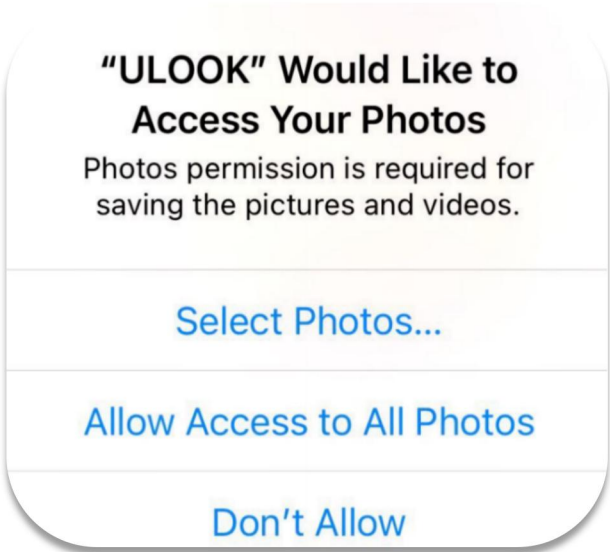

- The App will access the camera and photos.

# App Functions




- Enable or disable the access as required.

6.Close the App after use.


7.How to remove the App.

- Delete the [ULOOK] App.

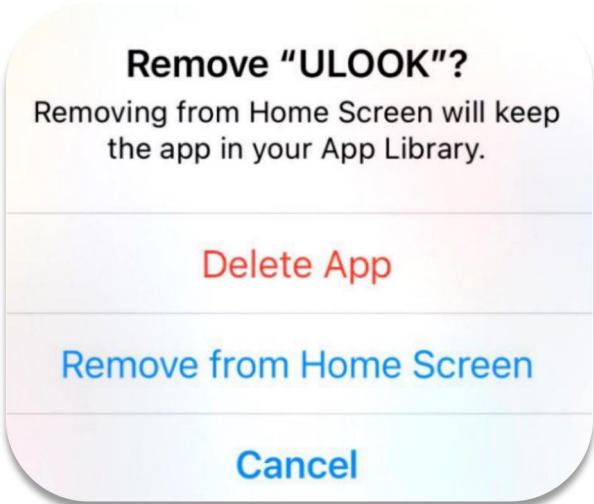



# App Functions

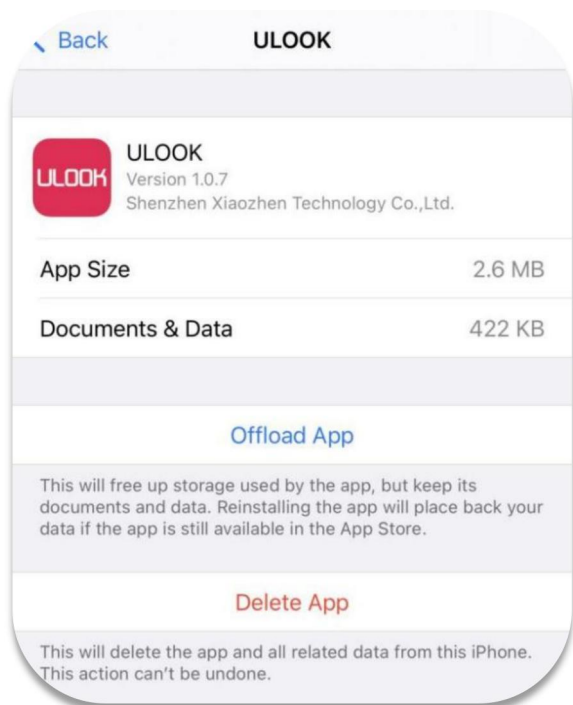

- iPhone: Settings → General → Storage → [ULOOK] App → Offload / Delete.

## Care Steps

Removing blackheads needs a combination of several steps and keeping regular pore cleansing.

- 1.Download and install the App. Check [App Download] for details.
- 2.Connect the App with the device. Check [App Functions] for details.

# Care Steps

## 3. Facial care before use.

- Remove makeup, then clean the face.



- Use the hot towel or facial steamer to steam the face for 5-8 mins to help open the pores.
- A hot water bath isn't the first choice because it will make moisture loss.

# Care Steps

- Applying blackhead softener to the nose can help soften blackheads and bring them to the surface.



4. Attach the appropriate suction head to the device as your needs. It is best to start with a small round suction head the first time.
5. Please test the suction power on your palm or forearm before use.
6. Correct cleaning areas.



- It will work well in the blue area.
- It should be careful to use on the gray area.

# Care Steps

7. Choose the appropriate suction power: Level 1 → Level 2 → Level

3. It is best to start with the lowest setting, keep it moving. Increase the suction level gradually according to skin feeling.



8. You can use pressure downwards, making the suction head tighter against the skin and slide it across your skin.


9. Two situations.

- The device sucks blackheads away, but just those on the surface.
- It will suck stubborn blackheads a little out at the first use, and they will be pulled out entirely after a few more sessions.
- Note: You can also use the included blackhead remover set to handle those blackheads. At this time, they can be pulled out more quickly than before.

# Care Steps


## 10.The correct ways for use.

- Move the probe in a horizontal direction instead of a vertical direction.
- Keep moving and do not stay in the same position for more than 2 seconds. You will see pores that don't pop, keep going over them to avoid suction marks.

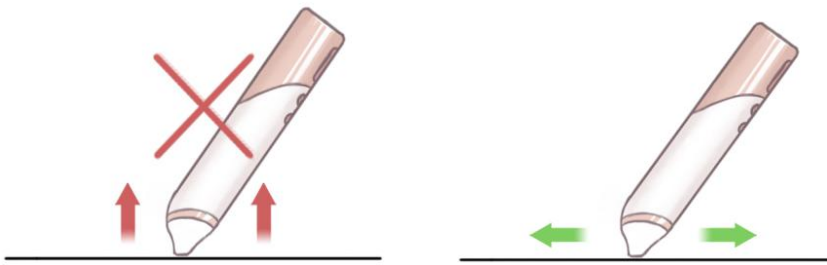

- Move the suction head in a circular or linear massage.
- Don't let the observation of the image distract you from the time when the probe stays on the skin.
- Tilt the side of the probe to release the pressure and take it off quickly.

# Care Steps

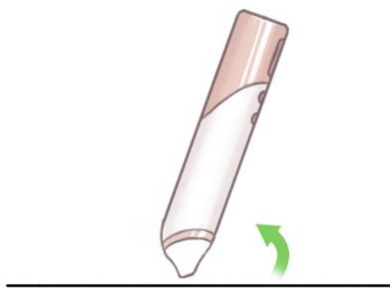

- Drag it while sucking the skin. After taking the suction head away, suck and drag like before again. Keep doing this repeatedly.

## 11. Facial care after use.

- After cleaning, apply pore minimizing toner or cold mist steam to help shrink pores.



- Then turn on the red LED light for circular massage for about 5 mins to promote better skin absorption of the essence.
- Switch to blue LED light to shrink pores further, as well as improve grease secretion.



# Care Steps

- Finally, apply a moisturizing mask, gel, or cream to hydrate the skin.
- It is normal for the skin to turn red slightly after suction. It will dissipate within 1 hour generally if you follow the manual.
- If the probe inadvertently stays on your skin for a long-time, resulting in redness or bruises. In that case, they will disappear in about 1-3 days.
- Note: This time is for reference only, depending on the skin damage. You can also apply ice to the skin.

12.Pore cleaning is recommended 1-2 times per week. Once a week for sensitive and dry skin, for oily and combination skin twice a week.

13.After two weeks: The pores are improved noticeably, and the facial skin becomes smoother.

# Video for Use

Scan the following QR code with your phone, watching the specific process of using the blackhead vacuum.

## Introduction


### Features

- The powerful suction can effectively and deeply suck blackheads grease, dirt, and exfoliation from pores.
- The camera inside allows you to directly observe the actual skin condition and the process of blackheads removing, avoiding blind cleaning, or having omissions.



# Introduction

- Three different suction heads.
- With both red LED light and blue LED light.
- Three levels of suction for various needs.

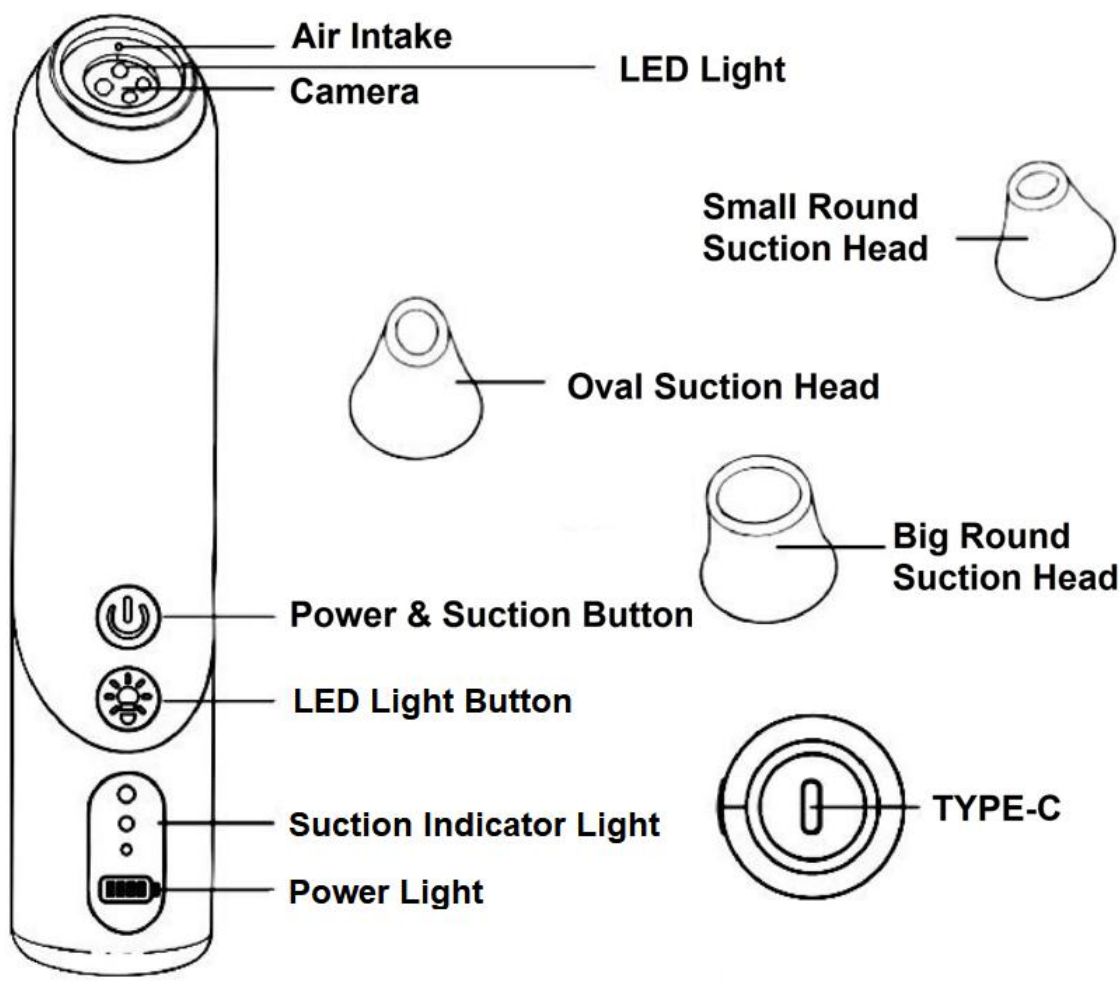

# Introduction

## Suction Heads

The suction heads are made up of soft silicone. Those are safe, harmless, and non-irritating. The edge of the suction head will not scratch the skin, more kind to sensitive skin.

### Big Round Suction Head

- Strong suction power and wide range.
- T-zone.
- It should be careful to use on thin skin.




### Small Round Suction Head

- Moderate suction power.
- The end and alae of the nose.
- Better for thin and sensitive skin.

### Oval Suction Head


- Moderate suction power.
- Nasal grooves and areas with wrinkles.

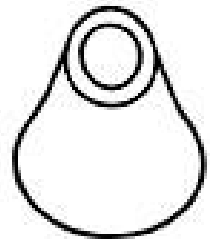

# Introduction

## Suction Heads Connection

Connect the suction head to the top of the device. Make sure that the suction head is fully connected to the device to avoid air leakage.

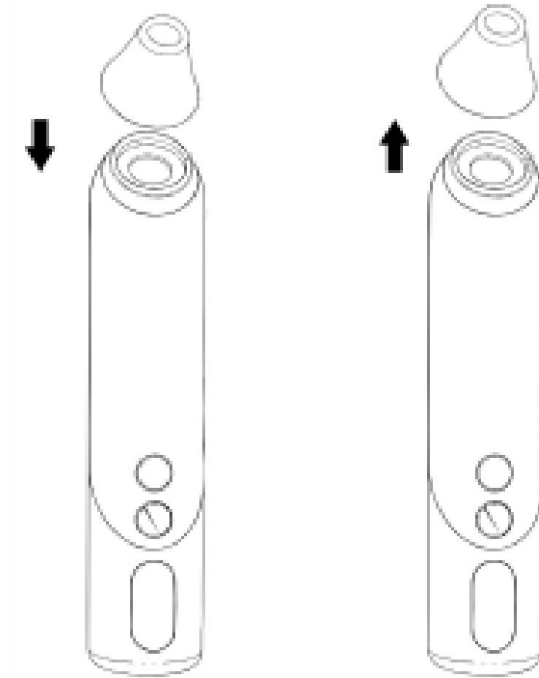

# Functions

## Power & Suction Button

- Press and hold for 2 seconds to turn on, and the device will be on standby. The LED white light lights up, and the display shows the green power light.

# Functions

- Short press the button on standby. In a cycle, the suction power will switch between [Standby on- Level 1 - Level 2 - Level 3 ]. At the same time, the corresponding round white light in the display will light up.
- Long press and hold the button for 2 seconds to turn it off on standby or working.

## LED Light Button

- Short press the button on standby or working.
- LED light colors will change once with each press. Colors [White - Red - Blue] will vary in turns by this separate button.

## Protection

- The device will automatically shut down if not operated within 15 mins, whether on standby or working.
- It will automatically shut down after five flashes of the power light if the power supply voltage is low to 3V during the use of the device.
- You cannot turn it on during charging.

# Attention

## **Read and follow the following safety terms carefully.**

- 1.Do not use it on the chin, corners of the eyes, around the eyes, and corners of the mouth.
- 2.Do not use on skin with wounds, rashes, performed plastic surgery, inflamed, itchy redness.
- 3.Consult a doctor before using this device with any skin problems.
- 4.Do not look at the running LED Light at close range.
- 5.Children do not use it, avoiding accidents and injuries.
- 6.It is prohibited to charge when wetted.
- 7.This product is not waterproof.
- 8.Keep the device and USB connector dry, but not be in direct sunlight or close to a heat source.
- 9.Do not knock or drop the device.
- 10.Stop using it immediately with uncomfortable feelings.
- 11.Stop using it with any malfunction and contact customer service.
- 12.Do not disassemble, modify or replace the battery inside.

# Charging

1. Charge the device fully before using it for the first time.
2. Use a standard qualified 5V/2A adapter for charging.
3. The red power light is on during charging.
4. The green power light is on when charged fully.
5. The red power light flashes when the power is low and needs charging.  
Note: Confirm whether the USB connector jack has poor contact or the power supply is standard if you cannot turn the power on after charging for 10 seconds.
6. Fully charge at least once every six months.
7. Charge, use, and keep the device at  $-10^{\circ}\text{C}$  to  $40^{\circ}\text{C}$ .
8. Do not charge it for more than 12 hours at the longest.

# Cleaning and Safety

1. Ensure the device is off when cleaning.
2. Do not wash it with water.
3. The device is not waterproof, so wipe it clean with a damp tissue or cloth after each use.

# Cleaning and Safety

- 4.Clean the lens with a cotton swab.
- 5.Do not use a corrosive solution to wipe the device.
- 6.Keep the device dry.

# Specification

Product Name	Blackhead Remover	Material	ABS & Silicone
Weight	82 g	Charge Input	5V/2A
Charging Interface	Type-C	Rated Power	3W
Charging Time	80 mins	Working Time	40-60 mins
Battery Capacity	500mAh	Working Temperature	At -10℃ to +40℃
Volume Level	≤55dB	Suction Level	58-65 KPa
Resolution	640*480/480*480	Waterproof	No
Adapted Systems	Available for iOS, Android 6.0, and above systems		
Package	Blackhead Remover * 1, Big Round Suction Head * 1, Small Round Suction Head * 1, Oval Suction Head * 1, USB Charging Cable * 1, User Manual * 1		

# Precautions of Usage

1. Please search [ULOOK] on App Store/Google Play to download if you cannot scan the QR code.
2. Charge fully and turn again if unable to turn it on for the first use.
3. Please turn on the device again and reconnect if connect is unsuccessful or unstable.
4. Stuttering or flickering of the screen in use may be due to the insufficient power of the device. Just charge and reconnect.
5. The device's temperature in use rises slightly, which is normal, and you do not need to worry.
6. Ensure steam your face and apply the blackhead softener before use. You will feel pores improved after several times instead of once use.
7. It is normal for the skin to turn red slightly after suction. It will dissipate within 1 hour generally if you follow the manual.





## FCC Statement

This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions: 1) this device may not cause harmful interference, and 2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

## FCC RF Radiation Exposure Statement

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.