



Dear Consumer: Thank you for purchasing a LOBOOR product. The LOBOOR brand focuses on the prevention and treatment of sub-health and chronic diseases...

Catalogue: 1. Catalogue, 2. Cautions, 3. Technical specifications, 4. Electrical schematics, 5. Product diagram, 6. Mode of operation, 7. Warm Tips, 8. Daily Maintenance, 9. Troubleshooting, 10. Description of product and parts recycling, 11. Precautions when disposing of products, 12. Product Warranty Card.

Cautions: Dangers: 1. Never use with the following medical electronic devices: 1. Cardiac pacemakers and other implantable medical electronic devices; 2. Life-sustaining medical electronic devices such as artificial heart lungs; 3. Electrocardiographs and other medical electronic devices...


When using the product, the electrodes should be closely fitted to the skin and evenly contacted. The electrodes of this product have a heating function and may cause low temperature burns when the same part is used for a long time (more than 15 minutes).

5. Pacemaker users: 6. Patients with metal built-ins: Do not use near the heart, above the neck, on the head, in the mouth, on the pubic area or on areas with skin disorders and do not massage the soles of both feet at the same time...

If the product does not work, please cut off the power immediately and take it off, otherwise high temperatures or short circuits may occur and accidents may occur. Also contact customer service for after-sales advice.

Recommendations: Please make sure that the mains button is working properly before using it. When using the product for the first time, or if the battery power has been depleted after a long period of time, please follow the instructions for charging to ensure that it works properly.

Technical specifications: Product name: Neck Massager; Rated charging voltage: 5V; Rated input power: 10W; Rated voltage range of the mainframe: DC3.2V-DC4.2V; Power rating of the main machine: 2W; Rated voltage range of the remote control: DC3.2V-DC4.2V; Remote control power rating: 0.5W; Mainframe battery capacity: 800mAh; Remote control battery capacity: 110mAh; Implementation standards: GB4706.1-2005/GB4706.10-2008.

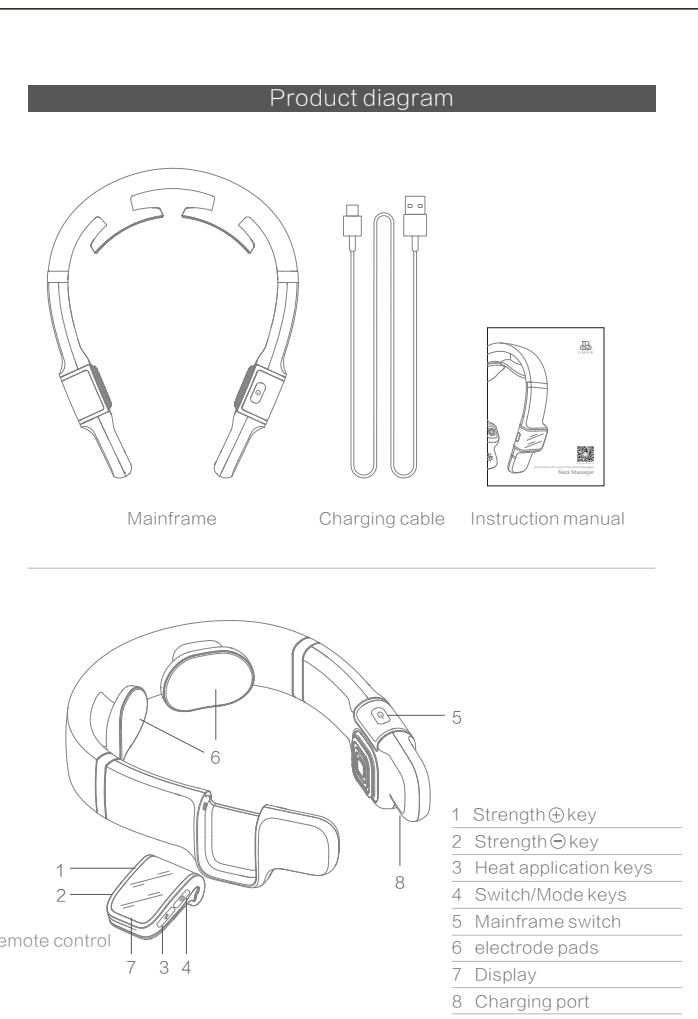

Wearing style: 1. When the red light of the main unit flashes, accompanied by the voice prompt "low battery", the main unit is low on battery, please charge it in time. When charging, use the USB charging cable to connect to a power adapter (5V/1A/2A) that meets the standard for home appliances or connect to the USB port of your computer for charging...

Instructions for use: 1. Before use, first moisten the skin of the neck with a wet towel, then wear the product on the neck, so that the electrode piece close to the skin of the neck, adjust the position of the product, open the heat and adjust the appropriate massage intensity and mode through the remote control.

5. The "Strength key" of the remote control is the strength reduction key, each short press will reduce the output strength by one grade, and accompanied by corresponding voice prompts, such as "strength 14", "strength 13"..... The product has a total of 16 output strengths (0-15) and can be adjusted to a minimum of 0 strength.

Troubleshooting: If you encounter problems when using the neck Massager, please refer to the following points for tips and solutions. If the neck massager is still faulty or not working, switch off the power.

Table with 2 columns: Status, Possible causes/processing. Rows include: Unable to start, No charging or no charging light indication, Shows normal but no sensation, Has a tingling sensation, Very weak stimulation.


Product and parts recycling instructions: Table of names and contents of harmful substances in products. Columns: Part name, Pb, Hg, Cd, Cr(VI), PBB, PBDE. Rows: Hardware, PCB Plate, Plastic parts.

Precautions when disposing of products: Please support the recycling of rechargeable batteries (lithium ion). The batteries used in this product are special rechargeable batteries and must not be used in other products.

The lithium-ion batteries used in this product are a valuable resource that can be recycled. Please dispose of them at the end of use to a professional organization that can recycle rechargeable batteries.

5. Please keep your warranty card in a safe place as it will not be replaced if lost. 6. In addition to the maintenance of the products specified in this warranty card, the user with a valid warranty card shall enjoy other rights within the scope of the relevant national laws and regulations during the repair period.

User warranty vouchers: Name, Model, ADD, Warranty Period, Invoice Number, Date of purchase. Maintenance records: Date of delivery, Date of calibration, Breakdowns, repairs, treatment. Quality certificate: Production models, Date of production, Date of inspection, Inspectors.

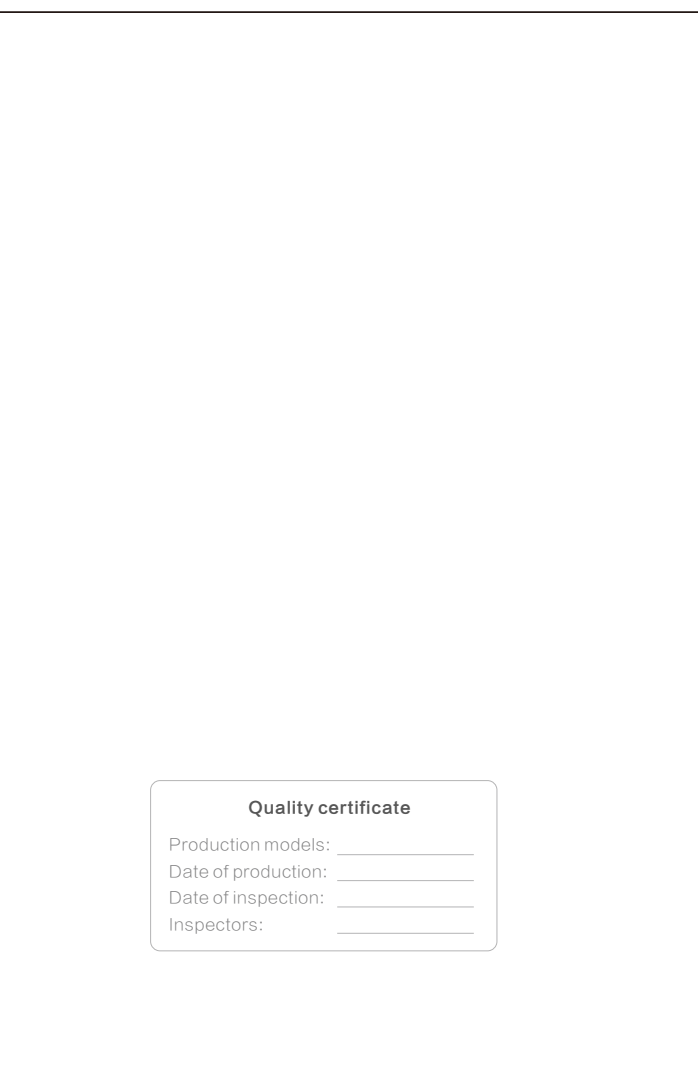

Marketing Centre: Guangzhou Warren Electronic Technology Co., Ltd. Address: Room 201, No. 111 Airport Road, Baiyun District, Guangzhou, China (Self-addressed 2F AD). Tel: 020-86342676. Manufacturer: Sichuan Yuka Medical Equipment Co., Ltd. Production Address: Second floor, North Wing, Building 102, Mianyang Export Processing Park, No. 291, East Feiyun Avenue, Mianyang. High-tech Zone. National Service Hotline: 400-8018-358.