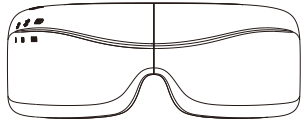
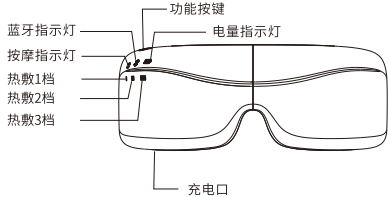
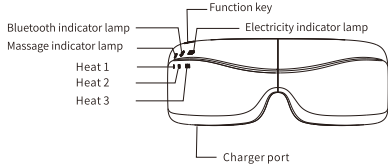



<div><div><div>D13 Heated eye mask</div><div>D13 护眼仪</div></div><div></div><div><div><div>使用前请仔细阅读本使用说明书,并妥善保管</div><div>Please read this instruction manual carefully before using the product and keep it in a safe place.</div></div></div></div>	<div><div>D13 护眼仪</div><div>产品型号: D13 工作电压: DC5V 充电电流: 0.9A 工作电流: 50mA-1100mA 使用功率: 3.7W 电池容量: 1200mAh 产品重量: 169g 执行标准: QB/T QB 1108.02-2016 质 检: 合格品 产 地: 中国制造</div><div><div>按键功能及充电指示</div><div></div></div></div>	<div><div>功能操作说明</div><div>1.长按3秒开关机，默认进入蓝牙模式，播放蓝牙模式提示音，连接后播放已连接提示音； 2.短按按键，开启1档加热模式，亮热敷1档灯/35度； 3.再次短按按键，开启2档加热模式，亮热敷2档灯/40度； 4.再次短按按键，开启3档加热模式，亮热敷3档灯/45度； 5.再次短按按键，关闭加热 6.双击按键，进入1档振动按摩模式，按摩指示灯亮绿灯； 7.再次双击按键，进入2档振动按摩模式，按摩指示灯亮蓝灯； 8.再次双击按键，关闭振动按摩模式</div><div><div>本机支持蓝牙连接；需要听音乐或其它声频时，在开机3分钟内连接。连接方法与通用蓝牙连接相同。蓝牙默认名称CFJG。</div></div></div>	<div><div>充电方法</div><div>1.充电时，请使用5V/1A或以上的正规充电适配器。 2.使用本产品前，请使用本工厂提供的充电线，给本产品充满电。 3.充电时电量指示灯长亮，充满亮蓝灯。 4.本产品一般充电时间约为2小时；不可长时间一直在充电，可能影响内部主板及电池使用寿命。 5.低电提示：红色提示灯常亮。</div><div><div>保养</div><div>1.使用前尽可能清理眼部周边皮肤；如果使用湿性眼贴同时使用本品，建议在本品绒布表面敷层保鲜膜； 2.本品可用牙刷粘少量清水或酒精清洁清洗；清洗后注意放置通风环境下自然晾干。 3.本产品不可长时间放阳光或高温环境。</div></div></div>	<div><div>注意事项</div><div>1.本产品采用锂电池，充电可反复使用；切记不可长时间保持充电状态，会严重影响电池及产品使用寿命。 2.本产品在不使用中，请保持在关机状态。（本品一个热敷疗程为15分钟；15分钟后自动关闭热敷功能；如果要继续使用时，请重新开启热敷）</div><div><div>警告：本品具有发热表面，对热度不敏感的人谨慎使用，避免长时间热敷引起灼伤。</div></div></div>
---	---	--	--	---

<div><div>D13 Heated eye mask</div><div>Model: D13 Working voltage: DC 5V Charging current: 0.9A Working current: 50mA-1100mA Power used: 3.7W Battery capacity: 1200mAh Product weight : 169g Executive standard: QB/T QB 1108.02-2016 Quality inspection: qualified products Place of origin: Made in China</div><div><div>Key function and charging indication</div><div></div></div></div>	<div><div>Functional Operation Description</div><div>1. Long press for 3 seconds to switch on and off the device, Bluetooth mode will be entered by default, and a Bluetooth mode prompt tone will be played. 2. Press the button to turn on the first gear heating mode, and light up the first gear heating lamp /35 degrees; 3. Press the button again to turn on the 2-gear heating mode, and light up the 2-gear heating lamp /40 degrees; 4. Press the button again to turn on the 3 grade heating mode, and light up the 3 grade heating lamp /45 degrees; 5. Press the button again to turn off the heating 6. Double-click the button to enter the vibration massage mode in the first gear, and the massage indicator light is green; 7. Double click the button again to enter the second vibration massage mode, and the massage indicator light is blue. 8. Double click the button again to turn off the vibration massage mode .</div><div><div>This machine supports Bluetooth connection; When you need to listen to music or other audio frequency, connect within 3 minutes after starting the machine. The connection method is the same as the universal Bluetooth connection.The default Bluetooth name is CFJG.</div></div></div>	<div><div>The charging method</div><div>1. When charging, please use a normal charging adapter of 5V/1A or above. 2. Before using the product, please use the charging cable provided by the factory to fully charge the product. 3. When charging, the battery indicator light is bright and blue light is full. 4. The charging time of this product is about 2 hours; Do not charge for a long time, which may affect the service life of the internal mainboard and battery. 5. Low power indicator: The red indicator is steady on.</div><div><div>Maintenance</div><div>1. Clean the skin around the eyes as much as possible before use; If the product is used with wet eye patch, it is recommended to apply the plastic wrap on the surface of the flannelette. 2. This product can be cleaned by toothbrush with a small amount of water or alcohol; After cleaning, place it in a ventilated environment to air dry. 3. The product should not be exposed to sunlight or high temperature for a long time.</div></div></div>	<div><div>Notes</div><div>1. This product adopts lithium battery, which can be used repeatedly after charging; Remember not to keep charging for a long time, it will seriously affect the battery and product life. 2. This product is not in use, please keep in the shutdown state. (A hot compress course of this product is 15 minutes; Automatically turn off the hot compress function after 15 minutes; If you want to continue using, please re-open the hot compress)</div><div><div>Warning: This product has a hot surface. Use with caution for people who are not sensitive to heat. Avoid burning caused by prolonged hot compress.</div></div></div>	<div><div>合格证</div><div>QUALIFIED CERTIFICATE</div><div><div>检验员 CHECKER</div><div>生产日期 DATE</div></div><div></div></div>
---	--	--	--	--

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.