


# 成品尺寸：92×125mm，127克哑粉纸，骑马钉

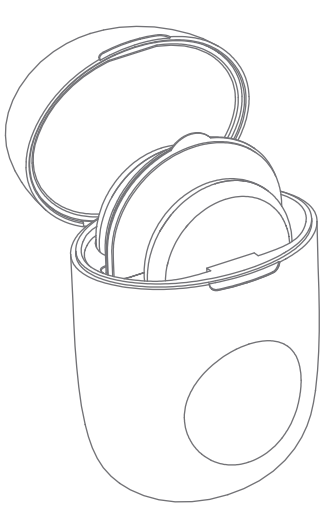

**Mini Pulse massager**  
Please read carefully before using the product

**• Catalogue**

- Safety Warning ----- 01
- Product Accessories ----- 03
- Product Specification And Introduction ----- 04
- Working Indicator Light And Charging Instructions ----- 05
- For Use ----- 06
- Storage And Maintenance ----- 10
- Massage Position And Simple Troubleshooting ----- 11
- Battery Removal And Disposal Methods ----- 12
- Product Warranty Card ----- 14

Please read this manual carefully before using.  
We are not responsible for using this product in the wrong way.


**⚠ Safety Warning**

- Please do not use with the following medical electronic devices to avoid accidents and injuries.
- Pacemakers and other implanted medical electronic devices, Artificial heart and lungs and other life-sustaining medical electronic devices, Electrocardiographs and other life-sustaining medical electronic devices.
- Do not use this product within 1 meter of short wave or microwave therapy equipment.
- Do not use this product while sleeping or bathing.

**Do not use the product for the following groups of people:**

- Infants, pregnant women, women in labor, people with skin allergies, people taking medication, drinking alcohol.
- People taking medication, drinking alcohol, or people with mobility problems.
- People with malignant tumors, heart nerve abnormalities, infectious diseases, acute diseases, people undergoing treatment, abnormalities in the body.
- People who have fever or are in menstruation or express their opinions.
- People who are in breastfeeding period.
- All kinds of pets.

- If you feel uncomfortable when using this product, please stop using it immediately and consulting a doctor.
- When using, please do not make metal substances (such as belts, bracelets, watches, necklaces, etc.) in contact with the pad.
- It is recommended that the same part of the using time not more than 30 minutes, please do not paste the pad on the body for a long time.
- When using, please keep the skin dry, and ensure there is no scar on the position that needs to be massaged.
- Please do not use this product during intense exercise.
- Do not put this product in high temperature, high humidity and dusty environment.
- Do not put the product in your mouth.
- Do not disassemble, repair or modify the product with out permission.
- Do not dispose of the product and waste batteries.
- If you need to move or remove the product, please be sure to turn off the controller.
- The elderly must use the product in the company of others.
- Recommended age: 16 years old.
- May appear phenomenon of rapid heartbeat in the first use, please remove the pad; wait for the heart beat normal and then use again, if still appear this phenomenon, please stopping.
- The product package does not contain a power adapter, charging cable of this product is TYPE-C.
- Charging box and the controller is not waterproof, please do not wash, if the surface is dirty, please use alcohol cotton to wipe it.

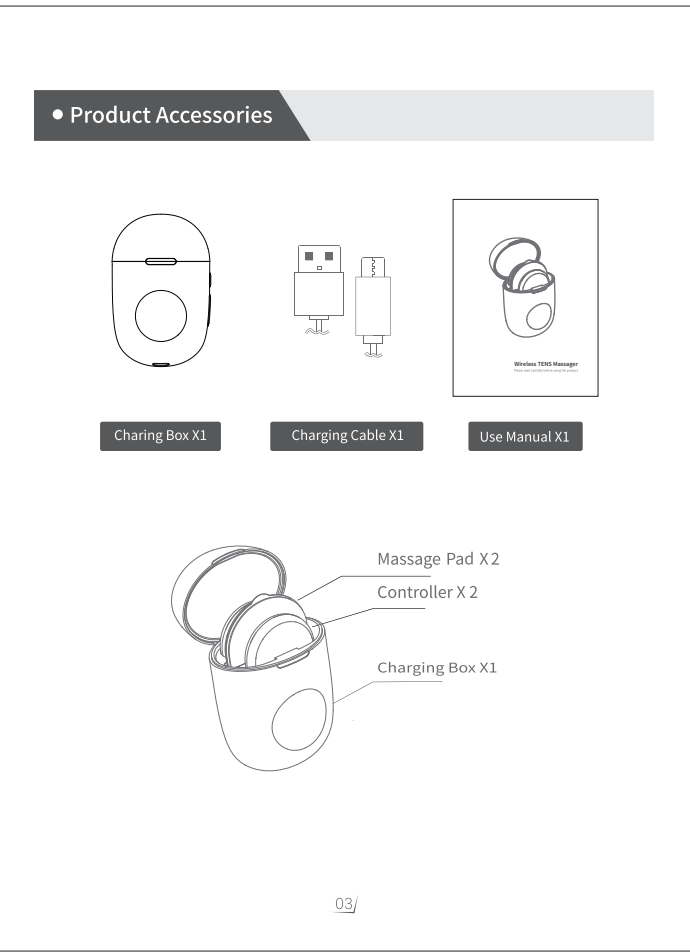


**• Specifications**


Model: U8PRO  
Weight: about 98g  
Controller timing: 15 minutes  
Charging box: Charging voltage: 5V±5%  
Rated power: 1.8W  
Working temperature: 10°C-45°C

Charging box battery capacity: 450mAh  
Controller battery capacity: 60mAh  
Size: 60.80\*46.10\*95mm  
Mainframe: Charging voltage: 5V±5%  
Rated power: 0.3W

**• Introduction of Product Key Diagram**

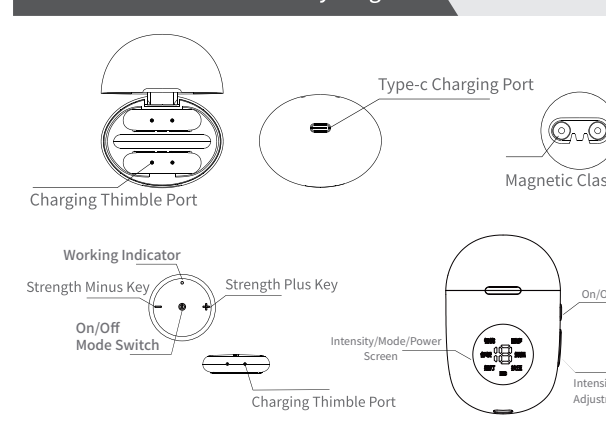

**• Working Indicator**

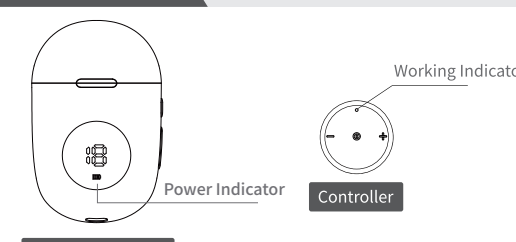


**• Charging Instructions**


- When the controller is connected to the pad, press the on/off button, the red light of the controller is on, start working.
- When the controller and pad separate, the red light slow flashes.
- When the machine capacity is not enough, the red light flashes rapidly, the machine sounds six times before automatically shutdown.
- When the controller is placed into the charging box, the machine sounds one time, the red light becomes brighter and then dimmed gradually, it starts charging. After Fully charged, the red light is always on for 20 seconds and then turns off automatically.
- When the controller is placed into the charging box, the red light lights up, displaying power. When the controller starts to charge, the red light of the charging box automatically turns off after about 3 seconds.
- When the charging box is charging, the red light flashes, and the red light is always on when it is full.

**• For Use**


**① Installation**  
Take out the pad from the charging box, make the pad fit the magnetic clasp of the controller.

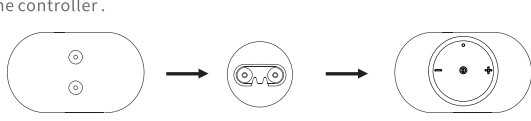


**② Stick**  
Tear off the protective film on the pad (put protective film back into the charging box to prevent loss), stick the pad on the parts that need to be massaged.

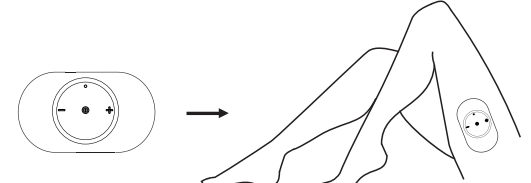


**③ Power**  
Press the On/Off key briefly to turn on the machine, the product will start working.

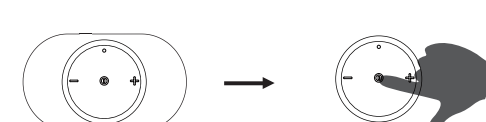


**④ Intensity Level Adjustment**  
Press "+" or "-" key to increase or decrease the intensity, there are 16 levels of intensities adjustment, the machine sounds twice when it reaches maximum level or minimum level, the other sounds once.

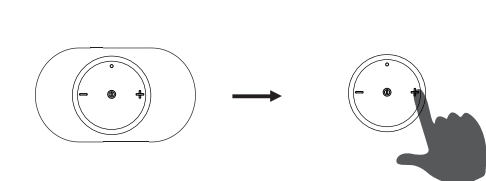


**⑤ Mode Adjustment**  
Short press the "O" key for mode switch, 5 modes can be cyclically switched (switch the first 4 mode the machine will sound once, switch to the last mode will sound twice), each mode switching, the intensity automatically turn to first level.

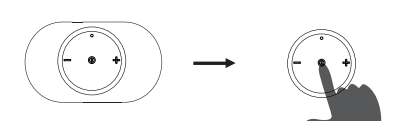


**⑥ Remote adjustment of the charging compartment**  
Long press the "O" key for 2 seconds to switch on, the remote display shows mode 1 intensity 0. The power indicator is always on. Short press the "O" key, 6 modes can be cyclically switched, short press the "+" or "-" key to strengthen or weaken the intensity, a total of 16 levels of intensity, long press the "O" key again for 2 seconds to turn off the machine.

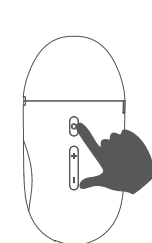

**⑦ Separation**  
Remove the controller when it under operation, the red light of the controller flashes for about 2 minutes after sound twice and then automatically shutdown, remove the controller and then assemble the pad to use again need to press the On/Off key on the controller to re-operation.

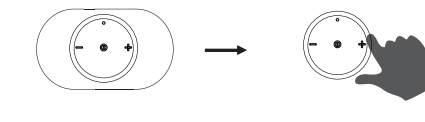

**⑧ End Of The Massage**  
The controller will automatically shutdown after 15 minutes when you are using or you can press the On/Off key to turn off the controller. And then stick protective film on the pad, and put them back to charging box.




**• Storage And Maintenance**

Pad storage: [Diagram] Controller Storage: [Diagram]


Separate the pad from the controller, stick the protective film on the pad and put both into the charging box.

- In order to extend the life of the built-in battery, if not used for a long time, please keep the controller fully charged, and charge it again in 2-3 months.
- If there is dust or stains on the massager surface, wipe it lightly with a wet towel (towel should not have water flowing out)
- If the pad stickiness is reduced, you can use a small amount of water to clean, stickiness will restore after getting dried.
- The pad can be used repeatedly about 30 times. If the pad stickiness is not enough, it can not achieve massage effect and it is easy to feel tingling in massaging. If pad stickiness cannot be restored, need to replace it with a new pad.
- Please do not use or storage the product in the following place:
  - Direct sunlight places
  - Dusty places
  - Places near smoke and fire
  - Strong electromagnetic field
  - High temperature and humidity places
  - Places where water is available

**Best Position For Massage**

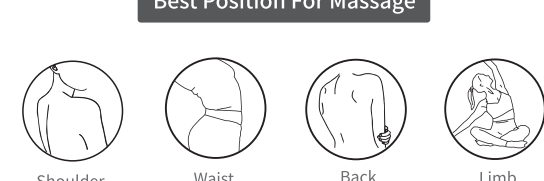

**Prohibited Massage Position**

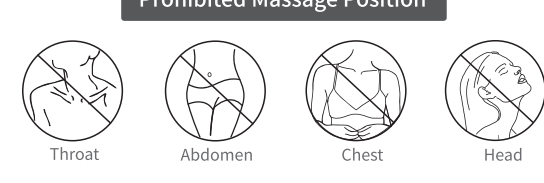


**• Simple Troubleshooting**

Can not be turned on or off	Try to charge the machine, if it still not be fixed, please contact after sale service
The machine cannot be charged	1. Please confirm the output voltage of the using power adapter is DC5V/1A and it is in normal working condition. 2. Replace the charging cable can use TYPE-C phone data cable, if it still not charge, please contact the after-sale service.
There is tingling sensation during use	Please make sure the pad best fits your skin, if the stickiness is not enough, please use water to clean it and use it again after getting dried. If every time you use it, there is tingling sensation in any position, you may belong to the sensitive body, and you are not suitable to this product.

**• Battery Removal And Disposal Methods**


- The built-in battery in this product contains substances that may pollute the environment, so remove the built-in lithium polymer battery before disposing of the product.
- Keep the power off before disassembling the product and removing the lithium polymer battery.
- The detailed steps for removing the lithium polymer battery are as follows. (1) Disassemble and remove the rechargeable battery in the order of ①④;

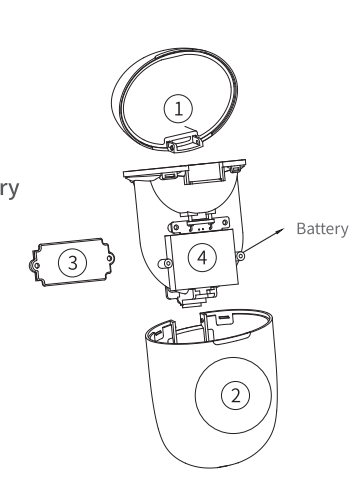

**Charging Box Battery Disassembly Order**

(2) The battery in the controller is disassembled and removed in the order of ①③;

**Controller Battery Disassembly Order**

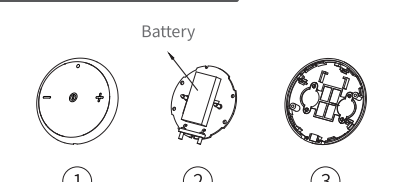

4. The removed lithium polymer batteries should be safely disposed of and discarded in the places specified by the government or recycling processing department

**• Product Warranty Card**

Product Model	Customer Name
Product Name	Contact Number
Date Of Purchase	Purchase Channel
Invoice Number	Contact Address

(Can be cut along here)

---

Product Model	Customer Name
Product Name	Contact Number
Date Of Purchase	Purchase Channel
Invoice Number	Contact Address

**Product Name:** Mini Pulse Massager  
**Product Model:** U8Pro  
**Production Date:** \_\_\_\_\_  
**Inspector:** \_\_\_\_\_  
**Inspection Conclusion:** \_\_\_\_\_

FCC Caution:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception,

which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

-- Reorient or relocate the receiving antenna.

-- Increase the separation between the equipment and receiver.

-- Connect the equipment into an outlet on a circuit different

from that to which the receiver is connected.

-- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

FCC ID: 2A9KD-U8PRO