

备注：封面封底反面留白

105mm

148mm

封底

封面

工艺说明:

图号：EM-724BT

名称：英文说明书

颜色：彩色

厚度：128克

材质：铜版纸

尺寸：210 mmx148 mm 对折 旗马钉装订

SC-EM724BT-IM REV 3

Contents

Before Use

Precautions	2
Scale Components	3
Installing (or Replacing) Batteries	4
Measuring Your Weight	6
Specifications	7
Maintaining the Scale	9
How the Wireless Connection Works	9
Before Requesting Service	10

Precautions

- ❖ Place your scale on a carpet or on a solid level surface, such as tile, wood, or laminate flooring.
- ❖ Do not place your scale on a soft, wet or slippery floor.
- ❖ Your scale contains delicate, high-precision parts. Avoid extremes in temperature, humidity, direct sunlight, shock, moisture and dust.
- ❖ Do not step on your scale with your feet wet. Be sure that your feet are dry to prevent slipping or serious injury.
- ❖ Step on the center portion of your scale gently.
- ❖ Do not step on the display window.

Indications for Use

This scale is designed for personal use.



1

2

Scale Components

Display

Weight Overload (Max + 11 lb)

EEEE

Low Battery Indicator

Lo

Front



Back



- ❖ The scale has an "Automatic power-on and off" function.

3

Installing (or Replacing) Batteries


Caution

- ❖ Insert each battery with its positive (+) and negative (-) terminals properly aligned with the corresponding symbols in the battery compartment.
- ❖ Improper installation of the batteries may cause battery leakage.
- ❖ Replace all four batteries with new ones, when the "Lo" symbol is displayed.
- ❖ Do not mix new and used batteries. It may cause damage to the scale or shorten battery life.
- ❖ The battery life is dependent upon the environmental temperature, battery type and usage.
- ❖ Remove the batteries from the battery compartment, when the scale is not to be used for a long time. They may leak and damage the scale.

4

Steps for Installing (or Replacing) Batteries

1. Remove the battery compartment cover on the underside of the scale by gently pushing the hook in the direction indicated by the arrow and lifting the cover backward.
2. Install four new batteries with the positive (+) and negative (-) terminals properly aligned with the corresponding symbols in the battery compartment.
3. Replace the cover by sliding it into the compartment and gently pressing into place.

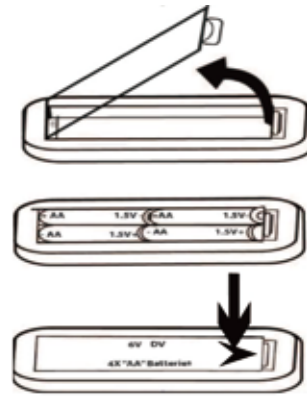

5


Measuring Your Weight

1. Place the scale on a solid level surface.
2. Step on the scale gently, placing both feet on the center portion of the scale.

The scale will sense when you step on it and automatically turn on.

"-----" Four flashing dashes will display until the measurement is taken.

Please stand still during the measurement.

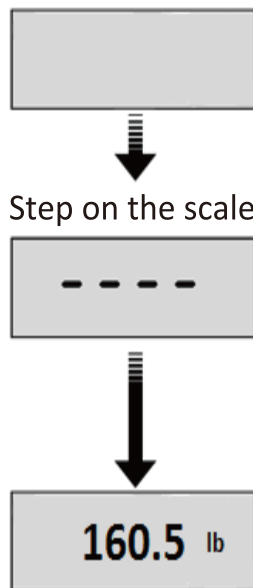

6

3. Your weight is displayed.

Once the reading is taken you will hear 3 loud beeps.

Simultaneously, your weight will be displayed, all at once, in 3 flashing instances.

Then your weight will be displayed for 5 seconds on the LCD.

4. Step off the scale.

The scale will automatically turn off after displaying your weight even if you are still standing on.

Specifications

Digital display	Liquid crystal display, 71.4 X 43.4 mm
Maximum capacity	200 kg / 450 lb
Minimum display	0.1 kg / 0.2 lb
Wireless communication	Bluetooth 4.0

7

Power source

Four type AA batteries (R6P or LR6, 1.5 volt x 4)
Auto on/off
Approximately 1000 measurements
Instruction manual and Batteries (4)

Battery life

Accessories

Compliance with European Directives

CE The scale complies with the statutory EMC (Electromagnetic Compatibility) directive 2004/108/EC.

Notes: Indoor use only. The displayed value may be adversely affected under extreme electromagnetic influences.

Compliance with FCC Rules

The scale complies with part 15 of the FCC rules and contains the FCC ID: 2ABGJEM724BT.

8

Maintaining the Scale

- ❖ Clean the scale with a soft dry cloth or a cloth moistened with water and a mild detergent. Never use alcohol or other harsh chemicals to clean the scale.

- ❖ Remove the batteries from the battery compartment, when the scale is not going to be used for a long time.

How the Wireless Connection Works

At the end of the measurement, your weight will be sent automatically to your tablet. You do not need to do anything.

9

Before Requesting Service

Cause/Problem	Corrective Action
Nothing is displayed by stepping on the scale.	Check the battery installation. Replace with four new batteries.
There is a large difference between estimated weight and data measurement.	Remove any foreign objects that might be touching the bottom of the scale.

If you cannot locate and fix the problem, contact your supplier.

SC-EM724BT-IM REV 3

10

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

11

12

深圳市乐福衡器有限公司

TITLE:

名称

说明书

VER:

版本号

V1.0

Page:

页码

DRAWING BY:

制图

DATE:

日期

CHECKED BY:

审核:

DATE:

日期

APPROVEDB:

批准:

DATE:

日期