




## SAFETY PRECAUTIONS

### CAUTION

 <b>COMPULSARY</b>	<ul style="list-style-type: none"> <li>◆ Please do not continuously use the massage chair over 20 minutes, and it is limited to massage a particular part of body in 5 minutes.</li> <li>◆ Please put the massage chair on flat ground for use.</li> <li>◆ When adjusting the legrest or backrest, please do not exert excessive force on it; in addition, keep children and pets well clear.</li> <li>◆ When removing power plug, please hold the insulation plug, do not pull power cord.</li> <li>◆ Please do not use the power supply that is beyond the specified electric voltage, as this may cause fire disaster or damage massage chair.</li> <li>◆ Before using, please uncover the back cushion to check whether backrest or fabric is damaged; if damaged, please stop using this massage chair immediately, and contact franchiser for repair immediately.</li> <li>◆ Please do not operate the controller with damp hands, and do not put heavy objects on controller.</li> <li>◆ Only people less than 120kgs are acceptable.</li> <li>◆ Please do not use other therapeutic equipments together with this massage chair.</li> <li>◆ During use, if you feel uncomfortable, please stop use immediately and consult your doctor.</li> <li>◆ This massage chair is designed and intended for non commercial household use only, do not use the massage chair for medical treatment.</li> <li>◆ The massage chair is equipped with heating function, so people who are insensitive to heat must use it carefully.</li> <li>◆ Always switch off power and remove the plug from the mains socket immediately after use and before cleaning.</li> <li>◆ Please do not smoke while using the massage chair.</li> </ul>
 <b>PROHIBITED</b>	<ul style="list-style-type: none"> <li>◆ To avoid causing discomfort, do not use the massage chair immediately after meals.</li> <li>◆ Do not fall into sleep on massage chair, when the massage chair is working.</li> <li>◆ Do not use the massage chair when you are drunk or feel bad.</li> <li>◆ If the upholstery becomes dirty, please use only known good quality leather cleaner, do not use chemical cleaners, petrol or detergents etc, and do not spray pesticide on it.</li> <li>◆ Do not sit on backrest, legrest, or armrest to avoid accidental injury or machine fault.</li> <li>◆ In case that PU cover, like back cushion or cloth cover is broken, please stop using the chair at once, and contact the distributor or repair man for help.</li> </ul>
 <b>REMOVE PLUG</b>	<ul style="list-style-type: none"> <li>◆ Before maintenance, please pull out power plug. In addition, please do not insert or unplug the power plug with damp hands, to avoid electric shock.</li> <li>◆ When the massage chair is not in use, please switch off power and remove the plug from socket. When power is off, please cut off power supply immediately and pull out the power plug from socket.</li> </ul>

## TROUBLE SHOOTING

No.	Problem	Possible Cause	Solution
1	The sound “Xu Xu”, “Ta Ta” is heard.	It is a normal sound made by the working motor or the sound is produced when the massage ball is rubbed against the cloth cover.	It is normal and no need for treatment.
2	The machine can not be started up.	The power supply plug is detached.	Insert the power supply plug into the socket again.
		The power supply switch is not turned on.	Turn on the switch
		The fuse of massage chair blew out.	Replace a new one with the same mode
3	No sound from the speaker or	The volume is turned to lowest.	Turn the volume up
4	Backrest or legrest can not ascend or descend.	Over pressure is put on backrest or legrest.	Reduce the pressure on backrest or legrest.
5	The airbag installed at seat cushion or legrest can not work.	To check whether there is air tube blockage.	To clear the air tube

If abnormal phenomenon still remains after above measures have been taken to the massage chair, please contact local distributor for repair and maintenance.

**Note:** Do not dismantle the product by yourself, and our company shall shoulder no responsibility to the damage of the product or the injury of people caused by improper use! Thanks your cooperation.

## ▶ PRODUCT STORAGE AND MAINTENANCE

### Maintenance

- ▶ Please twist the cotton cloth which has been dipped in neutral detergent before wiping the back cover, hand rest, legrest, air tube and some other plastic components, then use a dry cotton cloth to wipe them again.
- ▶ Use dry cloth to wipe the controller, DVD, and power supply box.
- ▶ Use a piece of cloth that has been dipped with neutral detergent to wipe the back cushion and pillow, then use a piece of wet cloth to wipe them, and finally put them to a well-ventilated place for natural drying.

### ⚠ Caution

- ▶ Make sure to pull out the power supply plug firstly before maintenance. Do not use wet hands to pull out the plug to avoid electric shock. When clean the back cover, please be cautious of detergent entering into the power supply box!
- ▶ Do not iron the upholstery.

### Storage method

- ▶ Wipe off dust and dirt from the product after use, if the product is not to be used for a long period of time, please cover the product with a dust cover (self provided).

### ⚠ Caution

- ▶ Do not store the product under direct sunlight or in places of high temperature. This may cause the color of the upholstery to fade.

### Kind reminding

During daily use, please check whether the following symptoms exist:


- ▶ Scorched smell exists.
- ▶ When touch the wire, the power is on, but sometimes off.
- ▶ The power supply cord is hot.
- ▶ Other abnormal symptom

### ⚠ Caution


- ▶ Please stop using the massage chair to avoid breakdowns and accidents occurring if any above symptom exists.
- ▶ Turn off the power supply switch and pull out the power supply plug, and then contact the local distributors or headquarters for repair and maintenance.
- ▶ Except all the normal operation and maintenance mentioned above, any other repair and maintenance can only be operated by our authorized dealer or service agent.

## ▶ NAME AND FUNCTION OF EACH COMPONENT

### External structure

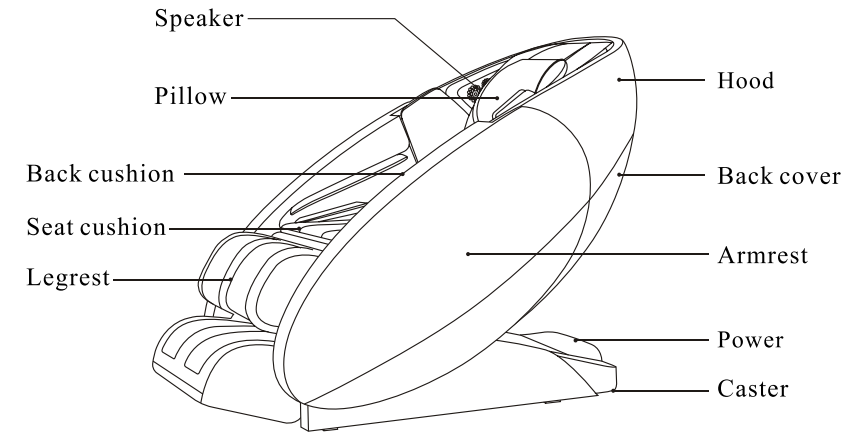


### Internal structure

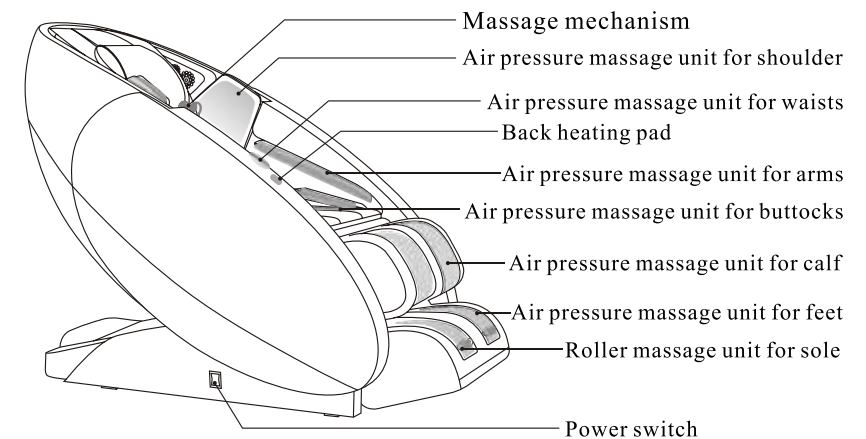

## NAME AND FUNCTION OF EACH COMPONENT


### Controller

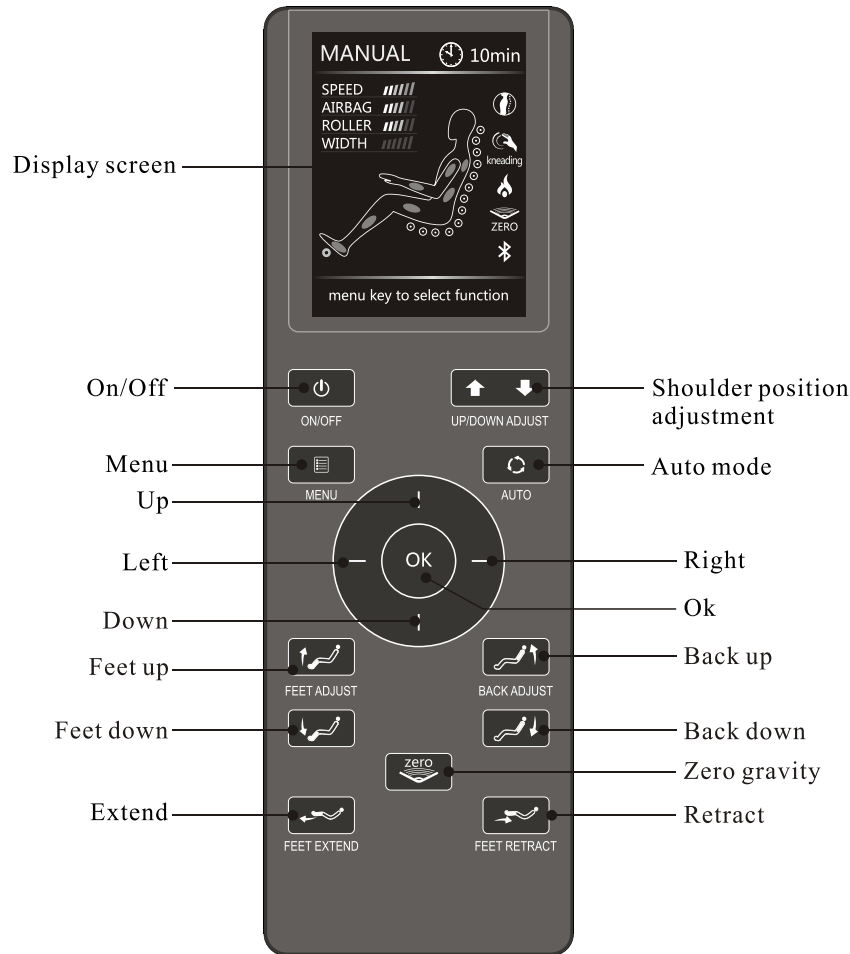

## PREPARATION BEFORE USING

### Move method

- ▶ There are moving rollers at the bottom of the chair.
- ▶ Raise the back to the highest point.
- ▶ Switch off and plug out the chair.
- ▶ Lift the front part of the chair as showed in the figure (one person raise the front part via holding moving holder, another person push the chair when holding the top of the chair), move the chair by the rollers.

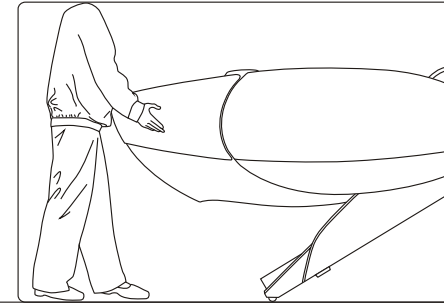


### ⚠ Caution

- ▶ Don't move it when there is people in it.
- ▶ Don't use the armrest while moving the chair, only hold the moving holder in the back and front

### Turn on power supply

- ▶ Plug the power cord into three core socket with well grounding.
- ▶ Switch on the power button locates at the bottom of the right (sitting direction) armrest.

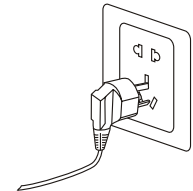


### ⚠ Caution

- ▶ Before switching on the chair, please make sure there is no damage on the power wire and other exposed wires.
- ▶ Make sure the power switch is off before the chair is connected to electricity.

## ▶ PREPARATION BEFORE USING

### Massage chair commissioning

- ▶ Before commissioning, please make sure there is enough space around the chair. Then use the controller to recline down the chair backrest to the limit position, rise up the legrest to the highest position and stretch out to the limit length.
- ▶ Testing all functions work normally, then power off the chair, return back to original position.  
(How to select the program, please refer to “ Controller instruction” and “display screen menu instruction” )

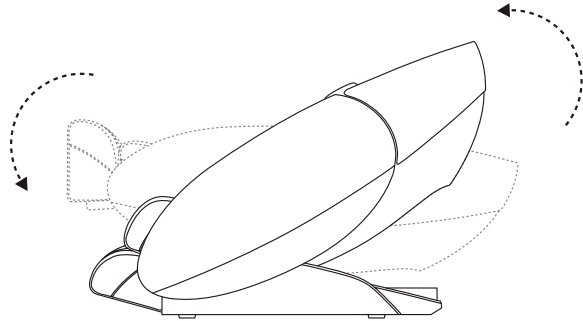

### ⚠ Caution









- ▶ If the legrest rises up, please do not seat on the chair.
- ▶ Do not make the legrest suffer from heavy loading suddenly, or it may damage the legrest mechanism.

### After using

- ▶ Please power off the chair and unplug the chair power cable from the power supply socket after using well the massage chair.
- ▶ Put well the controller inside of the pocket.
- ▶ Please power off the chair and unplug the chair from the power supply socket when the chair is not in use, to avoid children to start the chair accidentally.
- ▶ Please store well the chair if it is not in use for a long time.

## ▶ NAME AND FUNCTION OF EACH COMPONENT

### Instruction of Controller

-  **On/Off**
  - ▶ When power is on, click this button, display screen is lightening, then click up/down button to choose the auto mode, manual mode, air pressure massage, heat, roller massage, and set menu.
  - ▶ Click this button again to turn off all functions. The default working time is 20 minutes
-  **Menu**
  - ▶ Click this button to choose auto mode, manual mode, air pressure massage, heat, roller massage, and set menu.
-  **Shoulder position adjustment**
  - ▶ Under manual and point mode, press UP/DOWN to adjust the position of point massage.
-  **Auto mode**
  - ▶ Click this button to choose one of auto modes, such as “RECOVERY, EXTEND, RELAX, REFRESH, UPPER BODY AUTO, and LOWER BODY AUTO”
-  **Direction button**
  - ▶ In the menu and auto mode, click up, down, left, right button to choose the concerned function, then click “OK” to confirm.
-  **Zero gravity**
  - ▶ Click this button to turn On/Off zero gravity, then relevant icon in screen will be lightening. After adjusting the backrest and legrest position manually, zero gravity will be off.
-  **Feet up**
  - ▶ Press this button, the legrest will be moved up to the top position, releasing this button, it will stop.
-  **Feet down**
  - ▶ Press this button, the legrest will move down to the lowest position, releasing this button, it will stop.

## ▶ NAME AND FUNCTION OF EACH COMPONENT



### Back up

▶ Press this button, the backrest will be moved up to the original position, release this button, it will stop.



### Back down

▶ Press this button, the Backrest will be moved down to the lowest position, release this button, it will stop.



### Legrest extend

▶ Press this button, the legrest will be extended to the longest position, releasing this button, it will stop.



### Legrest retract

▶ Press this button, the legrest will be retracted to the original position, releasing this button, it will stop.




### WARNING

▶ Press power button after use it to make the chair in reset state. Please press any buttons in the remote controller to stop reset. Press the power button again to continue the operation or enter into massage state.

## Instruction of Display screen menu

### ▶ Auto mode

When power is on, click “ON/OFF” button, screen is lightening, then Click the “AUTO” button to start one of six auto modes, such as “RECOVERY, EXTEND, RELAX, REFRESH, UPPER BODY AUTO, and LOWER BODY AUTO”

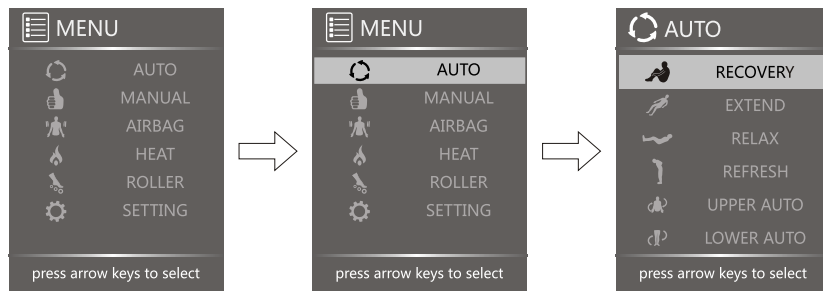

### Recovery

For someone after one day work and dinner, to relieve tiredness and recover vigour.


### Extend

For officers and drivers who sit for long time, to stretch full body and recover.

## ▶ PREPARATION BEFORE USING

### Floor protection

▶ Please put carpet or other pad under the massage chair, in order to avoid any unpredictable impression on the floor.

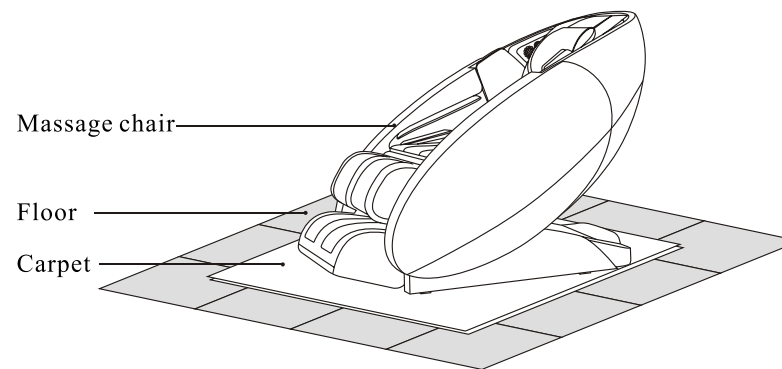

### Surrounding environment



Keep the massage chair away from damp places (for example: bathroom) to avoid electric leakage.



Keep the massage chair away from high temperature sources (for example: stove) to avoid fire or damage to leather material.

### Grounding


This massage chair belongs to Level One electrical appliance, make sure to use three-core power supply socket and to guarantee the ground wire connected to the power supply socket has obtained good ground connection to avoid electric leakage, electric shock and some other negative effect during use.

## ► PREPARATION BEFORE USING

### Installation position

Safety distance for the periphery of massage chair.

- Please make sure there is adequate space for installation this product:  
Back: minimum 50cm from the wall or any other blocks;  
Front: minimum 30cm from any blocks;
- In order to prevent from signal interference, please keep at least 1m away from the TV, radio and other video and audio devices.

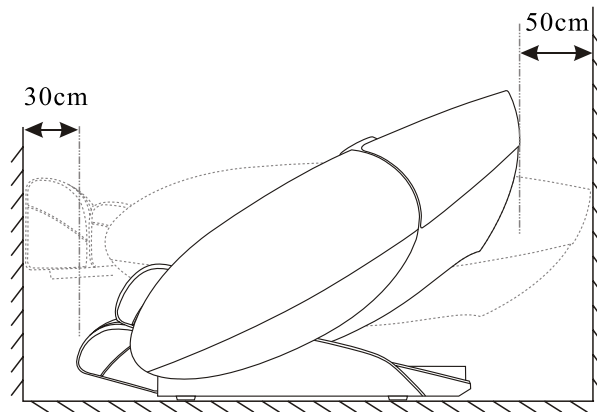

### ⚠ Caution

- Please check the power cord/cable in fine situation before using this masage chair to avoid any possible damage.
- Please confirm NO pet(s) or kid(s) staying under the chair backrest and footrest to avoid any possible accident.

## ► NAME AND FUNCTION OF EACH COMPONENT

### Relax

For early morning or noon leisure massage, to generate muscle energy.

### Refresh


For relieving the pains in neck, back and waist, combining with semi-auto mode, it can ease the pain effectively with focusing massage.

### Upper/Lower auto

For massaging different part of the body.


### ► Manual mode

Click “MENU” , choosing “MANUAL” to start manual program, then click up/down to choose: mode, part, width and speed.

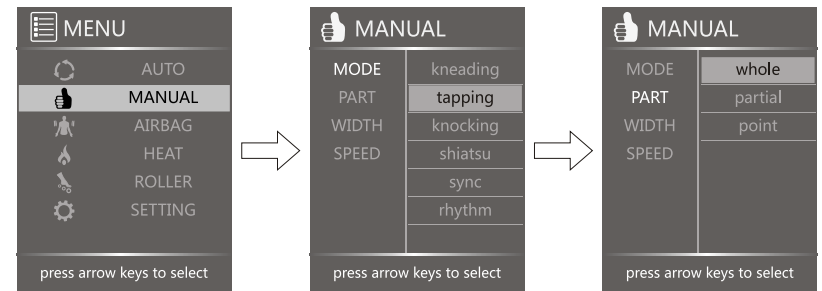

**Mode:** Click right button, then click up and down to circularly choose kneading, tapping, knocking, shiatsu, sync and rhythm massage.

**Part:** Click right button, then click up and down to circularly choose full, partial and point massage.

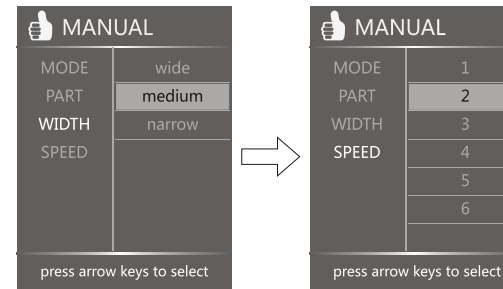


**Width:** Click right button, then click up and down to circularly choose wide, medium, narrow to adjust width of massage.

**Speed:** Click right button, then click up and down to circularly choose 1-6 levels of massage speed.

## ▶ NAME AND FUNCTION OF EACH COMPONENT

### ▶ Airbag

Click “MENU”, choosing “AIRBAG”, then click up/down key to choose: Part and Force.

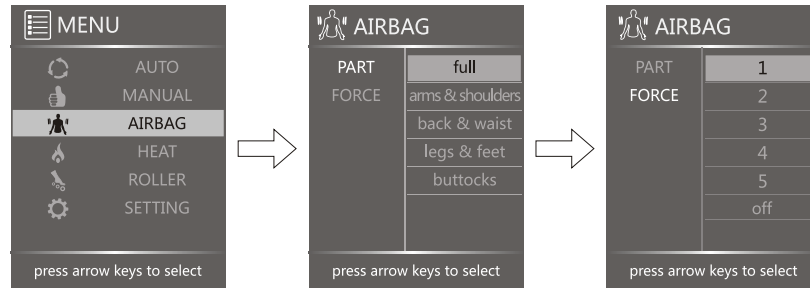


**Part:** Click right button, then click up and down to circularly choose full, arms & shoulders, back & waist, legs & feet, and buttocks

**Force:** Click right button, then click up and down to circularly choose 1-5 levels air strength.



### ▶ Heat

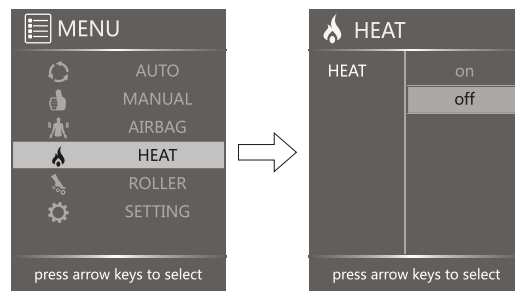
Click “MENU”, choosing “HEAT”, then click right to choose on or off.



## ▶ PREPARATION BEFORE USING

### Legrest Extend and Retract

▶ When the legrest in original state, in order to satisfy the most suitable massage position, press  or  to manually adjust the relative position between the footrest frame and leg.

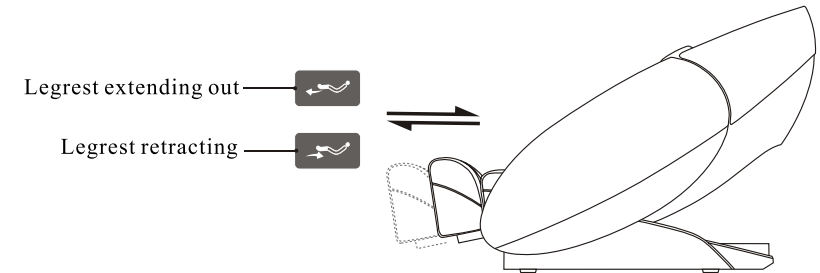

▶ Legrest can be extended or retracted automatically or manually according to different position and height, because of sensor inside

### WARNING

- ▶ Do not play in the chair and conduct improper operation.
- ▶ Before adjusting the angle between calf and backrest parts you need to make sure no obstacles is around.
- ▶ Do not insert leg, foot or other goods into chair's clearance when leg parts extend or reset.
- ▶ Do not get down from the chair when the chair is in the reset state.


### Blue tooth

▶ This masage chair is built-in Bluetooth module and speakers, we can just search it (the paring name is Novo) via a smartphone/tablet, then pair it with your smartphone/tablet to enjoy the musics.

## ► PREPARATION BEFORE USING

### Step 3 Other Accessories


- Plug-in the power cord to the chair power carbinet at back;
- Press the power button to swith on the chair,sit in the chair,find the chair remote control to start and enjoy the massage (attention please:if the chair is just unpack from the carton,we suggest you to use the chair remote to lift the chair backrest to upright position so you can sit in the chair)

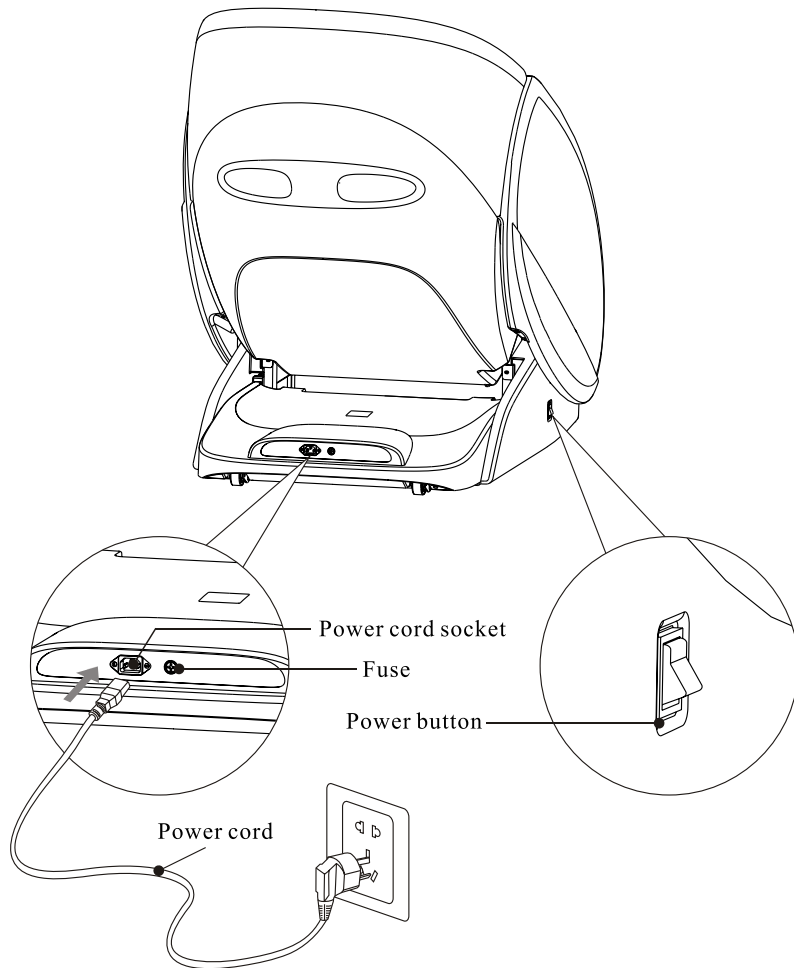

## ► NAME AND FUNCTION OF EACH COMPONENT

### ► Roller

Click “MENU” , choosing “ROLLER” , then click up and down to circularly choose 1-3 levels of speed and off.

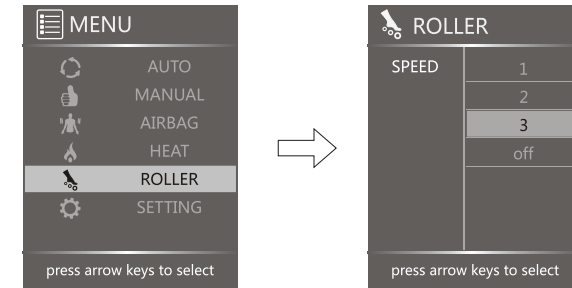


### ⚠ Caution

- Under the air massage of full body, the roller speed could not be adjusted.

### ► Set Menu

Press “Menu” button,select “Setting” to enter into Setting menu. Press “Up” or “Down” button to choose Time and Bluetooth setting.

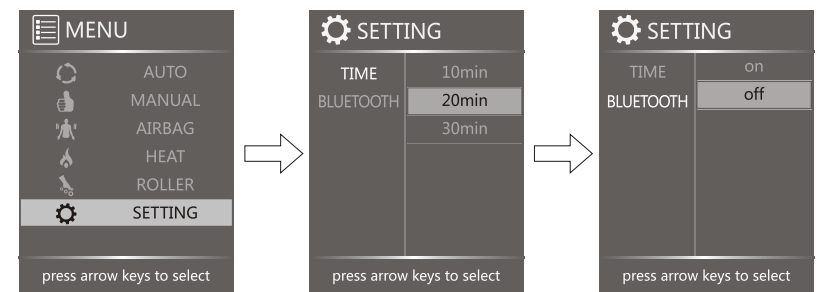

**Time:** Setting working time, firstly, press “Right” button, and then press the “Up” or “Down” button to adjust 10, 20, 30 minutes, and the default working time is 20 minutes.


**Bluetooth:** Setting “Bluetooth” , press the “Up” or “Down” button to open or close this function.




## NAME AND FUNCTION OF EACH COMPONENT

### Note:

- Under auto mode or manual mode, it will go back to home page as below drawing 10 seconds after without any operating.

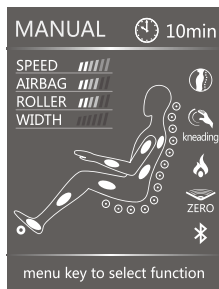

Home page in manual mode

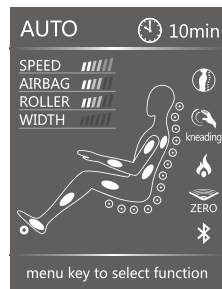


Home page in auto mode

### ⚠ Caution

- ▶ During working, click “MENU” twice, it will go back to home page immediately.
- ▶ In the home page, under manual mode or auto mode, click the left and right button to adjust air pressure intensity; under manual mode, click up and down, to adjust massage wheels speed.
- ▶ Click “AUTO” to choose auto program directly.

- Before running auto program, it will inspect shoulder position first. Except “Lower Auto” program

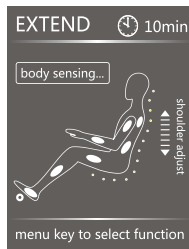

### ⚠ Caution

- ▶ Please do not disassemble the product to avoid any defective problems.
- ▶ Please do not split water and other liquid into the controller to avoid malfunction of keys and error code.
- ▶ Please do not put heavy things on the controller or power cable.

## PREPARATION BEFORE USING


### Installation Method

#### Step 1: Get out each part of the chair

- ▶ Unpack the cartons, carefully get out sewing units, accessories box and so on.
- ▶ Check all accessories with packing list.

#### Step 2 Installation of controller unit

- ▶ Connect the data line of the controller to the port on the inner side of right armrest (pay attention to the convex-concave between connector and port); then fix the screw cap. (see the figure below)

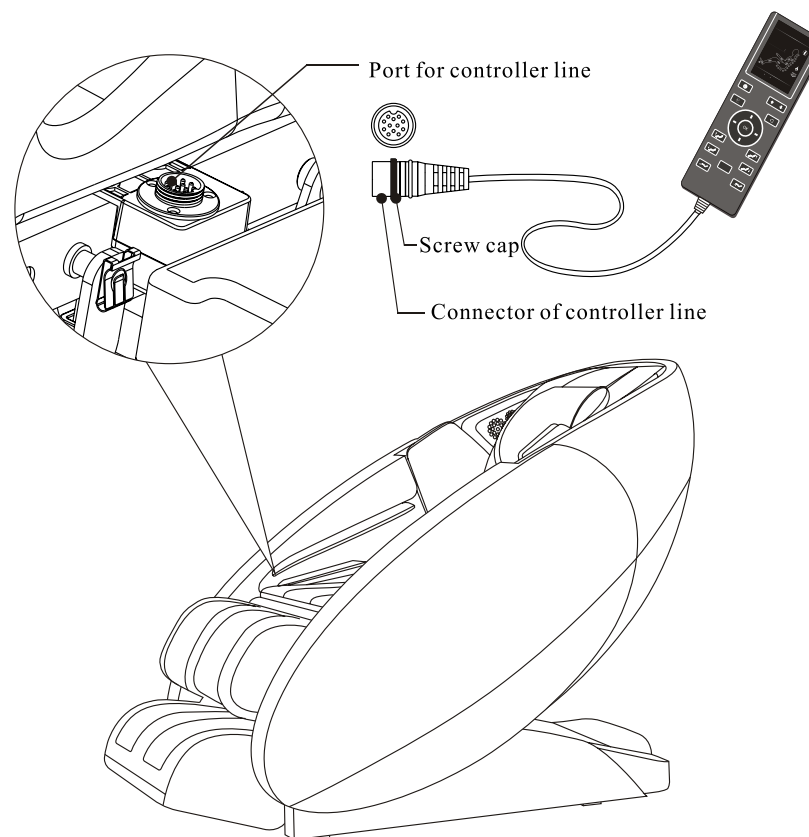

## **FCC WARNING**

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE 1: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.