PHILIPS


DiamondClean

9000


1

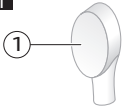2

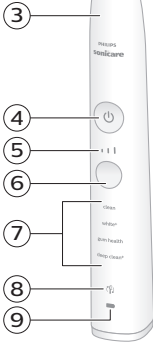4

5

6

7

8

9



10

11

12

13

14

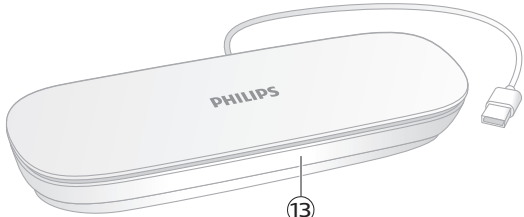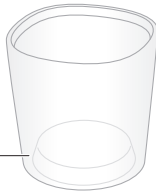— |

| —

— |

| —

English 6
Español 35
Français (Canada) 66

Introduction

Congratulations on your new Philips Sonicare power toothbrush! Superior plaque removal, whiter teeth and healthier gums are at your fingertips. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time. Further support and product registration are available to you at: <http://www.philips.com/welcome>.

Welcome to the Philips Sonicare family!

Intended Use

DiamondClean series power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. DiamondClean series power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGERS

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not immerse the charger in water or any other liquid.
- After cleaning, make sure the charger is completely dry before you connect it to a power outlet.
- Do not reach for a product that has fallen into water. Unplug immediately.

WARNINGS

To reduce the risk of burns, electrocution, fire or physical injury:

- Do not clean the brush head, the handle or the charger in the dishwasher.
- Do not use attachments other than those recommended by Philips Sonicare.
- Never use a charger with a damaged cord or plug.
- Always replace the charger with one of the original type in order to avoid a hazard.
- Do not plug into an outlet with a voltage other than specified on the bottom of the product. This product is designed to operate within a range of 100 to 240 volts. Voltage converters do not warrant voltage compatibility.
- Never force the plug into an outlet; if the plug does not easily fit into the outlet, discontinue use.
- Discontinue use if the product appears damaged in any way (brush head, handle or charger, etc.).

8 English

- The power cord cannot be replaced. If the power cord is damaged, discard the charger.
- Keep the power cord away from heated surfaces.
- This product contains no user-serviceable parts. Refer to 'Warranty and support' if the product no longer works properly or needs repair.
- Do not use the charger outdoors.
- Do not use the charger if dropped into water
- This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use as described in this booklet. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the product.
- Philips Sonicare is a personal care device and is not intended for use on multiple patients in a dental office or institution.
- This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see chapter 'Warranty and support').
- For connection to a supply not in the U.S.A., use a charger of the proper configuration for the power outlet.
- Use a UL-listed computer, USB Hub, or USB duplex receptacle when charging the handle in the travel case.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristle tufts dislodging.

- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle with soap and water after each use. Otherwise plastic cracking may occur.

MEDICAL WARNINGS

- 1 Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- 2 Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- 3 Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- 4 This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
- 5 Consult your physician if you have other medical concerns.

SAVE THESE INSTRUCTIONS

Your Philips Sonicare (Fig. 1)

- 1 Hygienic travel cap
- 2 Smart brush heads*
- 3 Handle
- 4 Power on/off button
- 5 Intensity light
- 6 Mode/intensity button
- 7 Modes
- 8 Brush head replacement reminder light
- 9 Battery indication

10 English

- 10 Charging glass
- 11 Charging stand
- 12 Charging base
- 13 (Charging)Travel case*
- 14 USB wall adapter*

* Note: The content of the box may vary based on the model purchased.

Sonicare app – Getting Started

The Sonicare app pairs with your toothbrush to provide you a connected experience. By connecting your toothbrush to your Sonicare app account, you will be able to:

- Easily review your progress report for ways to improve your brushing routine
- Receive personalized tips and recommendations for managing your oral health

To start using the Sonicare app:

- 1** Download the Sonicare app to your phone
- 2** Open the app and follow the guided steps
- 3** Pair your toothbrush with the app
- 4** Create your account
- 5** Complete firmware update, if prompted, to access the latest improvements and features.
- 6** Brush regularly

You are ready to start your connected experience. Be sure to brush regularly. You can receive updates to help improve your oral healthcare.

Note: Make sure your phone's Bluetooth is turned on when using the app. Your toothbrush uses the Bluetooth connection to transfer your brushing data to the app. If you have questions about why your brushing data is collected, be sure to review the

Sonicare Privacy Statement, available throughout the app setup process.

Sonicare app – Track and improve your oral health care

Your Sonicare power toothbrush is equipped with smart sensors to provide you with feedback about your brushing habits, including:

- Daily brushing sessions
- Time spent during each brushing session
- Applied brushing pressure
- Brush head replacement reminder (based on actual use)

The Sonicare app collects the data provided by the smart sensors and tracks your brushing habits over time.

Each of the features available in the app are designed to help you understand your brushing habits and provide you with simple recommendations for improving and maintaining your oral health. You will receive personalized brushing information in the Sonicare app.

Brush heads

To see what type of brush head you have, check the label at the bottom of the brush head. The symbol at the bottom of the brush head indicates that it is a smart brush head, equipped with the BrushSync technology.

12 English


| | | | |
|-------------------------|------------------------|--------------------------|------------------|
| Brush head name | Premium Plaque Control | Premium White | Premium Gum Care |
| Brush head label | C3 | W3 | G3 |
| Benefit | Plaque removal | Plaque and stain removal | Healthier gums |

To explore our full range of brush heads, go to www.philips.com/toothbrush-heads for more information.

Using your Philips Sonicare

Brushing instructions

- 1 Push the brush head firmly onto the handle. Wet the bristles and apply toothpaste.

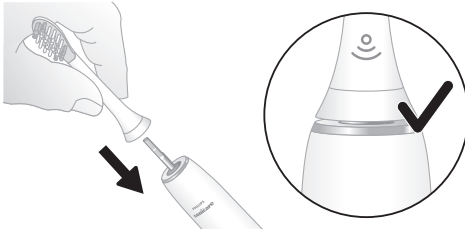
Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.


- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), apply light pressure to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush head in contact with the teeth at all times.

- 3 Press the power on/off button to turn on the Philips Sonicare.

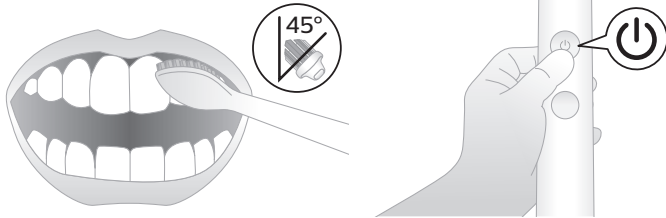
- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.

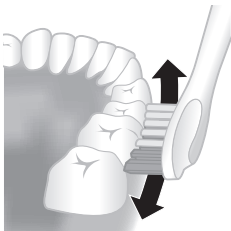Note: Different modes are available (see 'Brushing modes').

- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

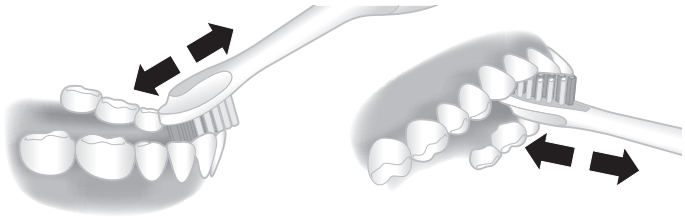Note:

- A change in vibration from the handle alerts you when you apply too much pressure (see 'Pressure sensor').
- Your Philips Sonicare toothbrush is safe to use on braces, dental modifications and dental restorations (fillings, crowns, veneers)
- Do not use your power toothbrush in the shower.

BrushSync technology

BrushSync technology enables your brush heads to communicate with your handle using a microchip. The symbol at the bottom of the brush head indicates that it is equipped with the BrushSync technology. This enables the following features:


- Brush head mode pairing
- Brush head replacement reminder

Brushing modes

Your power toothbrush is equipped with 4 different modes: clean, white+, gum health and deep clean+.

Based on the smart brush head you use, the optimal mode and intensity will be selected automatically by the BrushSync technology.

To manually change the mode, press the mode/intensity button while the toothbrush is off with the brush head attached.

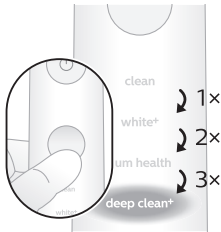*** Note:**

- All brush heads work with each mode.
- When the Philips Sonicare toothbrush is used in clinical studies, it must be used at high intensity, with the EasyStart feature turned off.

Brush head mode pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates to the handle and automatically pairs it with the recommended mode and intensity.







Note: If you change modes/intensity, the brush head will remember the last used setting. If you wish to brush with the recommended setting, you can find the recommended mode/intensity in 'Brushing Modes'.

Brush head mode pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates to the handle and automatically pairs it with the recommended mode and intensity.

Note: If you change modes/intensity, the brush head will remember the last used setting. If you wish to brush with the recommended setting, you can find the recommended mode/intensity in the following table.

| | clean and deep clean+ | white+ | gum health |
|------------------------------------|---|-----------------------------------|--------------------------------|
| Benefit | Plaque removal | Plaque and surface stains removal | Plaque removal and gum massage |
| Recommended brush head | Premium Plaque Control (C3) | Premium White (W3) | Premium Gum Care (G3) |
| Recommended intensity level | 3 | 3 | 3 |
| Total brushing time | clean: 2 minutes deep clean+: 3 minutes* | 2 minutes and 40 seconds | 3 minutes and 20 seconds |
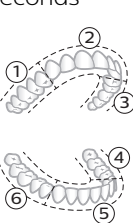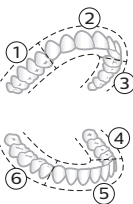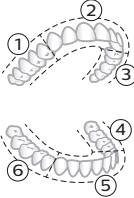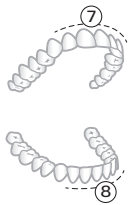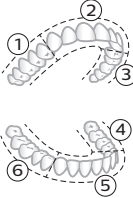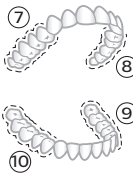
| | clean and deep clean+ | white+ | gum health |
|---------------------|--|---|--|
| How to brush | <p>clean: Brush each segment for 20 seconds</p>  <p>deep clean+: Brush each segment for 30 seconds*</p>  | <p>Brush each segment for 20 seconds,</p>  <p>then brush segments 7 and 8 for 20 seconds each.</p>  | <p>Brush each segment for 20 seconds,</p>  <p>then brush segments 7, 8, 9 and 10 for 20 seconds each.</p>  |

*Note: When brushing with the Sonicare app in 'deep clean+' mode, the brushing time will be 2 minutes.

Intensity settings

Your power toothbrush comes with 3 different intensity settings:


- low intensity (one bar)

18 English

- medium intensity (two bars)
- high intensity (three bars)

When attached to the handle, the smart brush head will automatically select the recommended intensity. To manually change the intensity, turn on the handle and press the mode/intensity button.

While toothbrush is ON, press this button to change the intensity setting

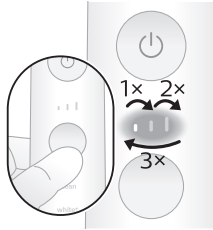Features

- BrushPacer
- Pressure sensor
- Brush head replacement reminder
- EasyStart

BrushPacer

The BrushPacer features indicates when you have to move to the next section of your teeth with a brief change in vibration. The toothbrush turns off automatically at the end of the brushing session.

Pressure sensor

Your Philips Sonicare measures the pressure you apply while brushing to protect your gums and teeth from damage.

If you apply too much pressure, the handle will change its vibration until you reduce the pressure.

Note: The Pressure Sensor comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

Brush head replacement reminder

Over time, the handle will track the wear of your smart brush head by measuring:

- The overall pressure you apply while brushing
- The total time you have brushed with your brush head

Once your brush head is no longer effective, the brush head replacement reminder will blink amber and the handle will make a series of beeps and tones.

Note: The brush head replacement reminder comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

EasyStart

The EasyStart feature gently increases the power of the handle over the first 14 brushings to help you get used to the brushing power of the Philips Sonicare.


The EasyStart feature works with all modes.

Note: The EasyStart feature comes deactivated with your product. To activate EasyStart, see 'Activating or deactivating features'.


Charging and battery status

Charging on charger base

- 1 Plug the cord of the charging base into an electrical outlet.
- 2 Place the charging stand or the charging glass on the charging base.

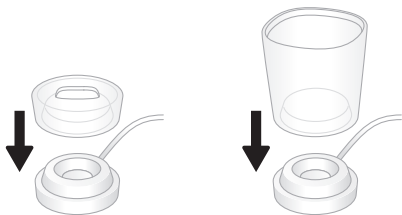- 3 Place the handle on the charging stand or in the charging glass.

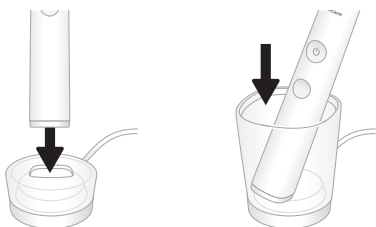- 4 The charger will charge your toothbrush until it is fully charged.

The handle is charging when the battery indication is blinking white. When the handle is fully charged, the battery indication illuminates white for 30 seconds and then turns off.

Note: Your handle comes pre-charged for first use. After first use, charge for at least 24 hours.

Charging with the travel case*

*Note: The charging travel case may not be included depending on the model purchased.

- 1 Take the USB cable from the bottom of the charging travel case.
- 2 Connect the charging travel case to 5V DC input power either from:
 - The USB adapter provided
 - a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
 - a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
 - a UL-listed or certified 1310 Class 2 power supply.

Important: The wall adapter must be a safe extra low voltage source, compatible with the electrical ratings marked on the bottom of the device.

- 3 The travel case charges your toothbrush until your toothbrush is fully charged.

Battery Status (handle on plugged in charger)

When the handle is placed on the charger or in the travel case, the battery indication will communicate the battery level.

Battery status

Battery indication

22 English

| | |
|----------|--|
| Charging | Blinks in white |
| Full | Illuminates white for 30 seconds and then turns off. |

Battery Status (when handle is not placed on Charger)

When you turn on the handle, the battery light at the bottom of the handle will indicate the status of the battery.

| Battery status | Battery indication | Sound |
|----------------|----------------------|------------------------------------|
| Full | Illuminated in white | - |
| Partially full | Blinks white | - |
| Low | Blinks amber | The handle emits 3 beeps |
| Empty | Blinks amber | The handle emits 2 sets of 5 beeps |

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder
- Pressure sensor

Step 1: Place the handle on the charging stand or in the charging glass.

Step 2: Press and hold power button for:

| | | |
|------------|---------------------------------|-----------------|
| Easy start | Brush head replacement reminder | Pressure sensor |
|------------|---------------------------------|-----------------|

| | | |
|-----------|-----------------|-----------------|
| 3 seconds | Up to 5 seconds | Up to 7 seconds |
|-----------|-----------------|-----------------|

Step 3: Let go power button when you hear

| | | |
|--------|-------------------------|----------------------------------|
| 1 beep | 1 beep and then 2 beeps | 1 beep, 2 beeps and then 3 beeps |
|--------|-------------------------|----------------------------------|

If you see the battery indication blink white 2 times and hear 3 tones from low to high, then the feature has been activated.

OR


If you see battery indication blink white once and hear 3 tones high to low, then the feature has been deactivated.

Cleaning


- All parts can be cleaned with a damp cloth or rinsed with warm water.


- Remove the brush head from the handle and rinse it thoroughly.

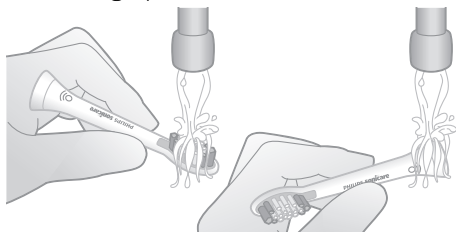- Rinse the entire handle, especially the brush head connection, at least once a week.

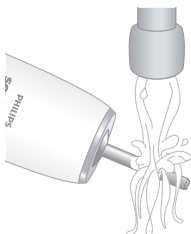Caution: Do not push on the rubber seal at the top of the handle. Gently clean around the rubber seal.

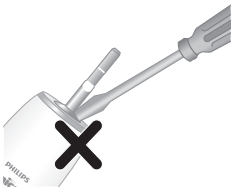Cautions:

- Do not clean the brush head, travel cap, handle, charger or travel case in the dishwasher.
- Do not use isopropyl alcohol, vinegar, bleach, or any other household cleaning products, to clean the handle as this may cause discoloration.
- Unplug the charger before you clean it.
- Do not use essential oils to clean the brush heads as this can cause damage.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Frequently asked questions

| Question | Answer |
|---|---|
| Can I leave my Sonicare on the charger between brushings? | Yes, you can always keep your Philips Sonicare toothbrush on a plugged-in charger. This does not affect the battery lifetime. |

Question

Answer

Will other brush heads fit my Philips Sonicare toothbrush?

All Philips Sonicare electric toothbrush brush heads that snap on will fit your toothbrush. However, not all are equipped with BrushSync technology. You can identify the brush heads with BrushSync technology by the icon at the base of the brush head.

How often do I need to replace my brush head?

We advise you to replace your Philips Sonicare brush head after 3 months of use (based on brushing twice a day for 2 minutes each session). The brush head replacement reminder will also notify you when it is time to replace your brush head.

Is my Philips Sonicare toothbrush waterproof?

Yes, your Philips Sonicare toothbrush is waterproof. However, never put your handle in a dishwasher to clean it. And do not use it in the shower.

Question**Answer**

What phones are compatible with the Philips Sonicare app?

The Philips Sonicare app support Android OS, and Apple iOS. For the most up to date version of the Philips Sonicare app check out the Android Play Store or Apple App Store. Mobile devices must have Bluetooth 4.0 or newer.

What other Philips Sonicare toothbrushes can I use with the Philips Sonicare app?

You can only use Philips Sonicare Flexcare Platinum Connected, DiamondClean Smart, DiamondClean 9000 and ExpertClean toothbrushes with the Sonicare app. Look for the Bluetooth symbol on the front panel of the package.

Disposal

- This product contains a rechargeable lithium-ion battery which must be disposed of properly.
- Contact your local town or city officials for battery disposal information. You can also call 1-800-822-8837 or visit **www.call2recycle.org** for battery drop-off locations.
- For assistance, visit our website **www.philips.com/support** or call 1-800-682-7664 toll free.
- Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit **www.recycle.philips.com**.

Removing the rechargeable battery at the end of lifetime

Warning

- Only remove the rechargeable battery when you discard the handle.
- Make sure the battery is completely empty when you remove it.
- Any attempt to remove the battery will void the warranty.

For assistance, visit our website at www.philips.com/support or call 1-800-243-3050 toll free.

Removing the rechargeable battery


Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it.

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.


- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.

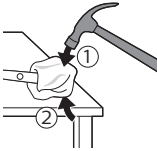
- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.


Note: You may have to hit on the end several times to break the internal snap connections.

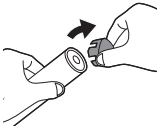
- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

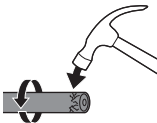
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

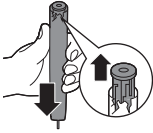- 6 Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

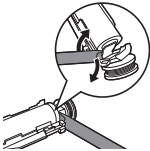
- 7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

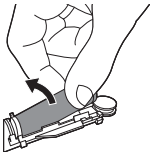- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

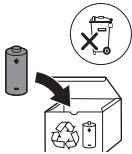- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

FCC & ISED Compliance

- Radio Equipment in this product operates at 13.56 MHz
- Maximum RF power transmitted by the Radio Equipment is 30.16dBm

Changes or modifications not expressly approved by Philips could void the user's FCC authorization to operate the toothbrush.

For the United States

This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any

32 English

interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment (handle) off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna (of the radio or television).
- Increase the separation between the equipment (handle) and receiver.
- Connect the equipment (charger) into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

For Canada

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's license-exempt RSS(s). Operation is subject to the following two conditions: 1. This device may not cause interference. 2. This device must accept any interference,

including interference that may cause undesired operation of the device.

Warranty and support

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided. Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection. Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center. Internet information: www.sonicare.com (North America) or www.philips.com/support (outside North America).

WARRANTY EXCLUSIONS

What is not covered under warranty:

- Brush heads.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES

34 English

LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

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©2019 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.
de Philips Oral Healthcare, LLC et (ou) KPNV.