TENS NECK MASSAGER

SHARPER IMAGE®

CONGRATULATIONS ON YOUR PURCHASE OF THE **TENS NECK MASSAGER**

Keep this instruction manual as it contains important information for future reference.

Due to continuous product improvements, the

picture on the manual may differ slightly from the actual product.

SHARPER IMAGE® NAME AND LOGOS (INCLUDING THE 77 LOGO AND TOMORROW'S TOMORROW) ARE TRADEMARKS OR REGISTERED TRADEMARKS. DISTRIBUTED BY AND © 2021 MERCHSOURCE, LLC. IRVINE, CA 92618. ALL RIGHTS RESERVED. FOR CUSTOMER SUPPORT, PLEASE CALL 1-800-374-2744.

Item # 1014706

Printed in China

RECHAREGABLE BATTERY

WARNINGS:

- \cdot I here is a risk of fire and personal injury if battery is uou-replaceable. · This product contains batteries that are
- bnuctaked, damaged, or misused.
- · Do not dismantle, crush, incinerate or short-circuit
- · Never expose battery to extreme temperatures or pattery.
- · Keep away from flammable materials and direct moisture.
- sunlight.
- agabrers. · Only charge using the provided charger, cable, and/or
- · Adult supervision is required while charging the battery. · Rechargeable battery is only to be charged by an adult.
- cyarding. · Allow battery to cool to room temperature before
- trom flammable products. · Cyarde ou a vou-tlammable surface and keep away
- smokes, or gives off a strong odor. · Discontinue charging it battery becomes hot, swells,
- enclosure, and other parts. In the event of such regularly examined for damage to the cord, plug, · paffery chargers used with the product are to be
- · Recycle or dispose of battery according to federal, yaz peeu tebalteg. damage, chargers must not be used until the damage
- · Keep these instructions for future reference. state, and local laws.

the battery cover.

alcohol wipe.

- \cdot to access the remote battery, use a coin to remove

· Do not submerge in water.

CARE AND TIPS

state, and local laws.

immediate medical attention.

explode or leak.

· Do not use detergent or soap when cleaning.

· Keep these instructions for future reference.

· Recycle or dispose of batteries according to federal,

swallowed or placed inside any part of the body, seek

immediately. Keep new and used batteries away from

hours and lead to death. Dispose of used batteries

can cause internal chemical burns in as little as two

Cell Battery. A swallowed Button or Coin Cell Battery

WARNING: This product contains a Button or Coin

WARNING: Contains button or coin cell battery.

· Do not dispose of batteries in fire, batteries may

· Do not short-circuit the battery supply terminals

children. If you think batteries might have been

Hazardous if swallowed-see instructions.

· Power down before cleaning. Clean with a 75%

FCC STATEMENT

This device complies with part 15 of the FCC Rules.

Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.


Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, this is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

CONTENTS

- · TENS Neck Massager
- · Wireless Remote

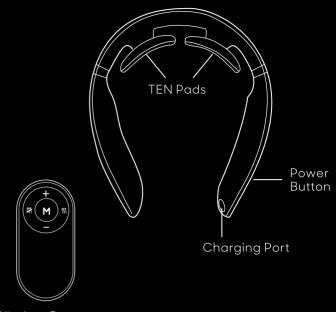

Wireless Remote



FEATURES


- · TENS massage with heat mode
- · 3 massage modes and 15 strength levels
- · Wireless remote with magnetic mounting


GETTING STARTED

- 1. Charge using the provided USB cable and unplug when fully charged.
- 2. Secure the rubber charging port cover.
- 3. Remove any electronic items or jewelry you are wearing.
- 4. Using a clean, damp towel, moisten the sides and back of your neck.
- 5. Place the massager around your neck. Make sure the TENS pads are in full contact with your skin.
- 6. Turn on the massager by pressing the power button.
 A tone will be heard when massager is turns on.
- Find your comfort level. Press the [+] button to increase massage intensity and the [-] button to decrease massage intensity.
 Note: Massager will turn off if left at the lowest setting.
- 8. Press the [M] button to toggle through the 3 massage modes.

 Note: changing modes resets massage intensity to lowest level.
- 9. Press the \$\infty\$ button to turn heat on and press the \$\infty\$ button to turn heat off.
- 10. Press and hold the [-] button turn off sound cues.

 Press and hold the [+] button to turn on sound cues.
- 11. Use for up to 15 minutes at a time but no more than 30 minutes per day.

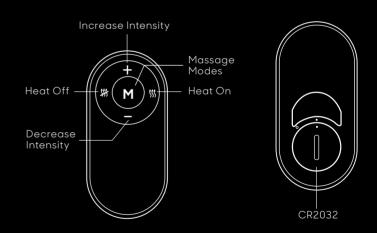

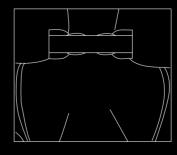



Warning: Use massager only as indicated. Do not use massager for any other body areas.