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

# QUICKSTART GUIDE

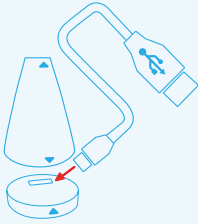
**1.** Charge your Playbrush until fully charged.

**2.** Insert any manual toothbrush into the Playbrush, aligning all 3 arrows with the toothbrush head.

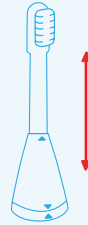**3.** Firmly shake your Playbrush twice vertically to turn it on.

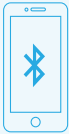**4.** Turn on Bluetooth on your smartphone/tablet.

**5.** Download one of our Playbrush games (search "Utoothia").

**6.** Launch the app and follow the instructions to connect your Playbrush. Once connected, Playbrush will flash green.

If you can't find the answers to your questions in this manual, you can also find some help here:

Visit our regularly updated FAQ  
[www.playbrush.com/faq](http://www.playbrush.com/faq)


Or send us a message  
[support@playbrush.com](mailto:support@playbrush.com)

# 1. CHARGING

Please charge your Playbrush fully before using it for the first time, usually taking a few hours.

1. Detach the white base from the coloured toothbrush holder.
2. Plug the micro-USB cable into the white base and recharge with any device that has a USB connection (e.g. computer or plug).
3. Look at the light colour:

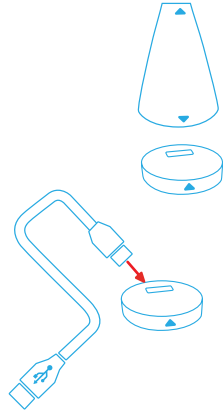**Blinking red and green:**

Playbrush needs charging.

**Blinking red:**

Playbrush is charging.

**Constant green light:**

Fully charged and ready to start playbrushing.

Once fully charged, re-attach the white base and coloured toothbrush holder, making sure all 3 arrows are aligned.

# 2. ASSEMBLING PLAYBRUSH

Insert your manual toothbrush, taking care to align the 3 arrows with the head of the brush. The bristles of the brush should be facing the same direction as the arrows.


### 3. TURNING IT ON

Playbrush uses its motion sensors to turn on and off.

**To turn Playbrush on**, firmly shake twice vertically as indicated, and the green light will begin to blink. Your Playbrush is now ready to connect. (If it blinks green and red then Playbrush needs charging.)

Your Playbrush will go into **sleep mode** automatically after being inactive for 1 minute. To reawaken, just move your Playbrush. Your device will not lose power whilst in sleep mode, so we recommend using sleep mode for day-to-day use.

**To turn Playbrush off**, shake twice again vertically, and the light will no longer blink.

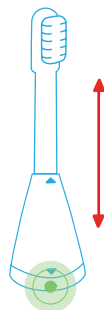
### 4. TURN BLUETOOTH ON

To use Playbrush, your smartphone/tablet needs to support Bluetooth Low Energy (Bluetooth Smart, Bluetooth 4.0, BLE). Please turn on Bluetooth now.

NOTE:

To connect Android phones to Playbrush via Bluetooth, please grant the permissions requested by the app.<sup>(1)</sup>
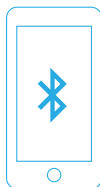
### 5. DOWNLOADING


Search 'Utoothia' on the App Store or Play Store and download one of our various games. Your smartphone or tablet must meet the following minimum requirements:

Apple: iPhone 4S / iPad 3 and iOS 8.

Android: Android 4.4 KitKat with Bluetooth Low Energy (BLE).

Kindle: 4th Generation Kindle Fire HD (2014) and Fire OS 4.


HELPFUL TIP: For younger children, we recommend starting off with our game 'Utoothia Sky'.
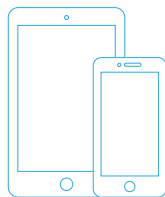
<sup>(1)</sup>In Android 6.0 (Marshmallow), Playbrush does need you to grant permission to access your location, but the permission is only for Bluetooth to work. Playbrush does not access your location.

## 6. FUN TOOTHBRUSHING

Launch a Playbrush game and follow the in-app instructions to connect your Playbrush. Once connected, Playbrush will start flashing green. We recommend brushing twice a day and we can't wait for you to begin your fun toothbrushing journey to great dental health!
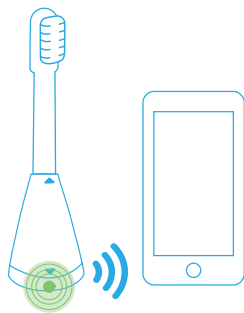
## IS YOUR PLAYBRUSH NOT WORKING?

### Playbrush Connection

- ✔ Please check that Bluetooth on your smartphone or tablet is turned on.
- ✔ Ensure location permissions are granted in the app for all Android 6.0 users.
- ✔ If you are still having connection issues, uninstall the app, turn Bluetooth off and on and reinstall the app.

### Brushing Accuracy

- ✔ The 3 arrows need to be correctly aligned with the toothbrush head.

### Charging

- ✔ Your Playbrush should be fully charged and the green light is visible. If it blinks red and green, Playbrush needs to be charged.
- ✔ Playbrush blinks red when you charge it. If it doesn't, you may need to reboot your Playbrush. To do this, charge for 1-2 minutes, remove the USB charging cable, and then plug it back in and charge for 12 hours.

If the above checklist has not solved your issue and your Playbrush is still not working, please contact us!

## CONTACT US

Visit our regularly updated FAQ [www.playbrush.com/faq](http://www.playbrush.com/faq)  
Or send us a message [support@playbrush.com](mailto:support@playbrush.com)

## SAFEGUARDS AND WARNINGS

The natural playfulness and temperament of children may result in unforeseeable situations and dangers for which the manufacturer cannot be held responsible. Children should be instructed on how to use Playbrush, draw their attention to potential dangers and keep children under supervision at all times. Playbrush GmbH doesn't accept any liability of any kind for accidents resulting from the use of Playbrush.

**Warning! Not suitable for children under 36 months due to small parts. Choking hazard.**

- When undergoing treatment of any oral care condition, consult your dental professional prior to use.
- Playbrush works with almost all manual toothbrushes but bigger toothbrush handles could expand the elasticated head.
- Use only for intended use and don't put Playbrush in the user's mouth.
- Playbrush is not waterproof when the coloured and the bottom white part are unattached, so don't clean with water or other liquids.
- Very strong magnetic fields should always be avoided as they can lead to malfunctions or damage.
- Under any circumstances, do not try to repair Playbrush on your own. The attempt to repair may cause damage and would void the warranty of your product.

### Certifications

Of course the safety of your children is our most important concern! Playbrush passed all relevant tests for Europe, the USA and Canada (EN71, BSEN 62115, BSEN 60825, REACH, ASTM F963, CPSIA, Wireless Regulatory RTTE Europe/FCC/IC, EMC for CE/FCC/IC, SAR, Bluetooth SIG, Radio Certificate USA/Canada )  
FCC ID: 2AF47PBR2 IC: 20694-PBR2 Model No.: PBR2

### Declaration of Conformity

Hereby, Playbrush GmbH declares that Playbrush complies with the essential requirements and other relevant provisions of the directives 1999/5/EC, 2009/84/EC and 2004/108/EC. The declaration of conformity may be consulted at [www.playbrush.com](http://www.playbrush.com).

### FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

No changes shall be made to the equipment without the manufacturer's permission as this may void the user's authority to operate the equipment.

This device has been designed and complies with the safety requirements for portable RF exposure in accordance with FCC rule part 2.1093 and KDB 447498 D01.

### RSS Compliance Statement

This device complies with Industry Canada's license-exempt RSS. Operation is subject to the following two conditions:

- (1) This device may not cause interference; and
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

This device has been designed and complies with the safety requirement for RF exposure in accordance with RSS-102, issue 5 for portable conditions.

### **Battery**

This product contains a non-serviceable, rechargeable lithium ion battery. DO NOT attempt to open or disassemble the unit to replace or remove the battery yourself. Please refer to original manufacturer for battery servicing and replacement. Follow safety instructions above carefully. If fluid should leak from the unit or rechargeable battery, do not touch. If fluid comes in contact with eyes, wash thoroughly with water and seek medical advice. Do not rub eyes. If you suspect that the battery is faulty through unexpected behavior such as not holding charge, please turn off and do not attempt to charge or use the product. Contact the original product manufacturer immediately.


### **Recycling Disposal**

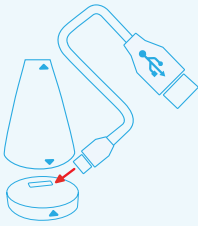This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product as household waste at the end of its useful life. This product was produced after Jan 1<sup>st</sup>, 2006.

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© 2016 Playbrush GmbH