# bewell<sup>®</sup> connect

Distribué par / Distributed by  
[www.visiomed-lab.com](http://www.visiomed-lab.com)

Visiomed SAS France  
Service clients / Customers service  
8, avenue Kléber  
75116 Paris - France  
Tel : +33 8 92 350 366 (0,34€/mn)  
Fax : +33 1 44 17 93 10  
[contact@visiomed-lab.fr](mailto:contact@visiomed-lab.fr)

Made in P. R. C.

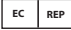
MyScale XL 200kg

BW-SC4

By Visiomed<sup>®</sup>

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BW-SC4

By **Visiomed**<sup>®</sup>

This product complies with the Reach (RoHS) standard.

## CONTENTS

<b>1. WARNINGS</b> .....	<b>15</b>
<b>2. INFORMATION</b> .....	<b>16</b>
<b>2-1. ABOUT THE BMI - BODY MASS INDEX</b> .....	<b>16</b>
<b>2-2. ABOUT DAILY CALORIE NEEDS</b> .....	<b>16</b>
<b>3. USE</b> .....	<b>16</b>
<b>3-1. PRECAUTIONS BEFORE USE</b> .....	<b>16</b>
<b>3-2. USING THE WEIGHING SCALES</b> .....	<b>17</b>
<b>3-3. CHANGING THE BATTERIES</b> .....	<b>18</b>
<b>4. CLEANING AND MAINTENANCE</b> .....	<b>18</b>
<b>5. AVOIDING MALFUNCTIONS</b> .....	<b>19</b>
<b>6. TECHNICAL SPECIFICATIONS</b> .....	<b>19</b>
<b>7. TROUBLESHOOTING GUIDE</b> .....	<b>20</b>

Learn how to install and use the BewellConnect application by connecting to our site  
[www.bewell-connect.com/install](http://www.bewell-connect.com/install)

*The manufacturer reserves the right to modify the technical characteristics of the product without prior notice.*

Dear customer,

You have just acquired this set of BewellConnect® connected and smart weighing scales and we thank you. We hope you get the best out of it and advise you to read this instruction manual carefully.

## 1. WARNINGS

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
- No modification of this equipment is allowed.
- Information regarding potential electromagnetic or other interference between the MyScale Analyzer and other devices together with advice on ways to avoid or minimize such interference.
- This device requires no calibration.
- The user must check that the equipment functions safely and see that it is in proper working condition before being used.
- Conditions that can unacceptably affect the ME equipment :
  - Effects of lint, dust, light
  - List of known devices or other sources that can potentially cause interference problems
  - Effects caused by pets, pests or children

- Use this appliance only for the use for which it is designed as described in this manual.
- This appliance may be used for personal use at home. Any other use could void the warranty and cause injury.
- Use this appliance in an ambient temperature range of 0 to 40°C.
- Do not expose this appliance to extreme temperature conditions > 60°C or < 20°C.
- Do not use this appliance at a relative humidity of more than 85%.
- This appliance must always be placed in a clean and dry place.
- Do not expose this appliance to sunlight or water.
- Do not expose this appliance to electric shocks.
- Do not use this appliance outdoors.
- Never drop the appliance.
- Follow the maintenance instructions specified in this manual.
- Do not attempt to open the appliance. In case of problems, contact your dealer.
- Do not leave this appliance within the reach of children.
- The measurement results are given for information. If in doubt about your results, please contact your doctor.
- Discontinue use of the device in case of anomalies or malfunction.
- Do not use any battery other than the mentioned batteries. Do not attempt to recharge non-rechargeable batteries. Do not throw them into a fire.
- Remove the batteries if the appliance is not used for a long period.
- This appliance is not designed to be used by persons (including children) whose physical, sensory or mental capabilities are reduced, or persons without experience or knowledge, unless they have been able to benefit, by the intermediary of a person responsible for their safety, from surveillance or prior instructions concerning use of the appliance.

**Important**

*Observe the polarity of the batteries. Incorrect polarity can cause damage and compromise the warranty of your appliance.*

**2. INFORMATION**

The human body is composed almost exclusively of the following elements: bone, muscle, organs, fat and water.

**2-1. ABOUT THE BMI - BODY MASS INDEX**

THE BMI (or BMI in English) is an indicator of your corpulence and potential medical risks, based on the ratio of weight to height.

<b>BMI DIAGNOSTIC (GIVEN AS A GUIDELINE ONLY)</b>	
<b>BMI</b>	<b>Diagnostic</b>
< 16.5	Malnutrition
16.5 - 18.5	Slim
18.5 - 25	Normal
25 - 30	Overweight
30 - 35	Moderate obesity
35 - 40	Severe obesity
> 40	Morbid obesity

**2-2. ABOUT DAILY CALORIE NEEDS**

All food provides calories and all calories provide energy. The human body uses this energy to live, breathe and perform all daily activities. A man consumes more energy than a woman. Usually, his energy needs are from 5 to 10% higher than those of a woman - with the exception of a pregnant or lactating woman, whose daily calorie needs ensure the good health of the mother and the baby.


**3. USE**

**3-1. PRECAUTIONS BEFORE USE**

- This device is not suitable for children under 10 years old. Do not leave children unsupervised when using the appliance.
- Place the scales on a sturdy, level and non-slip surface (avoid installing them on a rug, for example).
- The weight and fat mass vary during the day. We advise you to weigh yourself before meals to obtain comparable results. It is recommended that you should measure your weight and body fat mass at the same time each day, preferably on an empty stomach.

**CAUTION**  
**RISK OF EXPLOSION IF BATTERY IS REPLACED**  
**BY AN INCORRECT TYPE**  
**DISPOSE OF USED BATTERIES ACCORDING**  
**TO THE INSTRUCTIONS**

### 3-2. USING THE WEIGHING SCALES

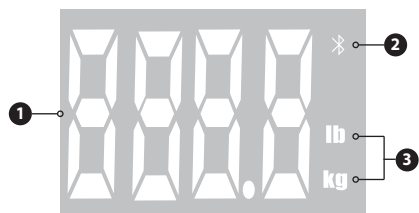1 Display

2 Bluetooth

3 Weight unit

**Overload warning :** Immediately remove anything on the scales if the screen displays Err. Irreversible damage could be caused to the appliance.

#### Choice of the weighing unit

Go to the App MyScale to select Kilos or Pounds.

#### Weighing mode

Step onto the scales. It lights up. The screen displays 0:0, and then your weight. The digits flash, and then stabilize and remain on continuously.

#### Programming of personal data

Schedule your personal data, such as age, size or sex directly in the MyScale application.

#### DOWNLOADING THE APPLICATION

Download the BewellConnect® application to your phone or tablet:

- by the App store or Google Play
  - or by flashing the QR Code on the side of the box.
- Then click the MyScale icon.

#### DATA TRANSFER

To follow the evolution of your weight and your proportions, connect your device (phone, tablet) by Bluetooth (connectivity Bluetooth 4.0 smart Low Energy, compatible from iPhone 4S).

Open the app and select your profile. Open MyScale and go to the measure screen before weighing yourself, in order to save your datas.

Switch on the scales by stepping on them.

The data transfer is automatic with the BewellConnect® application. Simply assign the measurement made to a profile. You can consult your history on the application.

### 3-3. CHANGING THE BATTERIES

The battery compartment is located on the back of the appliance.

- When using for the first time, insert 4 new AAA batteries in the scales.

- When the LCD screen displays the icon , the batteries are low or flat. Replace the flat batteries, taking care to observe the correct polarity. Incorrect polarity can cause damage and compromise the warranty of your appliance.

 Never insert rechargeable batteries. Use single-use batteries. Do not mix old and new batteries.

 Do not dispose of your old batteries with your household waste.

Take your flat batteries to a collecting agency.

If the appliance is not to be used for a long time, be sure to remove the batteries to avoid any leakage of liquid that may damage the unit.

### 4. CLEANING AND MAINTENANCE

- Do not wet the weighing plate. This could make it slippery and cause an accident.

- Store your scales in the horizontal position.

- Put your scales away in dry place away from humidity.

- Use a soft and dry cloth to clean the unit. If necessary, use a soft lightly moistened (tap water) cloth.

- Do not use alcohol, benzene, thinner or other chemicals to clean the appliance.

- This product is not waterproof. Avoid any contact with liquid.

- If the unit has been exposed to a temperature below 0°C, be sure to leave it for 30 minutes at room temperature before use. Without this waiting, the results could be distorted.

- If the unit is not used for a long period of time, be sure to remove the batteries (they could damage the unit). Observe the polarity of the batteries. Do not mix old and new batteries.
- Do not use rechargeable batteries. Do not dispose of your old batteries with your household waste.

## 5. AVOIDING MALFUNCTIONS

- Do not drop the appliance. Avoid shaking and shocks.
- Do not insert anything into the appliance, except batteries.
- Do not dismantle the unit.
- Do not put the appliance in water.
- Observe the temperature and humidity ranges stipulated in the specifications.

## 6. TECHNICAL SPECIFICATIONS

Product name	BewellConnect® MyScale XL 200kg		
Reference	BW-SC4		
Item	Units	Minimum	Maximum
Weight measurement range	kg	5	200
Accuracy of graduations			
Accuracy up to 50 kg	kg	0.5	
Accuracy up to 100 kg	kg	1	
Accuracy of 150 kg	kg	1.5	
Calories accuracy	Kcal	115	
Weight units	kg / lb		
Profiles	8 user profiles Age: 10-80 years Size: 100-200 cm		
Dimensions	350 x 350 x 28mm		
Weight	2kg		
Normal conditions of use	Temperature : 0 - 40°C Humidity : 85%		
Storage conditions	Temperature : -10 - 60°C Humidity : 85 %RH		
Bluetooth technology	4.0 smart Low Energy		
Power supply	4 batteries AAA (included)		

Battery life	3-4 months (for 4 daily measurements) Low battery indicator
Automatic shut-off	±10 seconds
Supplied accessories	A connected set of scales, instructions, a quick installation guide, 4 AAA batteries.

## 7. TROUBLESHOOTING GUIDE

If one of the following problems occurs when using your MyScale, refer to this troubleshooting guide. If the problem persists, please contact our customer support.

SCREEN ERROR	POSSIBLE CAUSES	SOLUTION
The icon  is displayed	Flat batteries	<i>Insert new batteries</i>
"O-LD" is displayed	Overload warning	<i>The maximum weight supported by the scales is 200 kg. Immediately remove the weight from the scales. Irreversible damage could be caused to the appliance</i>

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**EN:** Visiomed\* will repair or replace this product free of charge in the case of defective parts or manufacturing defects, in accordance with the conditions mentioned below as follows:

**DURATION:** 24 MONTHS RETURN TO WORKSHOP

**LIMITS AND EXCLUSIONS:** This guarantee concerns only the original final purchaser. A purchase invoice, or another proof of purchase, with this guarantee card will be required to obtain an after-sales service, in accordance with this guarantee. This guarantee card will not be extended to another person only the original final purchaser. This guarantee becomes void if the serial numbers on the product are modified, replaced, illegible, absent, or if repair has been carried out by a service not approved, including the user.

This guarantee covers only the defects of the material or parts, occurring during normal use of the product. It does not cover the damage caused during the transport of the apparatus, caused due to repairs being carried out by the distributor, by any modifications undertaken, any connection of equipment not approved by Visiomed\*, or causes contrary to those written in the user manual or notice. Moreover, the present guarantee does not cover damage due to falls, bad handling, bad installations, damage by fire, floods, lightning, or any other natural disaster. This guarantee does not cover the packing of the material, the accessories, the defects caused by commercial exposure of the product, show room, sale space, demonstration etc... Normal maintenance, cleaning and the replacement of parts where wear is normal, are not covered by the terms of this guarantee. Visiomed\* and its representatives and agents will not in any case be held responsible for any damage and consecutive damages due to the mishandling of this product. This guarantee is the only valid one at Visiomed\*, any other guarantee (commercial guarantee) except this one will not be taken into account.

**IMPORTANT:** During the guarantee period if you are dissatisfied with the repairs of this product, please contact the Visiomed\* customer service.

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**EN:** *Disposal of Waste Equipment by Users in Private Households in the European Union.*

*The symbol on the product or on its packaging indicates that this product must not be disposed of with your other household waste. Instead, it is your responsibility to dispose of your waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. The separate collection and recycling of your waste equipment at the time of disposal will help to conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For more information about where you can drop off your waste equipment for recycling, please contact your local city office or your household waste disposal service.*

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© 11/2015 BW-SC4
**EN:** *BewellConnect® is a brand registered by Visiomed Group SA. The constructor keeps the right to modify the photographs and characteristics without previous notice. Registered design. The logos and trademarks are the property of their respective owners.*