
STEP 1 Charging

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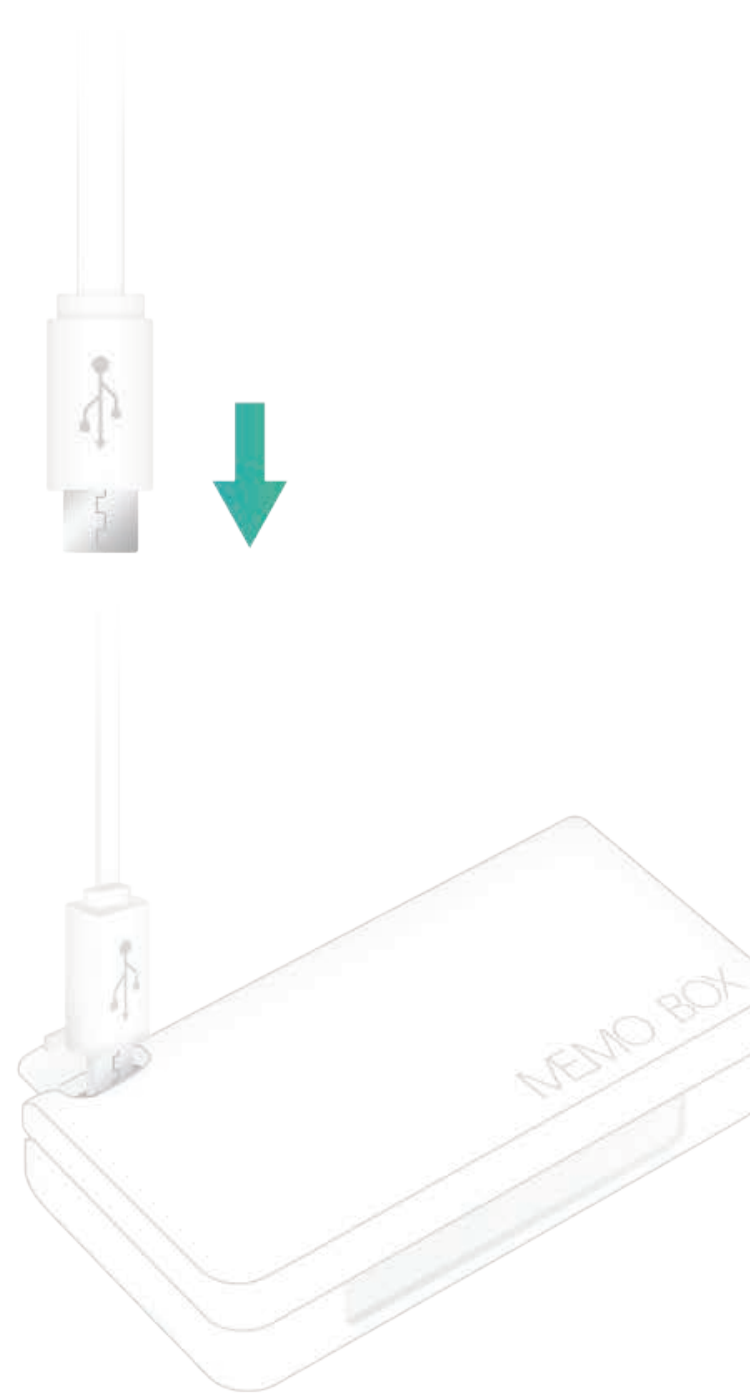
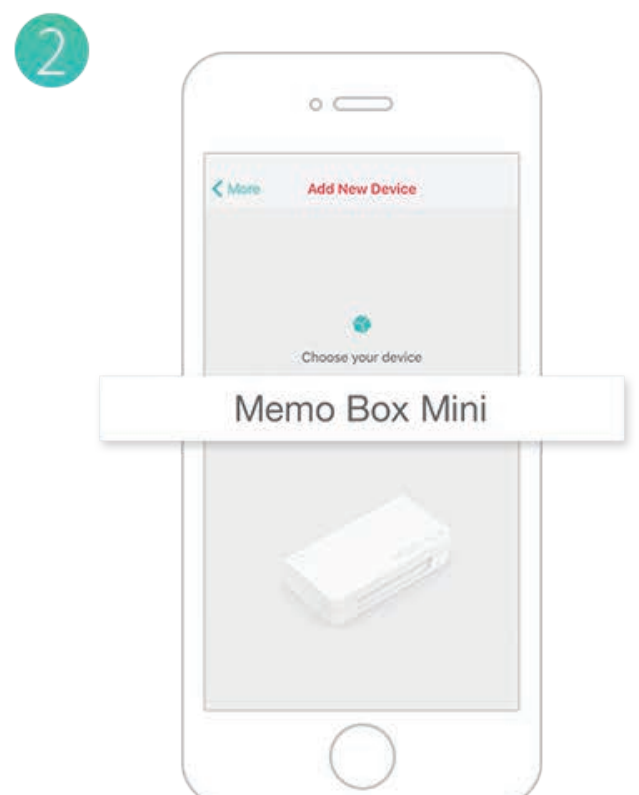
EN Choose your Memo Box and follow the on-screen prompts.

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STEP 4 Use Memo Box

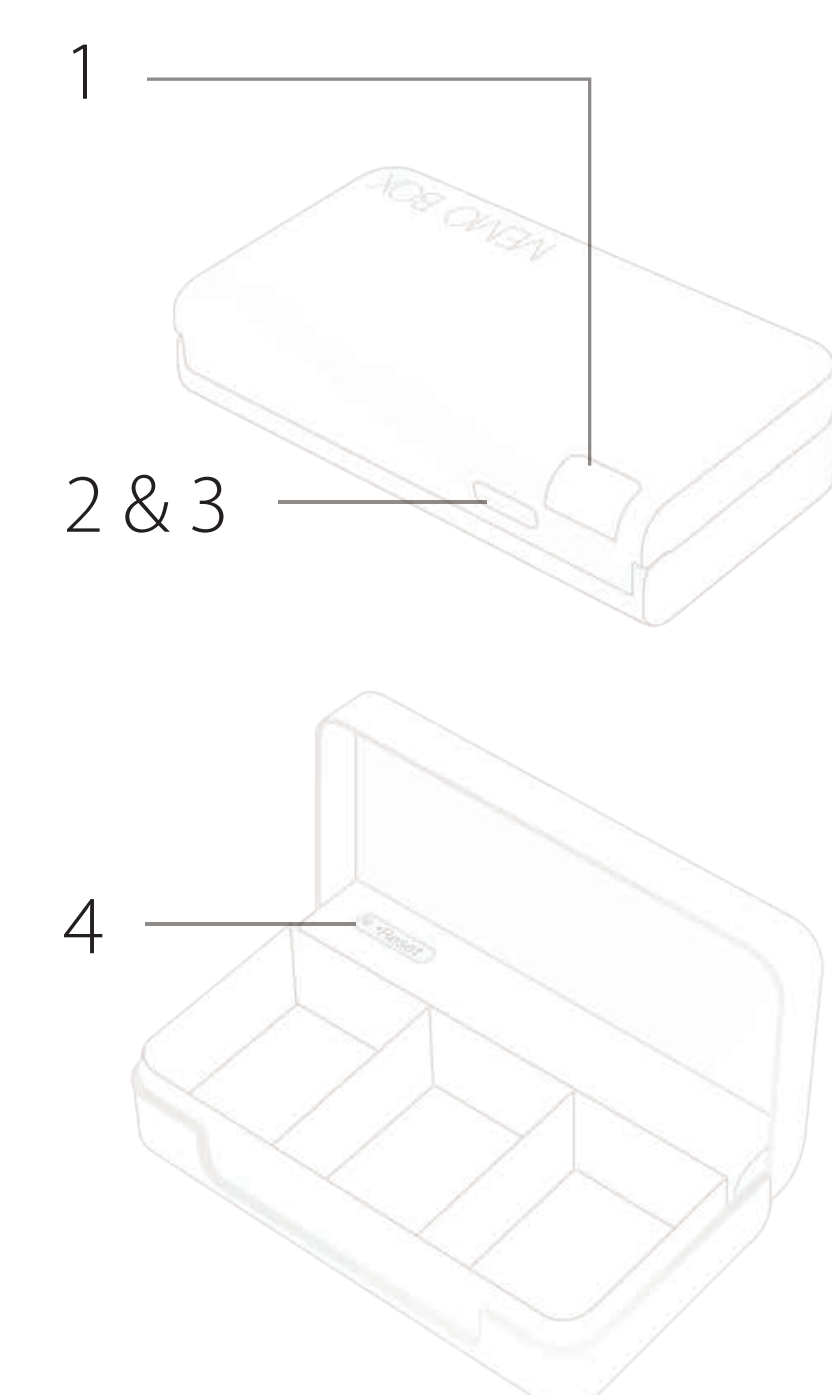EN 1. USB charging port
2. Indicator light

Red Potential double dosing
Green Time to take meds
Blue Charging

3. Snooze button
When Memo Box rings, press once to snooze.
* No 'beep' when pressed means battery is not connected / dead.

4. Reset button
Using a firm wire, such as a paper clip, push the pinhole reset button to reset Memo Box in case of device malfunctioning.

5. Clean Instruction
Wash the inner storage compartments with water, and wipe the other parts with a wet cloth. Do not spill water on the box shaft to avoid product damage.

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FCC ID: 2AH3P-M1603

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This modular complies with FCC RF radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

IC : 21387-M1603

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

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