


60*80mm 铜版纸128g

User Manual

Visible Blackhead Remover

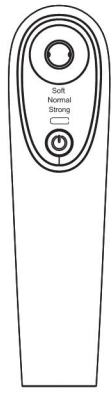

Dear customer, thanks for purchasing this product. Please read this user manual carefully and keep it for future reference.

Product Introduction

Visible Blackhead Remover is a multifunctional beauty instrument that can help clean pores, remove acne and blackheads. It uses wireless 1080P (300,000 pixel) microscope, which allows you to observe your skin status and every pore clearly, and thus help better remove the blackheads. The device can not only bring you shocking blackhead removing experience, but also prevent your skin from being damaged by excess use of the device.


Precautions

1. The battery may run out during long-distant transportation. It's recommended to fully charge the device before first use.
2. For better pore cleaning effect, it's recommended to use blackhead removing essence and/or facial steamer to open the pores before use the device. The pore

cleaning effect will be affected by closed pores.

3. It's recommended to test the suction with your forearm before using on your face.
4. To avoid injuries, do not allow the suction tip stay too long on one place, and do not press it too hard or pull it out vertically. Incorrect use may cause facial redness, which will fade away in a few days.
5. Do not apply the device to wounds or scars on your skin.
6. Please keep the body of the device away from water or other liquid.
7. Do not use chemical cleanser to clean the device, or the device may be corroded.
8. It's recommended to use the device once a week for dry skin and twice for oily skin. Do not use the device for more than 5 minutes every time. Excessive use may cause skin allergies.
9. Please clean the suction tips after each use to avoid the filter being blocked. Use alcohol with a concentration of less than 70% to clean the suction tips if needed.
10. Do not use the device on the delicate skin of children to avoid injuries.
11. The device is equipped with a Micro-USB charging cable, and a 5V/1A adapter is not included.

Product Overview

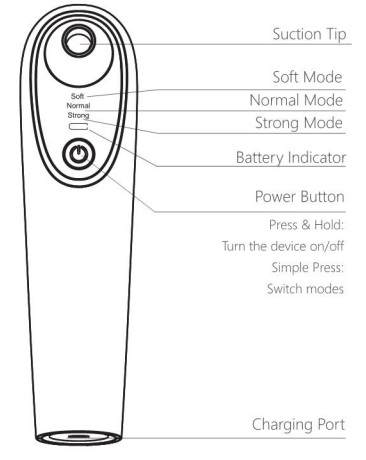

- Suction Tip
- Soft Mode
- Normal Mode
- Strong Mode
- Battery Indicator
- Power Button
- Press & Hold: Turn the device on/off
- Simple Press: Switch modes
- Charging Port

Charging Methods


1. The battery indicator flashes red when the battery is low, which reminds you to charge the device.
2. Connect the device with a 5V/1A adapter (not included) via the charging cable. Plug into the socket to start charging.

APP Instruction

I. Download the APP

Method 1 (iOS and Android System)
Search "X-Beauty" in the App Store or Google Play, click to download and install the APP.

Method 2 (iOS and Android System)
Scan the QR Code, then download and install the APP.




iOS System Version: 6.0 or higher
Android System

3. The battery indicator keeps red during charging and turns white once charging completed. Unplug the adapter then.

II. Connect the APP and the Device

1. Press and hold the Power Button on the device for 3 seconds to turn it on. The battery indicator will light on.
2. Turn on the WLAN of your phone and enter the WLAN setting interface. Find "XPRE150" and click to connect.
3. After connected, open the APP and click "V" to enter the camera working interface. The 6 LED fill lights on the head of the device will switch on.

III. APP Interface Instruction

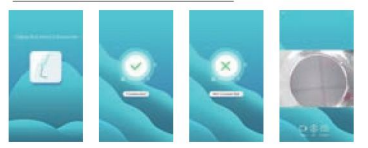



Application Methods

1. APP main interface
2. Connected: click "V" to start use the device with camera
3. Not Connected: please connect the APP and the device according to the method above
4. Camera working interface


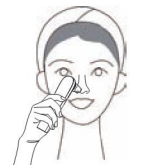
Icon Function

- [X] Icon: Click to disconnect the APP and the device
- [?] Icon: Click to enter the Help Center
- [AR] Icon: Click to zoom in or zoom out
- [Photo] Icon: Click to take pictures
- [Video] Icon: Click to record videos


2. Clean your face, use blackhead removing essence and facial steamer to help open the pores for better result.
Note: The effect may be affected if the pores are closed.


4. Short press the Power Button to choose a certain mode. The corresponding mode indicator will light on.

6. After clean, press and hold the Power Button to turn the device off.
Note: Do not use the device for more than 5 minutes each time.
Note: If you do not want to use the APP you can close the APP or reboot the device.



Note: The body of the device is not waterproof. Do not rinse under water or immerse it into water.



8. Then apply skin care products like facial mask, toner, cream, etc.

Store & Care

1. When cleaning, please make sure that the device is turned off.
2. Remove the suction tip from the device before washing it.
3. The main device is not waterproof. Clean it with wet cloth and not to wash it under tap or in the dish-washer.
4. If the device is damaged, stop use it.
5. Don't disassemble the device by yourself. If there's any problem or the device cannot be used normally, please stop using it and contact the customer service center.

Specifications

Material: ABS+PC
 Battery Capacity: 800 mAh
 Input Current: 500 mA(max)
 Output Power: 4W
 Product Size: 163*43*60mm
 Product Weight: 143g
 Package Size: 207*115*48mm

Package Contents

- 1 * Visible Blackhead Remover
- 2 * Sealing Ring
- 1 * Big Round Suction Tip
- 1 * Oval Suction Tip
- 1 * Skin Smooth Suction Tip
- 1 * Micro-USB Charging Cable
- 1 * User Manual

3. Select a suction tip and install it onto the device by pressing it down to the head of the device until it stopped.
Note: For first time of use, recommend to test the suction on your inner arm before use the device on your face.



5. Apply the suction tip to your face and slowly move it to remove blackheads.
Note: To avoid injuries, do not allow the suction tip stay on one place for more than 3 seconds, and do not press it too hard or pull it out vertically.



7. Remove the suction tip from the device by pulling it out with a slight angle. Wipe the dirt off with cotton swab and clean the suction tip with water. Then towel it dry or wait for natural air dry.



This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction