

## Compatibility

**Veta™ smart case** is a UV-protected, transparent case with sensors and electronics to monitor and communicate. Veta smart case connects to compatible iPhone®, iPad®, iPod touch® or Android™ mobile devices using Bluetooth® Smart. Veta smart case is compatible with:

- iPhone® 5 and newer
- iPod touch® 5th Gen & newer
- Samsung Galaxy® S5/S6
- Samsung Note® 4/5
- Google Nexus™ 5/6
- Android™ 4.4.4 devices and above with Bluetooth® Smart.

**Veta app** runs on iPhone® or Android™ smart phones and Wi-Fi enabled tablets and devices supporting Bluetooth® Smart, including iPad® and iPod touch®. These include recent versions of Android phones and tablets supporting Android version 4.4.4 or greater and Apple products supporting iOS version 8 or greater.

*Veta is not a medical device and is not intended, nor should be relied upon, to notify a caregiver, clinician or emergency services for the immediate delivery of care.*

Veta™ is a trademark of Aterica Inc. EpiPen® and EpiPen Jr® are registered trademarks of Mylan Inc. All other trademarks mentioned herein are trademarks of their respective owners. Mylan Inc. and Pfizer Canada Inc. do not sponsor or endorse the products referenced herein.

## Regulatory Compliance

This device complies with Part 15 of the FCC Rules and with Industry Canada's licence-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

This portable transmitter with its antenna complies with FCC/IC RF exposure limits for general population / uncontrolled exposure.

**Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.**

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes :

- (1) l'appareil ne doit pas produire de brouillage, et
- (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

# Getting Started with Veta

Congratulations on your Veta™ purchase!

Let's get it ready to use.



To locate the FCC/IC label, remove the battery cover from the battery compartment. The label is located on the side of the battery compartment.

**Model:** Z1604 **FCC ID:** 2AHBP-Z1604 **IC:** 21100-Z1604

VGS20160129




**www.aterica.com**

## Assemble Your Veta Smart Case

You have everything you need to get up and running:


- 1 empty Veta smart case with yellow cap attached
- 1 alternative green cap
- 2 AAA batteries
- The Veta 'smarts' — a cylinder that holds all the electronics, lights, sounds and batteries.

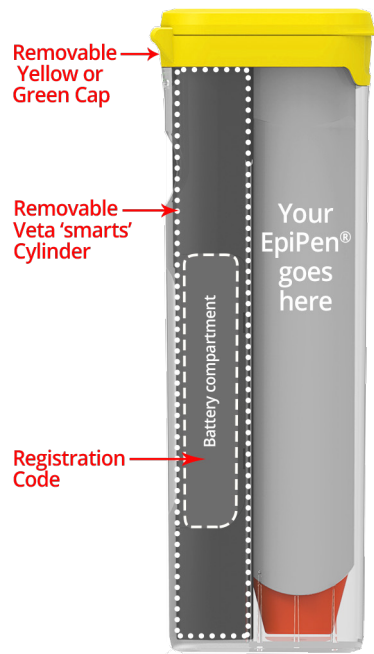

- 1. Choose a cap:** Choose a yellow or green cap (to match your EpiPen® dosage color). (To replace the yellow cap with the green cap, carefully remove the cap by pushing upwards with your thumb just under the cap's hinge. Then, snap on the green cap.)
- 2. Locate your registration code:** Each Veta smart case has a unique registration code located inside the battery compartment in the 'smarts'. Before continuing, remove the battery cover and record this unique number on the next page (see \*). (You'll record the other information during setup.)

- 3. Add batteries:** Unwrap and insert the two AAA batteries into the Veta smarts battery compartment as illustrated in the compartment. Veta flashes as it comes to life. Then, slide the battery cover into place.

- 4. Make the case smart:** Open the cap and slide the 'smarts' into the Veta smart case, with the rounded edge next to the cap hinge. Push the smarts down gently, until it snaps into place.

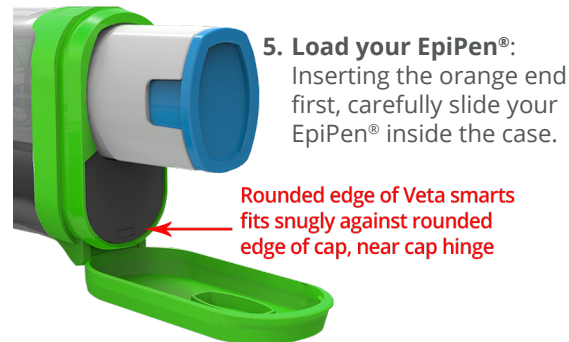

- 6. Ready:** Close the cap securely. Veta flashes again to confirm it's ready for setup.

## Download the Veta App

Download Veta app on as many mobile devices or computers as necessary. You will require an app for the Veta smart case owner and one for each person included in that person's support circle.



## Set up Your Veta Smart Case

The Veta app will guide you through registering, activating and setting up your Veta smart case. The app will help you define your custom privacy, temperature, zones, expiry date and support circle settings. Begin this process by opening the app on one of these devices.

Purchased On: \_\_\_\_\_  
Activated On: \_\_\_\_\_  
Model Number: \_\_\_\_\_  
\* Registration Code: \_\_\_\_\_  
Case Name: \_\_\_\_\_  
Smart Case Owner: \_\_\_\_\_



### Questions?

Call 1.855.476.1390 x202  
Visit [www.aterica.com/support](http://www.aterica.com/support)