

GROWP

Growp Scale | Model: BF1513-01

For safe use of the product, please read this User Manual first and save it for future reference. Please also keep safe the enclosed Warranty Card. Do not lose or discard it.

1 Product Use and Features

First, we would like to thank you for choosing Bongmi's product. This product offers a more accurate and more reliable way to help you understand your baby's growth.

Product Application: weight measurement

Product Features:

- Simple design — no button, no screen.
- Auto on by stepping on.
- Data displayed on and stored in the height ruler, no need to bend down to check.
- Auto uploading of height ruler data after syncing with BonBaby APP, permanent Cloud storage.
- Support baby-holding mode for baby weight measurement.

2 Precautions


⚠ Precautions for Use

- Do not jump on the product.
- Leave the infrared transmitter uncovered while measuring.
- Please place the scale on a flat and hard surface. Uneven surfaces or soft surfaces like carpet may result in inaccurate results.
- Please use it in a room lower than 3.5m with a light-colored ceiling; otherwise, weight data may not show on the height ruler.
- Please do not step on the scale with wet feet or when the scale surface is wet to avoid the risk of slipping.

⚠ Precautions for Storage

- Do not soak it directly in water.
- Do not leave it in high temperature, high humidity, dusty environment, direct sunlight or directly against the air conditioner.
- Do not leave it in places where a strong impact may occur, vibration is present or on an uneven surface.

3 Product Illustration

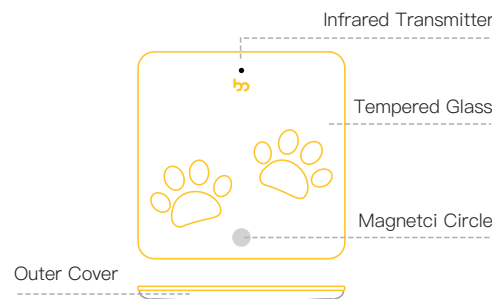


The scale consists of tempered glass, an infrared transmitter, a magnetic circle, an outer cover, etc.

4 Accessories List

Please check that the following items are included in the package: a Growp Scale, 4 1.5V (AAA) batteries, and a User Manual.

In case of anything missing, please contact our after-sales support.

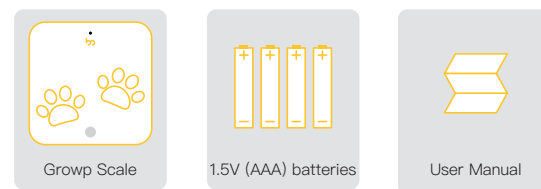

5 Battery Installation Instructions

Battery installation instructions for Growp Scale:

- Open the battery compartment cover.
- Insert 4 AAA batteries into the battery compartment of the scale with the positive and negative ends facing the correct directions as indicated.
- Put the battery cover back on in the reverse direction.

Friendly Reminders from Bongmi:

- If the Growp Scale is not going to be used for some time (3 months or above), please remove the batteries.
- Do not mix old batteries with new ones or different types of batteries to avoid affecting the life of the product.
- This product and used batteries shall be disposed in compliance with the electric product disposal regulations of the country of use.


6 How to Use

Step 1

Stand onto the scale with both feet and keep balance.




Step 2

When the Bongmi logo () is flashing, the scale is measuring the weight.



Step 3

Once measurement is finished, the Bongmi logo () will remain on and weight data will appear on the height ruler.



7 Download BonBaby App

Search "BonBaby" on App Store or Google Play to download



Data will be synced with your APP through Bluetooth, enabling scientific and quantitative management of your baby's growth and offering professional analysis and advice .

Apple devices: iPhone 5 or above with version at least iOS 8.0, iPhone 4/4s excluded.
Android devices: 4.3 or above and supports Bluetooth 4.0/4.1

8 Technical Specification

Product name:
Growp Scale

Division:
50g

Model:
BF1513-01

Unit:
Kg/Lbs

Size:
312mm*312mm*23.5mm

Weight:
1.91Kg(without battery)

Battery:
4 1.5V (AAA) batteries

Operating conditions:
Temperature: 0°C~40°C
Humidity: ≤85%RH

Battery endurance:
1 year

Storage and transportation conditions:
Temperature: -20°C~60°C
Humidity: ≤85%RH

Measurement:
3-150Kg

9 Troubleshooting

Symptoms	Reasons	Solutions
Bongmi logo does not flash when standing on the scale.	The batteries of the scale are dead or hardware damaged.	If it still doesn't work after replacing the batteries, please call our customer service.+86-400-1068-321
Bongmi logo is on, but the LED dot matrix screen doesn't show the measurement data.	The infrared transmitter is blocked.	Uncover the infrared transmitter.
	Batteries of the height ruler are running low.	Replace the batteries of the height ruler.
	The level ruler of the height ruler is unfolded.	Turn the level ruler till it's folded.
	Hardware of the height ruler is damaged.	Please call our customer service.+86-400-1068-321
	Room is higher than 3.5m.	Place it in a room lower than 3.5m

10 Maintenance and Care

- Keep the Growp Scale away from heat.
- This product is not waterproof. Do not let water drip or splash into the product. In case of exposure to water, please wipe with a dry cloth as soon as possible.
- Wipe the scale surface with a clean and soft cloth. If necessary, wipe it with a wrung-out cloth and then dry it with a dry cloth.
- Avoid wiping with a rough cloth in case of scratches on the scale.
- Do not wipe this product with chemical or corrosive solvents.

11 Warrantly and Service

Warrantly	
Product model:	Point of sales:
Purchased data:	Note:
Customer service: +86-400-1068-321	

Service:

- 1-year warranty (from date of purchase) is included under proper usage.
- Warranty doesn't cover parts other than Growp Scale itself.
- Under these circumstances, repair will be charged:
 - A: Artificial destruction, inappropriate use, or used for commercial purposes
 - B: Warranty expired
 - C: Failure to present the Warranty Card
 - D: Content of the Warranty Card is inconsistent with the actual object, or blurred and unrecognizable
 - E: Dirt and scratches on the product itself or accessories.
 - F: Repaired and modified by the user, damaged by natural disaster or external force

12 Disclaimer and Copyright Statement

Disclaimer:

This manual is written in accordance with available information and contents are subject to change without notice. Bongmi has made the best effort to guarantee the manual's accuracy and reliability. However, the company does not assume liability for loss or damage caused by omissions, inaccuracies or typographical errors in this manual. If there are discrepancies in this manual due to software upgrades, instructions and function dictated by the software.

Copyright Statement:

BONGMI INC. reserves all rights, and retains the final right to interpret and revise this handbook and the statements contained therein. Copy right of this manual belongs to BONGMI INC. Without prior written permission, no company or individual is allowed to reproduce, copy or translate this manual or any information contained therein.



GROWP

Growth Ruler | Model: BF1513-01

For safe use of the product, please read this User Manual first and save it for future reference. Please also keep safe the enclosed Warranty Card. Do not lose or discard it.

1 Product Use and Features

First, we would like to thank you for choosing Bongmi's product. This product offers a more accurate and more reliable way to help you understand your baby's growth.

Product Application: height measurement

Product Features:

1. Original measuring solution for easy and accurate height measurement
2. Cute shape, small and portable, suitable for use at home
3. Integrated Bluetooth, complementary with BonBaby APP. Auto transmission and permanent Cloud storage of data
4. Simple no-button design, turn and unfold the level ruler to activate height measuring
5. Off to sleep mode automatically, no need to turn off manually

2 Precautions


⚠ Precautions for Use

1. When measuring height, it is recommended to stick the ruler end to the middle of the metal circle on the scale.
2. Do not pull the measuring tape too fast so that the ruler end does not break away from the metal circle and hit the person on the scale.
3. When removing the ruler end from the magnetic circle, it is recommended to slide it out from the side of the metal circle.
4. After the measurement, please withdraw the measuring tape slowly. Withdrawing too fast may result in failure of placing it back to its original position.

⚠ Precautions for Storage

1. Do not soak it directly in water.
2. Do not leave it in high temperature, high humidity, dusty environment, direct sunlight or directly against the air conditioner.
3. Do not leave it in places where a strong impact may occur, vibration is present or on an uneven surface.

3 Product Illustration

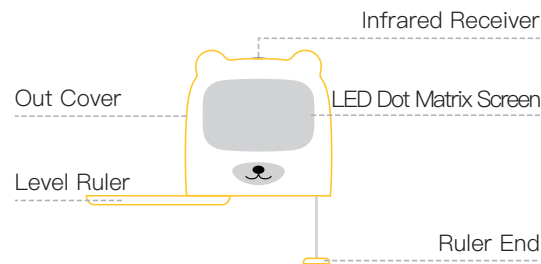


The height ruler consists of an outer cover, an infrared receiver, a LED dot matrix screen, a level ruler and a ruler end.

4 Accessories List

Please check that the following items are included in the package: a Growth Ruler, 2 1.5V (AAA) batteries, and a User Manual.

In case of anything missing, please contact our after-sales support.

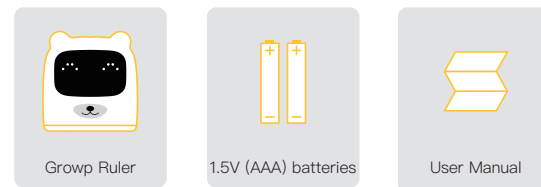

5 Battery Installation Instructions

Battery installation instructions for Growth Ruler:

1. Open the battery compartment cover.
2. Insert 2 AAA batteries into the battery compartment of the scale with the positive and negative ends facing the correct directions as indicated.
3. Put the battery cover back on in the reverse direction.

Low battery warning: if the battery icon flashes when turning the device on, it means that the batteries are running low. Please replace the batteries. Batteries with insufficient power may affect the accuracy of the measurement.

Friendly Reminders from Bongmi:

1. If the Growth Ruler is not going to be used for some time (3 months or above), please remove the batteries.
2. Do not mix old batteries with new ones or different types of batteries to avoid affecting the life of the product.
3. This product and used batteries shall be disposed in compliance with the electric product disposal regulations of the country of use.

6 How to Use

Step 1

Turn and unfold the level ruler, move the ruler end closer to the scale and stick it onto the metal circle.



Step 2

Let your child stand upright on the scale and do not move. Pull the height ruler up, attach the level ruler onto the top of your child's head and keep it balanced.



Step 3

The measurement is completed after 3 seconds; the height data will be displayed on the LED dot matrix screen and flash 3 times.



7 Download BonBaby App

Search "BonBaby" on App Store or Google Play to download



Data will be synced with your APP through Bluetooth, enabling scientific and quantitative management of your baby's growth and offering professional analysis and advice.

Apple devices: iPhone 5 or above with version at least iOS 8.0, iPhone 4/4s excluded.

Android devices: 4.3 or above and supports Bluetooth 4.0/4.1

8 Technical Specification

Product name: Growth Ruler

Model: BF1513-01

Size: 95mm*90mm*29mm

Battery: 2 1.5V (AAA) batteries

Measurement: 1-180cm

Wireless transmission: BLE 4.1

Transmission distance: 10m

System requirements:

iOS: iPhone 5 or above with version at least iOS 8.0, iPhone 4/4s excluded. Android: 4.3 or above and supports Bluetooth 4.0/4.1

Division: 0.2cm

Unit: cm/in

Weight: 120g(without battery)

Battery endurance: 1 year

Operating conditions: Temperature: 0°C-40°C Humidity: ≤85%RH

Storage and transportation conditions: Temperature: -20°C~60°C Humidity: ≤85%RH

10 Maintenance and Care

1. Keep the Growth Ruler away from heat.
2. This product is not waterproof. Do not let water drip or splash into the product. In case of exposure to water, please wipe with a dry cloth as soon as possible.
3. Wipe the scale surface with a clean and soft cloth. If necessary, wipe it with a wrung-out cloth and then dry it with a dry cloth.
4. Avoid wiping with a rough cloth in case of scratches on the scale.
5. Do not wipe this product with chemical or corrosive solvents.

11 Warranty and Service

Warranty	
Product model:	Point of sales:
Purchased data:	Note:
Customer service: +86-400-1068-321	

Service:

1. 1-year warranty (from date of purchase) is included under proper usage.
2. Warranty doesn't cover parts other than Growth Scale itself.
3. Under these circumstances, repair will be charged:
 - A: Artificial destruction, inappropriate use, or used for commercial purposes
 - B: Warranty expired
 - C: Failure to present the Warranty Card
 - D: Content of the Warranty Card is inconsistent with the actual object, or blurred and unrecognizable
 - E: Dirt and scratches on the product itself or accessories.
 - F: Repaired and modified by the user, damaged by natural disaster or external force

12 Disclaimer and Copyright Statement

Disclaimer:

This manual is written in accordance with available information and contents are subject to change without notice. Bongmi has made the best effort to guarantee the manual's accuracy and reliability. However, the company does not assume liability for loss or damage caused by omissions, inaccuracies or typographical errors in this manual. If there are discrepancies in this manual due to software upgrades, instructions and function dictated by the software.

Copyright Statement:

BONGMI INC. reserves all rights, and retains the final right to interpret and revise this handbook and the statements contained therein. Copy right of this manual belongs to BONGMI INC. Without prior written permission, no company or individual is allowed to reproduce, copy or translate this manual or any information contained therein.

FCC Caution:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.