# Hello

INNER BALANCE™  
BLUETOOTH®


## **Welcome to your new Inner Balance® tech product.**

This brief guide shows you how to pair your Inner Balance product with your device. Once you start the Inner Balance app, we'll be right there with you giving you a quick video tour of the app's features and what you can expect.

**Add Heart.™**

# Inner Balance®

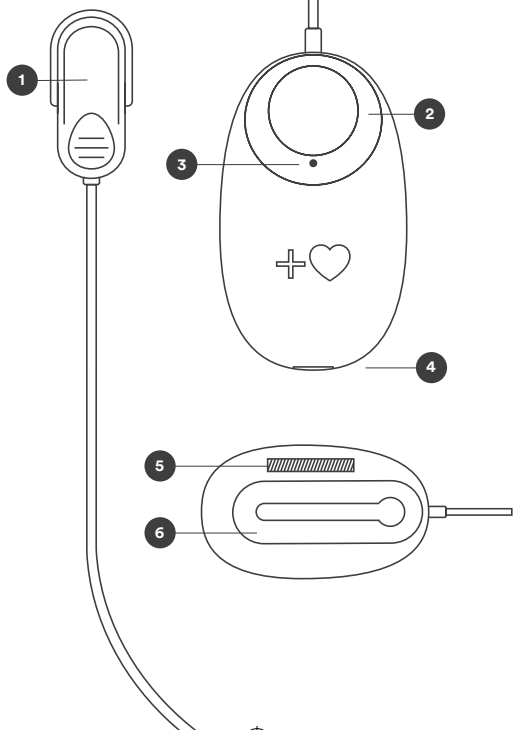- 1 **Ear clip** / Ohrclip / Oorclip /  
Sujeción para la oreja / 耳夾
- 2 **Power button** / Einschalttaste /  
Aan-uit knop / Botón de encendido /  
電源按鈕
- 3 **LED light** / LED-Licht / Led-lampje /  
Luz LED / LED 燈
- 4 **USB port** / USB-Anschluss / USB-poort /  
Puerto USB / USB 連接埠
- 5 **Identifier label** /  
Identifikationskennzeichnung /  
Identificatielabel /  
Etiqueta de identificador / 識別標籤
- 6 **Clothing clip** / Clip zur Befestigung an  
der Kleidung / Kledingclip / Sujeción para  
la ropa / 衣夾

## Pairing Instructions

1. Attach the USB cord to your Bluetooth® sensor and charge for 2 hours before use.
  2. Download and install the Inner Balance App from Google Play for your Android device or from the App Store for your iOS device.
  3. With your sensor charged, open the app.
  4. Once the Inner Balance app is open, press the power button on the sensor and attach the ear clip to your earlobe. Make sure the LED on the sensor is lit.
  5. Select the Inner Balance sensor from the list (ie: HeartMath-HRV-A123). If you are in the vicinity of other HeartMath sensors, make sure you select the one whose last 4 characters correspond to the identifier label located on the back of your sensor.
  6. In the app, if you do not see the scan dialog, tap the Start icon on the bottom toolbar to prompt it. If the Start icon is not visible, tap once on the display to bring the toolbar into view.
- 
- 

Now you're ready to begin your Inner Balance experience. Your Inner Balance will automatically pair with your phone in future sessions.

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# Warranty Statement

## Limited Warranty

The pulse sensor and peripheral accessories are warranted to the original owner for a period of One Year from the date of purchase against defects in manufacturing by HeartMath, and not by the dealer from whom the equipment was purchased.

This Limited Warranty is void if failure of the hardware has resulted from accident, abuse, or misapplication.

### What HeartMath will do:

If the pulse sensor develops a defect within the one year period, the sensor will be repaired or replaced with a new or thoroughly reconditioned and inspected unit. Proof of purchase is required.

If proof of purchase is missing or product was not purchased new, the warranty will be based on the manufacture date of the pulse sensor. For peripheral accessories, a One Year warranty only applies if purchased new and proof of purchase is provided.

For HeartMath's product knowledge base, technical or warranty assistance, please go to <http://www.heartmath.com/support>.

## FCC Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

HeartMath • 14700 West Park Avenue, Boulder Creek, CA 95006  
[www.heartmath.com](http://www.heartmath.com) • [info@heartmath.com](mailto:info@heartmath.com)

HeartMath is a registered trademark of Quantun Intech, Inc.