

1
QR CODE
<https://qra.colgate.com/27pwk-qr>
 currently redirecting to
<https://www.colgate.com.au/>

185mm

Colgate®
PULSE

Hi there, let's get to know each other

Scan to learn more

A brush that stays in tune with your smile

Replaceable head
 Amplify your smile with any Colgate Pulse expert designed brush heads

Keep it fresh by replacing your brush head every 3 months as recommended by dental professionals

Smart-Sonic technology
 For a superior clean*

Find your brushing rhythm with 2 modes
 Gentle Clean, Daily Clean

Smile Light Indicator
 Make your brush smile when you hit daily brushing goals or get into a healthy multi-day brushing rhythm

2 Minute timer and quad pacer
 Guides you to a better brushing routine by helping you reach your daily brushing goal

Rechargeable Brush
 One charge = 7 days when brushing twice daily

Power Source
 Slim charging base to style up your sink

Ergonomic Handle
 Your hand won't want to let go

*vs an ordinary manual toothbrush

An app that tracks your brushing routine

Sync up to the app for personalised tracking to guide you to better brushing

Brush Tracking
 Guides you to better brushing by tracking your brushing frequency and duration

Smile Rewards
 Get into a healthy brushing rhythm and get rewarded*

Stay connected - Online or Offline
 The app tracks your progress with or without the app open

Your way, everyday
 Choose to use the Colgate Connect app or save information to the Apple® Health app

*See colgate.com/apprends for more information on how you could get rewarded

GET STARTED

1. CHARGE YOUR TOOTHBRUSH
 Before your first use, we recommend fully charging your toothbrush. Plug the Pulse charging base into a certified AC adapter or other USB power source with a minimum power of 1.25W (5V = 0.25A).

Charge your Pulse electric toothbrush in a dry location away from water.

The Smile light pulses teal or white while charging
 The Smile light will be teal or white when fully charged
 The Smile light will flash orange when it needs to be recharged

Leave charger unplugged between charges.

2. SYNC YOUR TOOTHBRUSH
 Download the Colgate Connect app and make sure the Bluetooth® is turned on.

Create your profile and sync your toothbrush every 10 days to view your full brushing analysis and collect smile rewards.

The Smile light will flash blue when it recognises the brush during the connection process.

3. TURN IT ON AND FIND YOUR BRUSHING RHYTHM
 Before turning on the brush, press your Colgate Pulse brush head firmly into place. Press the on/off button to switch between brushing modes. The brush stays in tune with you and opens in the last mode used.

DAILY CLEAN
 Everyday clean
GENTLE CLEAN
 Low vibration soft clean

Pressing the button within a 3 second period will switch modes. If you press the button after 3 seconds the brush goes into standby mode.

4. SMILE COACHING & SMART TECHNOLOGY
 Be guided by the app or the smart sensors if brushing offline.

2 minute timer and quad pacer:
 Guides you to a complete clean by pulsing every 30 seconds and buzzing after 2 minutes

5. TURN YOUR BRUSH OFF
 After pressing the off button, the brush will remain in standby mode for 20 seconds.

You can take a quick break and if you resume brushing within the 20 second period your 2 minute brushing streak will not be impacted.

After 20 seconds in standby mode the brush will turn off.

6. MAKE YOUR BRUSH SMILE
 Get in the right flow and the brush will let you know.

Teal: You've hit your daily brushing goal of 2 minutes twice a day
Rainbow: You're in a healthy brushing rhythm with a multi-day streak
White: Default

Note: When charging the brush, the smile light will pulse in dark teal if you're in a healthy brushing rhythm and pulse white if not.

90mm 90mm 90mm 90mm

READ ALL INSTRUCTIONS BEFORE USING

IMPORTANT SAFEGUARDS
 When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

WARNING

- This toothbrush can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the toothbrush in a safe way and understand the hazards involved. Not suitable for children under 8. Keep out of reach of small children.
- Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Use only a certified AC Adapter or other USB power source with USB outlet and output voltage rating 5Vdc and output current rating of at least 250mA.
- USB power sources must only be used in dry locations away from water.

INSTRUCTIONS

CLEANING INSTRUCTIONS
 Clean your toothbrush head thoroughly after each use by gently pulling the brush head from the handle and rinsing under running water.

Regularly clean the handle and charging base of the toothbrush using a damp cloth. This is especially important if your toothpaste contains baking soda or other bicarbonates or peroxide. Otherwise, plastic cracking could occur.

STORAGE INSTRUCTIONS AND DEVICE MAINTENANCE
 The following maintenance instructions are important to ensure that your brush continues to work as it was designed. Failure to follow these instructions may cause your device to stop functioning altogether.

If you plan to stop using the toothbrush for an extended period, we recommend that you unplug the charging base, clean the toothbrush handle and heads, and store them in a cool dark place. If you start using your toothbrush again after a long period of non-use, place it on the charger for a latest 10 hours.

ELECTROMAGNETIC COMPATIBILITY
 The Colgate electric toothbrush complies with all standards regarding electromagnetic fields. It is safe to use if handled properly according to the instructions provided in this manual.

BLUETOOTH®
 The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and are used under license.

CHARGING BASE
 Input 5Vdc, 0.4A

HANDLE
 Model A2000C

NOTICE
 For energy savings purposes, charge only when the battery level is low. When not in use, keep the charger unplugged from the socket.

HOW TO RECYCLE YOUR ELECTRIC TOOTHBRUSH
 We're committed to reimagining a healthier future for all people, pets and planet through the commitments we make to sustainability and social impact - and the actions we take.

We've partnered with TerraCycle to make recycling your Electric Toothbrush easy! To learn how visit <https://www.colgate.com/electrical>

WARRANTY
 You can view the warranty of your purchased product by visiting our website at <https://shop.colgate.com/pages/warranty>.
 In the event the device fails to perform as intended, please contact Colgate-Palmolive Company at 1-823-862-9251 or www.colgate.com/pulse.

FCC STATEMENT
 WARNING: Changes or modifications to this unit not expressly approved by Colgate-Palmolive Company could void the user's authority to operate the equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.


NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 —Reorient or relocate the receiving antenna.
 —Increase the separation between the equipment and receiver.
 —Connect the equipment into an outlet on a circuit different from that

SAVE THESE INSTRUCTIONS



+61 02920 4418 email@tweak.com.au

FILE NAME PXXXXXXXXX Colgate Pulse Adopt Series 1 Manual 2022.ai
 PROJECT Colgate Pulse ETB 2022 PRINTER
 CLIENT Colgate-Palmolive CONTACT
 JOB NO CP3675 PRINT PROCESS
 DATE 02.02.22 MAX COLOURS
 REVISED 09.06.22 MIN DOT
 SOFTWARE Adobe Illustrator CC2020 SUBSTRATE XXXX
 OPERATOR AJC TRAPPING Not included
 NON PRINTING
 DIELINE

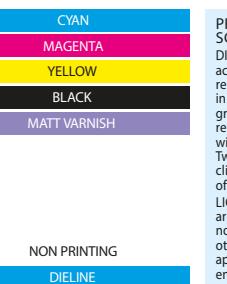

PRINTER TO TEST AND ENSURE THE SCANNABILITY OF ANY SUPPLIED BARCODE.
 DISCLAIMER: Tweak Design endeavours to ensure the quality of the artwork it creates. It is the client's responsibility to ensure that this digital artwork is correct in all areas. Please ensure to double-check spelling, grammar, layout and design. The client is responsible for final checking and approval and compliance with all appropriate regulations and client requirements. Tweak Design is not liable for any damage and/or loss to any client or third party, or any claims arising from, as a result of any errors or omissions in this digital artwork.
 LICENCE TERMS: Imagery and fonts used in this artwork are the property of their respective owners and must not be used to derive additional works or used in any way other than originally intended. For any new works the applicable licences must be acquired. Tweak assume the end user will adhere to applicable copyright laws.

WARNING OVERPRINT PREVIEW IS OFF.
 It is vital for accurate viewing to turn on overprint preview, for a guide to turning this on go to tweak.com.au/acrobat

0 10mm 50mm 100mm 150mm

XXXXXX