


Colgate®  
PLAQLESS PRO  
Smart Sensing Sonic Toothbrush

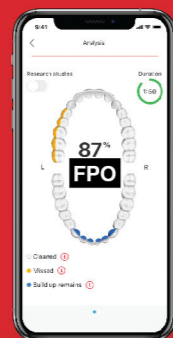

**PROVIDES REAL-TIME FEEDBACK TO HELP YOU IMPROVE YOUR ORAL HEALTH**

- Track your brushings over time with the Colgate® Connect app
- Brush on your own schedule and sync to the app later
- Coaches you to help improve your technique
- Quick convenient brushing guidance through Applewatch app

**KEY CONTENTS:**

2 brush head refills, travel case, & inductive charger

<http://colgate.com/connect/plaqlesspro>

\*vs. a manual toothbrush



The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and are used under license.



1-833-862-9251  
[www.colgate.com/plaqlesspro](http://www.colgate.com/plaqlesspro)



Dist. by:  
**COLGATE-PALMOLIVE COMPANY**  
New York, NY 10022 U.S.A.

PLAQLESS PRO  
Smart Sensing Sonic Toothbrush  
Makes the Invisible Visible

**HOW IT WORKS:**



**ADVANCED SMART SENSORS**  
Detect buildup on your teeth.

**INDICATOR LIGHT**  
Ring turns blue when buildup is present.

**TIMER**  
Quadrant timer illuminates every 30 seconds.

**INDICATOR LIGHT**  
Ring changes to white when area is done.

**MODE BUTTON**  
Press the mode button to switch to quiet mode.

PROVIDES REAL-TIME FEEDBACK TO HELP YOU IMPROVE YOUR ORAL HEALTH

Colgate®

PLAQLESS PRO  
Smart Sensing Sonic Toothbrush

Makes the Invisible Visible

Helps find & remove buildup in real-time. For a cleaner, healthier mouth.\*




\*vs. a manual toothbrush

Kit includes:




**IDENTIFY BUILDUP**  
Locates & identifies buildup in real-time, so you can brush it away.

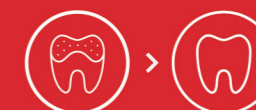

**KNOW WHEN YOU'RE CLEAN**  
Light ring changes from blue to white when area is clean.



**CHECK IN ANYTIME**  
Automatically syncs with the Colgate® Connect app.

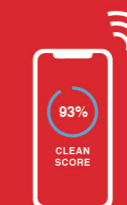

**TRACK & IMPROVE**  
Watch your Clean Score improve & earn rewards in the Colgate® Connect app.

Colgate®


# PLAQLESS PRO

SMART SENSING  
SONIC TOOTHBRUSH

PLAQ  
JUST MET  
ITS MATCH.

LET'S GET  
STARTED.


## LEARN ABOUT THE BRUSH

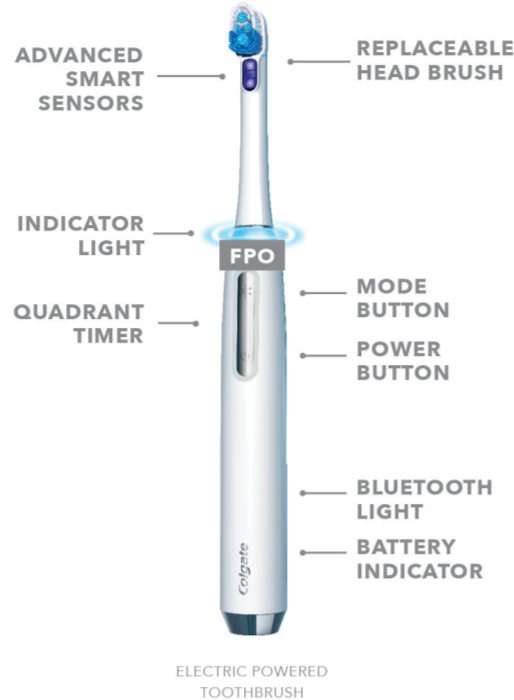

## KIT INCLUDES:



2 BRUSH  
HEAD REFILLS

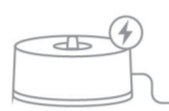


INDUCTIVE  
CHARGER



1 TRAVEL CASE

## INSTRUCTIONS

### FOR BEST RESULTS:

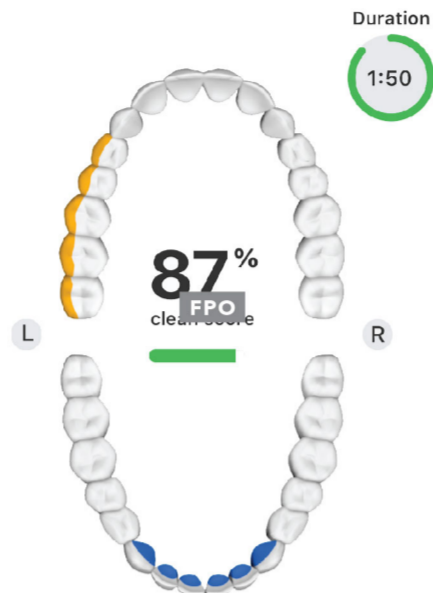

1. Start with your front teeth and work your way back towards your molars as you brush.



2. Brush in circular motion on each tooth. Continue brushing in that zone until the indicator light tells you to move on.



3. Hold the brush handle horizontally. When the handle is angled, the sensors may not be able to read your teeth.

## APP FUNCTIONALITY:

The Colgate® Connect app has additional features to enhance your experience with the Plaqless Pro toothbrush.

This includes:

- Video tutorials
- Personalized coaching in real-time
- Facts and history to help you improve your brushing over time

## CONNECTING YOUR BRUSH TO THE APP:

1. Launch the app.
2. Ensure that your phone or tablet has Bluetooth turned on.
3. Begin app introduction.
4. Turn on brush to pair, it is automatically in pairing mode. The Bluetooth® light will blink blue when paired.

## QUADRANT MODE:



Brushed for  
30 seconds



Brushed for  
1 minute



Brushed for  
90 seconds



Brushed for  
2 minutes

## KEY BRUSH INDICATORS:



Brush neck ring is green while initializing. Once initialized, the brush will be in feedback mode by default.

## FEEDBACK MODE: PROVIDES FEEDBACK ON BUILDUP

Full cleaning signal feedback: Brush motor lapses after 2 min and quadrant timer will be all white, Brush turns off after 4 min or when you press the power button.



**Blue light:** *Detecting*  
Continue brushing in zone.



**White light:** *Move on*  
Area is clean or you have brushed for enough time in zone.



**Red flashing light:** *Error: Correct angle*  
Sensor cannot see teeth. Rinse to clean sensor or hold the brush horizontally to resume detection.



**Light off:** *Brush not in mouth*  
When the brush detects it is outside the mouth, the light ring will turn off and the motor speed will reduce.

## QUIET MODE:

In quiet mode, the indicator light ring will not turn on to provide feedback as you brush. You will still have the 30 second quadrant timer to keep your pace.



The brush will automatically turn off after 2 minutes while in quiet mode. To change between feedback and quiet mode, press the mode button on the brush.

## REPLACING THE BRUSH HEAD:



1. To remove, pull the brush head straight up.



2. Position the new brush head so that the bristles are aligned with the quadrant timer. Push the brush head all the way down.



3. Your toothbrush is now ready to use again! After every use, rinse toothbrush head thoroughly, shake off excess water and store upright.

## BATTERY INDICATORS:

Battery indicator is green when battery is charged, orange when battery is low. Charge when indicator turns red. Battery indicator will flash green while charging and solid green when complete.

## BLUETOOTH INDICATORS:

Bluetooth indicator is green when connected and blue during firmware update.

## CHARGER:

To coil or uncoil the cord from the charger base, push on the post to pop the cord spool from the outer housing. Adjust cord length as desired and reinsert spool into housing lining up cord exit with groove. Plug USB connector into provided AC adapter & place toothbrush on charger to charge.

## IMPORTANT SAFEGUARDS

### READ ALL INSTRUCTIONS BEFORE USE

When using the electrical products, especially when children are present, basic safety precautions should always be followed including the following:

## DANGER

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not touch the AC adapter with wet hands.
- AC adaptors and chargers should always be completely dry when connecting to power.

## WARNING

- Adult supervision required for persons with disabilities. Replace refill if brush head becomes loose due to damage.
- Not intended for children. Keep out of reach of children 3 years of age and younger.
- The supply cord cannot be replaced. If the cord is damaged, the appliance must be scrapped.

## SAVE THESE INSTRUCTIONS

- When charging in a bathroom, lavatory, or other wet location, use only the AC adapter and charger provided. Other USB power sources must only be used in dry locations away from water.
- Stop using this toothbrush if it is damaged in any way (cord, plug, brush head, toothbrush handle, charger and/or AC adaptor). Return the product to a service center for examination and repair. Go to [www.colgate.com/connect](http://www.colgate.com/connect) to review how to handle a product for examination and repair.

- Use this product only for its intended use as described in this manual. Do not use brushheads, charger, or attachments not recommended by Colgate.
- Keep the cord away from heated surfaces.
- This product is only to be used with the power supply unit provided with it.
- Use the charger only within rated capacity of outlets or voltages specified.
- No user serviceable parts. Do not disassemble the handle as it may cause personal injury or device failure.
- If after using this toothbrush you experience excessive bleeding or if bleeding continues to occur after a week of use, please contact your dentist as it may be an indication of a more serious oral care issue.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

## NOTICE:

For energy savings purposes, charge only when the battery level is low. When not in use, keep the charger unplugged from the socket.

## ELECTROMAGNETIC COMPATIBILITY:

The Colgate® Plaqless electric toothbrush complies with all standards regarding electromagnetic fields. It is safe to use if handled properly according to the instructions provided in this manual.

## STORAGE INSTRUCTIONS AND DEVICE MAINTENANCE:

The following maintenance instructions are important to ensure that your brush continues to work as it was designed. Failure to follow these instructions may cause your brush to stop working as desired or to stop functioning altogether.

- It is recommended that you change your brush head every three months to ensure that the effectiveness of your brushing is maintained. If you do plan to stop using your toothbrush for an extended period, we recommend that you unplug the charging base, clean the toothbrush handle and heads, and store them in a cool, dark place.
- If you start using your toothbrush again after a long period of not using it, place it on the charger and, as soon as the level of battery charge allows, pair it with your phone.
- Clean your toothbrush thoroughly after each use by rinsing it under running water. Regularly clean the handle and the charging base of the Colgate® Plaqless Pro toothbrush using a damp cloth. This is especially important if your toothpaste contains baking soda or other bicarbonates or peroxide, otherwise, plastic cracking could occur.

## TROUBLESHOOTING:

- In the event the device fails to perform as intended, please contact Colgate-Palmolive Company at 1-833-862-9251 or [www.colgate.com/plaqlesspro](http://www.colgate.com/plaqlesspro)

## PRODUCT DISPOSAL:

- Please dispose of the device properly in accordance with local regulations.
- Do not dispose of battery in fire.

## WARRANTY:

- You can view the warranty of your purchased product by visiting our website at <https://shop.colgate.com/pages/warranty>.

## SPECIFICATIONS:

- Battery: NiMH
- Use and store at room temperature
- Power Supply Unit Input: 100-240VAC, 50/60HZ
- Handle Model PD1: Bluetooth® BT 5.0

## FCC STATEMENT:

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device must not be co-located or operating in conjunction with any other antenna or transmitter. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Distributed by  
Colgate-Palmolive Company  
New York, NY 10022  
U.S.A.

Call us at 1-833-862-9251  
Website: [www.colgate.com/plaqlesspro](http://www.colgate.com/plaqlesspro)  
Model Number: PD1

P93939390