


0. Preparations for use -Power on




1. Press the push button on the right side of the cushion to start the machine.

0. Preparations for use - Install apps

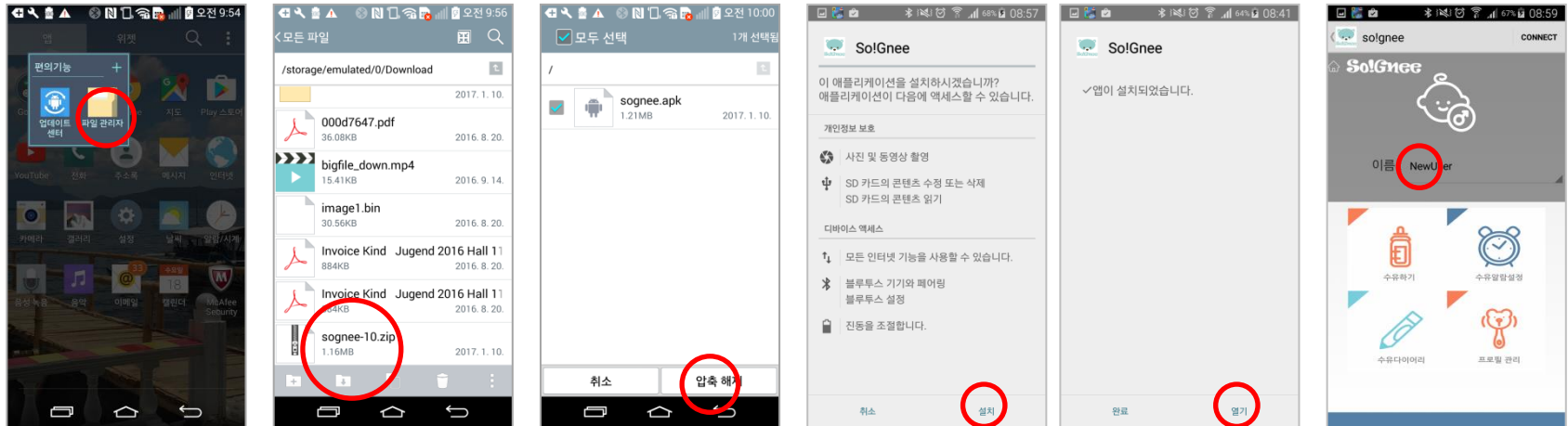


1. Select [My Files] → [Download] → [sognee] to download
cf: Download "Alzip" and "Uncompress"
2. Install downloaded apps
3. When you open an installed app, Default baby name will be "NewUser"

0. Preparations for use - Connect apps

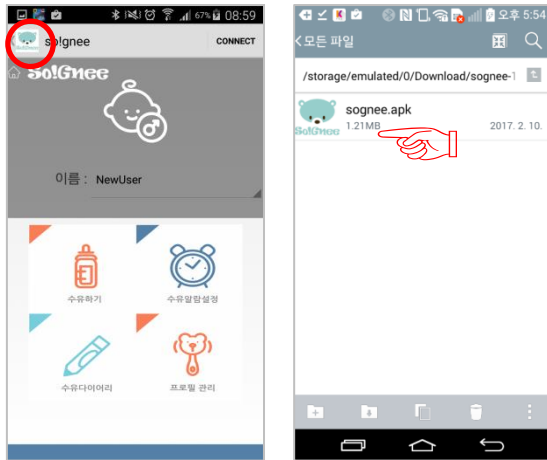

1. Touch the "Bear" logo to select the app you want to connect to.
2. Touch the app to connect.

1. Baby Registration

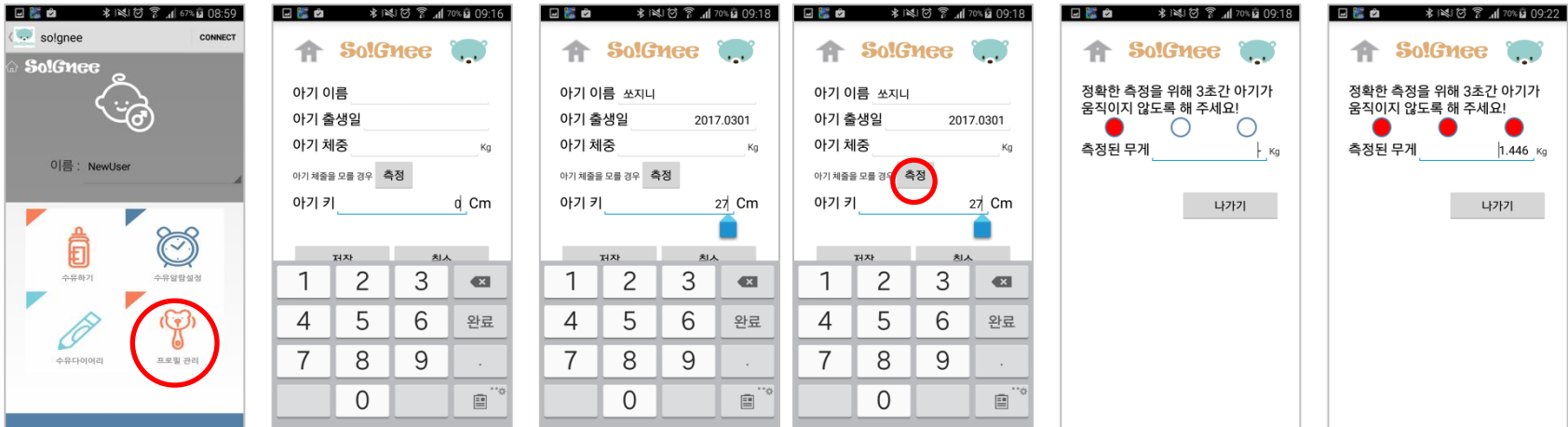

1. Baby Information Registration

Select [Profile], enter baby name, birthday, and height.


2. Weighing

Place the baby on the "Scale" and select the [Measure] button to measure baby's weight.

If you know your baby's weight, just fill it in.

3. Save the measured information. Be sure to click 'Save'.

2. Setting (automatic calculation mode / manual calculation mode)

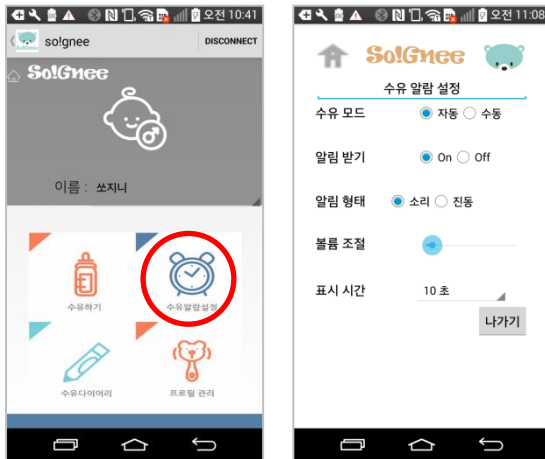



1. Set calculation mode

In [Settings], select between 'Automatic calculation' and 'Manual calculation'.

The default setting is 'Automatic calculation'


2-1. Setting (automatic calculation mode / manual calculation mode)

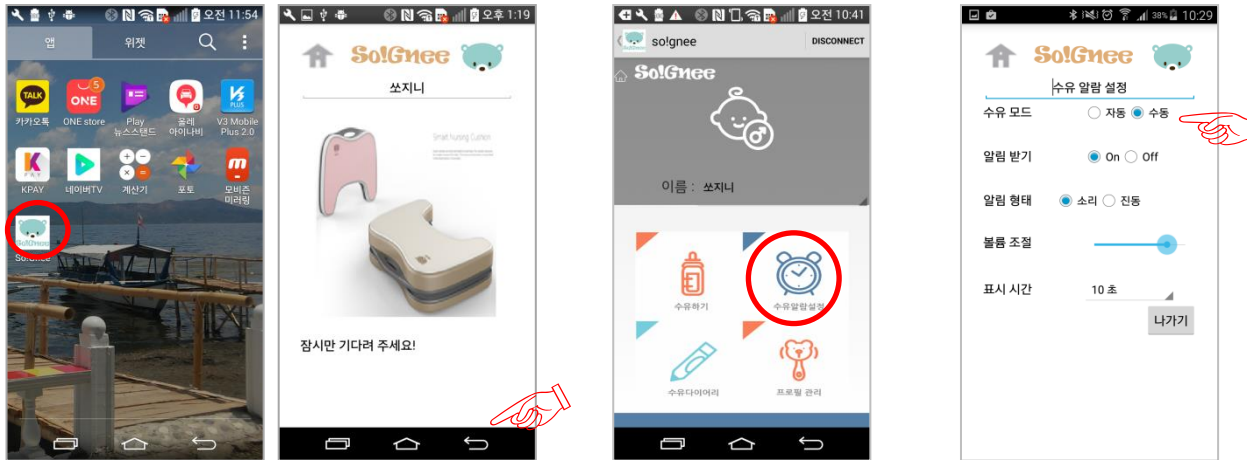


1. When you select the [So! Gnee] app, 'Auto' mode is set by default
2. Press 'BACK' () to go to [HOME].
3. Under Settings, select Manual.

3. Feeding measurement (Automatic calculation mode)




1. Select the [So! Gnee] app
2. Wait for a while, and if you see "Please do not let the baby move," put the baby on the feeding cushion.
3. If you see 'Breastfeeding', start breastfeeding.
4. When feeding is complete, press [PAUSE].
5. Check the results.
6. Pressing [Feeding complete] automatically saves the measurement results.
7. If you press  (return) without pressing [Feeding complete], the amount of feed will not be saved.


3-1. Additional feeding measurement (automatic calculation mode)

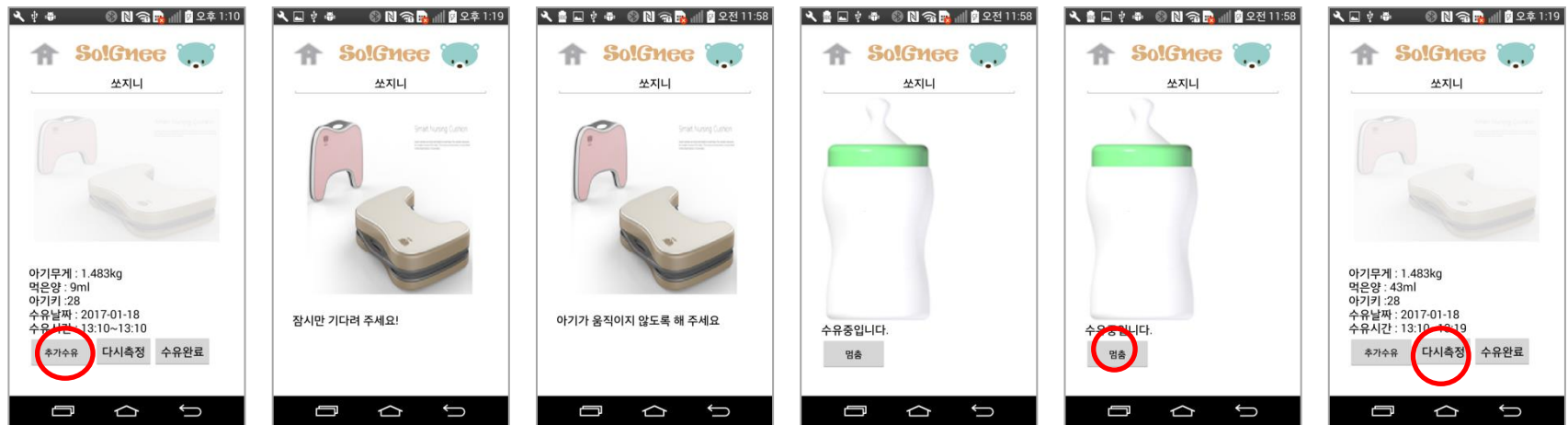


1. Get your baby off the cushion.
2. Select [Additional Feeding].
3. When feeding is complete, press [Pause] to see the results.
4. If the result is not correct, press [Re-measure] to re-measure.
5. Pressing [Feeding complete] saves the measurement result

* Note: If you click 'Return ()' without pressing [Feeding completed], the result of feeding measurement will not be saved.

3-1. Feeding measurement (manual calculation mode)

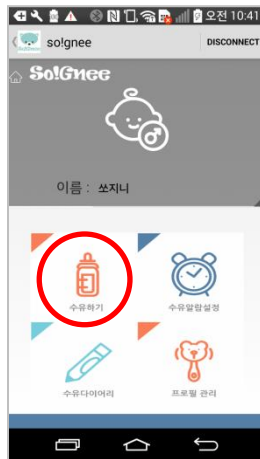

1. Select the [So! Gnee] app
2. Wait for a while, and if you see "Please do not let the baby move," put the baby on the feeding cushion
3. If you see 'Breastfeeding', start breastfeeding.
4. When feeding is complete, press [PAUSE].
5. Check the results.
6. Pressing [Feeding complete] automatically saves the measurement results.
7. If you press  (return) without pressing [Feeding complete], the amount of feed will not be saved.

FCC Information

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions :

- (1) This Device may not cause harmful interface, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for CLASS B digital device, pursuant to Part 15 of FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try correct the interference by one or more of the following measures:

- 1.1. Reorient or relocate the receiving antenna.
- 1.2. Increase the separation between the equipment and receiver.
- 1.3. Connect the equipment into an outlet on a circuit different from that to which receiver is connected.
- 1.4. Consult the dealer or experienced radio/TV technician for help.

WARNING

Changes or modifications not expressly approved by the manufacturer could void the user's authority to operate the equipment.

"CAUTION : Exposure to Radio Frequency Radiation.