

hapee®  
Easycare

Hapee Sonic Vibration Toothbrush  
Model: B017B/02

## User Manual

Hapee, keep you smile everyday!



### CAUTIONS

#### Prohibition


- This product is not suitable for physical, sensory or mental disability people, or people lack of experience and knowledge (include children) using alone.
- Do not modify, disassemble or repair this product by yourself. Please contact customer service if you have any problem with our products.
- In case of abnormality or malfunction, please stop using immediately. In this case, the use of this product may lead to accidental injury. If you have any question, please contact the customer service.
- Do not immerse the toothbrush in water or other liquids for a long time.
- Do not put the toothbrush and the charging cradle under the sun light. It should be stored in dry and ventilated environment.
- Do not place any metal objects such as coins, etc., 10cm around the charging cradle.

#### Warning

- Please do not keep it within the reach of children or infants, or it may cause accidents and injuries.
- This product is only suitable for cleaning teeth, gums and tongue. Please do not use it for other purposes. If you feel any discomfort after use, please discontinue using it.
- If you have a bleeding condition after use of this toothbrush or bleed for more than a week, please consult your dentist.
- People with severe periodontitis, who are being treated, or who are suffering from oral discomfort, should consult a dentist before use, or they may cause damage to the teeth and gums
- This product complies with the safety standard for electromagnetic devices, and any consultation with a pacemaker or other implantable device. Please consult a doctor before use.
- Do not use toothbrush refill head which recommended by other manufacturer.
- When charging, please keep away from strong electromagnetic field, such as induction cooker.

### SPECIFICATIONS

Charging	Wireless charger
Toothbrush rated voltage	3.7V
Toothbrush rated power	1.8W
Charging cradle rated input	DC5V/500mA
Fully charged	About 6 hours
Battery life	About 30 days

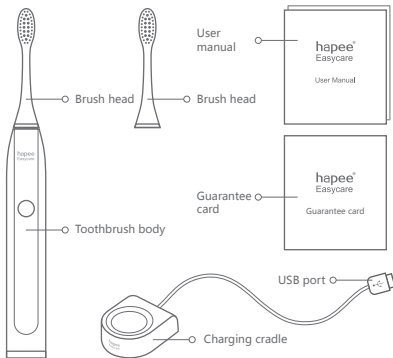

#### Packing list

Item	Quantity
Toothbrush body	1
Brush head	2
Charging cradle	1
User manual	1
Guarantee card	1

### FCC Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.  
This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:  
-- Reorient or relocate the receiving antenna.  
-- Increase the separation between the equipment and receiver.  
-- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.  
-- Consult the dealer or an experienced radio/TV technician for help.  
This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.



Zhejiang Hapee Oral Health Co., Corp.

Address:  
No. 1412 Zhouan Road, Jiaxing, Zhejiang Province  
Customer service number:  
400 088 5112 (only available in China)  
Website:  
www.hapee.com.cn

### USER INSTRUCTIONS

- For the first time use, remove the motor protection cover from the main body of the toothbrush.
- Put the brush head on the main body of the toothbrush carefully (there is some space between the brush head and the main body of the toothbrush, which is necessary for the normal shaking of the toothbrush head).
- Apply proper amount of toothpaste on the brush. (Tip: in order to prevent toothpaste from splashing, please put the brush in the mouth and then start the toothbrush)
- Press the start button to start brushing, and you can choose the suitable cleaning mode as needed.
- After the toothbrush is started, there is a short pause every 30 seconds to remind the user to switch the oral cleaning area. The toothbrush starts for 2 minutes and the toothbrush automatically stops working (user can restart the toothbrush to brush teeth for longer time)

### FUNCTION INSTRUCTIONS

#### Start-up button

Short-press the start button only once, and start the toothbrush. Change the brush cleaning mode every time you press the start button after started. Hold the button for 1 second, the toothbrush stops working.

#### Clean

Remove the food residue and tartar which can not be clean by ordinary toothbrush. Suitable for daily use of the general population.

#### White

Effectively remove stains, tobacco stains, coffee stains and other stubborn stains, make your teeth stay healthy and bright white. Suitable for use by a particular population.

#### Polish

Effectively remove stubborn stains, polishing teeth. Suitable for use by a particular population.

#### Soft

Gently removes dirt from the teeth and protects the gums. It is suitable for people with mild gingivitis or slight gum bleeding.

#### Sensitive


Gently cleaning teeth, while relieving the teeth sensitive. Suitable for sensitive teeth, such as elderly and pregnant women for daily use.

#### Bluetooth Indicator light

When a Bluetooth link is established with the mobile terminal app, the Bluetooth indicator lights up.

#### Power Indicator light

Enough power - not bright, low power - flashing red light

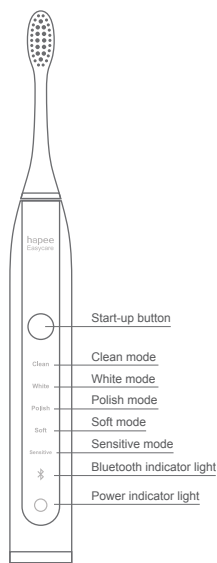

### CHARGING INSTRUCTIONS

- Breathing blinking: low battery
- Short blinking: in charging
- Always ON: fully charged

#### Low battery


During use, when the power indicator light blinks like breathing continuously, it indicates that the toothbrush is lack of electricity and needs recharging. When the battery is used up, the toothbrush cannot start and the power indicator will not turned on.

#### Battery charging

Place the main body of the toothbrush on the charging cradle and make sure the charging cradle is plugged in. The power indicator blinks shortly and continuously as charging.

#### Fully Charged

When the electric indicator continuously lights on, the charging is finished. (For better of usage the battery, please charge the battery for more than 12 hours at the first time charge)

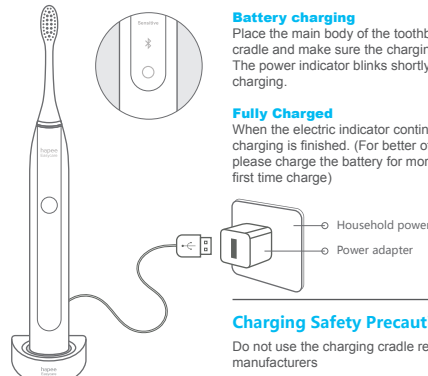

#### Charging Safety Precautions

Do not use the charging cradle recommended by other manufacturers

Do not use power adapter or other power supply that does not conform to the rated input specification of the charging cradle, otherwise it may cause fire or electric shock due to short circuit.

Do not use the wet hand to connect the power adapter to the home outlet, which may cause injury.

### TOOTHBRUSH MAINTENANCE



#### Cleaning brush head

Rinse the brush head with clean water every time after brushing teeth



#### Rinse the toothbrush body

Rinse with clean water and dry naturally. Do not immerse the toothbrush body in water for a long time. Do not use high temperature water or other corrosive liquids to clean the main body of the toothbrush.



#### Clean charging cradle

Clean the dirt with a soft cloth and clean the charging cradle before pulling the USB port of the charging cradle out of the power supply.



#### Change brush head

The dentist recommends changing the brush every three months. Please use the exclusive brush head specified in this product.