

Portable Smart Pillbox

USER MANUL



* Please, carefully reading the user manul before using your device;

* Intactly, keep your warranty card.

■ Admire Your Attention

Thank you for using our automatic pill dispenser Leboxmade by XiamenZayata Technology Inc. The product is only for medicine usage reminding. The medicine, medicine administration time and dosage must be according to the prescription.


■ Advanced Features

- Exquisite portable
- Intelligence operation
- Double alerts
- Automatic record

■ Contents

Name	Quantity
Pill Box	1 piece
Battery	2-AAA
User Manul	1 piece
Warranty Card	1 piece
Portable Bag	1 piece

■ Control & Parts

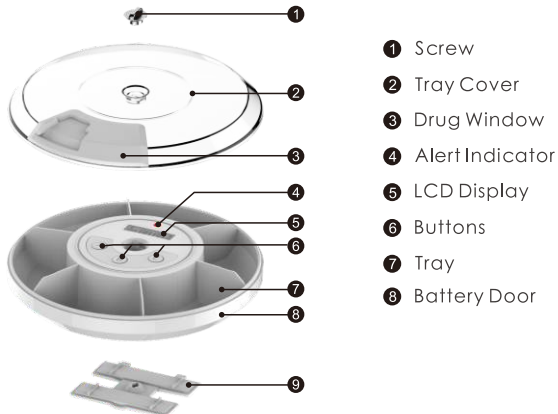

■ Insert the Batteries

1. Open the battery door, insert the batteries(2-AAA) as illustrated .



2. Close the battery door.

* Clean tray with a soft damp cloth or clean sponge only.

■ Connect Pillbox With APP

1 · Scan the QR code below through a browser on phone device.

· Visit the website on phone device:

“<http://www.zayata.com/dl/med10.php?> ” .

· Iphone user can search “Lebox” at App Store.





Scan & Download APP




Lebox

2. Click the software icon to enter the home page after installation.

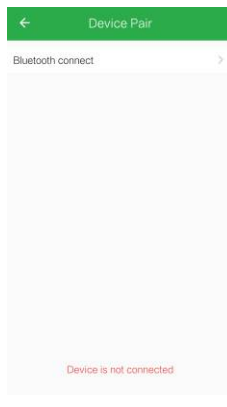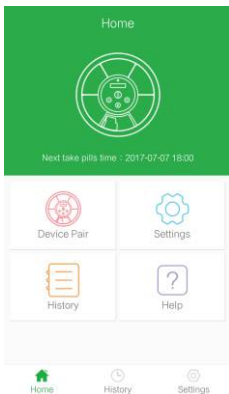

3. Click the item **【device matching】** to enter the next interface, choose one of the following operations.

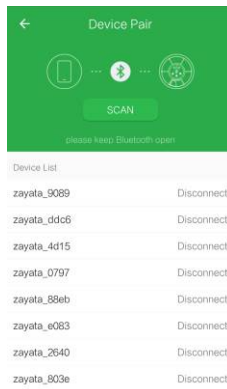


a . Chose the Bluetooth , click the **【scan】** option , chose your pill-box type to connect. Bluetooth icon shows on the pillbox screen, it means the App connected. Bluetooth icon disappears, it means the App disconnected.

b . Easy connection


Chose the **【easy connection】** , then scan the QR code to connect. Bluetooth icon shows on the pillbox screen, it means the App connected. Bluetooth icon disappears, it means the App disconnected.




■ Setting Alarm on Your Application

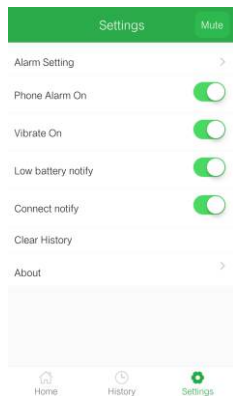


1. Click the return arrow to the home page, then click the item **【settings】** .

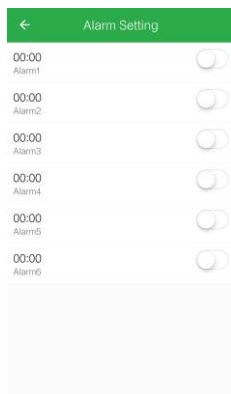

2. Click the **【alarm setting】** , setting the alarm by your requirement and conserve.

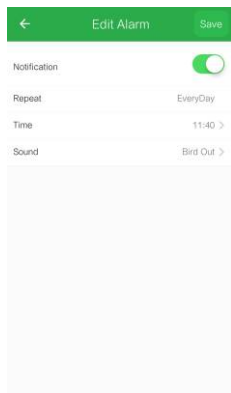

3. After setting all alarms, then back to set other parts.

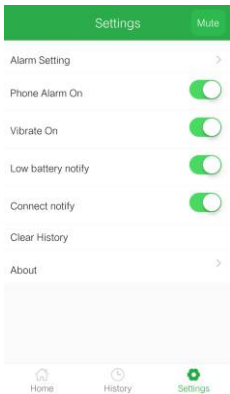

4. Without beep alarm
To stop the alerting noise
without purging your medi-
-cine out, click the item
【mute】.

■ Manual Adjustment

■ Displayed

● Real time

Can be settled both on
your application and
your device.



● Pill reminding

The alarm mark will be shown on the LCD display
after the pill reminding settled down. When the pill
time is up, the device will "beep" and lighting to
remind the user.

● Battery low power remind

When the mark "⊗" showed on the display, please
change the new batteries.

■ Open The Cover



Rotate the screw on top,
then open the cover.

■ Checking Medicine-taken Records



1. Click the item 【medica-
tion records】, it will show
your medication records
in real time and mark the
time when you missed the
dosage (the medicine
information should be fill-
ed when you setting the
alarm)

* the yellow circle means
you have filled the medi-
cine name on that date
on the application.

* the red dot means you
forgot to take your medi-
cine on that date;

* take your medicine 10
minutes later after the
alarm will be marked
【forget taking】.

■ Setting The Clock



Step 1:

Press and hold button A over
2 seconds until "—" stop fla-
shing.



Step 2:

Press button B to adjust hour
digits. Press button C to
adjust minute digits.



Step 3:

Confirm clock settings by pres-
sing Button A. The setting page
will log out if no operation in 10
seconds.

■ Setting & Checking Alarm Times



Step1:

Press button A, when the clock symbol turns up, press Button B(hour) and Button C (minute) setting the alarm clock.



If you wish to cancel the alarm, press Button A until "--:--" is displayed.



Step2:

Confirm your desired alarm time by pressing the Button A. Then "alarm2" will show up.



Step3:

Repeat steps 1-2 to set the Alarm 2,3,4,5 and 6.



Step4:

After confirming Alarm 6, press Button A back to the clock mode. (it will also automatic confirm and return to the clock mode if you have no operation within 15 seconds).

■ Testing Your Dispenser



1. Press and hold Button B for about 5 seconds: screen all displayed, the indicator light on ,the buzzer rings.
2. Hold your dispenser do the down medicine gesture: indicator light on and 1 will be displayed.
3. Press Button A: indicator light on and 2 will be displayed
4. Shake your dispenser: indicator light on and 3 will be displayed.
5. Rotate the cover by trays ,the LED light on by times and 4 will be displayed.
6. Press Button C to exit test mode and return to clock mode.

■ Locking Your Dispenser



Place the cover on bottom of unit, rotate the cover until totally fit with the tray, then insert the screw and turn counterclockwise.

■ Dispense Pills & Stop The Alarm



To stop the alarm and dispense the pills, rotate the drug window to the tray, without this act the alarm clock and the indicator will keep alerting continually 5 minutes.

■ Low Battery Alert & Standby Wakeup



A battery mark will show up to indicate the batteries are low and should be replaced to ensure proper function. To save power, Lebox will stand by and shut down if no operation in 30 seconds. Touch it to wake up.

Warranty Card

Dear customers:

Thank you for trust our product, please carefully read and retain our warranty card and show it when your dispenser need maintenance.

Repair date	Repair Content	Support staff

■ About warranty:

- 1 .The products are guaranteed to be free from defects in materials and workmanship under normal and intended use for a period of 1 year from the date of purchase.
- 2 .Under the following situation, we will ask a charge for helping you fixing the products: in case of damage in transit, inadequate care or neglect, abuse, abnormal use, misuse, accidents, normal wear and tear, damage due to environmental or natural elements, failure to follow product instructions, immersion in water, battery leakage or improper installation, improper storage or maintenance of the products, or improper service or cause not arising out of defects in materials and workmanship.

■ FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Warning:

§ 15.21 Information to user.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

§ 15.105 Information to the user.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions,

may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

* RF warning for Portable device:

The device has been evaluated to meet general RF exposure requirement.
The device can be used in portable exposure condition without restriction.



Scan & download App