



User Manual

YAKOOK

SMART MEDICINE CHECKER

Quick Guide

01. Contents of the Package

Contents of the Package



Yakook




Battery
(CR2477)



Quick Guide


How to Attach it

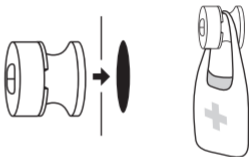

- ▶ Use double-sided tape on the bottom of the product to attach the product to a place where it is easy to hang
- ▶ Hang the medicine bag in the hollow in the center of the product.


02. Power Control & Battery Change

How to turn on the power

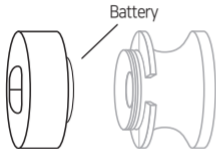

The battery is already inserted when purchasing the product, so you can use it right after assembling the product.

How to change the battery

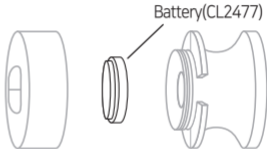

When replacing the battery of Yakook, disassemble the top and bottom plates, remove the used one, and insert a new one.

03. App Installation

Medication checker only works with its exclusive app.
Take a picture of the QR code below or check out the app store
or play store to find the exclusive Yakook app.

YAKOOK

<http://yakook.co.kr/en/app>



- ▶ Follow the instructions provided by the app for the connection between Yakook app and Yakook button.

04. Cautions

Refer to the safety precautions and the product guide for using on Yakook website to use this product safely.

Safety Precautions.

<http://yakook.co.kr/en/guide>



05. Customer Support / Product Specifications

If you experience any problems using this product, contact customer support by emailing appinfo@medisolution.co.kr

The limited warranty period for this product is one year from the date of purchase. If the product is defective or a functional error occurs during this period, you can be exchanged the product for free or get its firmware updated.

| | | | |
|----------------------|------------------|-------------------------|-----------------|
| Connectivity | BLE 4.2 | Working Temp | 20 to 120°F |
| Optimal Range | Up to 30 ft. | Battery. | CR2477 (3v) |
| Dimensions | 2.2 x 2.0" | Certification | R-R-Mdp-YAK-001 |
| Weight | 8.8 oz | FCC ID | XXXXXX-XXXXX |
| Frequency | 2402 to 2480 MHz | Frequency output | -2 dBm ± 1 dB |

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

CAUTION : Any Changes or modifications not expressly approved by the manufacturer could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- **Reorient or relocate the receiving antenna.**
- **Increase the separation between the equipment and receiver.**
- **Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.**
- **Consult the dealer or an experienced radio/TV technician for help.**

A minimum separation distance of 20 cm must be maintained between the antenna and the person for this appliance to satisfy the RF exposure requirements.



YAKOOK

SMART MEDICINE CHECKER

Main Features

05. Main Features



Button Test

If you press the button when it is not the time to notify you of the medication, the green indicator flashes once.

Flash once per click



No Device Added

If you press the button of Yakook that is not added on the app, the red indicator flashes once.

Flash once per click



No Notification Scheduled

Enter the medication time on the app. The green indicator flashes every 10 seconds for 5 minutes from the medication time.

Flash every 10 min.
With no schedules



※ Synchronizing the app and the device is recommended.

05. Main Features



Medication Notification

Enter the medication time on the app. The green indicator flashes every 10 seconds for 5 minutes from the medication time.

Flash every 10 sec.



Medication Input

Take medicine at the notification and press the button to record the medication. The green indicator flashes three times.

Flash three times with medication input



Duplicate Medication Warning

If you press the button within five minutes again after the medication, the red indicator flashes three times.

Flash three times for duplicated warning



05. Main Features



Medication Reminding

If you do not press the medication button despite the medication reminder, the red indicator flashes every 10 seconds for an hour.

Flash every 10 sec.



※ If a reminder for another medication is scheduled within an hour, it only flashes until the next medication time.

Power turned on

The green indicator flashes for three seconds when the battery is inserted.


Lights on for 3 sec.



Low Battery

When the battery remains at 20%, the green and red indicators flash alternately once every 10 minutes.

Green & red flashes

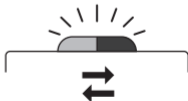



**MEDI PLUS
SOLUTION**

MediPlus Solution Co., Ltd.

2nd Fl, 39 Seongbuk-ro, Seongbuk-gu, Seoul 02880, Korea