210 mm

Intelligent Wireless Fitness Apparatus

Thank you for purchasing Intelligent wireless fitness apparatus.


This product is designed for training your muscles via electrical stimulation(it


To avoid causing harm and loss due to user's improper use, please refer to the "safety precaution "and then use the unit after fully understand.



is not a medical equipment).

Please keep this "instruction manual". Content Safety Precautions P. 2 Parts Identifications P. 6 Precautions for Use P. 7 Installing the Battery(Remote control) P. 8 Charging(Main machine) P. 8 Usage of Electrode Gel Pads P. 9 P. 11 Basic Usage P. 15 Storage Maintenance P. 16 Q&A P. 17 **Trouble Shooting** P. 19 Product Specifications and Usage Conditions P. 21 Consumables(Sold Separately) P. 22 Warranty Card P. 23

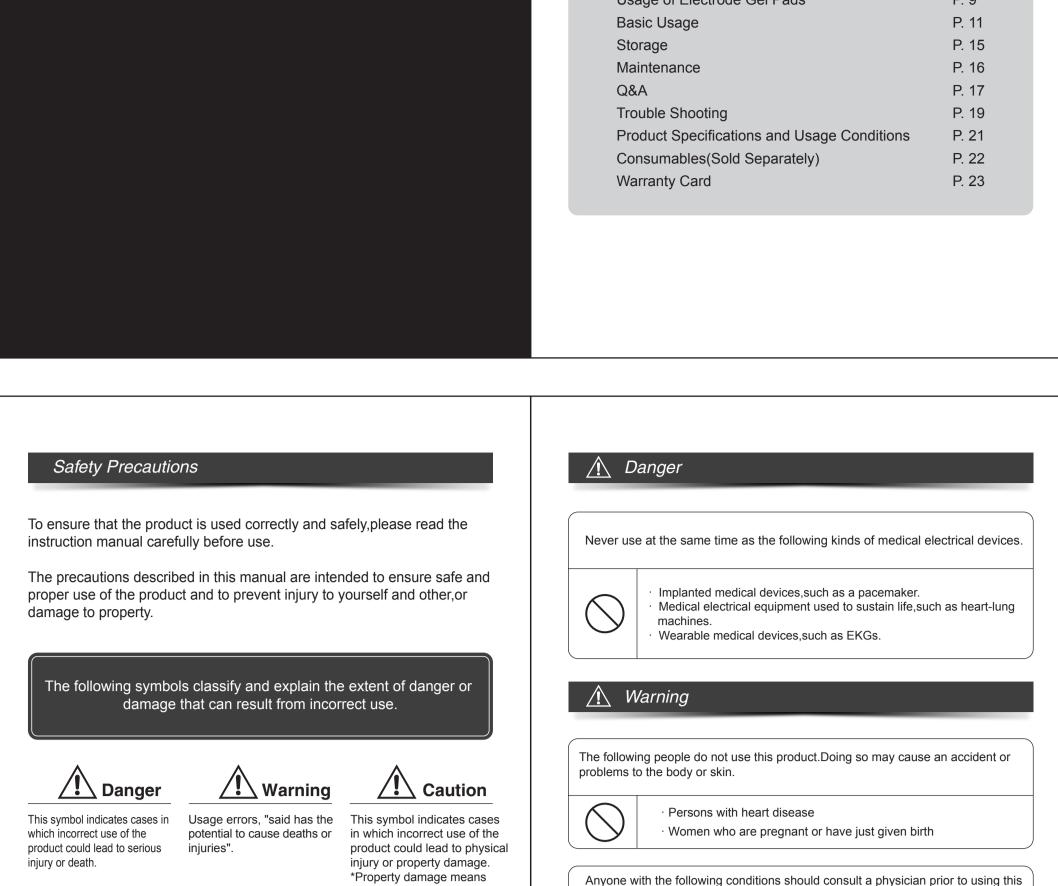

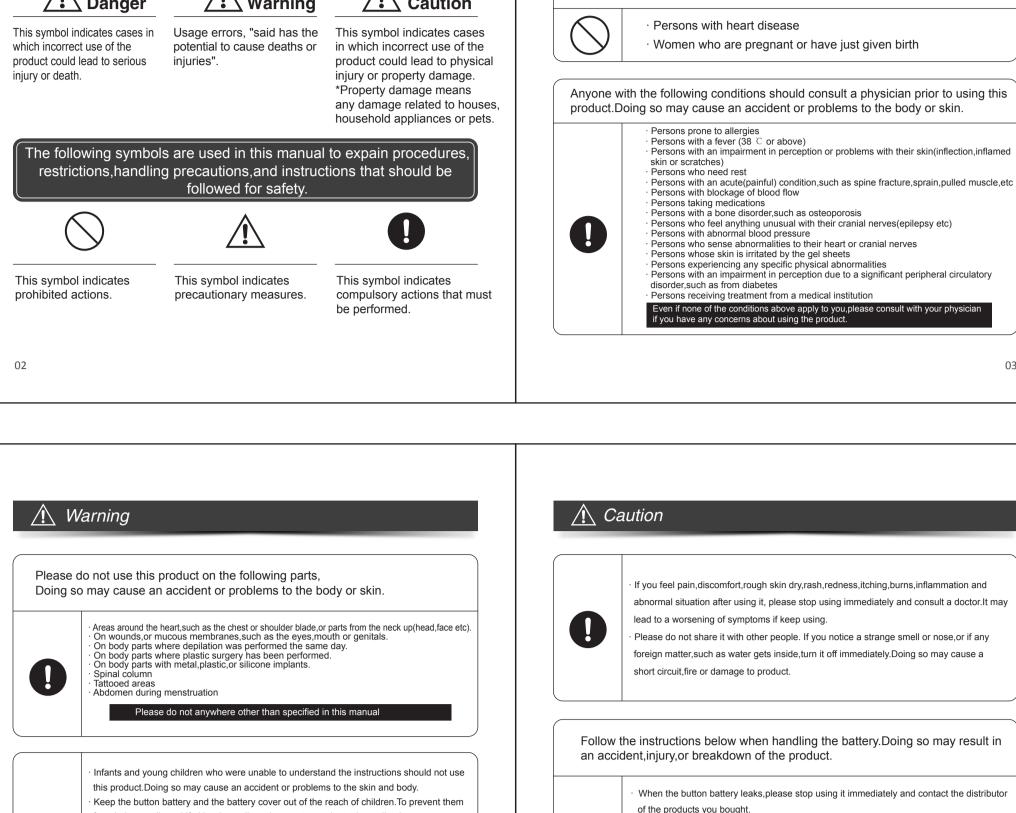

If liquid from the button battery gets in your eyes,it may cause damage or blindness .Do

not rub them immediately, rinse them with plenty of clean water for at least 1 minute and

If you get liquid from the battery on your skin or clothes,it may harm your skin.Rinse your

When installing the button battery, be sure to put it in facing the right way. When the button battery runs out, please take it out immediately.

Do not bend or fold the gel sheet, otherwise the gel sheet may be damaged.

Stick the gel sheet to your skin,if there is a gap where it is stuck on,electricity will not flow

uniformly across the entire surface of the gel sheet and you may feel a strong stimulus.

· Always adhere the gel sheet to the electrodes when using the product .Do not touch the

skin immediately with plenty of clean water.

Follow the instructions below when handling gel sheets.

· Please do not use damaged, deteriorated, dirty gel sheet.

Turn off the product before removing the gel sheet.

electrodes directly when the power is on.

Doing so may result in burns.

from being swallowed. If either is swallowed, contact your doctor immediately.

Follow the precautions below when using the product. Doing so may result in an

In such case, continuing to use it may result in an accident or injury.

Do not touch the metal electrodes. Doing so may result in shock.

Do not heat the product itself or the gel sheet in a microwave.

may make you feel ill or nauseous.

accident, electric shock, injury, damage or getting it dirty.

Do not use the product on young children or pets.

· Do not use it while exercising or driving an automobile.

Do not use in a humid location, such as in the bathroom.

Do not let children or pets play with the product.

Do not spill any liquids on the product.

The electrode button is aligned with the electrode hole of main machine.

4. Attach the pad with main machine to your body.

Peeling Method of Electrode Pad

When removing the gel pads, please turn off it, it may cause pain or burns.

2. Put it back to the storage sheet, then put it into the storage bag for proper storage for

1. Please slowly take the gel pad off after using.

the next use.

10


3. Remove the film on the gel pad.


Do not disassemble, modify or attempt repairs. Doing so may result in an accident or injury.

Do not subject to shocks, such as by dropping or hitting the product. Doing so may damage it

Do not use in combination with any other electrical device or beauty equipment. Doing so

Do not use for any purpose other than as intended and according to the instructions in the

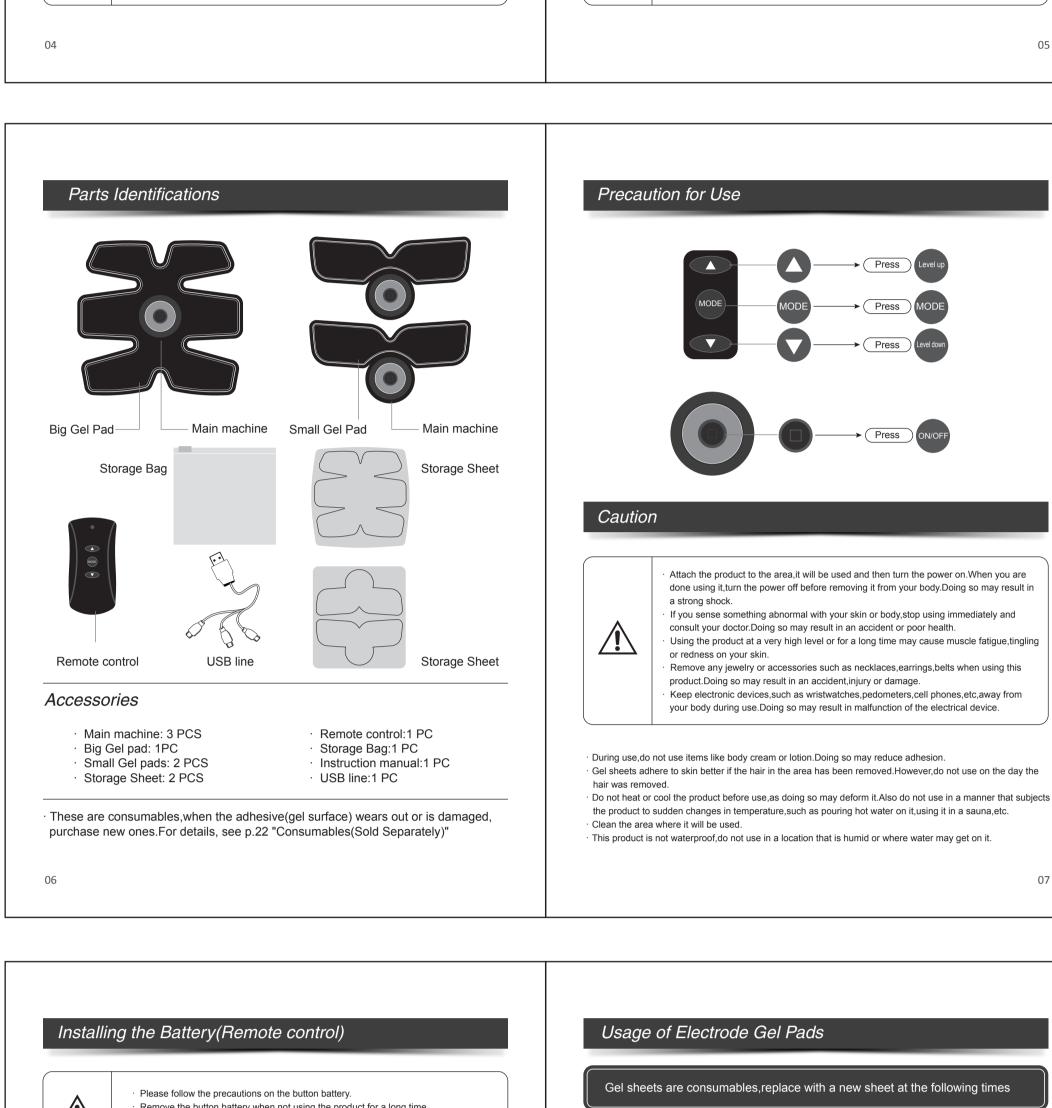

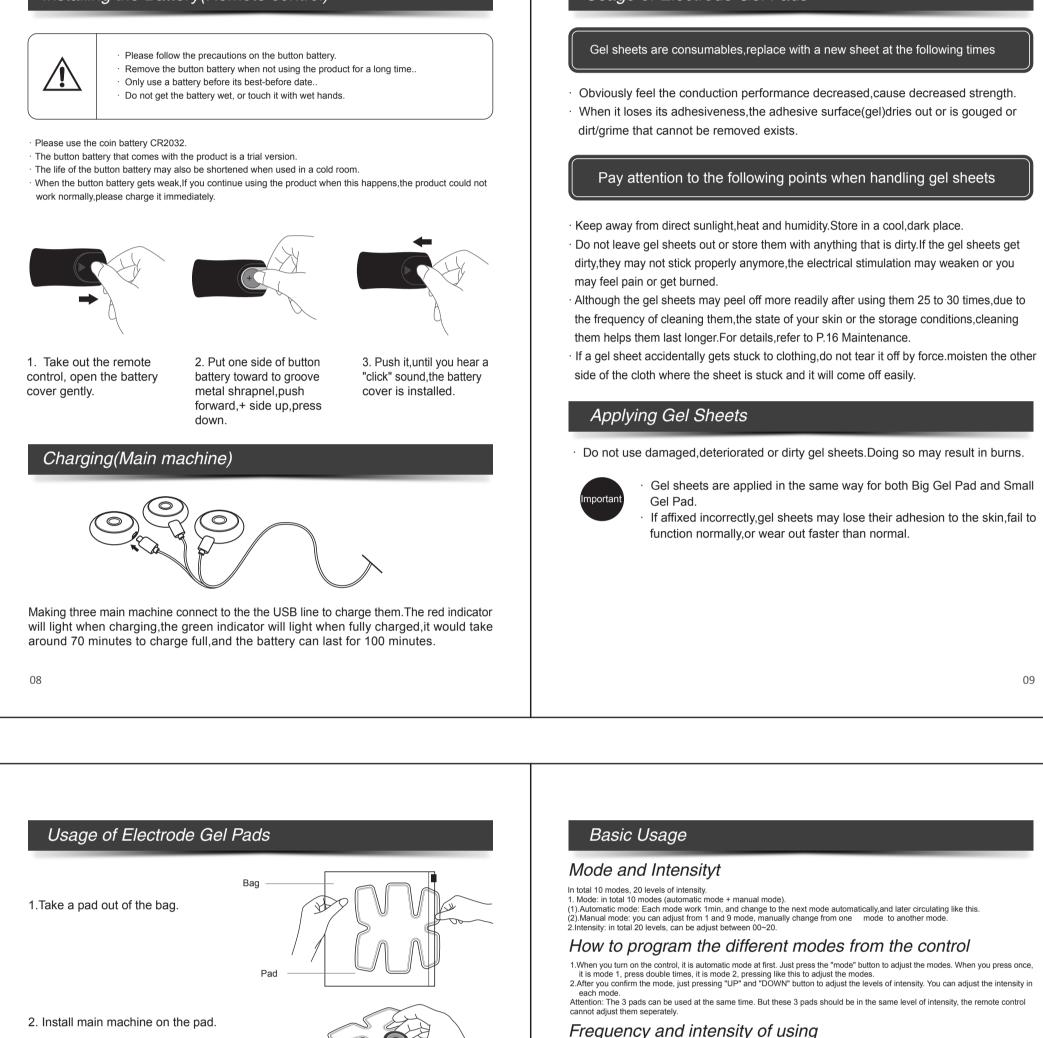

Use once per day per area used, use a weak level when using it for the first time, once you are

Stick the gel sheet to your skin, if there is a gap where it is stuck on, electricity will not flow uniformly across the entire surface of the gel sheet and you may feel a strong stimulus or pain

If you can't feel the stimulation readily, reposition the product 1 to 2 cm at a time up

or down. When repositioning the product, turn off the power.

Waist

11


used to a weak level, adjust the level as desired.


Wipe any sweat off the area to be used before attaching. Touch the adhesive (gel surface) as little as possible.

1. Attach the product to the area where you wish to use it.

Abdomen

The hydrogel is attached to the following sketch map:

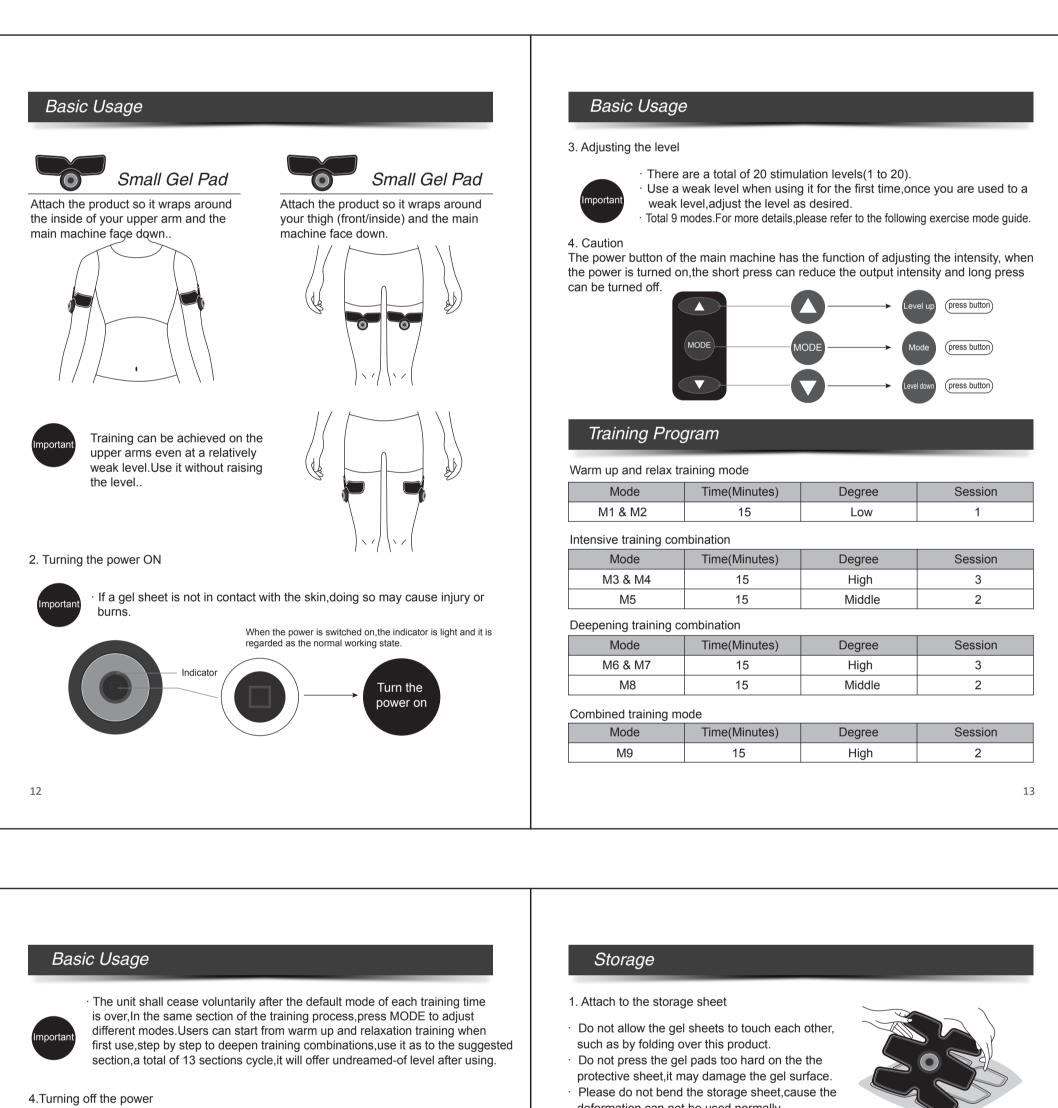

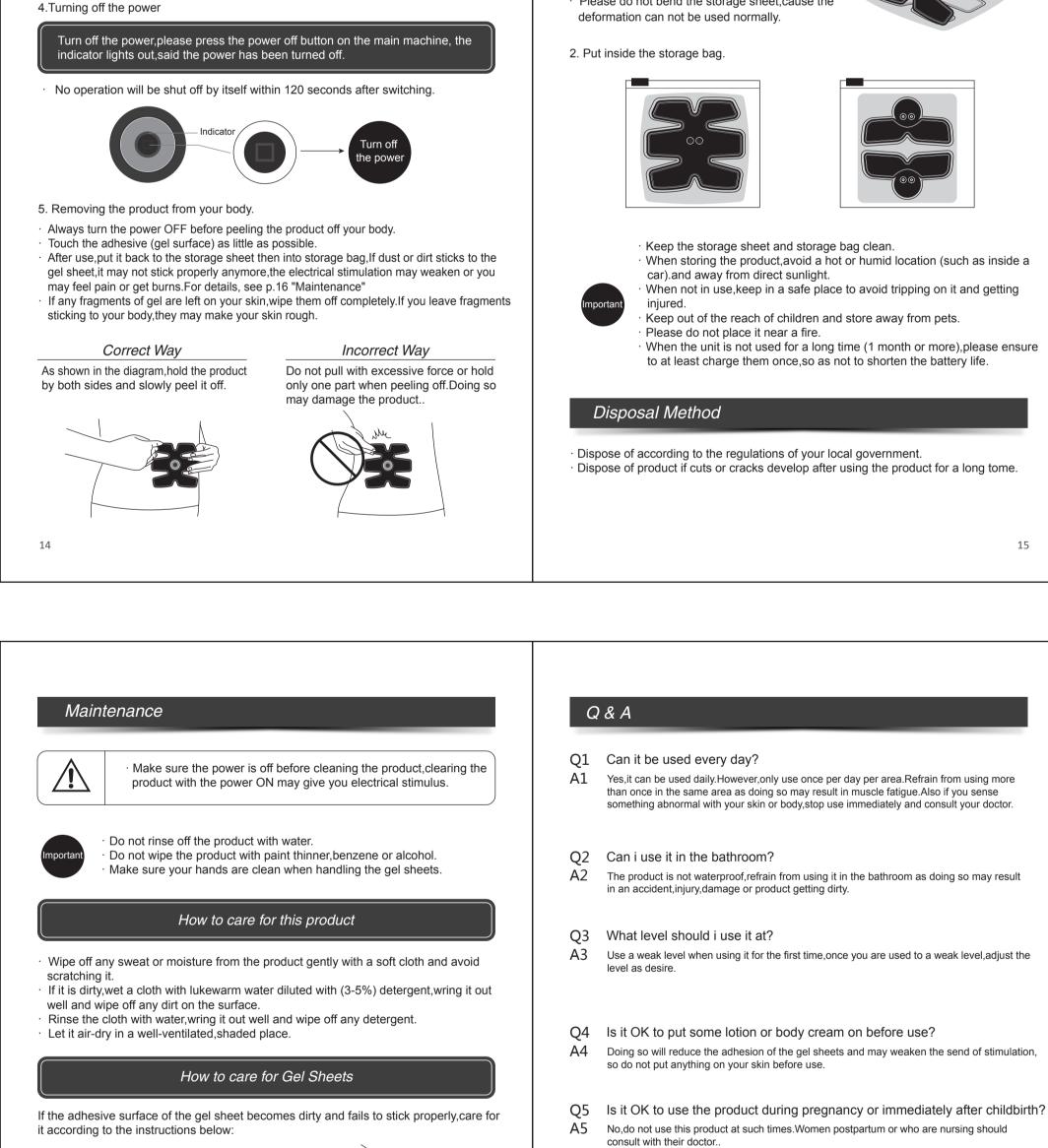

Q6

Trouble Shooting

What part of the body and circumstances can not be used?

Do not use this product for parts of the neck (head, face, etc.), spine, chest and shoulder

blades and other parts around the heart, In addition, Please do not use this product who

17

those with heart disease. For more details, please refer to P.2 Safety precautions.

In the interest of product improvement, these specifications are subject to change without

· Get some drops of water on your fingertips and rub the adhesive surface (gel surface) and remove any dirt. · Shake off as much water as you can and leave it

shaded location.

16

Q & A


adhesive (gel side) face up to air-dry in a well-ventilated,

adhesive surface.(gel side) and shorten its life.

Do not clean it for a long time or too frequently. Doing so may damage the

· Allow to air-dry until there is no moisture on the surface

and then affix the product to the storage sheet.

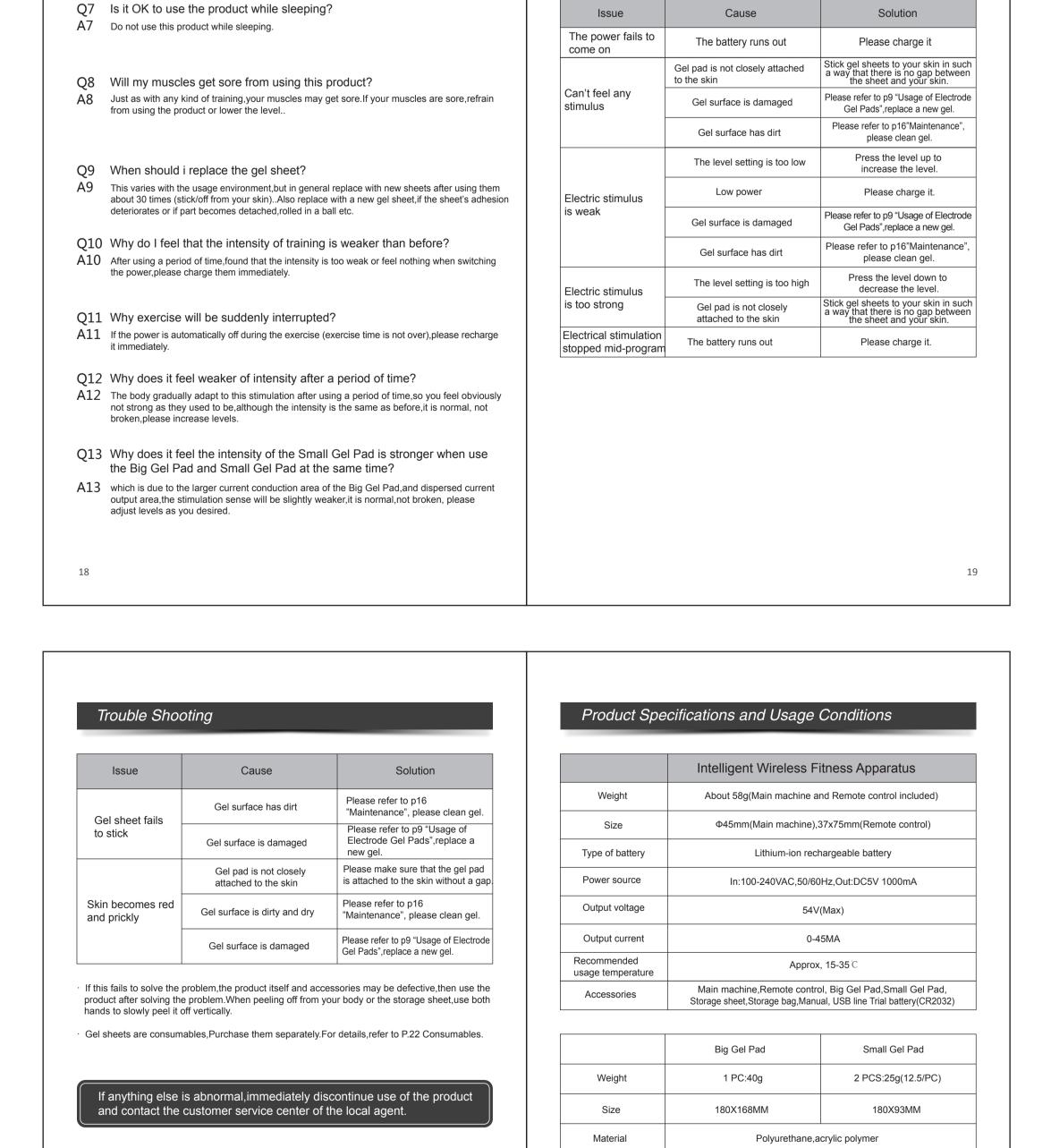

Consumables (Sold Separately)	Warranty Card	
Gel sheets are consumables,when the adhesive (gel surface) wears out or is damaged, purchase new ones.	Valid for warranty 1 year after the date of pure	hase
The following problems may arise if anything other than genuine gel sheets are used:	Name	Model No.
Failure of the product to perform as intended. Partial damage or breakdown of the product	Purchase Date:	
· Sensation of a strong shock along with burns or rough skin. Please note that we shall bear no responsibility whatsoever for any problems that result	Consumer's Add:	
from using anything other than genuine gel sheets. To purchase genuine gel sheets,consult the customer service station of the local agent.	Consumer's Name:	
	Consumer's Tel:	
Bin not need Constituted need	Dealer Name:	
Big gel pad Small gel pad	Dealer Add:	
	Dealer Tel:	
Gel sheets		
22		
FCC Statement		

22	23
FCC Statement	
This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.	
Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.	
NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: 1) Reorient or relocate the receiving antenna. 2) Reorient or relocate the receiving antenna. 3) Reorient or relocate the receiving antenna. 4) Consult the dealer or an experienced radio/TV technician for help important announcement.	
The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.	