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BP-003RX / BP-003TX

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Daiper Alarm Devices

BP-003RX / BP-003TX

**Alarm Devices for Dedicated Daiper
Dedicated Diaper; Main Device & Sub-Device**

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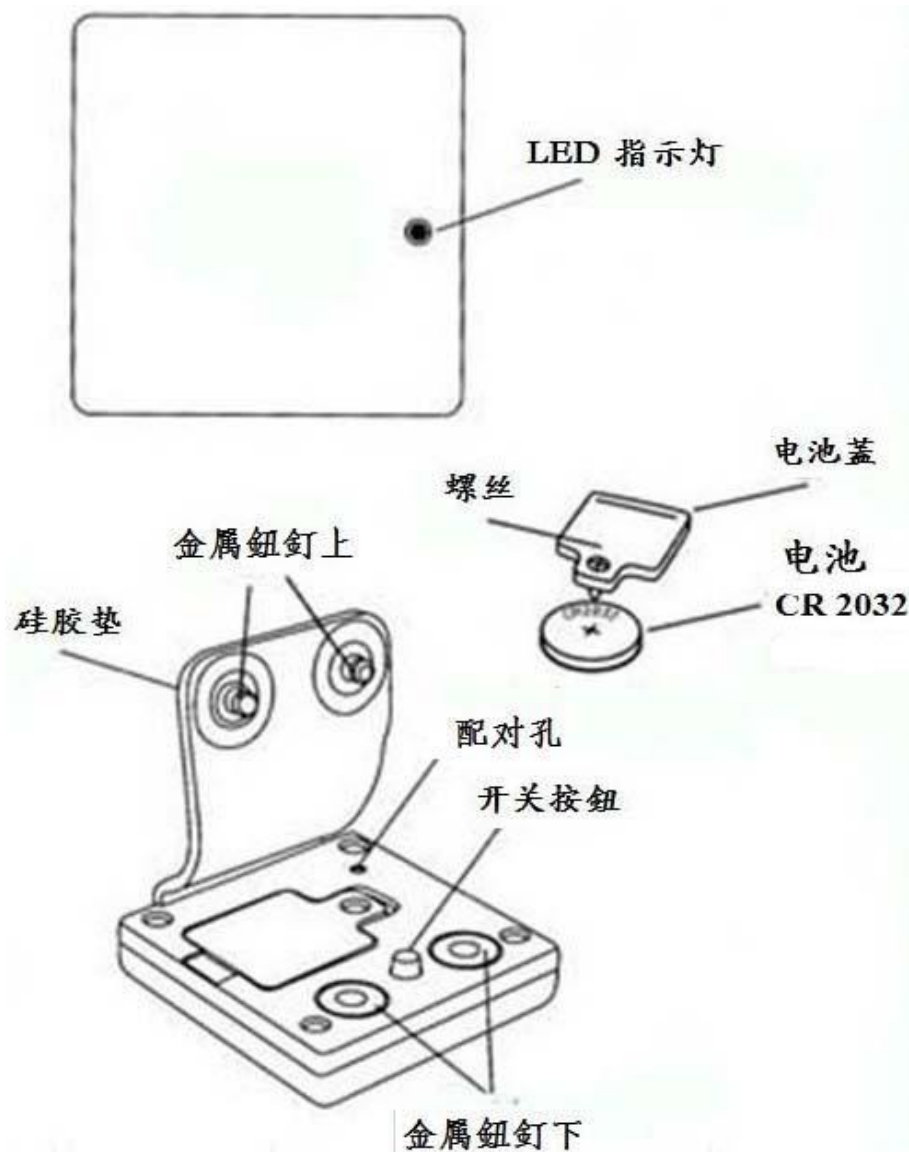
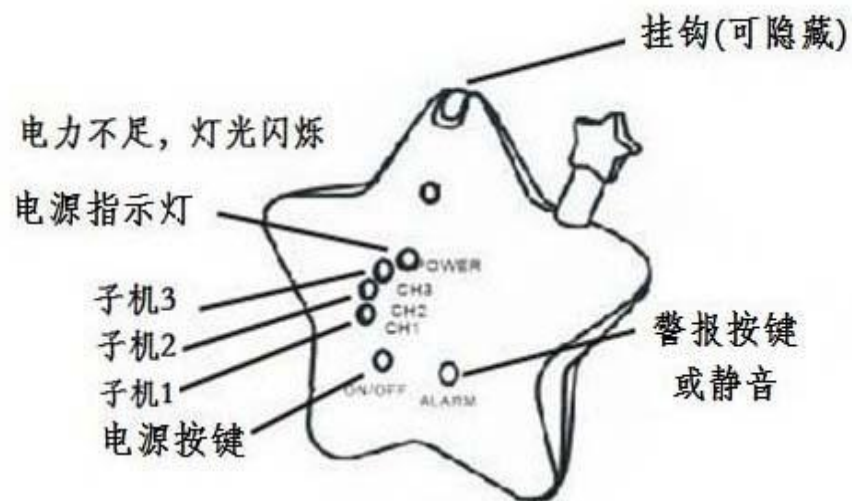
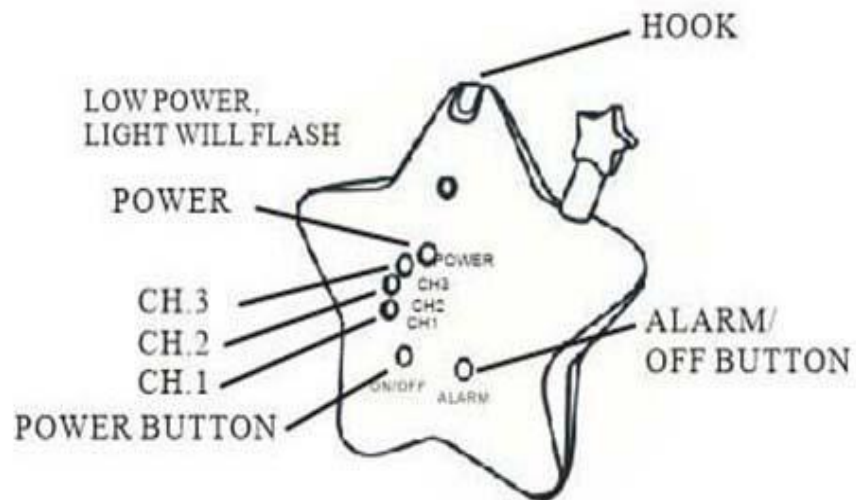
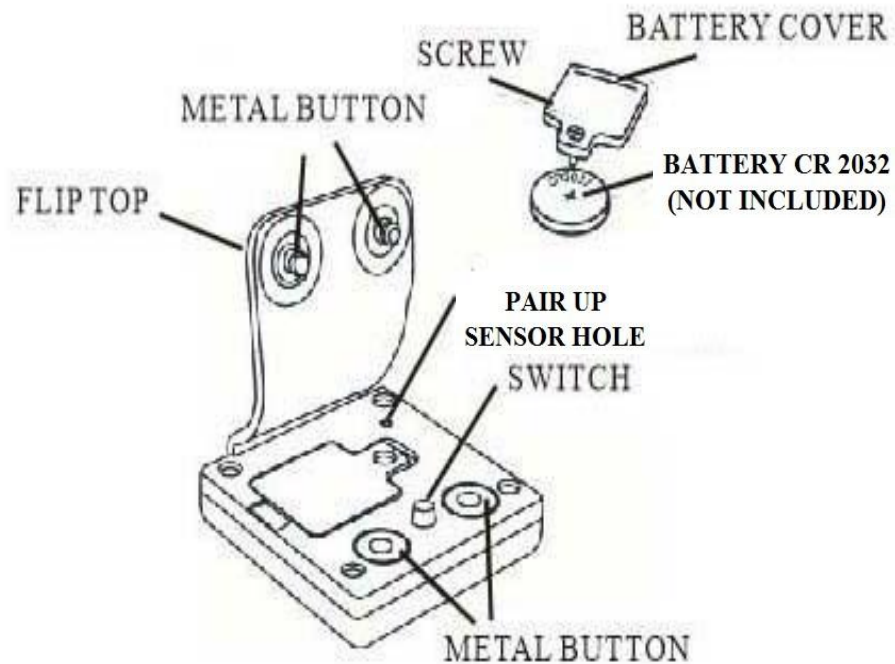
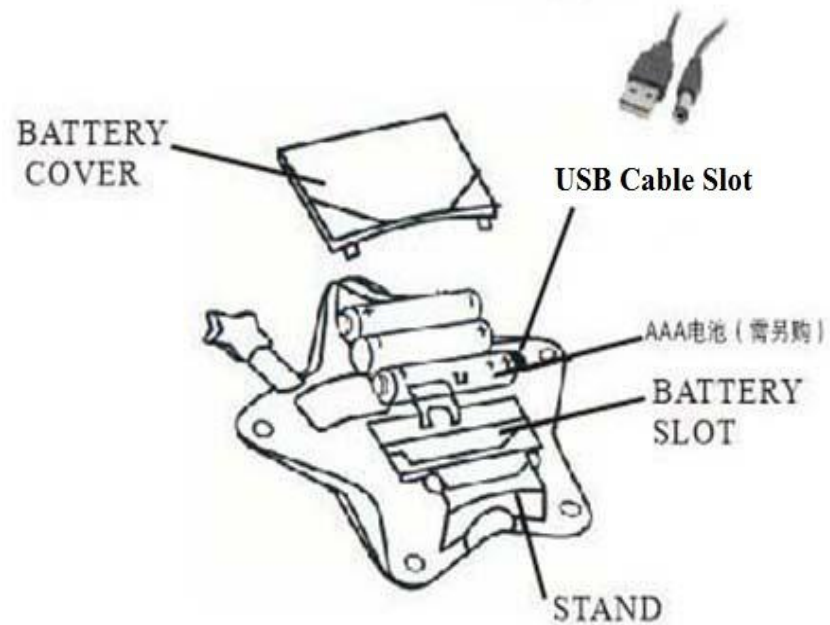
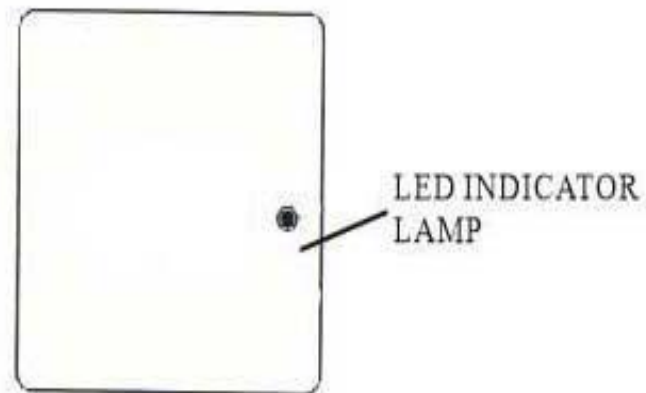
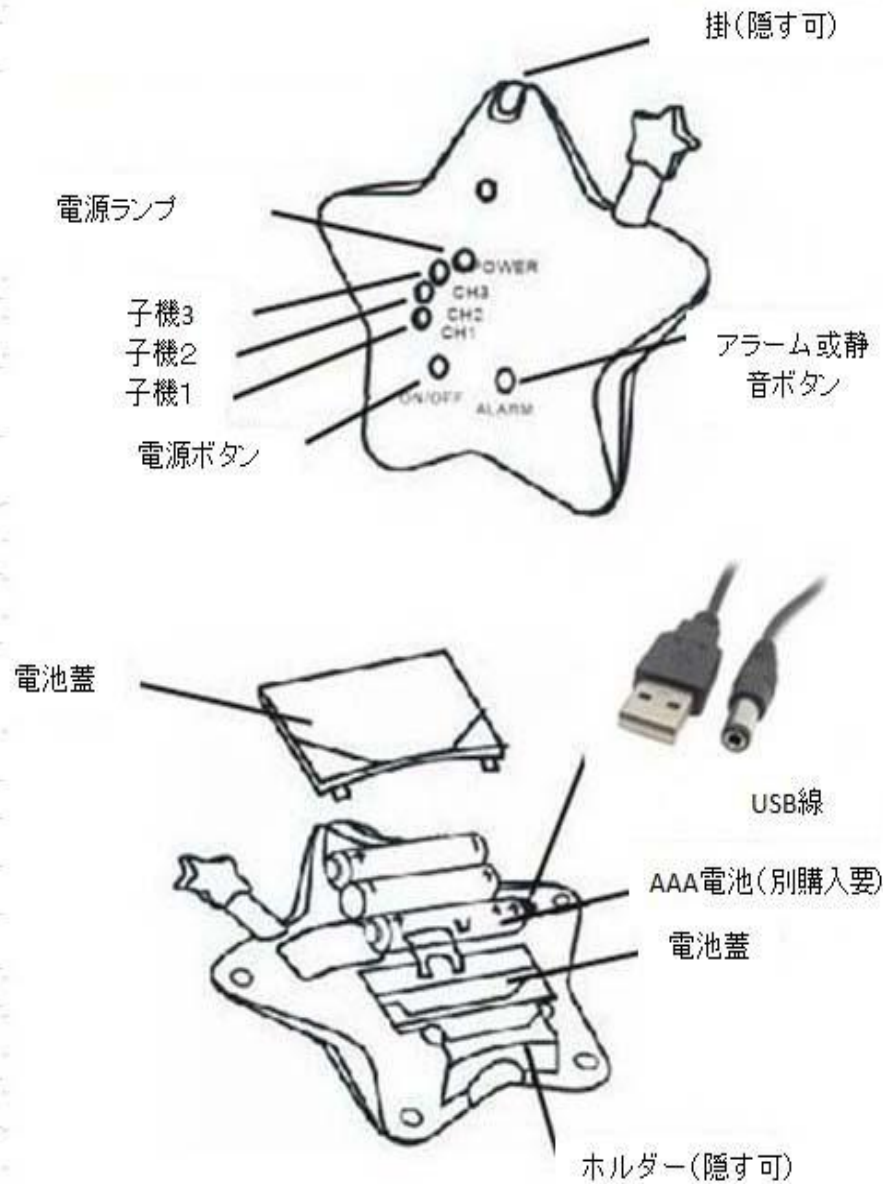
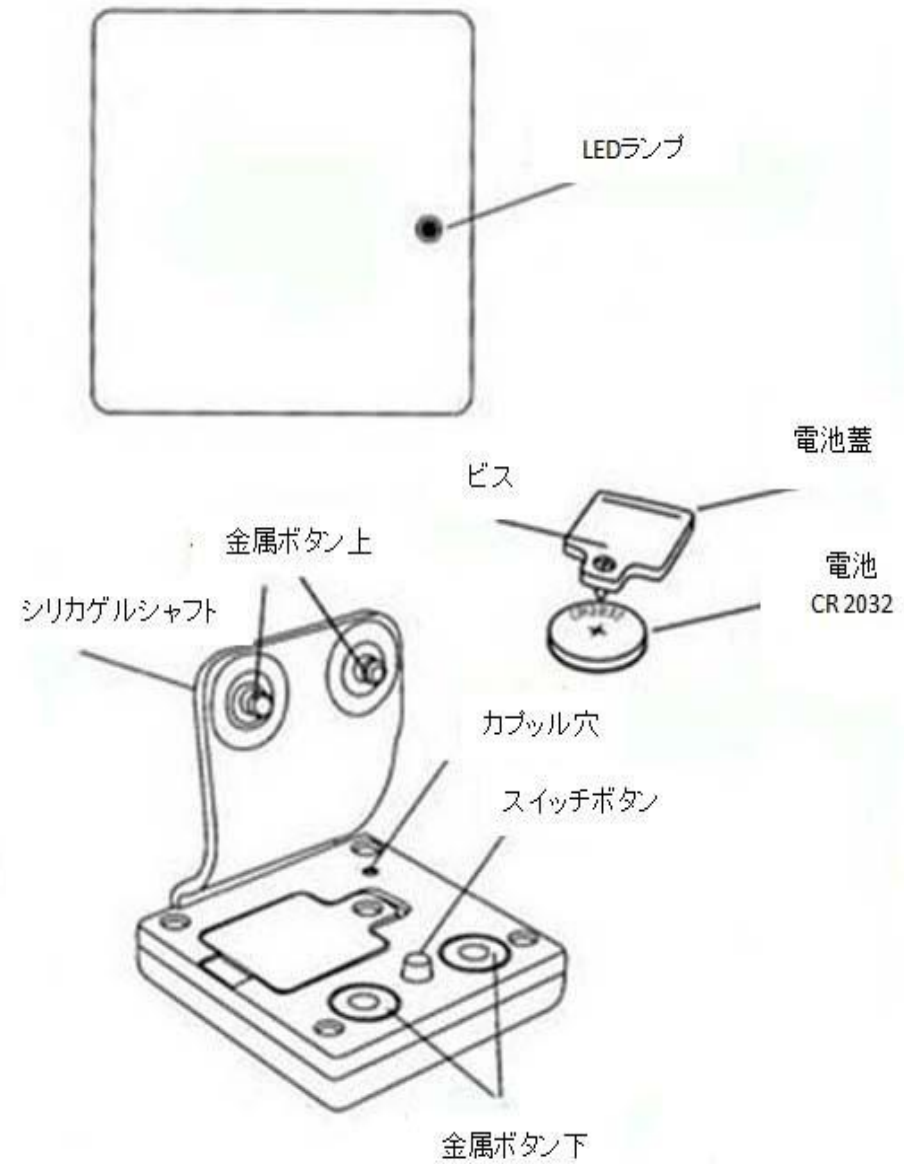
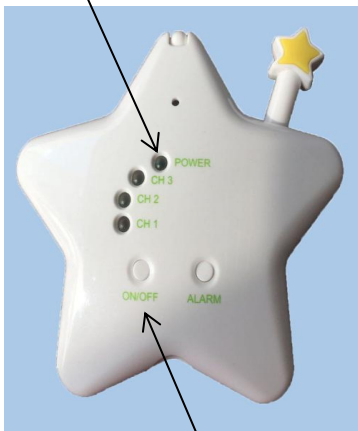
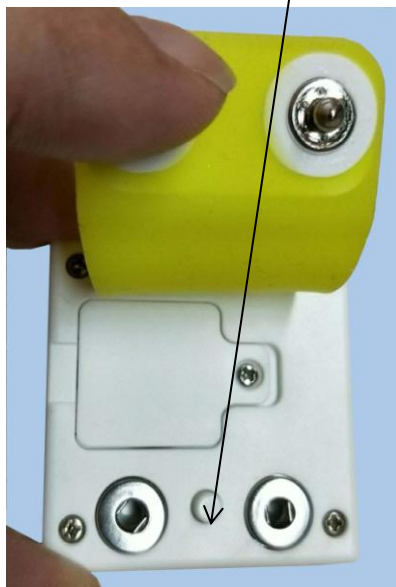
Step 1 : Switch On Main Device

Switch on the button and the red light will be on.

Press ON/OFF button

Step 2 : Reset the Main Device



How to Reset

Press On/Off and Alarm button at the same time, wait for 2-3 second, 3 green lights will be flash once, reset is finish.

Step 3 : Operate the Sub-Device

Make sure the sub-device with power supply, press the back side button and check the light at the front side, if the green light is flashing, that means the device is ready.

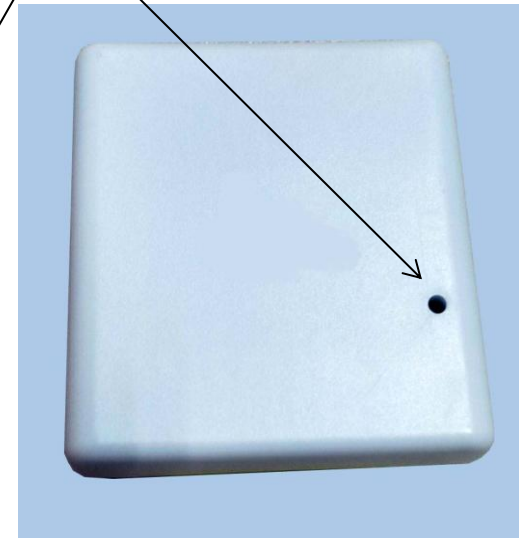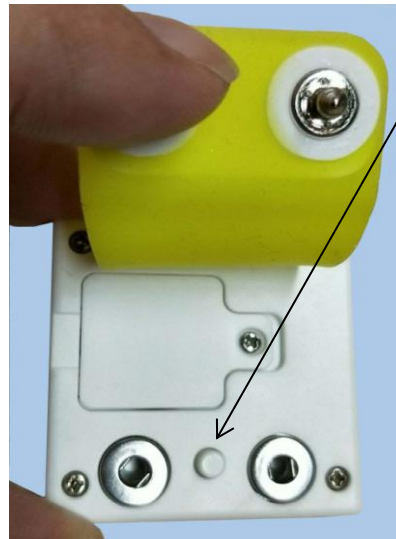Step 4 : How to Pairing

When start to pairing, please switch on the main

Insert the pin (included) into the pair up sensor hole, and main device the light "CH-1" will be flashing, pairing is finish. If you need more than one, than the 2nd and 3rd sub-device are same procedure (as CH-1), than the CH-2 & CH-3 will flashing respectively.

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
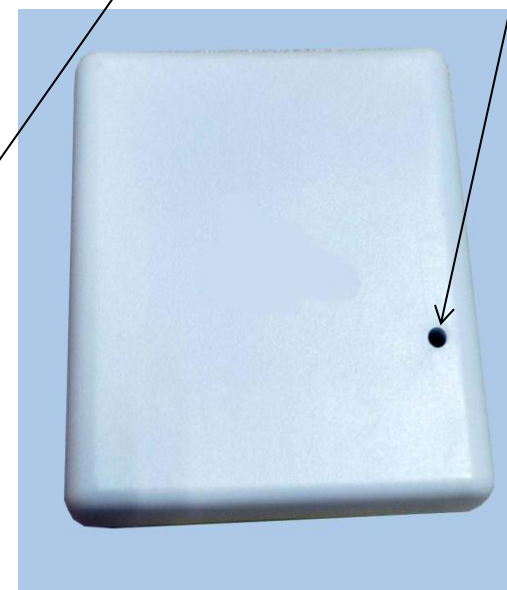
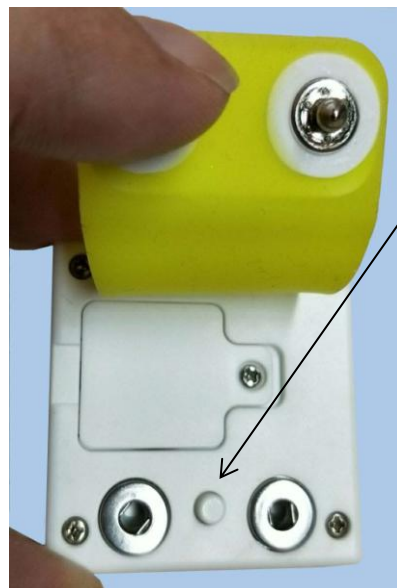
Step 5 : Set-up for Dedicate Daiper

After pairing, main device the CH-1 green light will non-stop flashing, just press the button at the sub-device (as photo), CH-1 light will be off (but still in stand by position), take the sub-device and set up on the daiper's carbon stripes, and the set up is finish.


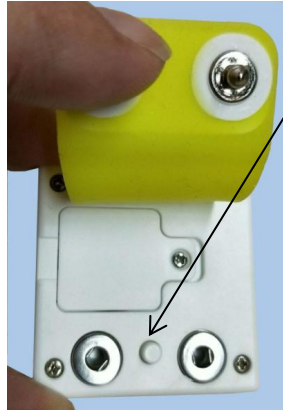
Step 6 : Final Check

After set-up, the main device can put within 20 meters distance (maximum 20 meters for in-house and 60 meters for out-door), wear-up this diaper for necessary people, after piss and poo, main device will alarm it, that means the user need to renew the daiper, if only press the alarm button, the green light will continue to flash, after 10 second the alarm will on again, until the sub-device separate with used daiper, than green light will off. If set-up the sub-device on the new dedicate daiper, than both device will working again.

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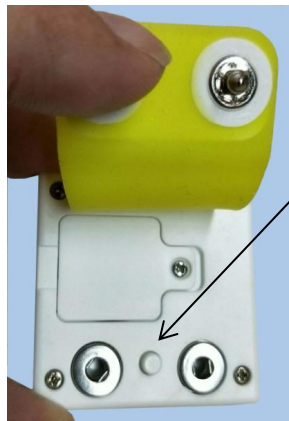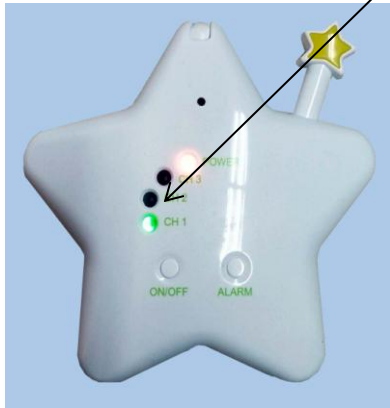Warning

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.**
- Increase the separation between the equipment and receiver.**
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.**
- Consult the dealer or an experienced radio/TV technician for help.**

--Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

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