SOOCAS

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SOOCAS

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USER MANUAL

SOOCAS sonic electric toothbrush X3 USB Edition

## IMPORTANT

### ⚠️\*Note for Users

The brushing efficiency of the product is 7 times higher than that of a manual toothbrush. Gums may feel itchy or bleed slightly when you use the product for the first time because gums need time to adapt to the high-efficiency cleaning. This will generally last for one or two weeks. If it remains after two weeks, please stop using it and consult your dentist.

### \*WARNING

1. Users who have received tooth orthodontic, filling or plantation should consult dentists prior to use;
2. SOOCAS sonic electric toothbrush complies with regulations for electric magnetic compatibility. Users of cardiac pacemakers or other implanted devices should consult doctors or device manufacturers prior to use;
3. The product is designed for tooth, gum and tongue cleaning. Do not use it for any other purpose. If you feel unwell or painful, stop using the product and consult a doctor;
4. Replace the brush head that deforms or whose bristles are curving. You should replace the brush head once every three months. Brush head that is abraded should be replaced immediately. Do not use brush heads other than those recommended by SOOCAS;
5. If your toothpaste contains superoxide, sodium bicarbonate or bicarbonate chemicals (commonly used in whitening toothpaste), the brush head must be cleaned using soap and water after use to prevent plastic from cracking.

### \*PROHIBITION

1. Do not immerse the charging dock in water or use it in any liquid;
2. Make sure that the plug on the charging dock is completely dry before connecting the dock to the power supply;
3. The product is equipped with rechargeable battery. Do not throw it into fire or put, use or recharge it in a high-temperature environment;
4. Do not hit or drop the product. If it is damaged because of collision or falling, stop using it;
5. Do not modify, disassemble or repair the product;
6. Keep the product away from a strong electromagnetic field such as an induction cooker during recharging.

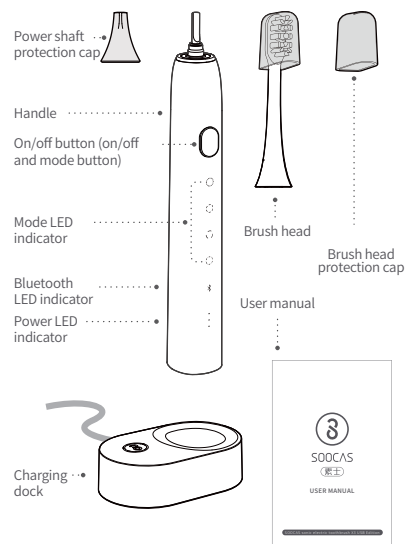
### \*IMPORTANT

1. Make sure that the product is not damaged prior to use. Stop using it if it is damaged or falls heavily;
2. Make sure the voltage of the power supply is consistent with the rated voltage of the product before connecting it to the power supply;
3. Do not use a damaged charging dock;
4. Do not use the product if the handle or brush head is damaged;
5. Do not use the charging dock in high temperatures, such as under direct sunlight.
6. Do not use the toothbrush during a bath or shower.

### \*REMINING

1. Do not share one brush head with others;
2. Do not twine the power lead on the charging dock when the toothbrush is being recharged;
3. Do not put metal objects such as coins and clips on the charging dock;
4. Disconnect the dock from the power supply if the product is not to be used in a long period of time to prevent a fire or other dangers.

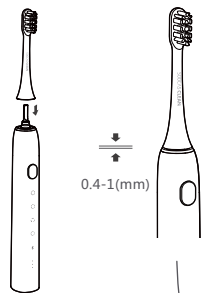
## PRODUCT DESCRIPTION

## USING YOUR TOOTHBRUSH

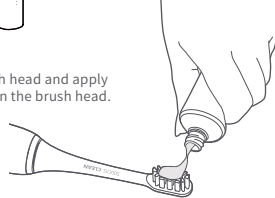
### STEP 1:

- Install the brush head on the handle and make sure that there is a space between the two parts. (The space can prevent the head from colliding with the handle when it vibrates quickly and can improve the motor's power. Because stains can easily hide in the space, we suggest that you disassemble and clean the brush head regularly.)

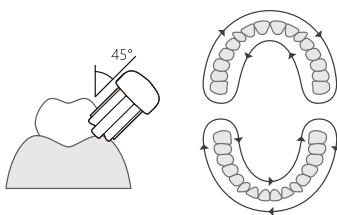
### Step 2:

- Wet the brush head and apply toothpaste on the brush head.

### STEP 3:

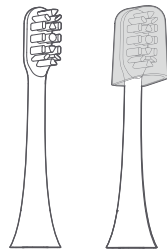
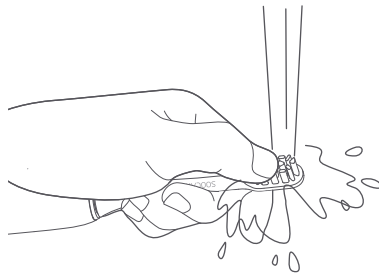
- Guide the brush head slowly from tooth to tooth with a 45 degree angle between bristles and gums line.
- You can press the on/off button to change the brushing mode.
- The toothbrush may vibrate frequently during use to remind you that you should brush another place.

	On/off/mode switch
	• Switch on: Press and release the button
	• Mode: Press the button when the toothbrush is on to change brushing mode.
	• Poweroff: During starting-up, switch to Mode IV to shut down/Long press to shut down.
	• If any exception occurs, press and hold the button for 12 seconds to restore default setting.

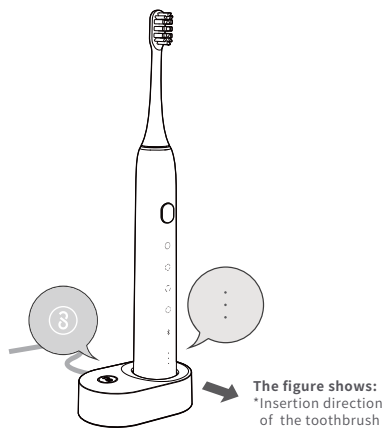
### CLEAN BRUSH HEAD

- After brushing, rinse your brush head under running water.
- After cleaning, cover the protective cap to keep the toothbrush clean.

## CHARGING

- After the charging dock is connected to power supply, the LOGO on the dock lights up;
- Clean the toothbrush and put it on the dock as shown in the figure. The LOGO light becomes stronger;
- During the charging process, the energy LED indicator will flicker. After it is fully charged, three LED indicators will remain on for about 2 minutes before they light off;
- The charging process takes more than 16 hours when the energy LED indicator is red.

## FOUR BRUSHING MODES PROFESSIONAL AND EFFICIENT

### Daily Clean

The brush head can vibrate up to 31,000 times per minute to create effective liquid flows in your mouth and remove dental plaque and stains quickly. This mode is developed for daily clean.

### Sensitive

The sensitive mode can effectively brush your teeth and protect sensitive gums. This mode is suitable for users who have gums swelling and pain or received tooth implantation or orthodontic treatment.

### Gum Care/First Use

This mode is designed to massage sensitive gums and create a more comfortable brushing experience. It is developed for who have troubles such as gums bleeding and oral ulcer or who just received tooth extraction. We recommend using this mode for one week if you are using SOOCAS sonic toothbrush for the first time so that your teeth can adapt to the change from a manual toothbrush to electric toothbrush.

### White/Custom

The toothbrush vibrates 34,800 times in one minute and 37,200 times in the next time. This mode is designed to scrub off smoking, tea, coffee and other stubborn stains. You can use Custom function in the APP to customize your own brushing mode.

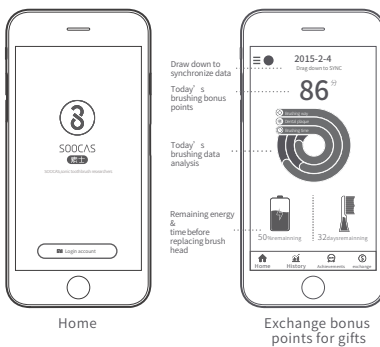
### SOOCAS APP

Download SOOCAS APP. You can get bonus points after toothbrushing, which can be exchanged for gifts.

#### Download APP

Search SOOCAS in the app store or scan the following QR code to download the APP.

\*Note: SOOCAS has the power of interpretation for the above-mentioned.

## DATA SYNCHRONIZATION

- Download the APP. Register a MI account and log in. Click the avatar to enter My Device.

Active Bluetooth on mobile phone and search SOOCAS.

When the Bluetooth LED indicator flickers, press the on/off button.

Step 1

Step 2

Step 3

## ENVIRONMENTAL NOTICE

Product contains rechargeable lithium battery, which may be contaminating. Please take the battery out before you dispose it or send it to a recycling center appointed by local government and send the battery to a designated waste collection station. If you cannot take the battery out, you can send the product to our after-sales service center. We will remove the battery and dispose it in an environmental-friendly manner.

## Specifications

### Handle

Name of Product:  
SOOCAS Sonic Electric  
Toothbrush  
Model: X3 USB  
Rated voltage/power: DC 3.7 V  
Rated power: 2 W  
Charging time: about 16 hours  
Waterproof grade: IPX7  
Country of production: China  
Standards : GB-4706.1-2005,  
GB -4706.59-2008

### Charging dock

Charging method:  
inductive charging  
Type: sealed charging dock  
Waterproof grade: IPX7  
Input: 5 V  $\equiv$  0.5 A

## HAZARDOUS SUBSTANCES

Hazardous substances and content in product

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This form is prepared in accordance with SJ/T11364.  
 ○ : It indicates that the content of the hazardous substance in this part is lower than the required level of GB/T26572.  
 X : It indicates that the content of the hazardous substance in this part is higher than the required level of GB/T26572.

\*Brush head has passed the FDA standard test

\* This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**CE FC**  **RoHS**

 **IPX7**  

**\* Tips:**

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

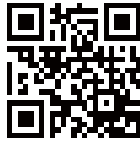
If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

You are welcomed to visit [WWW.SOOCAS.COM](http://WWW.SOOCAS.COM)  
Service hotline: 400-788-4008  
SOOCAS (Shenzhen) Technology Co., Ltd.  
Room 201, Building A, No.1, No.1 Qianwan Road,  
Qianhai Shenzhen-Hong Kong Cooperative zone,  
Shenzhen City, China

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