
碧波庭
BIBOTING

BIBOTING BB-333
Healthcare Massager

Product
Manual

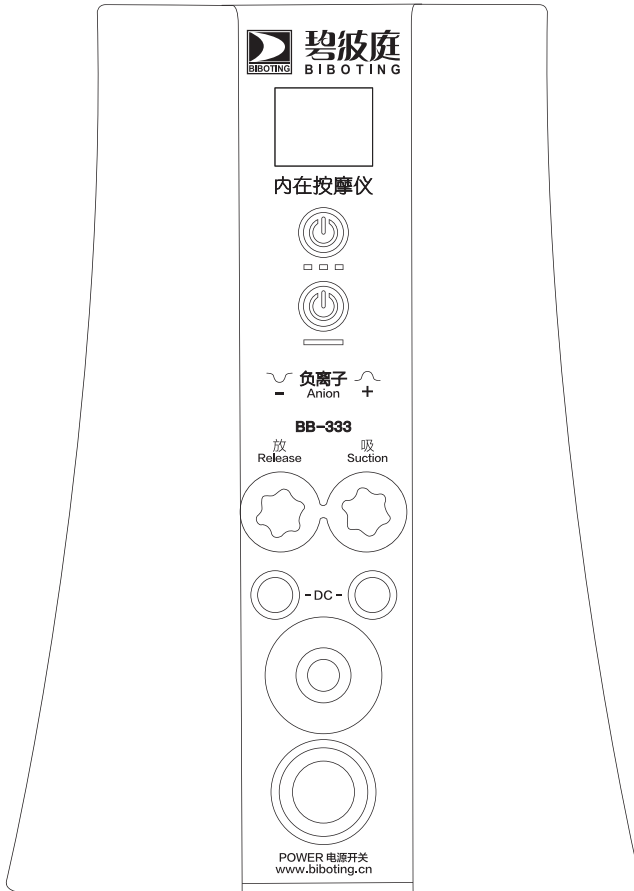This instrument is used for massage and it is not a replacement of any treatment. It is suitable for home use.

BIBOTING Authenticity Statement

1. 碧波庭, BIBOTING and H34D are trademarks of BIBOTING International CO., LTD.
2. There is a warranty code for product with regular purchase registration that you can check with BIBOTING Company.
3. Please read the precautions at the side and the bottom of the instrument carefully before use.
4. Please check carefully in order to avoid purchasing counterfeit products to protect your consumer rights and for your safety.

BIBOTING Healthcare Massager

- Thank you for purchasing BIBOTING Healthcare Massager.
- Please read the operation manual carefully before use and pay special attention to safety matters in order to operate and use the product correctly.

Contents

Precautions	01
Description of Functional Instruction	05
Product Operation	07
Product Accessories Manual	08
Product Maintenance	09
Troubleshooting	10
After-Sales Service	11
Specifications	13

BIBOTING International CO., LTD

Contact No.: (+886) 3-250-4567

Address: Floor 8, 66-7 Sec. 2, Nankan Rd, Luzhu District, Taoyuan City, Taiwan.

Factory Registration Certificate No. (TAIWAN): 99-687589-00

ISO9001 Registration Certificate No: 83165191534

BIBOTING Official LINE

BIBOTING Official FB

Precautions

- Please consult your doctor before use if you are receiving drugs or have medical condition.
- This instrument is used for massage and it is not a replacement of any treatment.
- Using the product on the following parts are strictly prohibited:
 - ① Artery position on both sides of the neck
 - ② Trachea position around the neck
 - ③ Area around eyes and eyeball
 - ④ Area around ears and ear canals
 - ⑤ Temple on both sides of the head
 - ⑥ Varicose veins or wounds
 - ⑦ Nostril
 - ⑧ Mouth
 - ⑨ Anus
 - ⑩ Urethra
- DO NOT USE under the following physical conditions:
 - ① Breast cancer or other cancer patients
 - ② Those who use Ogilvy injection
 - ③ Those with cardiac catheter or pacemaker installed
 - ④ Pregnant women
 - ⑤ Those with difficulty repairing diabetic wounds
 - ⑥ Those with major illness
 - ⑦ Those with severe hypertension
 - ⑧ Heart disease patients
 - ⑨ Epileptic patients
 - ⑩ Patients with skin disease such skin ulcer, infectious inflammatory skin, cellulitis and etc.
 - ⑪ Those with severe anemia (usage may cause dizziness)
 - ⑫ Patient with surgical procedure

- Do NOT place the product near the heater, cooker or direct contact with the sun.
- Do NOT use the product in the environment that is too humid or over 104° F (40° C) in order to avoid electric shock or mechanical failure.
- Moisture accumulation in the internal electronic parts can easily affect the normal operation if the product is moved from a low temperature environment to a warmer place. Please leave the product at normal room temperature for 1 hour before using it.
- This product is NOT suitable for use on children. Please supervise and ensure children do NOT use this product as toy to ensure their safety.
- Please stop using the instrument immediately if you feel unwell in the process.
- Please be sure to use the designated essential oils and skin care products from special ingredients and original supplies to avoid eroding the cup material or the phenomenon such as blisters, broken skin, bruises, pore blockage and other adverse effects from happening. The company shall not bear any responsibility for the above consequences.
- The suction level is recommended to start from level 4 to 5, gradually strengthen the suction when you are feeling more comfortable.
- Do NOT use more than 30 mins on each breast, belly and hip, rest for at least 6 hours before repeat the use. Do NOT use more than 2 mins on spots other than breast, belly and hip, move to other body part before continue using the device.
- Healthy Cup must be disinfected before and after use, do ensure vibrating part does not contact with water during cleaning.
- Do not use the product within 2 hours after consuming food.

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using (this appliance).

DANGER

To reduce the risk of electric shock:

- (1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons:

- (1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking on or taking off parts.
- (2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- (3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- (4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- (5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- (6) Do not carry this appliance by supply cord or use cord as a handle.
- (7) Keep the cord away from heated surfaces.
- (8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- (9) Never drop or insert any object into any opening.
- (10) Do not use outdoors.
- (11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- (12) To disconnect, turn all controls to the off position, then remove plug from outlet.

Grounded Products

(1) The following item shall be included following the heading “WARNING” Connect this appliance to a properly grounded outlet only. See Grounding Instructions.

WARNING

(a) For all grounded, cord-connected products:


GROUNDING INSTRUCTIONS

This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.


DANGER

Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product – if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

- (b) For a grounded, cord-connected product rated less than 15 A and intended for use on a nominal 120-V supply circuit, the instructions in either (1) or (2):
 - (2) This product is for use on a nominal 120-V circuit and has a grounding plug that looks like the plug illustrated in sketch A in Figure 86.1. Make sure that the product is connected to an outlet having the same configuration as the plug. No adapter should be used with this product.

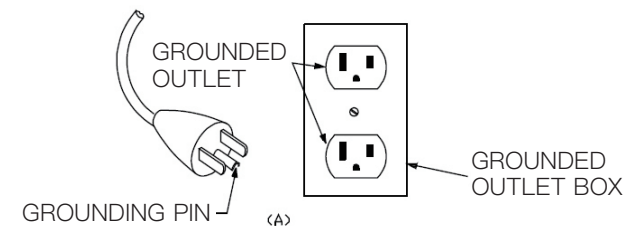Description of Functional Instruction

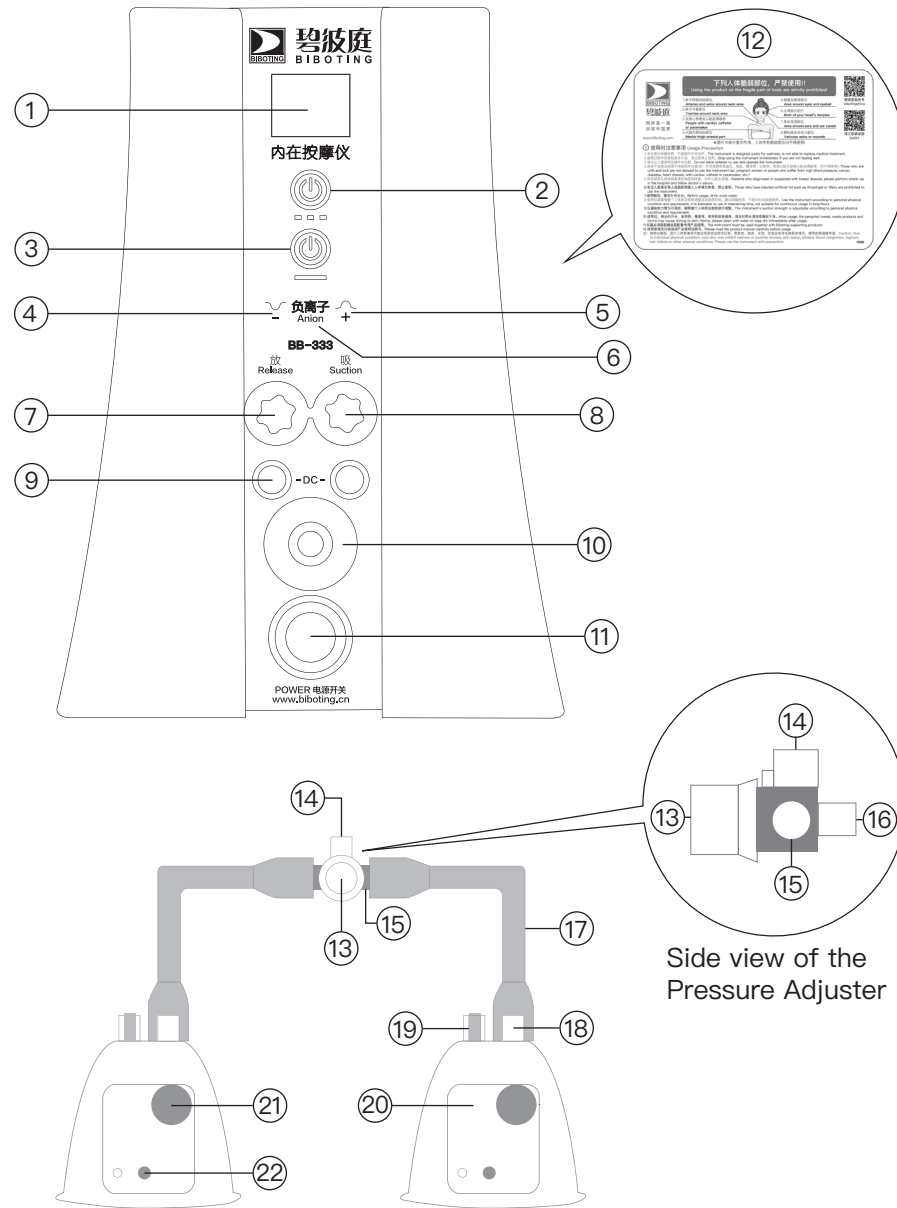- ① Count down Timer
- ② Suction and Release button: Press once to start, press again to stop.
- ③ Suction button: Press once to start, press again to stop.
- ④ Deflating indicator
- ⑤ Suction indicator
- ⑥ Negative ions indicator
- ⑦ Release timer button: the minimum time setting is 0.0 second.
- ⑧ Suction timer button: the minimum time setting is 0.0 second.
- ⑨ Healthcare Massager DC power point
- ⑩ Suction connector: Connect with long tube.
- ⑪ Power: ON / OFF power button.
- ⑫ Precautions for using the instrument
- ⑬ Pressure Adjuster: This knob has settings from 0–10, 0 is minimum and 10 is maximum which represent the suction strength. It is highly recommended to start with level 5.
- ⑭ The upper end interface of Pressure Adjuster: Connect interface with short tube for the use of Healthy Cup.
- ⑮ The interface for both ends of Pressure Adjuster: Both interfaces are used to connect with short tubes.
- ⑯ The back interface of Pressure Adjuster: This interface is used to connect with long tube.
- ⑰ Short tube: Connect to Healthy Cup.
- ⑱ Healthy Cup interface: Connect with short tube.
- ⑲ Safety relief valve: Pull up to release the Healthy Cup.
- ⑳ Rhythm cube
- ㉑ Frequency adjusting knob: Adjust vibration frequency of the Health Cup.
- ㉒ Healthy Cup DC power point: Power supply for Healthy Cup operation.

After switch on the power button, the default settings are:


- Complete cycle time: 30 minutes
- Deflating frequency: 0.6 second
- Suction frequency: 2.4 seconds
- Negative ion activation

Healthy Cup manual

- Pull up safety relief valve or adjust Pressure Adjuster to reduce suction power if the suction power is too strong.
- There is a rhythm cube on Healthy Cup, turn the switch clockwise to increase the frequency or anti-clockwise to reduce the frequency.
- Healthy Cup sizes are 8, 10, 12 and 15 (according to the scale on the cup), refer to table below for a suitable Healthy Cup size.


Suitable Bra Size	Suitable Healthy Cup
A or B	No.8
B or C	No.10
C or D	No.12
D or E	No.15

Product Operation

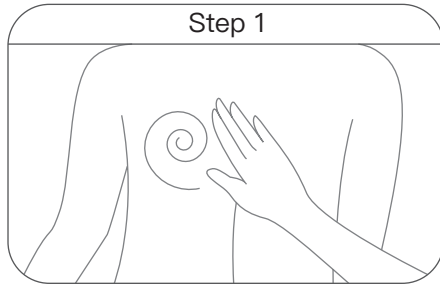
Step 1

1. First wipe clean around the body part.
 2. Apply essential oils.
 3. Apply nourishing cream.
- Apply evenly around the body part.

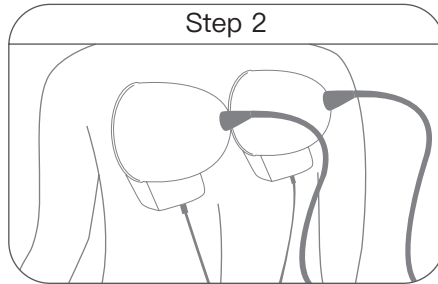
Step 2

1. Switch on the power button.
2. Place Healthy Cup on the body part.

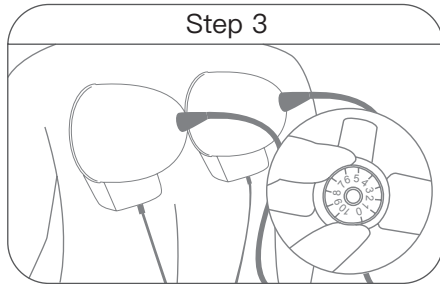
Step 3

1. Use Pressure Adjuster to adjust the suction level.
 2. Please start from a low pressure, then adjust to a comfortable level.
- Adjust the level according to individual capacity.

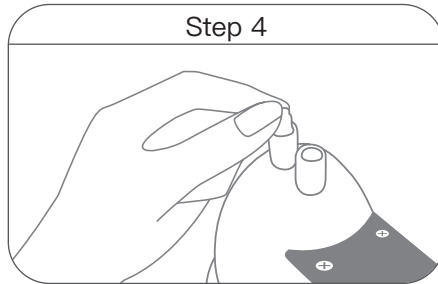Step 4

1. Pull up the safety relief valve to remove Healthy Cup.

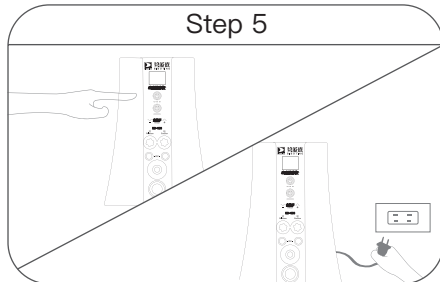Step 5

1. The instrument will stop automatically after the preset operation time.
2. Press the blue or red button to continue using after the instrument stopped.
3. Please remove the power plug when not in use.

【Simple Troubleshooting】


- Please ensure the power plug is on.
- Verify if the long and short tubes are connected correctly and turn Pressure Adjuster to a comfortable level (or a scale of about 5).
- Check if the DC cable is properly connected with Healthy Cup.

Product Accessories Manual

- Flat Cup / 1 pcs (for face)


- Trump cup / 1 pcs (for private parts)

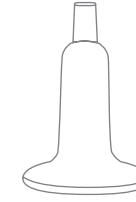
- DC cable (180cm) / 2 pcs

- Long tube (180cm) / 1 pcs

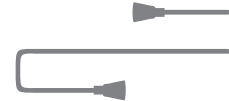
- Pressure Adjuster and 2 short tubes (40cm) / 1 set

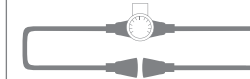- Car charger DC cable (200cm) / 1 pcs

- Healthy Cup No.8, No.10, No.12, No.15 / 2 pcs

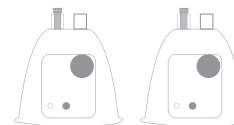- Listed Cup / 2 group (8 pcs)

Product Maintenance

Accessories	Cleaning	Attention
A. Healthy Cup	Please use standard cleaning solution or disinfection alcohol then use fibrous towel (do not use tissue or cotton towel) to wipe the cup in and out, if the black deflated rubber fastened in the cup falls, please use cotton swab with alcohol to clean the part directly, after that restore the rubber to the original position.	Cups product with physical contact must be cleaned and disinfected after using to avoid cross-infection of skin, and do not share this product with others for personal hygiene.
B. Listed Cup		
C. DC Power Cable	If the power cable is stained by essential oil or cream product, remove cable from instrument before cleaning, then use warm water on fibrous towel to wipe clean.	To avoid reducing the life cycle of the power cable as a result of short circuit due to excessive entanglement, please do not fold the power cable after cleaning.
D. Pressure Adjuster	Set the level to 0 and place into hot water for at least 5 mins, then clean Pressure Adjuster with cold water.	Place at room temperature to dry.
E. Long tube (180cm) / Short tube (40cm)	Clean inside and outside of the tube with warm water or cleaning solution.	Ensure hands are clean from any essential oil or cream products before switching on power button and other on/off buttons on the panel to avoid affecting the normal operation and life cycle of the instrument.

Troubleshooting

- When experience any problem using BIBOTING Healthcare Massager, please refer to the following for possible solutions.
- If the solution is not helpful, please turn off the power and remove the plug, and contact BIBOTING Service Centre immediately.
- Users should not disassemble or repair BIBOTING Healthcare Massager. Damage caused by self-disassembly or self-repair will invalidate the product warranty, and the company shall not bear any responsibility for the consequence.

Situations	Possible reasons / Solutions
Rhythm box on Healthy Cup is unable to vibrate.	<ol style="list-style-type: none"> ① Please ensure the DC cable is connected properly. ② Test by turning frequency adjusting knob on the rhythm box clockwise and anti-clockwise. ③ Replace with another DC cable or Healthy Cup, it must be sent for repair if it cannot be started. (This product is a consumable item. Do not crash it during the 1-year warranty period to avoid affecting its normal function.)
Switch on the power button and nothing appears on screen.	<ol style="list-style-type: none"> ① Restart the power button. ② Relocate to another power outlet and test, it must be sent for repair if it cannot be started.
Pressure Adjuster does not work.	<ol style="list-style-type: none"> ① There is foreign matter or grease blockage in the Pressure Adjuster. Please refer to cleaning instruction. (P.9)
Suction becomes weaker or could only feel the strength around suction hole of the main processor.	<ol style="list-style-type: none"> ① One of the tubes is not connected correctly or properly. Please check the tubes and reconnect properly. ② Tube might be broken especially at both ends, please replace the tube if it is broken. ③ There is foreign matter or grease blockage in the tube, please refer to cleaning instruction. (P.9) ④ There is foreign matter or grease blockage in suction hole of main processor. Please refer to cleaning instruction. ⑤ The stopper in Healthy Cup or Listed Cup is detached or broken, which causes the suction becomes weaker.

After-Sales Service

Warranty

These clauses apply only to our products under normal operating condition, the clauses do not apply if the damage is caused by the act of god or it is a human-induced damage.

- We provide 3-year free warranty repair for instrument main processor. A repair fee will be charged according to the condition of instrument if the 3-year quality assurance is over, and it is still within the 20-year instrument warranty period.
- The rhythm box on Healthy Cup is warranted for 1 year.
- If there are glitches of the instrument and cable accessories that affect the use within 30 days of purchase, they can be replaced after testing and verification.
- Healthy Cup and Listed Cups Set – Flat Cup, Feminine Cup and Trump Cup, are not covered under warranty, please examine the product upon receipt, any damage unrelated to production quality is not liable for free replacement hence must to be re-purchased.
- In the event the product is repaired or exchanged, the subsequent warranty period will be based on the original product hence there is no warranty period extension.
- Malfunction of instrument main processor due to natural disaster or man-made incident will be charged a service fee according to the faulty condition as a standard.
- This clause does not apply to instrument that has been dismantled by unauthorized service personnel.

Disclaimer

- The physical reaction caused by products or skin care products not by our company has nothing to do with us.
- Our products do not have any efficacy, please consult a professional medical institution if you have physical discomfort.
- Claims that products are from Biboting, our company reserves the power to legal recourse. If you use products not from the company but claim they are Biboting products, the company reserves the right to pursue legal prosecution.

Service Centre

- RESGREEN BIBOTING Showroom
Address: (2A-06 Showroom, Liancheng Garden Lobby) No. 9, Xiangfu Middle Road, Yuhua District, Changsha City, Hunan Province, China.
Trader Address: Floor 8, No. 66-7, Section 2, Nankan Road, Luzhu District, Taoyuan City, Taiwan.
- Service Centre in Shenzhen
(For product maintenance and repair only, not open for public)
Customer Service No.: (+86) 171-8748-8849 ∙ (+86) 170-9142-2298
Service Centre Address: No. 171, Guziyi Road, Xixiang, Baoan District, Shenzhen City, China. (A District, First Floor, Building B, Gate 3, Weicheng Industrial Zone)
- BIBOTING International CO., LTD
Customer Service No.: (+886) 3-250-4567
Service Centre Address: Floor 8, 66-7 Sec. 2, Nankan Rd, Luzhu District, Taoyuan City, Taiwan.
- Service Centre in Malaysia BBT BEAUTYCARE (Asia) SDN. BHD.
Customer Service No.: (+60) 12-267-6723
Service Centre Address: Unit 53-7, Level 7, The Boulevard, Mid Valley City, Lingkaran Syed Putra, 59200, Kuala Lumpur, Malaysia.

Repair Service Process

- ① For repair, please visit after-sale service website (<http://cnsrv.biboting.com/login.php>), login to register for the repair service.
- ② Please pack the product properly and send to the service address.
- ③ In China, product can be sent to the Shenzhen service base, while in other regions product should be sent to regional office in Taiwan.

Scan the code to login to the after-sales website.

Specification

Product name: BIBOTING Healthcare Massager

Product model: BB-333

Rated Voltage: 100V-240V \sim or 12Vdc

Rated frequency: 50-60Hz

Power: 35W

Dimension: Approx.270 (L) x130 (W) x180 (H) mm

Net weight: Approx.2.5kg

Origin: Taiwan, China

Executive standard number: GB4706.1-2005

Bluetooth system description

Version: V4.1

Operating Frequency Band: 2.4GHz ~ 2.48GHz unlicensed ISM band

Power: 7.5dBm (MAX)

Voltage: 1.8V-3.6V

NCC Warning

1. Without permission, any company, firm or user shall not alter the frequency, increase the power, or change the features and functions of the original design of the certified low-power RF motors.
2. The use of low-power RF motors shall neither affect flight safety nor interfere with legitimate communication. If interference is found, it should be suspended immediately until interference is no longer exists after improvement is made. Legitimate communication in the preceding paragraph refers to radio communication operating in accordance with the provisions of the Telecommunications Act. Low-power RF motors must withstand interference from legitimate communication or radiological, radiative electrical equipment for industrial, scientific, and medical applications.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/ TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.


For License-Exempt Radio Apparatus (RSS-GEN)

This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

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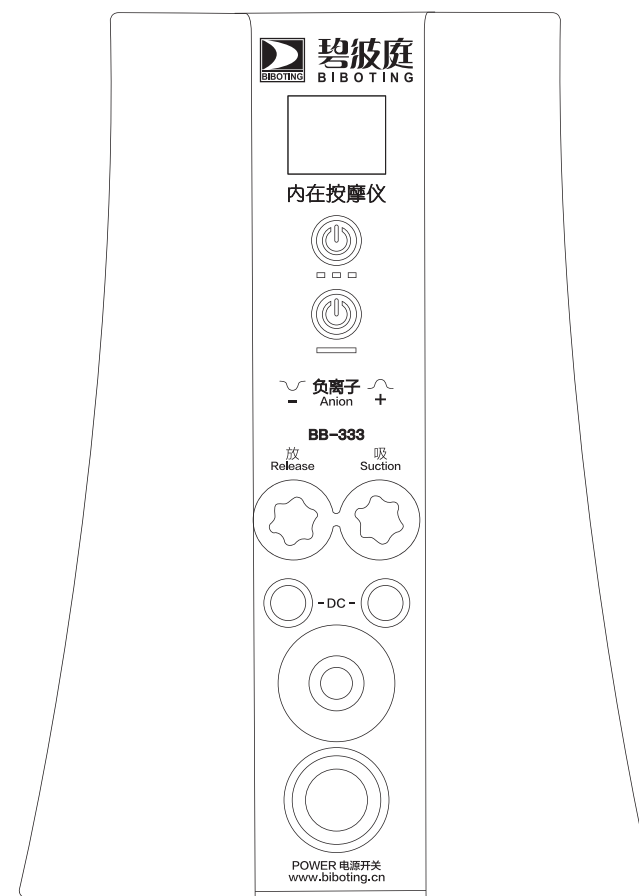BIBOTING International CO., LTD

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Adresse : Floor 8, 66-7 Sec. 2, Nankan Rd, Luzhu District, Taoyuan City, Taiwan.

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