

Eye Massager Smart Eye

OE-0999



catalogue

Thanks for purchasing our product, pls read carefully before use.specially for attention part,and make sure of using correctly. pls keep it to check in futher time.

Brief introduction	— — — — — — — — — —	02
Graphene introduction	— — — — — — — — — —	03
Infrared sensing	— — — — — — — — — —	04
application	— — — — — — — — — —	05
packing list	— — — — — — — — — —	06
use method	— — — — — — — — — —	07
Function	— — — — — — — — — —	08
working mode &attentions	— — — — — — — — — —	09
basic fault & troubleshooting method	— — — — — — — — — —	10
maintenance	— — — — — — — — — —	11
Ship&storage&specification	— — — — — — — — — —	12
QC certificate&Warranty card	— — — — — — — — — —	13

If the design and changes of product functions are involved, without prior notice, please refer to the actual product.

introduction

Product principle

the intelligent air pressure eye massage adopts the principle of traditional Chinese meridian and beauty, combined with the micro-chip control technology,improve eye blood circulation by soft warm, vibration and simulation of acupressure on the eye Built-in speaker and soothing music and relax you even more.

- 1.Intelligent air pressure : vibration to acupuncture point of eye and eyelid, Dredges the eye and head nervous system
- 2.Constant temperature soft warm: it takes by carbon fiber heating element,portable &constant temperature for eye massage, activates eye blood circulation, increases eye nutrient supply ability, effectively relieves eye dryness.
- 3.Vibration massage : vibrating the eye and eyelid, relieving blood circulation to reduce the elasticity of the ciliary muscle, enhancing the ability of the eye muscle to adjust.
- 4.Music: Relax with soothing music in 15 minutes.

Graphene

Graphene A planar film composed of carbon atoms. It is currently the world's lightest and most tough new nanomaterial due to its special nanostructure and high quality physicochemical properties. Graphene has important application prospects in electronics, optics, magnetism, biomedicine and drug delivery, and has been recognized as a "future revolution material"

Graphene Feature

portable convenient



stable performance





long life



infrared portable physiotherapy:

most close to body of wave physiotherapist:

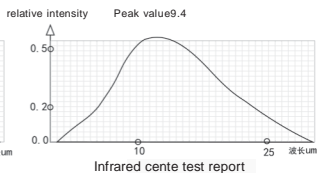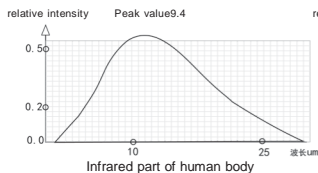

When graphene is heated, the far-infrared wavelength emitted is very close to the far-infrared wavelength of the human body, both of which are 6-14um. When the two wavelengths interact equally, resonance occurs. When resonance occurs with cells, The far-infrared heat energy penetrates into the subcutaneous tissue, uses capillary expansion, promotes blood circulation, strengthens tissue metabolism, increases regenerative capacity, and improves immunity of the body, thereby achieving medical care.

light of lif --Infrared sensing:

light of life power of health:

Far-infrared rays have been discovered as early as 1800 years ago. It is an electromagnetic wave that is invisible to the naked eye. Part of the sun is a dizzying effect of life and growth. It is inextricably linked to the life of the Earth. Aerospace scientists have been in existence. The human living conditions in the spacecraft under vacuum, low temperature and overload conditions are investigated. It is known that far infrared rays with a wavelength of 6-12 microns in sunlight are essential for human survival. This wavelength of light is related to the human body. The far-infrared rays emitted are of the same wavelength, so that the human body can absorb the light energy and exchange and transmit the human body. Therefore, people call this long-range far infrared ray "the light of life."

Introduction of 10 types'far infrared physiotherapy:



application

- ▶ Large, middle and primary school students, myopia, hyperopia, over-eye or long-wearing glasses.
- ▶ computer operators, accountants, designers, secretaries, teachers, long-time eye workers.
- ▶ eye bags dark circles insomnia.



Large, middle and primary school students, myopia



hyperopia



over-eye or long-wearing glasses



Designer secretary Teacher



Driver



long-time eye workers



Computer operator



Insomnia



Long-term mobile phone player

packing list



User manual
(including warranty/certificate)



Storage bag

--

Prepare before using

- ▶ Please confirm if the product has electricity. If the battery is too low and the machine function is turned off or cannot be turned on, please charge it before using it. Be sure to unplug the power adaptor after charging, to avoid charging for a long time.

charge

1. Connect the power adapter to the product USB port.
2. The charging LED light is red during charging, and the indicator light is blue when fully charged (or the LGD display battery power symbol is full).
3. During the charging process, the product does not support the booting work (charging in the working state, the product starts the self-protection program, and automatically stops working).
4. The main unit is powered by a built-in rechargeable battery, and the charging time is 3-4 hours.
5. Prompt for low battery charge during use. Please pause the use and recharge the battery.

Multiple charging methods

charging methods	computer USB port
	mobile
	standard power adaptor(recommend)


working time

The default working time for booting is 15 minutes. If you need to extend the usage time, you need to boot it again.

Function :

- 1.ON/OFF button: long press for 3 seconds to boot, long press for 3 seconds to shut down,Connect to Bluetooth,Double click to play music,Double-click ON/OFF again to pause playing music.Short press to switch modes.
- 2.Play music with Bluetooth, speaker Reminder under connect status.
- 3.3D vibration mode-left side -right side- left&right side (upgrade motor)
- 4.This product cannot adjust the volume, so set the volume.80% of the volume, suitable for most people, other people Adjust the volume level with your phone.
- 5.This Bluetooth name is HongWai-999.Please connect Bluetooth before playing music for the first time.
- 6.With low battery voice reminder function,When the battery voltage is low, the speaker will prompt "Battery low"


7.Infrared operation:

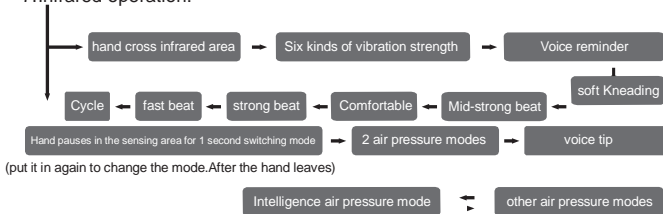

Attentions:

The six force function switches can only be switched in the vibration mode, otherwise it is invalid. Switching between the two air pressure modes only switches in the air pressure mode, otherwise it is invalid.

Function mode

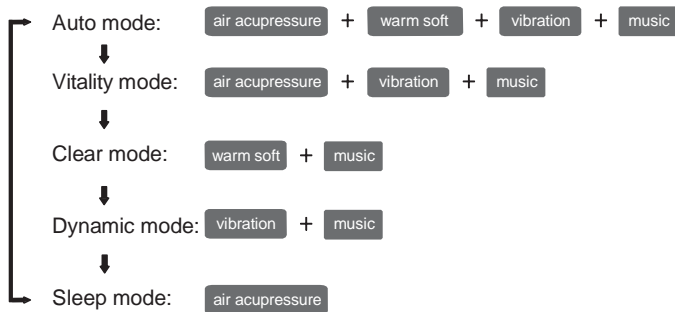

Attentions

To ensure that you use the massager correctly, please refer to the precautions before use:

1. Do not use the eye massager while walking outdoors or driving.
2. Do not use the eye massager in a hot or humid place.
3. There are symptoms of eye diseases that cannot be used. Major eye patients such as eye surgery, glaucoma, cataract, and retinal detachment have been banned.
4. Patients with deep stroke, hypertension and hypotension should be used with caution. Please consult a professional doctor. If there is any discomfort during use, please stop using it immediately.

5. Persons with mental and physical discomfort or mental retardation should be used with caution. Please use it under the supervision of the accompanying staff. If there is any discomfort during use, please stop using it immediately.
6. Please avoid using the eye massage device when drinking alcohol or physical discomfort.
7. Eye massager has a hot compress function, which must be taken care of when using temperature sensitive people.
8. Put the massage device on the eye area then turn it on again to avoid damage to the air bag.
9. The contact lenses, lenses, hairpins, hair accessories, etc. worn before use should be removed first.
10. This massage device is a personal health care product and cannot replace medicines and medical equipment.
11. If you feel pain or serious discomfort during use, please stop using it immediately and consult a physician.
12. This massager is designed as a home electronic health care product and is not intended for use in diagnostic, therapeutic and other commercial applications.
13. If cleaning is required, wipe it with a cotton cloth soaked with medical alcohol. Do not wash the product with water. Do not soak it in liquid.
14. Do not disassemble this product by yourself.
15. Please use the power adapter configured by our company and pay attention to electricity safety.
16. Please use this product for closed eyes. After use, the eyes will be blurred. After closing eyes for 1-3 minutes, the eyes will slowly clear.

basic fault and troubleshooting method

basic fault	possible cause	troubleshooting method
failure to start up	too low power	full charge in time
	auto-protect power off	hold it 5 seconds for charge,let it to be unlocked
	unable to charge	damage adapter
low massage strength	too low power	charge in time

maintenance

- 1.clean with a dry soft cloth if there is dirt on the surface.
- 2.clean with a soft cloth dampened with medical alcohol
if particularly dirty.
- 3.Make the product in a safe place after using.
- 4.take measures to prevent moisture, anti-rat,anti
-pollution.
- 5.please ask a professional technician to repair it if
troubleshooting . Do not disassemble it at will.

Atmospheric pressure: 86 kpa~106 kpa

Ambient temperature: -10°C~65°C

Relative humidity: ≤80% RH, no corrosive gas and well ventilated room.

Transportation: This product is equipped with simple anti-shock facilities, suitable for air, rail, road and ship transportation. It should avoid rain, snow splash, inversion and collision. When the storage period is more than 6 months, it should be checked out of the box and energized for half an hour. It cannot be close to the ground, walls and roof. The room should be well ventilated to avoid strong sunlight and corrosive gas erosion.

specification

rating voltage: 3.7V

rating power: 3.75W

input voltage:5V/1A

- -

certificate of approval

production date: _____

QC: _____

Warranty card (warranty record)

Name	
Tel	
Address	
E-mail	
Product name	
purchase Date	

Warranty date	Warranty content	Remarks

- 1.The company guarantees that the product purchased by the user will be free of charge for parts and maintenance within one year from the date of purchase if the product is used normally according to the instructions.
- 2.For more than one year, if you need to replace parts, the company only charges the cost of parts and repairs free of charge.
- 3.Please return the faulty machine to our company or send it to your local distributor. The company will handle it immediately.

FCC Warning

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE 1: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE 2: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Company: TIC AUDIO INC

Address: 15224 E STAFFORD ST ,CITY OF INDUSTRY,CA 91744

Name: RONGBIN CHEN

Tel.: +1-626-9680211

E-mail: ys33l@yiletec.cn