


正面

反面

ST-2060



FCC ID: 2AS4J-ST-2060

USER MANUAL

86-0755-96093099

ADD: Rm 602, building 4B, Software Industry Base, Haitian Road,  
Nanshan district, Shenzhen,China

SHENZHEN SLEEPTAILOR TECHNOLOGY CO.,LTD

I. Component Name

Exterior schematic drawing:

Front view

Back view

Function keys

Mica

Magnetic lap joint

Hook&loop

Mica

Heating Film

Mica


Main accessories:


User Manual

USB cable temperature control line

II.Model and Specification

Brand	SLEEP TAILOR™
Item Name	GRAPHENE PHYSIOTHERAPY MISC EYE MASK
Temperature control mode	Line Control
Power input	5V 2A
Power connection	Plug Use Cable into Charging Point to Connect the Power Supply
Dimension	690*88mm
Operating power	5~8W
Service environment	Temperature: -20℃ To 80℃; Humidity: <85%RH


2. Long press  button to enable/disable heating function.

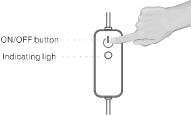

Note: The product is set to work for 1 hour each time. Please restart it if you want to continue

III.Use Instruction

1.Connect the Eye Mask with the power bank (self-prepared. Make sure the power bank and computer opening) or computer via USB cable temperature control line.



① heating

② light

③ adjust

④ power switch

⑤ power switch

⑥ power switch

⑦ power switch

⑧ power switch

⑨ power switch

⑩ power switch

⑪ power switch

⑫ power switch

⑬ power switch

⑭ power switch

⑮ power switch

⑯ power switch

⑰ power switch

⑱ power switch

⑲ power switch

⑳ power switch

㉑ power switch

㉒ power switch

㉓ power switch

㉔ power switch

㉕ power switch

㉖ power switch

㉗ power switch

㉘ power switch

㉙ power switch

㉚ power switch

㉛ power switch

㉜ power switch

㉝ power switch

㉞ power switch

㉟ power switch

㊱ power switch

㊲ power switch

㊳ power switch

㊴ power switch

㊵ power switch

㊶ power switch

㊷ power switch

㊸ power switch

㊹ power switch

㊺ power switch



㊻ power switch


㊼ power switch

㊽ power switch

㊾ power switch

㊿ power switch

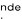3.Under heating state (  light is on),slightly press  button to switch different temperature levels (RED, YELLOW GREEN, BLUE).

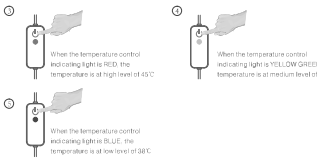

IV.Product Features

Product Features

1.High-tech new material graphene is used as the heating element, with advantages of light, convenient, stable performance, and long service life;  
2.The far-infrared wavelength emitted by the graphene heating film is very close to the human body far-infrared wavelength, providing physical therapy;  
3.Provide even sunbath heating, without producing electromagnetic radiation, safe and comfortable;  
4.Low driving voltage required for product working, high safety factor;  
5.Three-level intelligent temperature control, convenient to use;  
6.One-piece hot-pressing cutting technology, zero pressure, 360° shading, comfortable feel.

V.Precautions

Precautions

1.When the product is not in use, please timely turn off the power to avoid long-term heating;  
2.At the end of use, please remove graphene eye mask after turn off the power for 5 minutes, so that the eyes are adapted to external light;  
3.Recommended duration for each use shall not exceed 45 minutes;  
4.It is not recommended to use during long time sleeping;  
5.Please cut off the power while cleaning the product. Product can be cleaned but rubbing is forbidden. Ensure to use it again after drying off;  
6>Please store in a dry and cool place, avoid exposure and moisture;  
7.There is circuit inside, please do not disassemble and bend the heating part to avoid damage.

VI.Inapplicable for:

1. Infants;  
2. People without ability of self-care.

VII. Warranty Clause

1. Based on the national guidelines on warranties for repairs, replacements and returns, from the date of purchase, we will accept returns within 7 days, and offer replacements within 15 days as well as 1 year of complimentary maintenance services (for product wear and tear only).  
2. The following is outside the scope of this warranty (including, but not limited to): improper use, using the device beyond the recommended operating conditions; product has been tampered with by non-authorized repairs or damage to product is man-made (including

What is Graphene?

Graphene is a flat film consists of carbon atoms, the lightest and thinnest new nano-material with the best hardness and toughness ever tested in the world. Due to its special nano-structure and superior physical-chemical property, Graphene has a brilliant prospect to be applied in such fields as electronics, optics, magnetics, biomedicine, catalysis, energy storage and sensor. It is recognized as "future material" and "revolutionary material" of 21st century.

Light of Life - Far Infrared

Light of Life - Energy Source of Health  
Far infrared, discovered as early as 1,800 years ago, is an electromagnetic wave

Portable Far Infrared Physical Therapy

Ideal Light Physiotherapist for the Human Body

Ten Key Effects of Far-Infrared Physical Therapy

Visualizes molecular activation

Anaesthetic effect

Boosts metabolism

Treats pain in the joints

Strengthens immunity function

Regulates autonomic nervous system

Skincare and beauty

Improves circulation systems

Anti-inflammatory effect

Improves blood circulation

FCC Warning:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.  
This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.  
However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:  
• Reorient or relocate the receiving antenna.  
• Increase the separation between the equipment and receiver.  
• Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.  
• Consult the dealer or an experienced radio/TV technician for help.

Caution:

Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.  
  
This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:  
(1) This device may not cause harmful interference, and  
(2) this device must accept any interference received, including interference that may cause undesired operation.  
  
The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

www.sleeptailor.com