

User Guide

USB bluetooth body fat scale

FCC ID: 2ATTY-CF370BT



Mobile&Device



IOS 7.0



Android 4.3



Bluetooth 4.0

For First Use

1. Activate on

for first use, please press the UNIT button on the bottom side of the scale to be activated. Later on to use, stand on the scale and it will be auto on.

2. App Installation

Scan the QR code in the first page, or download "OKOK International" from App Store or Andriod App store.

3. Registration

Turn on the bluetooth, turn on the OKOK, follow the App steps to register the user profile.

Twelve health index checking

1. Default Eight



Body fat



Body hydration



Body muscle



Body bone



Metabolism



Weight



BMI



Visceral Fat

2. Additional four index, release them or not from the setting



Protein




Obesity Degree



Body Age



Weight without Fat

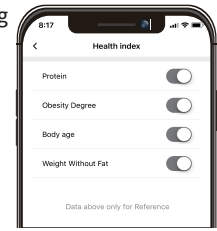

Multi Users

Total 8 users share one user account in one device;

Auto indentify users according to different weight memory.

Target weight

Users can set their own target weight to make the health plan

Attention

1. Please make sure good contact between bare feet and the metal electrodes, otherwise it will display your weight only.
2. Ensure the scale is on a flat, even and firm surface. Avoid carpet or soft surfaces.
3. The surface of the scale will be slippery if wet. Please make sure both the glass and your feet are dry before use.
4. This scale is not water proof. Do not expose to moisture.
5. Clean the scale with a slightly damp cloth. DO not use any chemical/abrasive cleaning agents to clean the scale.
6. Please remove the batteries, if the scale is not to being used for a long time.
7. The appliance is intended solely for domestic use within the home.
- 8.Children should be supervised to ensure that they do not play with the appliance.

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help. The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction