# ***RENPHO***

Live Healthy, Live Smart

**Model:SD-002**

Electric Eye Massager  
User Manual

## Attention



---

Please read and comprehend all the contents of this manual and operate the product according to the instructions.

1. **DO NOT** use if you have had an eye operation, or if you have any eye diseases such as glaucoma, cataract, and amotio retinae.
2. **DO NOT** clean the product with water; clean it with cotton soaked in alcohol
3. This product has a hot surface; people sensitive to heat should take extra care when using this.
4. Juveniles, disabled people, and the elderly should only use this product when properly supervised.
5. The product is intended for personal use; it is not intended to replace medical care and medical devices.
6. **DO NOT** use if you have any swelling in the eye area.
7. **DO NOT** disassemble the product, **DO NOT** impact or squeeze its interior lithium battery or expose it to fire;
8. **DO NOT** use if the battery has been immersed in water
9. **DO NOT** use in high temperatures or humid environments. Avoid direct exposure to the sun.
10. Consult your doctor if you have sensitive skin, allergic dermatitis, or orthopedic bio materials.
11. Remove contact lenses, hairpins, and hair accessories before using.
12. Start the eye massager after putting it on in order to avoid damage to the air bag due to over pressure.
13. Only use the charger (5V 1A, not included) to recharge the eye massager.
14. You must have your eyes closed while massaging.
15. **DO NOT** use sharp objects to scratch this product

# Product Structure

---
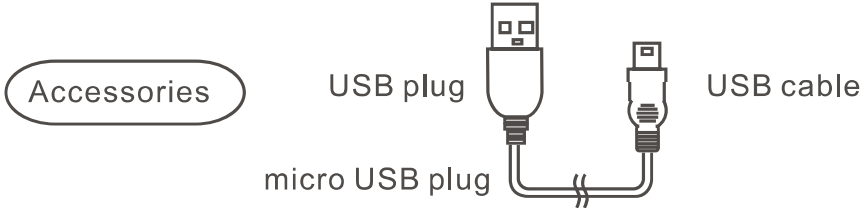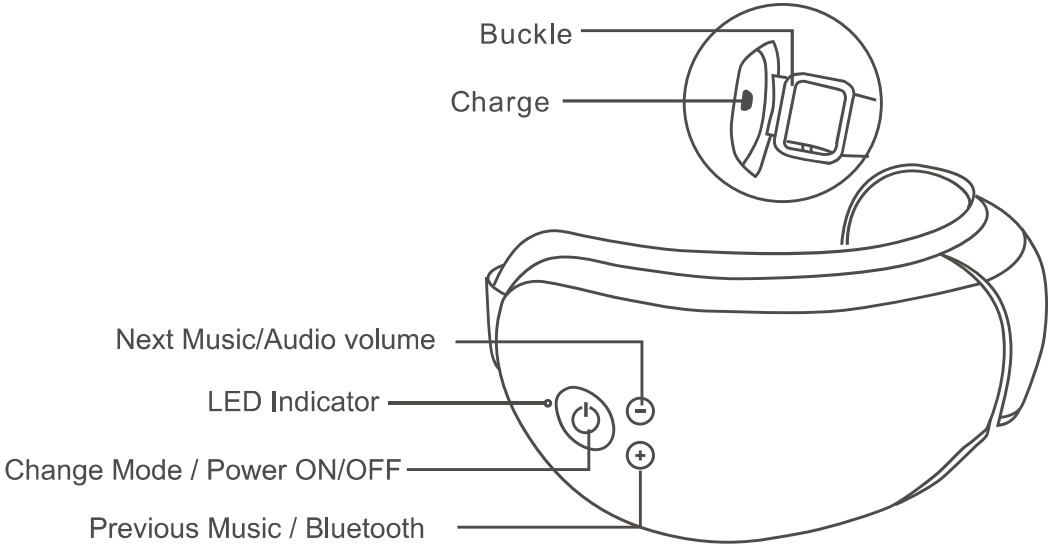
# Technical parameters

---

Product name: RENPHO Eye Massager  
Model: SD-002  
Input: 5V  $\equiv$  1A  
Rated power: 5W


# What's included

---

- 1 x RENPHO Eye Massager
- 1 x USB Cable
- 1 x User Manual
- 1 x Warranty Card

## Wearing the device

Align the cloth cover to your eyes, put on the device, adjust the elastic band, adjust the device to the most comfortable position. Hold the power button for 3-5 seconds to turn on the eye massager

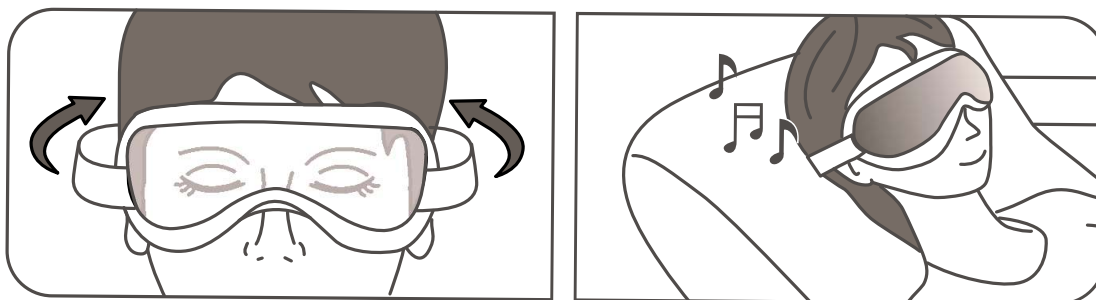Attention


Remove any glasses, contact lenses, hairpins and hair accessories before use to avoid discomfort. **Do not** operate the device with wet hands.

## During operation

Default setting: Clear mode = air pressure + heat compression + music

### Press power button to change mode:


Hold power button 3-5 seconds to turn on/off the eye massager. Gently press power button to change modes. Audio voice will confirm mode change.
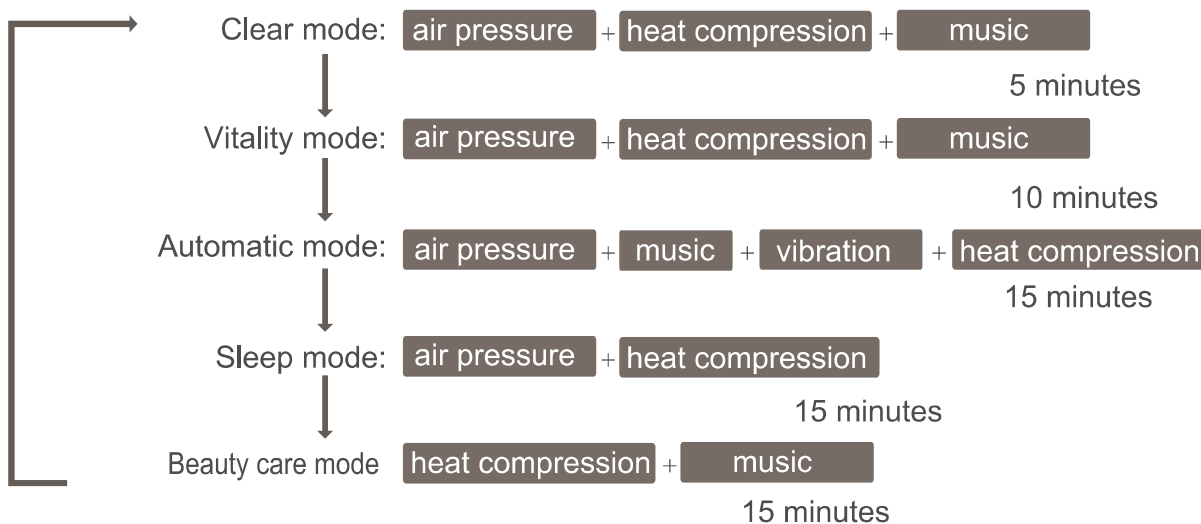
# Charging

---

1. Use the 5V mobile charger(not included) for charging your eye massager.
2. Insert USB cable into the socket of the mobile charger.
3. Softly insert the micro USB plug into the socket of the device, making sure the micro USB is inserted in the right direction.
4. Plug the charger into power.
5. While charging, the red LED indicator will show. When done charging, the LED will be green. The device will not work while charging is in progress.
6. The charge duration will be 2-3 hours. Disconnect the power adapter after charging is complete.
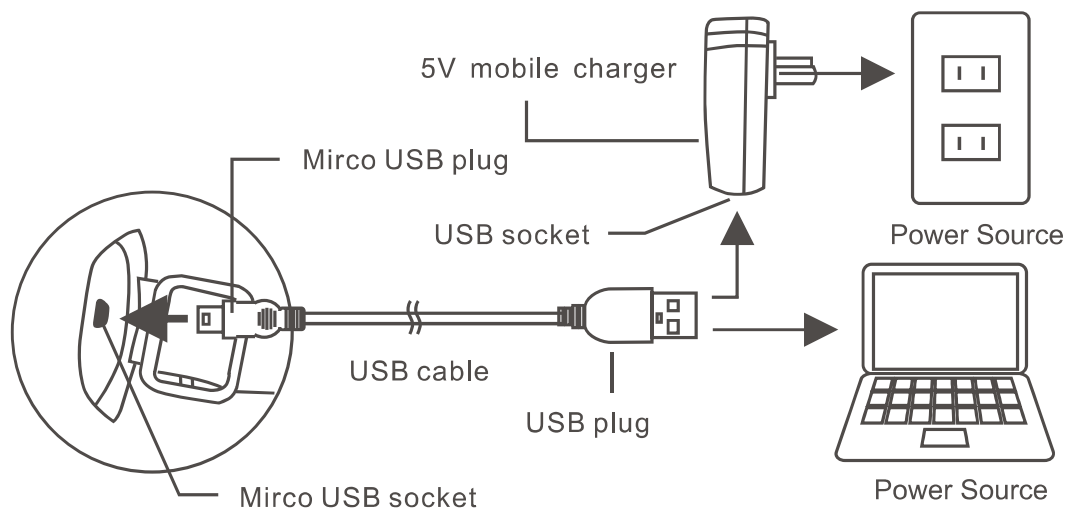
## Music adjustment

---

Press “-” button: play next track

Hold “-” button: adjust volume

Press “+” button: play previous track

Hold “+” button: enter Bluetooth pairing mode.

## Bluetooth pairing

---

Hold “+” button for 3 seconds to enter Bluetooth pairing mode, on your phone or tablet enter your Bluetooth settings and select “eye massager” to connect.

To quit the Bluetooth pairing mode, hold “+” button till the LED indicator stop flashing.

## FAQ and Troubleshooting

---

<b>Problem</b>	<b>Possible causes</b>	<b>Actions</b>
There is a sound mechanical or air sound coming from device	Normal sound of air pump inflating, motor working, or air pressure	No need to do anything
Device turns off automatically. Device does not turn on. Buttons are not responsive. Air pressure is weak.	Battery is low	Charge the device

If there are other problems, turn off the device and contact us at [support@renpho.com](mailto:support@renpho.com). RENPHO is not responsible for any damage or injuries that occur if you attempt to disassemble or repair the product yourself.

## Cleaning and Maintenance

---

Please turn off the power when cleaning, do not use corrosive detergents such as benzene or diluents to clean the product.  
**Do not** expose the device to excessive sun or hot temperature.

## Guarantee

---

This product comes with one-year warranty from the date of purchase.

Warranty will be voided due to any of the following:

Damaged due to the improper usage or storage.

Damaged by attempting to repair or disassemble the device.

Period of warranty exceeded.

## Need more help

---

Contact us via e-mail or phone

E-mail: [support@renpho.com](mailto:support@renpho.com)

Phone: 1-(844)-417-0149

Service time: PST 9.00am-16.30pm Monday-Friday

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NHJMG02

## **FCC Warning Statement**

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

## **JOICOM CORPORATION**

165 E. Freedom Ave., Anaheim, CA 92801

E-mail: [support@renpho.com](mailto:support@renpho.com)