


THERAGUN Elite

Wireless Charging Stand

Wireless Charger

Theragun Elite Wireless Charging Stand

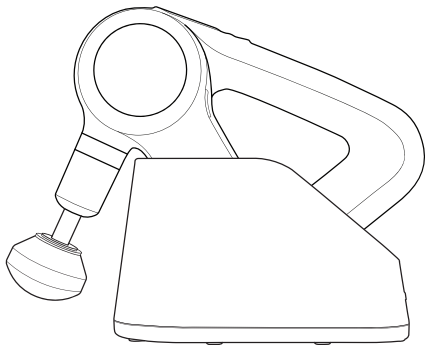

Theragun Elite not included with the Wireless Charging Stand

Charging your Theragun Elite



1. Plug the power adapter into the Wireless Charging Stand's charging port.
2. Turn off your Theragun Elite.
3. Place your Theragun Elite in the cradle of the Wireless Charging Stand and push forward until fully inserted.
4. The LED light of the wireless charging stand will light up when properly connected. The Theragun Elite is fully charged when the battery icon on the device's OLED screen is full.

The Theragun Elite can remain in the Wireless Charging Stand even when the device is fully charged.

THERAGUN

Born in Los Angeles, CA.
Change the way you *move*.