



CONTENTS

- What's in the box
- Welcome
- Important notices
- Warnings and important safeguards
- Main features
- Preparing for use
- Charging the brush
- Operating the brush
- Programming motor speed
- Cleaning and storage
- Maintenance
- Troubleshooting
- Warranty
- Replaceable brush heads

WHAT'S IN THE BOX

- Oralucent toothbrush handle with built-in red and blue LED light
- Replaceable brush head
- USB charging cable
- Wall charging plug with USB port
- Charging base
- User manual
- 1 oz • 30 ml Oralucent light-activated amplifying O₂Gel

Scan this code for the Quick Start Video
and to access our support page



WELCOME

Congratulations on taking charge of your dental health. Get ready to experience a whole new level of clean & bright.

The Oralucent toothbrush harnesses advanced light technology that outshines conventional toothbrushes. Brushing with the Oralucent phototherapy toothbrush is different from any other toothbrush you may have used, so please **read all instructions prior to using for the first time**. In just minutes, you'll be on your way to enhanced oral health.

IMPORTANT NOTICES

The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.

Unsupervised children should not play with the appliance. This toothbrush and O₂Gel, like all whitening products, is not intended for use by children under 12 years of age.

The product is a personal oral device and is not intended for sharing among multiple users. Utilize this product only for its intended use as described in this user manual. Discontinue the use of product if it appears damaged in any way.

CONSULT YOUR DENTAL PROFESSIONAL BEFORE USING THE PRODUCT IF:

- You have had oral surgery within the last 60 days
- You are currently receiving gum treatment or treatment for any oral condition
- You are taking light sensitive medications such as aminolevulinic acid or porphyrin-containing medications or have sensitivity to light
- If pregnant or nursing, before using any product containing peroxide, such as Oralucent O₂Gel
- You have medical concerns

WARNINGS AND IMPORTANT SAFEGUARDS

Do not

- Submerge the unit or charging dock in water or use while in bath or shower.
- Position the unit or charging dock in a wet area or where it is in danger of being submerged in water.
- Touch the USB port or plug with wet or damp hands.
- Attempt to open the unit, repair or replace the battery. The unit is not user serviceable.
- Drop the unit or charging base. Do not attempt to use if unit has been accidentally dropped as it may have been damaged and subject to malfunction.
- Operate the unit if it appears visibly damaged, cracked, or warped in any way, or if you suspect damage or malfunction.
- Use brush heads or accessories not recommended by manufacturer.
- Use for any other purpose than cleaning teeth and mouth.
- Point any emitting light directly into your eyes.

MAIN FEATURES

EYE COMFORT TECHNOLOGY

Patented technology turns lights on in your mouth while you are brushing your teeth, eliminating eye discomfort from bright light.

When you begin brushing, the lights in brush head will turn on. Lights shut off when removed from mouth. Note: Lights will stay on with continued brushing past 2 minutes to encourage longer brushing time.


BLUE LIGHT TECHNOLOGY

Uses the same type of blue light used in dental offices for whitening treatments.

Note: When combined with optional Oralucent light-activated O₂Gel enclosed in this kit, the antimicrobial and whitening benefits are accelerated.

RED LIGHT TECHNOLOGY

Therapeutic red light has been shown to stimulate circulation, reduce inflammation and promote cellular renewal.

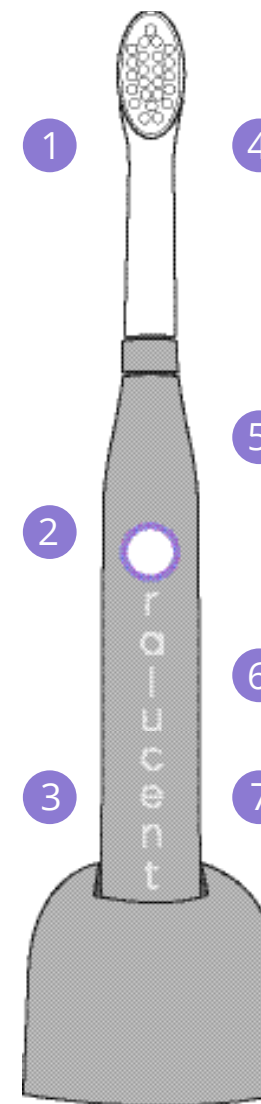

1

4

CUSTOMIZABLE MOTOR SPEED

Your brush is preset with a resettable medium motor speed. To select from the 6 motor speed levels, press and continue to hold down the button to activate the looping speed ramp. You will hear and feel the motor cycling through the speed options. When you reach your preferred speed, release the button. Your selection will lock in and can be reset at any time with the same steps.

5

MANUAL MODE OPTION

To change to manual mode (no motor), press and continue to hold down the button. You will hear and feel the speed ramp increasing, before looping back to the start at motor level 0. When the setting with no motor (level 0) is reached, release the button. Manual mode will lock in for future use.

6

GENTLE BRISTLES

The non-abrasive, soft premium bristles provide a gentle – yet thorough – clean.

7

DETACHABLE CHARGING CORD

Once toothbrush is charged, the cord can be detached from the bottom of charging base to clear counter space.

PREPARING FOR USE

1. Turn the base so that bottom is visible then insert the L-shaped end of the USB cable into the bottom of the charging base.
2. Insert the other end of the USB cable into the USB port of the round wall charger provided, and then insert the round wall charger into a standard US outlet.

Connect to power source

Position the charging base on a level, non-slippery surface in a convenient location near the desired power source. Be careful to situate both pieces so they are not in danger of falling or sliding into your sink, as neither the charging base nor the unit should ever be submerged in water.

CHARGING THE BRUSH

Before first use, fully charge the battery for at least 4 hours. After connecting the charging base to a power source with the previous 2 steps, place the brush into the charging dock.

Power button will pulse purple when charging. A solid purple button indicates that the brush is fully charged.

A full charge will last at least one week, assuming 2 minutes of brushing, twice a day. While it is not necessary to place the brush into the charging dock after each use, doing so will help it dry properly and maintain a secure, upright position. Note: It is normal for the toothbrush to be slightly warm when charging.

OPERATING THE BRUSH

Your Oralucient toothbrush is easy to use. When ready to brush, wet the bristles thoroughly and add your favorite toothpaste and optional light-activated O₂Gel. Simply press the button to start. Once you start brushing, the lights will come on in your mouth and intensify to therapeutic levels. The lights will turn off when you stop brushing. Note: A red ring around button indicates the brush needs recharging.

FOR OPTIMAL RESULTS:

2 DROPS (TOOTHPASTE + O₂GEL)

2 MINUTES

2 TIMES PER DAY

PROGRAMMING MOTOR SPEED

To customize your speed, press the button and continue to hold down to activate the looping speed ramp. You will hear and feel the motor cycling through the 6 speed options. When you reach your preferred level, release the button to lock it in.

Manual mode/No motor can be set by cycling through speed options until speed level 0 is reached. Release the button to lock it in to memory for future use. Speed can be reset at any time by following the same steps.

CLEANING AND STORAGE

After each use

Rinse the brush head thoroughly under running water until it is free of all remaining toothpaste and debris. If necessary, wipe clean your unit with a damp, soft cloth. Allow to air dry upright in its base.

When replacing brush head

Remove the brush head from the unit by grasping with one hand while holding the handle in the other hand. Pull the brush head straight off, without twisting. Wipe away any accumulated toothpaste and debris on the inner shaft with built-in light.

MAINTENANCE:

Keep the unit, charging base and all small parts in a dry place and away from children.

For the purposes of recharging the battery, only use the detachable supply unit provided with the appliance.

The appliance contains batteries that are non-replaceable.

EASY TROUBLESHOOTING TIPS

Disconnect the charging dock from power source and wipe clean with a soft, damp cloth. Allow all parts to dry before reconnecting.

Problem: The toothbrush will not operate

Possible cause: Toothbrush may require charging.

Solution: Check that the toothbrush is charged. If not, charge the brush fully.

Problem: The toothbrush is not charging

Possible cause: Toothbrush may not be properly inserted into the charging dock. The charging base may not be connected to a power source. Cables may be damaged and require replacement.

Solution: Confirm that the toothbrush is fully inserted into the charging base. Check that the charging base is connected to a power source. Safely disconnect from power source. Inspect cables for damage. If damaged, do not operate the unit. Contact customer care for replacement parts at **(877) 379-8781**.

Problem: Light is coming on before I have the toothbrush in my mouth

Possible cause: Inner shaft under brush head may have accumulated water, debris or residual toothpaste and/or gel.

Solution: Remove brush head and wipe inner shaft clean and dry. Also make sure that the plastic base of brush head where the light shines through is clear of build-up. If lights continue to come on outside of mouth, contact **customercare@oralucent.com**.

Problem: The light is not staying on while I'm brushing

Possible cause: Bristles may not be saturated with tap water and toothpaste.

Solution: Check that bristles are thoroughly wet and maintain a firm grip on the handle while brushing. If lights are still not coming on after brushing for 10 seconds, contact Customer Care.

Problem: The brush is not allowing me to program the speed I want

Possible cause: Button may not have been pressed and continuously held down to reach preferred setting in cycle.

Solution: Check that you are pressing and continuing to hold down the button, to activate the looping speed ramp before releasing the button at your preferred level.

Problem: Increased pressure is required when brushing for sufficient cleaning

Possible cause: Brush head may be worn or damaged.

Solution: Replace the brush head at least every three months, as recommended by the American Dental Association, or when worn or damaged. Visit **oralucent.com** to order replacement brush heads.

WARRANTY

Your Oralucient toothbrush has a 1-year warranty from the date of purchase.

If malfunction cannot be remedied by following the troubleshooting instruction, discontinue use immediately and contact Oralucient Customer Care at **customercare@oralucent.com** or **877-379-8781**.

Within the warranty period and with proof of purchase, we will replace the complete appliance, free of charge, for any defects in the appliance resulting from faults in materials or workmanship. If the unit is damaged, or if the product fails to function, do not dismantle, exchange, or manipulate any part of it, as there is risk of electric shock. Please contact Customer Care.

EXCLUSIONS:

Brush heads

- Damage due to improper use
- Defects that have a negligible effect on the value or operation of the appliance
- Damage due to unauthorized repair
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading
- Damage caused disassembling the product, damage caused by using an attachment not made by Oralucient, or damages caused by negligence

Changes or modifications not expressly approved by Oralucient, Inc. could void the user's authority to operate the equipment.

This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. The equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

To obtain service within the warranty period, contact Customer Care at **customercare@oralucent.com** or **(877) 379-8781** for a Return Merchandise Authorization (RMA) number and return shipping address and instructions. Do not return to the address on the original packaging.

REPLACEABLE BRUSH HEADS

Replace brush heads when worn or every three months, whichever comes first.

Oralucient replaceable brush heads are designed exclusively for use with the Oralucient toothbrush. To maintain maximum effectiveness, it is important to replace them regularly – before they become worn.

Detaching the brush head:

- Grip the brush handle firmly in one hand
- Use two fingers on your other hand to pull brush head straight off of the brush handle – do not twist

Attaching the brush head:

Slide the open end of the brush head onto the handle shaft, aligning the brush head so the bristles face the same direction as the button on the handle, and push to secure. You should hear a click when the brush head is properly attached.

To protect brush heads against premature wear or damage, follow these tips


- Rinse brush heads well after each use.
- Never rinse the brush heads with hot water. Doing so may damage the fibers or warp the entire brush head.
- Do not exert unnecessary pressure while brushing. Use only the amount of pressure required to maintain contact with your teeth and gums.
- Allow brush heads to air dry thoroughly between uses.
- To reorder brush heads or any Oralucient products, visit **Oralucent.com** or call **(877) 379-8781**.

DISCLAIMER

This product is not intended to diagnose, cure or prevent any disease. Oralucient, Inc. makes no claims, representations or warranties regarding the ability of this product to cure any physical ailments. A qualified health professional should be consulted with regards to any condition requiring medical attention.

This product contains a lithium-ion battery that must be disposed of properly according to local laws and regulations, and not as household waste. When the product reaches its end of life, check into your state and municipal disposal and recycling options.



 **oralucent**
CUSTOMER CARE

 **oralucent.com**

 **customercare@oralucent.com**

 **877- 379-8781**

Distributed by Oralucient, Inc.
Long Beach, California 90803