Venom Go

EN Operating Instructions BG Инструкции за експлоатация CS Návod k použití
DA Betjeningsvejledning DE Bedienungsanleitung EL Οδηγίες λειτουργίας
ES Instrucciones de uso ET Kasutusjuhised FI Käyttöohjeet FR Mode d'emploi HR Upute za rad
HU Használati útmutató IT Istruzioni per l'uso LT Naudojimo instrukcijos MS Arahan operasi
NL Gebruiksaanwijzing PL Instrukcje użytkowania PT Instruções de operação
RO Instrucțiuni de utilizare RU Инструкции по эксплуатации SK Návod na použitie
SL Navodila za uporabo SV Driftanvisningar TH คู่มือการใช้งาน ZH-TW 繁體中文 操作說明

EN

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

READ ALL INSTRUCTIONS BEFORE USING THE VENOM GO.

DANGER

To reduce the risk of electric shock:

- Always unplug this device from the electrical outlet immediately after charging and before cleaning or using.
- Do not reach for a device that has fallen into water. Unplug immediately.
- Do not use while bathing or in a shower.
- Do not place or store device where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
- Keep dry. Do not operate in a wet or moist condition.

WARNING

Read all safety warnings and instructions before using this device. Failure to follow the warnings and instructions may result in electric shock, fire and/or serious injury.

WARNING

TO REDUCE RISK OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- Use this device only for its intended use as described in this manual.
- This is not a toy. For adult use only. Do not use if injured.
- This device should never be left unattended when plugged in. Always unplug this device from the electrical outlet immediately after charging, prior to use and before cleaning. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Never operate this device if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or has gotten wet. Return the device to a service center for examination and repair.
- Do not carry this device by the cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Never drop or insert any object into any opening.
- Suitable for dry location ONLY. Do not use outdoors or on wet surfaces and do not expose to water or rain.
- Use product only according to these instructional materials. In addition, instructional use materials are available at hyperice.com.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

- Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the device frequently to reduce the risk of blistering. Immediately remove from body and turn off at first sign of discomfort.
- This device must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- If the device is worn while sleeping with the controls set to any temperature the user may suffer skin burns or heat stroke.
- Do not crush. Avoid folding with sharp creases. Do not fold Venom Go pad while in use.
- Do not modify or attempt to repair the device. Contains no user-serviceable parts. Please contact Hyperice for service.
- Prevent unintentional starting. Ensure the power button is in the off position before picking up, carrying, or storing the device.

CAUTION

When using an electrical device, basic precautions should always be followed, including the following:

- For adult use only.
- This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge of how to use this device, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- Children under the age of three are not to use this device due to their inability to react to overheating.
- The device has a heated surface. Persons sensitive to heat must be careful when using the device.
- This device must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- Children should be supervised to ensure that they do not play with the device.
- Unplug the charger from battery after charging or prior to use. Venom Go will not function while charger is plugged into battery.
- Never leave the Venom Go operating or charging unattended.
- Do not immerse device in water. Clean only with a damp cloth.
- Keep liquids away from device.
- Use only on a dry clean surface of the body as shown in this instructional material.
- Do not use Venom Go on wet skin or over topical ointment, especially ones that contain heat producing ingredients as skin burns could result.
- Burns can occur regardless of control setting. Check skin under product frequently and immediately remove from body and turn off at first sign of discomfort.
- Do not place any metallic objects between Venom Go and body, including pins or other metallic fasteners.
- Do not attempt to disassemble.
- Do not tamper or alter the Venom Go in any way.
- Carefully examine the device before each use. Contact manufacturer if it shows any sign of deterioration, such as checking, blistering, or cracking.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not sit on, press against, or crush Venom Go.
- Do not crush the device by placing items on top of it during storage.
- Cleaning and user maintenance shall not be made by children without supervision.

PLEASE DO NOT USE THE VENOM GO WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY:

- Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device.

These contra indications do not mean that you are not able to use a vibration or exercise device but we do advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

SAVE THESE INSTRUCTIONS

Symbol	Description	Location
IP21	Degree of protection against ingress of water	On pod
	Read instructions before use	On pod
	Direct current	On pod and in manual
	Manufacturer's name and address	On pod
	Separate collection for waste electrical and electronic equipment. This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmentally safe recycling.	On pod
	Caution symbol to indicate the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself.	On pod and in manual
	European Conformity	On pod
	Australia RCM	On pod
	Made in China	On pod
	UK Mark	On pod
	North American NRTL Mark	On pod
	Recycle symbol	On pod
	Heat control	On pod
	Power on/off	On pod
	Vibration control	On pod

BLUETOOTH® WIRELESS TECHNOLOGY

The Bluetooth® word mark and logos are owned by Bluetooth SIG, Inc., and any use of such marks by Hyperice is under license. In the unlikely event of loss of a stable Bluetooth connection, the system will attempt to re-establish its connection automatically. The Venom Go is completely autonomous, and will continue operating normally, even during a loss of connectivity. If this device does cause

interference, which can be determined by turning the device off and on, the user is encouraged to try to correct the interference by reorienting or relocating the device, increasing the separation between equipment and the device, or connecting the device to a different outlet on a circuit if it is plugged in.

The Venom Go device uses Bluetooth 5.0 wireless technology with the following radio specifications: FCC ID: 2AWQY-VENOMGO, IC: 26929-VENOMGO, Frequency: 2402 to 2480 MHz

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

This equipment complies with FCC/ISED radiation exposure limits set forth for an uncontrolled environment and meets the FCC radio frequency (RF) Exposure Guidelines and RSS-102 of the ISED radio frequency (RF) Exposure rules. This equipment has very low levels of RF energy that are deemed to comply without testing of specific absorption rate (SAR).

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer for help.

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INTERNAL BATTERY INFORMATION

The Venom Go device is equipped with a rechargeable lithium polymer battery. The internal battery is designed to allow use of the Venom Go anywhere—even when power outlets aren't available. The Venom Go may need to be charged before first use. A fully charged battery will provide power for 35+ minutes of continuous use. It takes approximately 2.5 hours to fully charge the battery when the device is plugged in and not in use.

- This device contains a lithium-polymer battery. Batteries are only replaceable by a skilled worker with a special tool. Contact Hyperice for service or repair. When the battery is at end of life, the device shall be properly disposed of. The battery must be complied with safely at an appropriate e-waste disposal or recycling facility.

WARNING

- Do not charge the Venom Go overnight or leave the Venom Go unattended when charging or in use.
- Follow all charging instructions and do not charge the battery pack or device outside of the temperature range specified in the instructions. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
- Nominal temperature range for charging the battery is 0°C ~ 45°C (32°F ~ 113°F).

- Under abusive conditions, liquid may be ejected from the battery; avoid contact. If contact accidentally occurs, flush with water. If liquid contacts eyes, seek medical help. Liquid ejected from the battery may cause irritation or burns.
- Do not use a battery pack or device that is damaged or modified. Damaged or modified batteries may exhibit unpredictable behavior resulting in fire, explosion or risk of injury.
- The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, do not expose a battery pack or device to fire or excessive temperature. Exposure to fire or temperature above 130°C (265°F) may cause explosion.
- In the event of the battery pack leaking, do not allow any liquid to come into contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical attention.
- In case of battery rupture, major leakage or explosion, evacuate the area and quarantine the contaminated area. Provide good ventilation to clear out any fumes, gases or pungent odors.
- Eyes – Rinse eyes with plenty of water for 15 minutes; Seek immediate medical attention.
- Skin – Rinse affected area with plenty of water and soap or take a shower for 15 minutes.
- Inhalation – Expose the person to fresh air and use artificial respiration if needed; Seek medical attention if necessary.
- Ingestion – Consult a physician or local poison control center immediately
- Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

CHARGING

To charge, connect the USB-C cord into the battery's charging port and plug the cord into a USB-A adapter. Then plug the adapter into the wall. Fully charge the battery for up to three hours before first use. A LED light pulses to show the battery level and indicate active charging. Full charge is indicated when the LED light stops pulsing and remains fully illuminated green. Pressing the power button will turn on the device and will display the battery level. Device may be recharged at any time and at any battery level. It is not recommended to fully deplete the battery to the orange LED level. Recharge the device monthly to best prolong battery life.

OPERATING

Ensure battery is charged. Peel the Venom Go pad from the liner and place the pad on your desired body part. On the pod, press the power button to activate the device. The Hyperice logo will illuminate red to show Venom Go is activated. Choose the desired heat and vibration options. Then connect the pod to the pad using the magnetic connection. The device will automatically start when the pod connects to the pad. You can adjust the heat or vibration to the new levels at any time. The device will turn off automatically after 10 minutes without further input. The device will vibrate twice to signify the end of a session or will vibrate three times when it has a low battery. When you are done using the Venom Go, adhere the pad to the liner for storage. The average run-time is around 35 minutes at level three heat, and longer at lower levels. Each pad lasts up to 20 applications.


CONNECT TO THE HYPERICE APP

To connect your device to the Hyperice App via Bluetooth® make sure that your device is turned on, Bluetooth® is turned on in your phone, and your device is within range of your phone. The Bluetooth LED on the pod will pulse blue, which means it is waiting to pair. If prompted, tap "Scan for Devices." Select your device when it pops up on the screen. Once you pair your Venom Go, HyperSmart™ will allow you to automatically start your session and control the heat and vibration settings from within the Hyperice App.

MAINTENANCE, CLEANING, AND STORAGE

Reapply pad to liner after use and store in a cool, dry location. Use a towel or soft cloth to wipe any moisture from device. Keep device away from heat sources, liquids, children, and pets.

SPECIFICATIONS

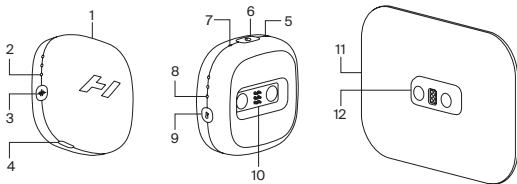1. Venom Go pod
2. Vibration pattern indicators
3. Vibration pattern button
4. Charging port
5. Battery life indicator
6. Power button (ON/OFF)
7. Bluetooth® indicator
8. Heat level indicators
9. Heat level button
10. Pad attachment magnetic connection
11. Venom Go pad
12. Device attachment magnetic connection

Vibration frequency: 88 Hz

Heat: Level 1—105°F (45°C), Level 2—109°F (42.8°C), Level 3—113°F (45°C)

Weight: Pod—2.5 oz., Single pad—0.6 oz., USB-C cord—0.5 oz.

CAUTION

Not intended for medical use in hospitals. Intended for household or commercial use.

When not in use, it is recommended to store the device between -10°C ~ 25°C.

WARRANTY

This product is covered by a limited warranty from Hyperice. Please visit hyperice.com/warranty to review the warranty in your country.

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hyperice.com

Hyperice, Inc., 525 Technology Drive, Suite 100
Irvine, California 92618 USA

Copyright 2022 © Hyperice, Inc. All rights reserved.



170-00167 Rev 2

**Let's
get
started**


Meet your new Venom Go

Front view

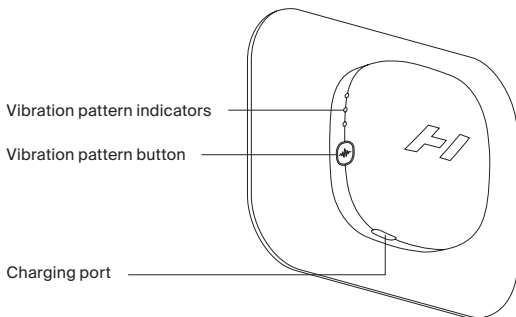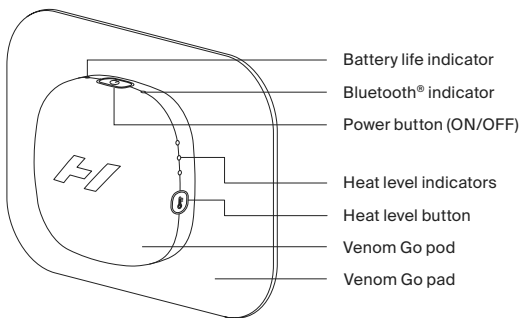Back view

Magnetic connection

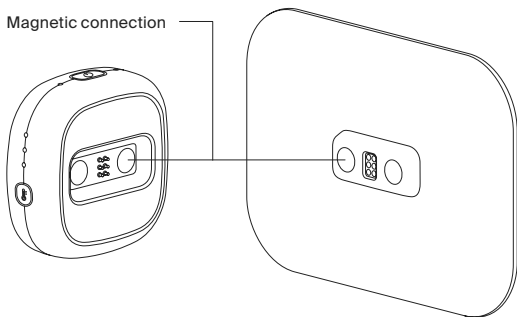
Please read the Operating Instructions before first use.

Before you start

Download the Hyperice App

Get the most out of your Venom Go, or any Hyperice device with the Hyperice App. Automatically stop and start your sessions and remotely control the heat and vibration.

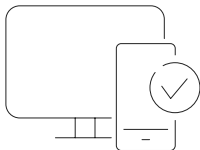
Register your device

Activate your warranty and ensure easy returns, repairs, or refunds at: hyperice.com/register-product

Charge your new Venom Go


Plug in your new Venom Go with the provided Hyperice USB-C until fully charged.

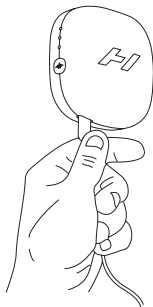Start your session

Power up


Turn your Venom Go ON by pressing the power button (top) once.

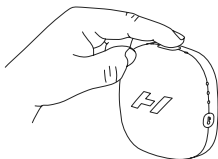Apply pad and connect pod

Peel the Venom Go pad from liner and apply to skin. Select heat level and vibration pattern. Connect pod by using the magnetic connection (shown below / right).


Note: Heat and vibration will not be felt until pod is connected.

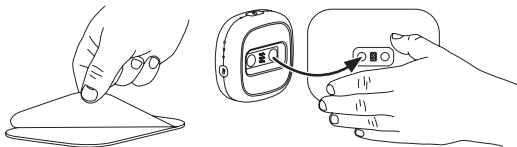WARNING: Only use on clean, dry skin. Do not use while perspiring. Do not dispose of pad liner. Reapply pad to liner after use and store in a cool, dry location.

Adjust temperature

Choose your desired temperature by pressing the heat level button (right side) once per level. The white lights will indicate the selected temperature.


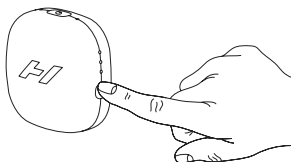
Add vibration

Choose your vibration by pressing the vibration pattern button (left side) once per pattern. The white lights will indicate the selected vibration pattern.

Enjoy your session

Feel the soothing relief of the Venom Go. Your session will automatically shut off after 10 minutes. Restart as needed.

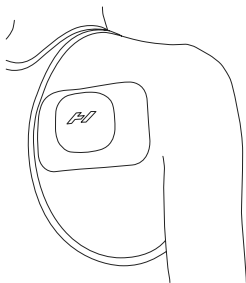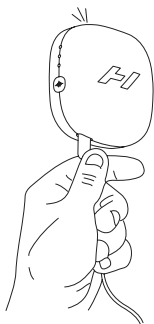Recharging

Full charge is indicated when battery indicator light on the pod remains green and fully illuminated.

Caring for your Venom Go

Make sure the power is OFF and the battery charger is not attached.
Use a damp, clean cloth to gently wipe pod and pad (if necessary).
Reapply pad to liner after use.

Store in a cool, dry area, out of direct sunlight when not in use.

Discard pad when it no longer adheres to skin. Additional pads can be purchased at hyperice.com.

We're here for you

HyperCare®

Get award-winning support from our HyperCare team – a group of specialists dedicated to your overall well-being and expert guidance on Hyperice products.

1.855.734.7224

hyperice.com

customersupport@hyperice.com

If outside of the United States, please visit hyperice.com/contact