

After-sales Service

**Thank you for purchasing Nicefeel oral irrigator.
Please read the following precautions:**

1. Please ask the seller for a valid sales document as the basis for the warranty when purchasing the oral irrigator.
2. The product is covered by free maintenance for failure due to non-human factors within one year according to relevant state regulations (please register by our support email with your order number for the warranty);
3. Costs and service fees will be charged as appropriate during the warranty period in the following cases:
 - (1). Damage caused by improper use, maintenance and storage by the consumers;
 - (2). Damage caused by dismantling of unauthorized repairer;
 - (3). No relevant sales document and it has been over one year after the production date;
 - (4). Damage caused by commercial use;
 - (5). Damage caused by force majeure (earthquakes, fires, etc.).
4. Please fill in the form below with care, and send it together with the product to be repaired to the service site designated by the manufacturer.

Warranty Record


Customer name		Tel	
Product model		Product color	
Date of purchase		Where the product is purchased (seller)	
Maintenance date			
Mailing address			

Fly Cat Electrical Co.,Ltd.

Bldg 28th, Rm. 101- 401 and Bldg 30th, Rm. 101, 301&401,
Longwangmiao Industrial Park, Baishixia East Community,
Fuyong Street, Bao'an District, Shenzhen, Guangdong, China
Tel: (86)0755-2729-8495
Fax: 0755-2727-3346
Website: www.fly-cat.com.cn

Nicefeel[®] Oral Irrigator

FC2520 Series Instruction Manual

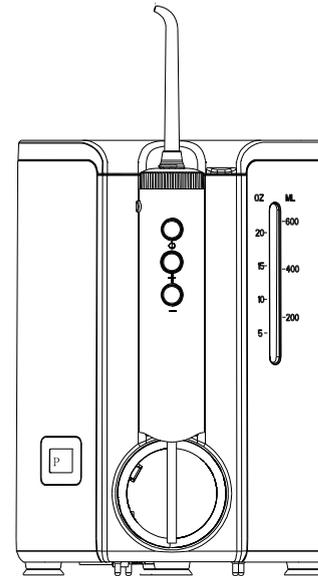

Important Safety Precautions

Basic safety precautions should be followed when using electrical products, especially when there are children in the house.

Please read all the outlined instructions before using

Dangers

To reduce the risk of electric shock, please pay attention to the following items:

- DO NOT handle the charger with wet hands.
- DO NOT immerse the product in water or any other liquid.
- DO NOT use the product while taking a bath.
- DO NOT place or store the product in a place where it may fall into the bathtub, shower or sink.
- DO NOT pick up the product falling into the water immediately, make sure it is unplugged immediately before taking any other action.
- Check whether the adaptor is damaged before using the product for the first time after purchase.

In order to reduce the risk of ambustion, electric shock, fire, serious injury or damage:

- DO NOT plug the product into an adapter that does not meet the specified voltage.
- DO NOT use the product if the cable or plug is broken. If it does not work normally or is immersed in any liquid, please contact your local dealer.
- DO NOT lead the water flow under your tongue, or into your ear, nose or other minor organs. This product is capable of producing pressure that may cause serious damage to these organs. Please refer to the operating instructions for proper use.
- Please follow the instructions or recommendations from dental professionals when using the product.
- A small amount of bleeding is normal when using the product for the first time. If the bleeding is serious, stop using and ask your dentist for advice.
- Fill the reservoir with water or other solutions recommended by the dental professionals.
- DO NOT drop or insert any foreign objects into the product.
- Keep chargers and wires away from hot surfaces.
- DO NOT operate the product in places where oxygen or aerosol sprays are used.
- DO NOT use mouthwash or rinse solution containing iodine or disinfectant.
- DO NOT use the product when wearing dental accessories, which should be remove before use.
- DO NOT use the product if there is an open wound on the tongue or in the mouth.
- If your doctor recommends you to take antibiotics before surgery, please consult a dentist before using the product or any other oral hygiene adjuvant.
- For household use only, DO NOT use the product outdoor.
- Children and individuals with special needs should be closely guided and supervised for using the product properly.
- Children should be supervised to ensure proper use of the product.

Please keep this User Manual!

Thank you for purchasing our Oral Irrigator!!

Thank you for purchasing our oral irrigator!

Welcome to use Nicefeel oral irrigator, which will give you healthier teeth and brighter smile!


Cleaning teeth with pulse water flow, the oral irrigator can remove harmful bacteria hidden deep in the crevices between teeth and below the gum line, where cannot be cleaned by brushing teeth only.

Traditional oral care:

Traditional toothbrushes cannot remove food debris and harmful bacteria hiding deeply.



They can only reach a limited part below the gum line.

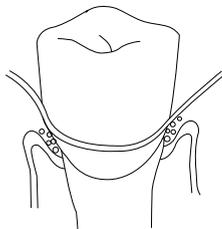

Nicefeel Oral Irrigator

Remove harmful bacteria hiding deep in the crevices between teeth and below the gum line.


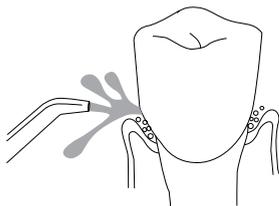


Table of Content

Important Safety Precautions-----	1
Contents-----	4
Product Structure-----	5
Function nozzle-----	6
Using Instruction for Product-----	7
Using Instruction for Professional Nozzles-----	11
Cleaning and Maintenance-----	13
Common Failures and Handling Methods-----	14


The Company reserves the right of final interpretation for any error in or misunderstanding of the content of this User Manual; any product or data update will be incorporated into the new version without prior notice; and in case of any change in the appearance and color of the product, the product inside the package shall prevail.

Product Structure

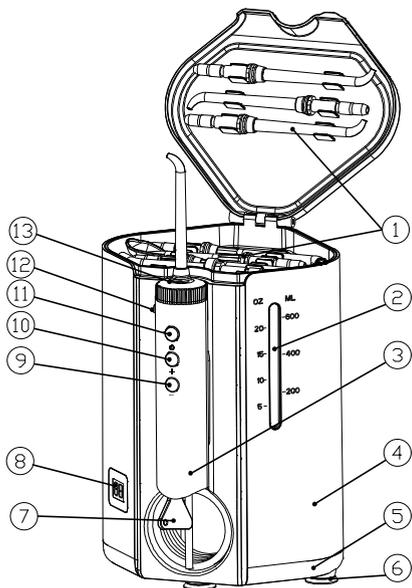

- | | | |
|--------------------------|-----------------------------|-------------------------------|
| 1. Nozzle combination | 2. Water tank (hidden) | 3. Remotely controlled handle |
| 4. Body | 5. Base | 6. Base suction cup |
| 7. key | 8. Gear display | 9. Turn-down button |
| 10. Turn-up button | 11. On/off and pause button | 12. Nozzle removal button |
| 13. 360° rotating nozzle | | |

5

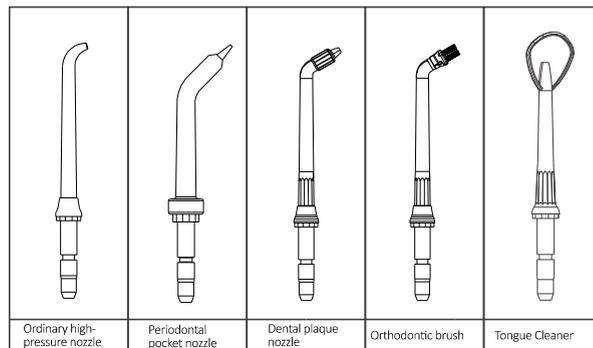

* Ordinary high-pressure nozzle: comprising three colored rings for three persons. (Family Pack)

* For the using instructions of other four special nozzles (tongue scraper, periodontal pocket nozzle, orthodontic brush and dental plaque nozzle), please refer to pages 11-12 of the Manual.

Gum bleeding is normal due to different tolerance when using the water flosser for the first time. Please feel free to use!


6

Using Instruction

Plug the wire into the socket. If the socket is controlled by a wall switch, make sure the switch is turned on. The on/off and pause buttons and two gear buttons are in the front of the product handle, and the latter (pause button) can remotely control the gear so as to control the water pressure. After the power is connected, a small display screen on the left side of the device's body will display 'P (Pause)'. Press the On/Off button to turn on, and gears 01-10 will be displayed on the screen, from which you can choose one to control the water pressure during use. A timing of 3 minutes is set between turning-on and off (Note: You need to press the On/off button again to operate the product after it is turned off by timing).

Preparing water tank


Remove the water tank cap from the top of the device. After filling the tank with warm water, put the cap back and press it firmly. The water tank is hidden inside the body, on the right side of which there is a long window where you can clearly see the water level inside the water tank. The surface of the large body is marked with capacity in ounces and milliliters.

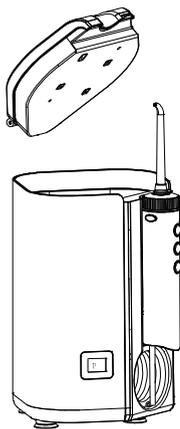

7


Insert and remove the nozzle

Insert the nozzle into the middle of the knob at the top of the handle of the water flosser. If the nozzle is locked in place, the colored ring shall be flush with the end of the knob. To remove the brush from the handle, please press the nozzle eject button and pull the nozzle out of the handle.

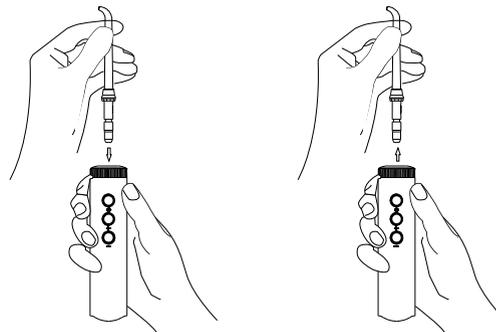

Adjust pressure settings



When using the product for the first time, please adjust the pressure gear on the handle of the product to the minimum value (the gear displayed on the body can be adjusted in gears 01-10). Start from the lowest gear, turn on the switch on the handle, and gradually increase the pressure after the oral cavity adapts to it. You can also increase the pressure under the guidance of a professional dentist.

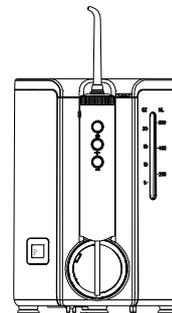

8


Using Oral Irrigator

Lean over the sink and place the injector head into the mouth. Align the nozzle to your teeth. Turn on the pushing switch on the handle to adjust the water pressure to the desired value. Control the direction of water flow by turning the knob on the top of the handle.

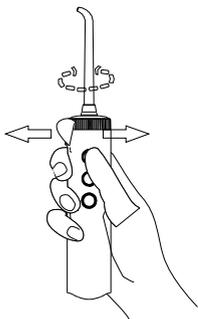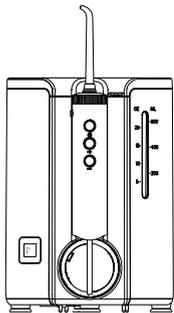

Use Tips


Adjust the water jet flow to 90 degree to the gum line. Close your lips to prevent the water from flowing sideways while allowing the water in your mouth to flow freely into the sink. For best effect, start cleaning from the molar teeth (back teeth) towards the front teeth. Slide the injector head along the gum line with short pauses between teeth. DO NOT stop until the areas around the gums and among the teeth are cleaned.

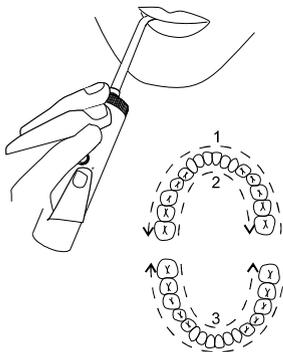

9


Pause control

You can turn off the ON/OFF button on the handle of the oral irrigator at any time and put the handle back on the body to temporarily stop the water flow.

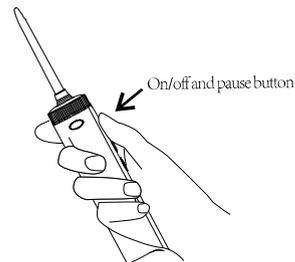

After using

Turn off the product, pour off the liquid left in the tank, and pull off the plug. When the device is turned on again, it will restore the gear before turning off based on the memory of the main body.

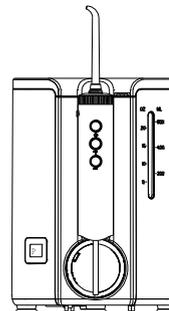


10

Using Instruction for Nozzles


High-pressure nozzle

Rinses deep into the mouth thoroughly to clean the parts hard to be cleaned with toothbrush. More effective than floss and making the gums healthier.

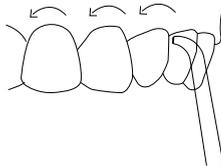


Periodontal pocket brush

The periodontal pocket brush can inject water or antibacterial fluid into specific parts below the gum line. When using the periodontal pocket brush, please lean over the wash basin, place the soft brush head against the teeth at an angle of 45 degrees, and gently place the brush head into the teeth socket below the gum line. Turn on the device and minimize the pressure, gently place the brush on the crevice between teeth, and then clean it along the direction of the gum line. Continue to move along the gum line until the flushing is completed.

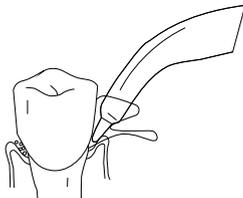

Tongue Cleaner


Lean over the wash basin. Place the tongue cleaner on the tongue and start the oral irrigator: use the lowest pressure first and then increase the pressure based on the degree of comfort. Gently place the tongue cleaner in the middle of the back part of the tongue and gently pull it to the front against the tongue.

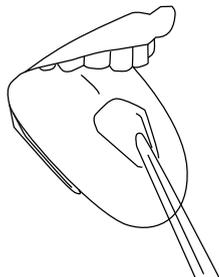

11


Orthodontic nozzle

The special orthodontic nozzle is designed to wipe and rinse the areas hard for braces and other dental appliances to reach. When using the orthodontic nozzle, lean over the wash basin to wash your teeth with the nozzle, and water flow will take the debris away.

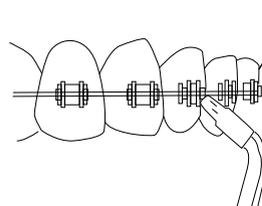

Dental plaque nozzle

Align the water with the gums at an angle of 90 degrees and slightly close your lips to prevent water from splashing while allowing the water to flow out from the corners of the mouth to the wash basin. Start flushing from the molar teeth (back teeth) and gradually to the front teeth. The nozzle shall be close enough to the teeth so that the bristles can gently touch the teeth. Gently move the nozzle along the gums with a short pause between two teeth to allow the water to pass between the teeth. Repeat the same procedure on each side of the teeth until all the part between the teeth as well as the front, back, top and bottom of the teeth are cleaned.

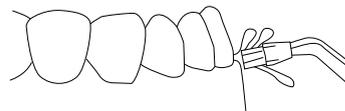

12

Cleaning:

If any special solution is used in the water tank, please rinse the device properly. Fill some warm water into the water tank.

Point the nozzle toward the wash basin. Turn on the device and do not stop it until the water tank is empty. This will help to prevent blockage of the nozzle and extend the service life of the device.

Clean the product with a soft cloth and a mild non-abrasive cleaner should be used when needed.

Be sure to pull off the plug from the power socket before cleaning.

Remove the water tank and run the oral irrigator till there is no residue left.

Maintenance

1. Each time after using, be sure to pull off the power plug and place the water tank upside down on a flat surface. After the water has been completely drained and wiped dry, place the water tank in a cool and dry place together with the cleaned main body. When not in use for a long time, the product shall be stored in a packing box or a carton with a similar storage function.

2. In case of machine failure, you can consult with the seller or call the Company to inquire about the repair method.


3. Do not disassemble or attempt to dissolve the product by yourself, as the damage caused by such behavior is not covered by the Company's warranty.

4. When the power cord is damaged, please go to the maintenance place designated by the Company or the repair shop with relevant qualification for repair. DO NOT replace the power cord or plug by yourself.

5. If the warranty period expires, the company will still provide relevant service, with appropriate material cost and labor fees charged.

Replace handle battery

If the battery requires replacement, use the key to align the locking hole on the bottom of the handle, and then turn it. The battery in the handle pops out. Each time you change the battery, the pressure is the same as the last time. If it is not used for a long time, remove the battery if the battery fluid damages the handle.

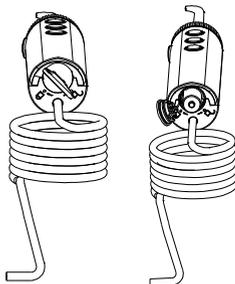

13

Common Failures and Handling Methods

Phenomenon	Reason	Handling Method
Leakage between the nozzle or the tongue cleaner and the handle knob	The brush does not fit well with the handle	Remove the brush and insert it well
Cracking of periodontal pocket brush	The brush is seriously worn	Replace a new brush
Leakage on the hose connecting the base to the handle of the oral irrigator	Hose damage	Replace it by yourself or contact dealer for after sale service
Water tank leaking	The spring is not installed or is installed at incorrect position	Install the spring in the correct position, or contact the dealer for repair
Pressure is not enough	Improper installation of the water tank	Reinstall the water tank and press on the water tank firmly
Device cannot start	Motor does not run	Confirm whether there is power input, or contact the seller for repair
Buttons on the handle cannot be used	Check whether the battery fails or is assembled incorrectly	Assemble or replace a new battery

We guarantee that this product is free of defects in materials and processes, and the warranty period is one year from the date of purchase. Please keep your receipt as proof of the date of purchase. If the product is not abused, misused, altered or damaged after purchase but used in accordance with the User Instruction, we will replace any defective part(s) of the product at our sole discretion. This Limited Warranty does not cover the attachments such as standard nozzles, tongue scrapers, periodontal pocket brushes, orthodontic nozzles or dental plaque nozzles. This warranty gives the user a clear legal rights, and you may also have other appropriate rights depending on country of residence.

Thank you again for purchasing Nicefeel Oral Irrigator!

14

FCC Warning

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE 1: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE 2: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.