


## VISUAL BLACKHEAD VACUUM DEVICE

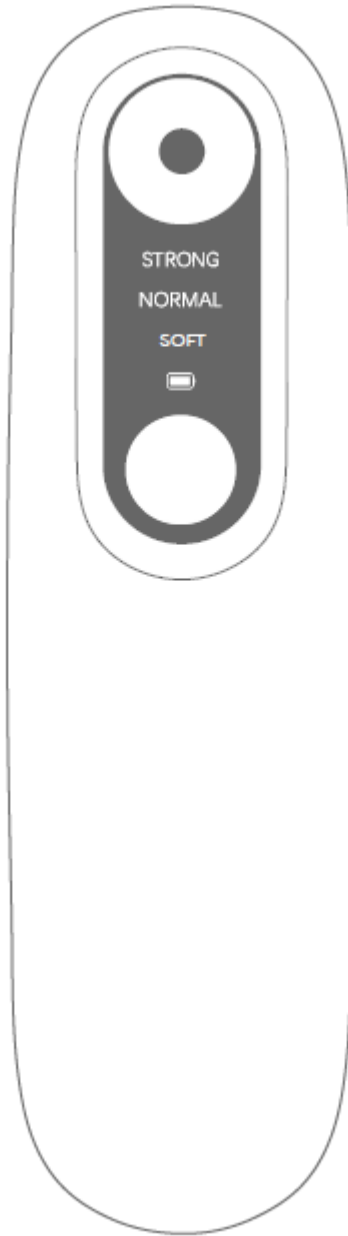

**Read the manual carefully before use**

# Manual

Thank you very much for purchasing this product. In order to ensure correct and safe use, please read carefully and keep this manual carefully before use.

## Attention

In order for you to use this product correctly, please read and follow the following safety guidelines carefully.

1. This product is not waterproof. Please do not immerse the product in any form of liquid. Do not use the product near water-bathed bathtubs, shower trays, washbasins or other water-containing utensils. Please keep the product and power adapter dry.
2. Before using this product, please dry the water or other liquids on your hands to prevent the water or other liquids from entering the interior and damaging the machine.
3. Do not overuse in an area, it may cause the skin to become red and swelling. It is recommended that the stay time in the same area should not exceed 3 seconds, and each nursing time should not exceed 5 minutes.
4. Do not use on injured or fragile skin, and on the chin and mouth, Use in fragile areas such as corners and around the eyes.
5. Do not use on skin with open wounds, acne or pimples.
6. Do not use it on cuts, areas that have undergone cosmetic surgery, areas where the skin is inflamed, itchy and red.
7. If you feel uncomfortable while using it, please stop using it immediately.
8. Continuous use of the body for a long time will generate heat, this phenomenon is not a product failure, Please rest assured to use.
9. Continuous use of the body for a long time will generate heat, this phenomenon is not a product failure, Please rest assured to use.
10. If you have any skin problems, please consult a doctor before using this product
11. Please put the product out of the reach of children, do not let children use, it will cause accidents and injuries.

## Introduction

The visual blackhead instrument uses non-invasive vacuum negative pressure technology, with 3 special suction heads, which can help effectively remove blackheads, pimples, dirt and aging keratin, clear pores and improve muscle skin. Built-in high-definition camera, which can directly monitor the skin condition and clearly observe the blackhead quilt. During the removal process, avoid blindly cleaning or leaving clean skin.



标准吸头




精准吸头




区域吸头



密封圈\*2

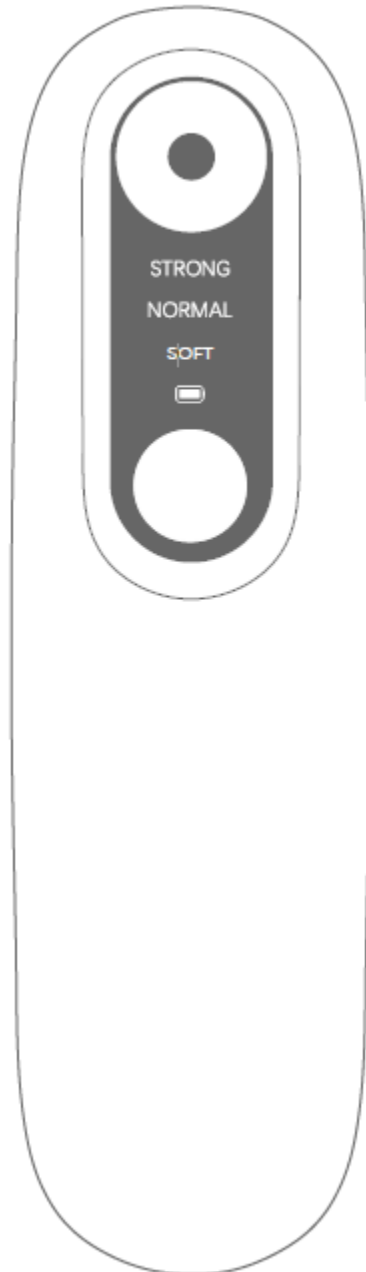


## Structure

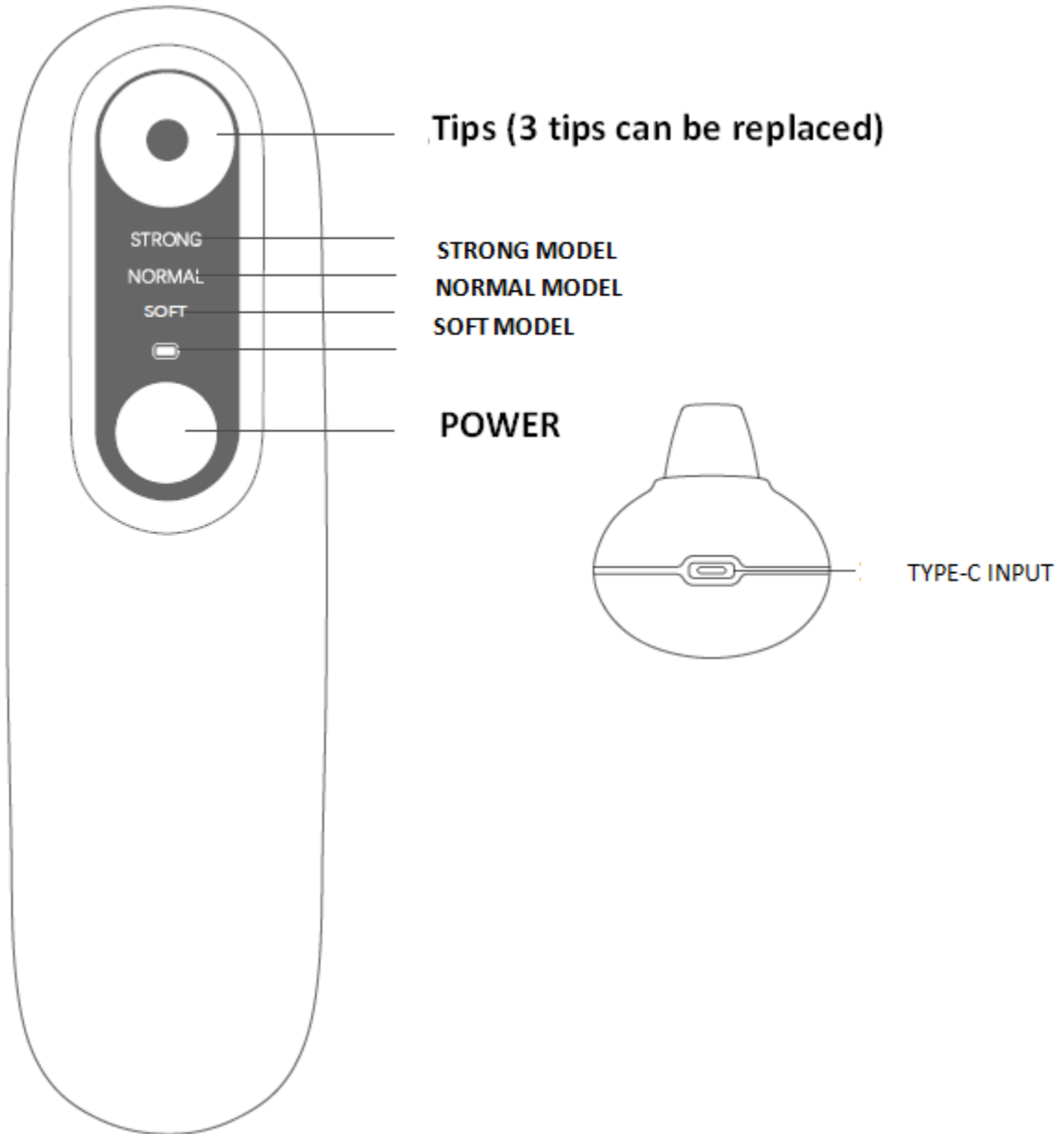

## Function

Press and hold the power button for 2 seconds, the instrument is turned on, mode & power indicator Light.

- Shortly press the button after switching on to switch the mode: Soft (soft)/ Normal / Strong/ Standby.
  - When the battery is low, the battery indicator light flashes slowly.
  - The white indicator of the battery indicator flashes slowly during charging, and the indicator of the battery indicator when fully charged
- The white light is always on.
- Each gear works for about 5 minutes. Standby mode after 5 minutes, wait
- The machine will automatically shut down if there is no operation within 10 minutes.
- Press and hold the power button for 2 seconds to shut down.

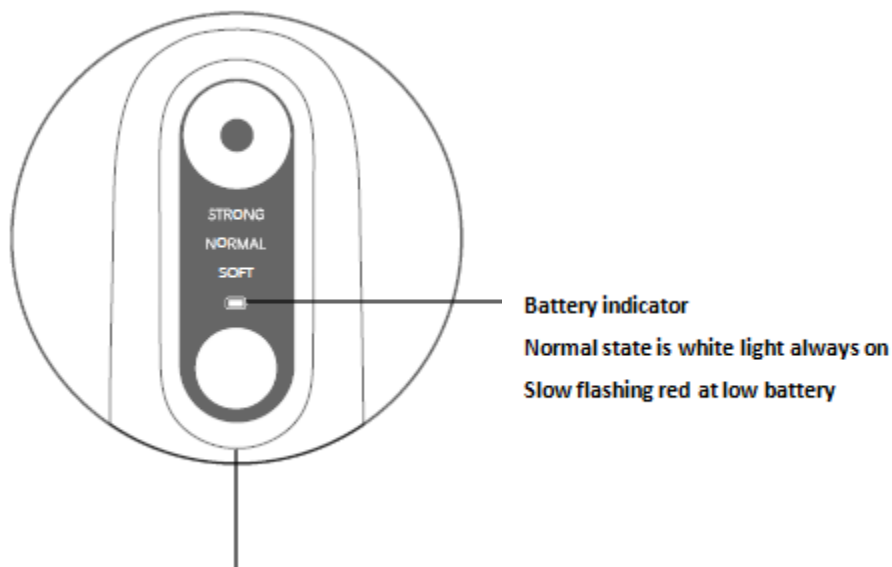

Press and hold for 2 seconds to turn the machine on or off  
Short press cycle to switch working mode

## Procedure

### APP download and installation

After scanning the QR code in the user manual with your mobile phone, click to download and install.

(If you use WeChat to scan, you need to manually jump to the browser to open)



## APP function

The camera connected to the instrument through the APP is more straight during use  
Easy to clean pores.

- ① After turning on the instrument, set the WIFI connection ID on the mobile phone, such as: UDU Beauty-1234.
- ② Open the "UDU" APP and enter the home page. If the device is not connected, the APP will The operation guide prompt appears.
- ③ After clicking "Open", you can see the transmission of the instrument camera in the APP window image.



- ④ When using, you can adjust the work of the blackhead according to the skin condition of the face position.
- ⑤ After clicking "Video", the image transmitted by the camera of the instrument can be saved to the mobile phone.
- ⑥ After clicking "Take Photo", the image transmitted by the instrument camera can be taken as a photo Save to phone.
- ⑦ If you need to switch windows, you can click "AR" to switch windows.

## Nursing

Download and install the APP (see APP download and installation steps for details)

② Connect the instrument through APP (see APP function description for details)

③ Clean the face before use, apply a hot compress or use a steamer to help pores

Open, and then use the blackhead export solution to soften the blackhead, you can easily suck out the blackhead.

④ Choose the tip suitable for your needs, install it on the top of the instrument and on the palm or forearm


Test the suction power.

⑤ When in use, you can press the button to select the appropriate mode according to your skin condition.


Soft: Soft mode; Normal: Normal mode; Strong: Strong mode.

⑥ Make sure that the suction head is close to the skin when using it, and there must be no gaps, otherwise it will not form a vacuum

Negative pressure.

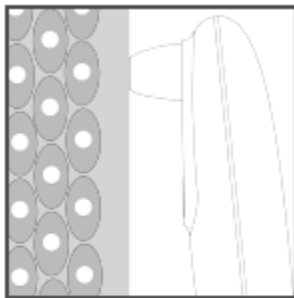

✓

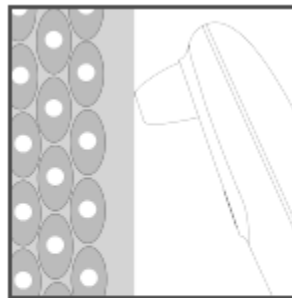

✗

⑦ Keep moving when using the instrument, the stay time in the same part should not exceed 3 seconds,

Each use time does not exceed 5 minutes, enter the standby mode after 5 minutes, standby 10

minute

No operation will automatically shut down.

⑧ After using, long press the power button to turn off the instrument, you can apply ice to the skin,

Or use with pore astringent to calm the skin.

⑨ Facial care is recommended 1-2 times a week.



## Charging

- ① Please charge the instrument before using it for the first time.
  - ② Please use standard qualified DC5V/1A adapter for charging.
  - ③ During charging, the power indicator white light flashes breathing.
  - ④ When fully charged, the power indicator white light is always on.
  - ⑤ When the power is low, the red indicator light flashes slowly. Please use it after fully charging.
- ※ If the connection has not started after 10 seconds of charging, please confirm whether the USB cable is light Poor micro-contact.

## Maintenance and protection

- ① When cleaning, please make sure that the instrument is turned off, the instrument can not be washed to avoid electric shock.
- ② When cleaning the tip, please remove the tip from the machine and rinse it with warm water.
- ③ The instrument is not waterproof, you can wipe the body clean with a wet paper towel or a damp cloth after each use, cut  
Do not rinse with water.
- ④ The lens of the product camera is cleaned with a cotton swab. It is forbidden to touch the lens with hard objects.
- ⑤ It is forbidden to use corrosive solution to wipe the body.
- ⑥ Please store the instrument in a normal dry indoor environment.
- ⑦ It is strictly forbidden to disassemble this product by yourself. If any malfunction occurs or cannot be used normally, please stop  
Stop using and contact online customer service.

## Specification

Name: VISUAL BLACKHEAD VACUUM DEVICE

Model: DZ02

Rating: DC5V/ 1A ( Built-in lithium battery 3.7Vdc, 1000mAh)

Size: 158.5\*45\*48mm

## Warranty

The warranty period of this product may vary according to the country of sale

There are different lengths of time, we will allow the law in the country

Try to extend it as much as possible.

● The warranty period of this product is 1 year from the date of purchase. The accessory warranty is from the date of purchase

From 1 year.

● If the instrument has quality problems during the warranty period, it will be repaired or replaced free of charge. But by you

Machine damage caused by improper use or non-standard operation is not covered by the warranty.

The cost price is for you to replace the damaged parts or replace them with new ones. The ownership of the old instrument belongs to the manufacturer.

● Please show your warranty card when you need warranty. The warranty date of this product is as The date is the start date. If you do not have a purchase date, it will be calculated according to the date of manufacture.

In the following cases, warranty or return may not be provided

Do not undertake to reissue the purchase invoice and warranty card for you, please keep your purchase invoice and warranty

card.

- ① For product damage caused by unpredictable events such as fire, earthquake, flood, etc.
- ② Instrument failure and damage caused by your negligence, misuse, or use under other abnormal conditions.
- ③ Instrument damage caused by improper disassembly and improper storage (such as rat damage or liquid penetration).
- ④ Damage caused by compulsory use of products in non-product use environments, including but not limited to

Work load, etc.

- ⑤ Product damage caused by circumstances beyond the control of other manufacturers.
- ⑥ The incidental damage to the function or use of this product (such as damage to your business interests, business

Interruption, other losses caused by using the instrument, the use of this product requires replacement of equipment and devices

Or the cost incurred by the service, your customers' claims, etc.), the manufacturer does not Responsibility.



**FCC Caution:**

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment.