


说明书尺寸:60*100mm 材质:105g双铜 风琴折(五折六页)



Visible Acne Remover

USER MANUAL

Thank you for your purchasing. For accurate and safe use, please read the manual carefully before use the product and keep it well.




The illustrations used in the manual only for reference, please refer to the actual product.

Product Introduction

This product is a personal care cleaning instrument, suitable for blackheads, acne, fat granules, pustules and pimples. Built-in high-definition camera, via APP connects to smart phone or tablet, directly watch the real-time removal process and your skin condition clearly. The cleaning head adopts professional 304 stainless steel, which is rust-proof and easy to disinfect and wash clean, simple design, compact and easy to store.


Package List:

- 1. Visible Acne Remover
- 2. User Manual
- 3. Charging Cable

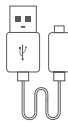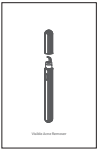


Visible Acne Remover User Manual Charging Cable

Product Structure

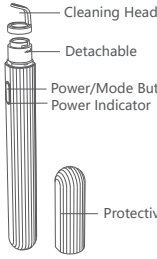

Cleaning Head
Detachable
Power/Mode Button
Power Indicator

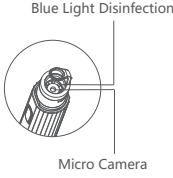

Blue Light Disinfection
Micro Camera
Type-C Charging Port


Functions

- 1. Long Press "Power/Mode" button 1s, turned on the device, the white indicator light flashes and enter into "deep cleaning mode".
- 2. Short press the "Power/Mode" button to switch to the "Blue Light Disinfection Mode", the blue light on the head is on.
- 3. Short press the "Power/Mode" button again, the product will cyclically switched between the two modes.
- 4. In any working mode, long press the "mode/open key" for 1s, the product will turn off, and working indicator will be off.

Operation

1. Two Ways to download APP(Support Android 5.1 and IOS 9.0 version or above)



- ① Download the APP from AppStore or Google Play;
- ② Scan the QR code to download the APP. (If use WeChat to scan, need to manually jump to the browser to open)



Compatible for both IOS and Android system



2. Connect Product WiFi

- ① Turn on the product and APP
- ② Enter the WIFI setting
- ③ Connect phone with product WIFI "UDU Facial-xxxx".



3. Open App


Open App and click "Start up", then watch the live video.



Using Steps

1. Deep Cleaning Mode

- ① After basic facial cleansing, apply a hot towel to compress and steam face to open the pores.
- ② Put the cleaning head against the skin, cleaning in downward direction and always keeping the whole cleaning head closely against the skin, refer to below picture. That performs better.

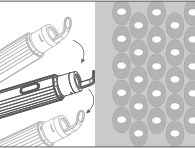

2. Warm Tips:

- ① Use alcohol cotton pads to disinfect the cleaning head before using the product, and clean face and hands to prevent bacteria invading the skin;
- ② After use, remove the cleaning head and clean it in time;
- ③ Not recommended for inflamed skin.

2. Blue Light Disinfection Mode

After finish cleaning, short press the "Power /Mode" button to switch to the blue light mode, to disinfect the cleaning head, and the blue light will automatically go out after disinfection finish.

✳ **Tips: Please cover the protective cap to avoid harm to eyes during disinfecting the cleaning head.**

About charging

- ① Please fully charge the product before first time use.
- ② Please use standard qualified DC5V adapter for charging.
- ③ The red indicator light is always on during charging.
- ④ The red indicator light off when it is fully charged.
- ⑤ The red indicator light flashes slowly when it's in low battery power, and please fully recharge before next time use.

✳ **When charging doesn't start after 10 seconds, please check if the Type-C cable port is fully plugged in or any poor connection issue.**

Device Cleaning and Maintenance


- ① When cleaning, make sure the device is turned off avoiding electrical shock.
- ② Please take off the cleaning head from the instrument and rinse by warm water.
- ③ The instrument body is not waterproof, please wipe it clean with a wet tissue or cloth after each use. Do not rinse the body with water.
- ④ Clean the camera lens with a cotton swab, and avoid hard objects collide with the lens.
- ⑤ It is forbidden to wipe the body with corrosive liquid.
- ⑥ Please store the device in a dry indoor environment.
- ⑦ It is strictly forbidden to disassemble the instrument by yourself. If there occurs any failure or malfunction, please stop using and contact online customer service immediately.

Specifications

Product Name: Visible Acne Remover
Model Number: DZ2103
Input: DC5V ~~~ 1A
Power: about 1W
Battery capacity: 250mAh
Working time: about 70 mins
Product size: 145x14.5x13.5 mm
Executive standard: GB 4706.1-2005 GB 4706.15-2008

Hazardous Substance

Part Name	Hazardous Substance					
	Lead (pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent chromium (Cr(VI))	Polybrominated biphenyls (PBB)	Polybrominated diphenyl ethers (PBDE)
Shell	○	○	○	○	○	○
PCBA	✕	○	○	○	○	○
Battery	○	○	○	○	○	○
LED	○	○	○	○	○	○
Cleaning Head	○	○	○	○	○	○
Type-C Charging Cable	○	○	○	○	○	○
Camera Module	○	○	○	○	○	○
Packaging Material	○	○	○	○	○	○



This table is compiled with the regulations of SJ/T11364.
O: Indicates the content of the hazardous substance is below the limitation specified in GB/T26572.
X: Indicates that the content of the hazardous substance exceeds the limitation specified in GB/T 26572.
Note: Components that contain toxic, harmful or elemental elements have no substitutes due to the limitation of the global development level.

Visible Acne Remover
DZ2103

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction

FCC ID: 2AYJD-DZ2103