

90MM

90MM

WOLONOW Smart Body Fat Scale

User Guide

WOLONOW

Please read through the user guide before using the product and keep it in good condition

Warning: This product must not be used by:

• Patients with a cardiac pacemaker or special population who have other electronic product to maintain their vital signs in their body. When this product is used to measure BMI, weak current flows through the human body, which current may affect the proper operation of the medical apparatus inside the body and lead to a severe consequence.

• Pregnant women. Those pregnant women who need to use this product should seek for their doctor's advice first.

• This product should not be used on an unstable or slippery surface to avoid falling.

Note: The device may not be able to obtain precise body data from the following people due to their special nature. Only the body weight information can be accurate.

• Children in growth stage

• Osteoporosis patients

• Elderly people

• Professional athletes

• Edema patients

• Patients with cold or fever

• Pregnant women

• Patients with dialysis

Please try to avoid measuring your weight if:

• The environment is humid, the scale surface is slippery or the product is not steady

• You have just finish sauna and shower and your body surface is not completely dry

• You have just finish intense exercises and your body is not completely quiet down

• You have just drunk a large volume of water or you have just eaten

• You have drunk a large volume of alcohol

Before measuring your weight:

• Turn on or restart the scale and wait for the " " icon stops flashing on the scale surface.

• Keep the scale static for 5 to 10 seconds for the first time use or after moving the scale to allow calibration of the sensor with the ground.

Maintenance and Care


• Please use clean wet cloth to wipe to clean the scale surface.

• Please do not use any detergent to clean the scale surface. Corrosive detergents may cause damage to the ITO conductive coating and affect the precision in body fat measurement. "ITO coating is used with the body fat measurement chip to measure your body fat"

Product Description

Thank you for choosing WOLONOW Smart Body Fat Scale


Hidden High-Resolution LED Screen

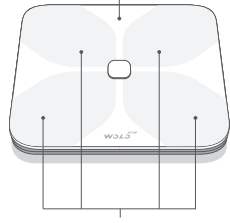

Imperceptible ITO Coating

Battery Cover


Scale Leg with High-Precision Sensor

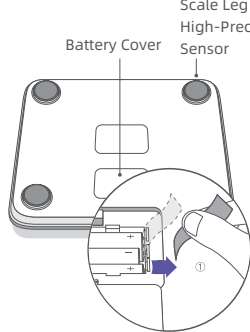


① Note: Please remove the insulating strip so that the product can be switched on.

Installation of WOLONOW App

Search for "WOLO" in Android App Market or App Store, or scan the QR code below to download and install WOLO App for connection with the device.





Scan QR Code to download the official app

The device or the app may not work properly if the app is downloaded through other channels. Before communication with your device, make sure your device is either:

• iPhone / iPod Touch with the system of iOS 9.0 or above; or

• Mobile phone with the system of Android 4.3 or above and with Bluetooth 4.0.


Enable Bluetooth on the phone and open WOLO App

After the Bluetooth on your phone is enabled and WOLO App is open, stand on the scale and measure your weight, and the data can be synched, displayed and recorded in your phone app.

Measuring Weight and Body Fat


Please stand on the scale barefooted in the posture as shown below for measuring your weight and body fat. Stand on the scale with your feet paralleled and aligned by the center of the scale surface.

Stand Upright

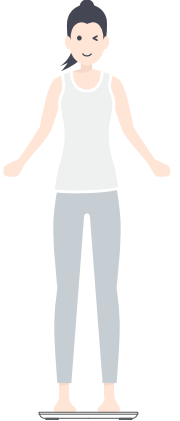


✓

Blocking the Screen

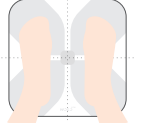


✗

Bend your Knee on the Scale



✗


Squat on the Scale




✗

Information on Display

Body Fat — BF



Weight Unit — kg

Low Battery — 

Weight Number

Note: When the red light flashes, the battery is low.

Product Performance Fault List

Name	Performance Fault
WOLONOW Smart Body Fat Scale	Failure of function listed on the user guide No display / missing display on the screen Unable to switch on the device or unable to access to data on WOLO APP Cracking on the enclosure due to structural defect or material factor

After Sale

The device is supplied with lifetime technical support and 1-year warranty for hardware.

FCC Requirement

FCC ID: 2A25Z-ITO260

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and

(2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.

- Increase the separation between the equipment and receiver.

- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

- Consult the dealer or an experienced radio/TV technician for help.

Basic Parameters

Product Name: WOLONOW Smart Body Fat Scale

Model: BFSITO - 260

Unit of Measurement: lb

Communication Mode: Bluetooth 4.0

Dimensions: 260 x 260 x25 (mm)

Execution Standard: Q/XF 001-2021

Net Weight: 1.13 kg

Weight Range: 5kg - 150kg

Measurement Precision: ±50g

Operating Ambient Temperature: - 5°C to 40°C

Enclosure Material: ABS engineering plastics

Battery Type: 1.5v AAA Battery x 3

WOLONOW

More detailed information is available from www.wolonow.cn. Or, call WOLONOW's customer service hotline 400-999-0509 for more inquiries.

Manufacturer: Star Atmosphere Technology (Shenzhen) Co., LTD

Address: Room 2706A, Block A, Sunshine TDI Building, Nanshan District, Shenzhen 518000, China

User Guide V1.0

风琴折页，128g哑粉字体印刷/单色 k 80 或冷灰11C /字体：阿里巴巴普惠体

切刀 压线 齿刀 开槽 半穿 粘合

Scale 1:1 Unit : mm

材料： 工艺： 尺寸：

印刷颜色： CMYK

设计： XXX

版本： XX