


Product instructions



Smart Sterilization Ultrasonic Teeth Cleaner

Product Description

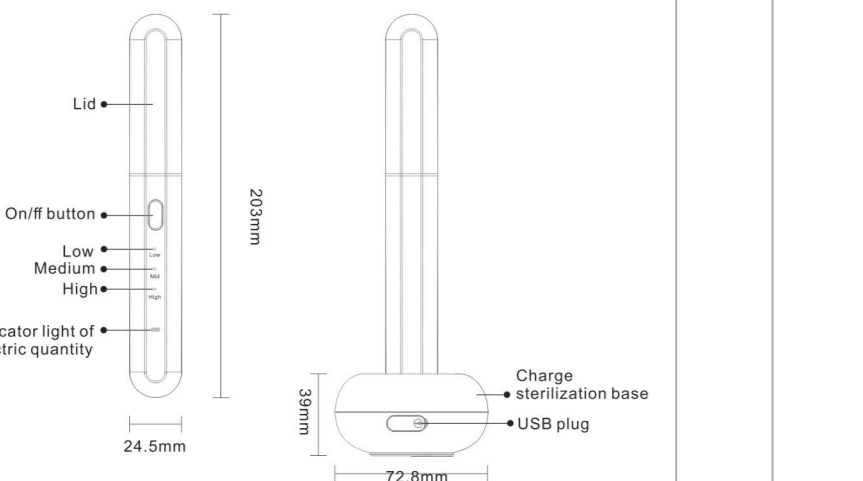

- On/off button**
On/off: Hold for 3s
Intensity adjustment: Low is the default mode, press once to switch to medium and press again to convert to high
- Charging base**
Plug to charge and disinfect
- Indicator light of electric quantity**
In the event that the battery is sufficient, the light is in green
The red light will flicker during the charging and it will turn to green when the battery is fully charged.

Turn on JKK tooth cleaner

- Scan QR code to download the mobile APP and install it
Open "Cvideo" in your browser and download and install it




Turn on JKK tooth cleaner

- Hold the On/Off button for 3s to turn the cleaner on.



- Open app

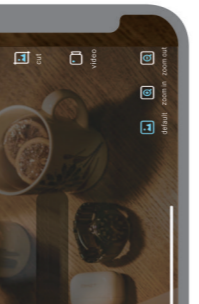


Connecting product WiFi

- Connectivity to WiFi "JKK-M20"



- Return to APP to use it

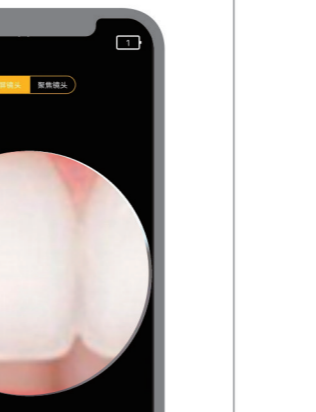

Replace the working head

- Use special tool to screw the working head tight and make sure to mount in place so as to avoid falling off





- The curved part of the working head shall aim at the camera, otherwise, the camera will not discover the tooth



Charging/Sterilization base

Magnet base for charging and plugging
Forward plug is for charging and sterilization and reverse plug is for fixing

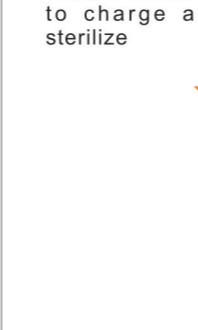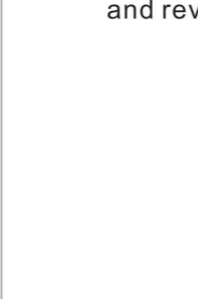

Storage bin for sterilization base

- Rotate the storage bin to replace the working head of the tooth cleaner



- External TYPE-C port can sterilize the working head while electrifying and the sterilization time is 5min and it will automatically stop the sterilization afterwards.

- External TYPE-C port can sterilize while electrifying



Open the USB plug → Inset USB cable to start the charging


Troubleshooting

Fault	Resolvent
No response to power on	several times according to the correct starting method Please make sure that there is electricity and please use the matched charging base
The camera doesn't light up	Re-start the tooth cleaner Please make sure that there is electricity
App crash	Unload the current APP, re-scan the QR code to install again
The indicator light flickers in red	Low battery, please charge the battery.
Product can correctly connect to APP but fail to display the real-time image	WiFi Disconnect the mobile phone from the WiFi of tooth cleaner and re-start it for connection again
The WiFi connection is instable	Re-start the tooth cleaner and re-connect APP

If the problems fail to be solved in accordance with the above operations, please contact the customer service and we will provide you with a feasible solution accordingly

Product parameters

Product name	Intelligent visual ultrasonic tooth cleaner
Weight/packing box:	145g/400g
Network standard	IEEE 802.11b/g/n
Antenna	Built-in FPC antenna
Working frequency	2.4GHz
Image transmission rate	30fps
Image sensor	CMOS
Working temperature	-10~50°C
Battery	520MAH polymer battery
Battery life	180 minu
Charging time	About 3h
Input parameter	5v-0.5a
Working voltage	4.2V
Rated power	5W
Lens diameter	4mm
Optimal focal length	2.5-4mm
Pixels	5 million

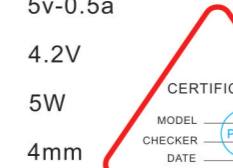

Notice:

- Please make sure that the lens are clean, otherwise, please use swab to make it clean.
- Please observe the ambient environment and do not use it while someone is running around so as to avoid crashing
- The product is not suitable for minors under the age of 16.
- It is prohibited for children to use this product so as to avoid injury.
- Do not wash the main machine with water or use it in liquid so as to avoid damage caused by soaking.
- The built-in battery of the product must be charged once a month if it is not used frequently.
- Protect the product from direct sun exposure, in particular, the main machine, so as to avoid ageing.
- During operating, the temperature of the product will raise slightly (it will not exceed 35°C maximumly), please feel free to use.
- Please try several times before use so as to ensure the safety and best use experience.

Warranty

User information:	
User name:	Contact number:
Consignee address:	
Seller:	Contact number:

Answers to frequently asked questions about tooth cleaner

- Will it hurt when clean the tooth?
No, it will not. The ultrasonic will not work unless it contacts hard objects, such as tooth or glass. When it contacts the soft tissue, such as, gum, it will not work, therefore, you will not feel any pain.
- Does the operation complex? Can the ordinary people use it?
The household tooth-cleaning (non-medical use) has been simplified in design. 99% ordinary consumers can use it easily. For any questions, please consult the online customer service about how to operate it.
- Who can't use it?
It is not recommended for pregnant woman, children and people with a cardiac peacemaker to use it; Do not use it if you have sensitive teeth, artificial teeth, porcelain teeth, periodontitis, dental caries, pathological teeth, thin natural tooth enamel (about 0.1% seen black spots on the teeth after use).
- Can it be used for brushing the tooth? No.
It has the completely different principle to the electric toothbrush. The toothbrush is mechanical motor vibration, whereas, this product is ultrasonic energy exchange, which cannot be used as a toothbrush.
- Does the product have warranty?
Yes, it does. You can enjoy 1-year free maintenance and 7-days unconditional return for quality issues that are not caused by people. You have nothing to worry about it.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.