


AquaSonic Clinical AQ

Ultra-Whitening Toothbrush
with Touchscreen and AI-Powered Brushing Coach

User Manual

Model: AQLCD

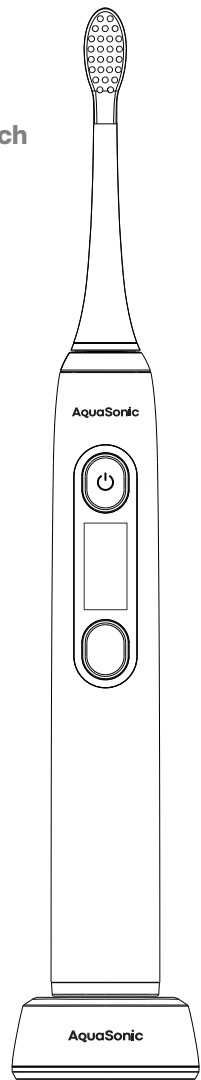

BLUETOOTH CONNECTION APP

Before using this product, please read these instructions completely and save this manual for future use

IMPORTANT SAFEGUARDS

This power toothbrush is safe and effective when used according to directions. Like any electrical device, however, basic safety precautions should be followed, especially when children are present.

DANGER-WATER

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not place or store the charger base where it can fall or be pulled into a bathtub, sink, or toilet. Do not submerge it into water or any other liquid.
- Do not reach for this product if it has fallen into water or any other liquid. Unplug it immediately.
- Do not use this product while bathing.
- Unplug the charger base before cleaning it, and make sure it is dry before plugging it into an electrical outlet.

DANGER-VOLTAGE

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug the charger base into an outlet with a voltage other than that which is specified on the bottom of the charger.
- Voltage converters do not guarantee compatibility.
- This device has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in an outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not bypass this safety feature.
- If the plug does not fit easily into an outlet, discontinue its use.

MEDICAL WARNINGS


- If you have had oral surgery or gum surgery within the past two months, consult your dentist before using this product.
- Consult your dentist if excessive bleeding occurs while using this product or if minor bleeding continues to occur after one week of use.
- Consult your physician or implant device manufacturer prior to using this product if you have any questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- Never use while sleeping or drowsy.

SAFETY WARNINGS

- This toothbrush should not be used by children, or by persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning the use of this product by an adult responsible for their safety. Cleaning and user maintenance should not be done by children without supervision.
- Children should be supervised to ensure they do not play with this toothbrush.
- This toothbrush should only be used for its intended purpose as described in this user manual. It is designed to clean teeth, gums, and tongue only. Do not use any accessories or brush heads other than those recommended by the manufacturer.
- This toothbrush is a personal care device and is not intended to be shared or used by more than one person in a dental office or institution.
- If your toothpaste includes peroxide, baking soda, or other bicarbonates which are common in whitening toothpastes, please ensure that you rinse your brush head and handle thoroughly between uses, to avoid cracking the plastic.
- Never operate this product if it has a damaged cord, plug, or charger base.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Do not use this toothbrush if it is not functioning properly. If it has been dropped, damaged, or dropped into water, discontinue use immediately. Discontinue use immediately if the product appears damaged in any way.
- Keep the toothbrush, charger base and power cord away from heated surfaces.
- If you experience any oral discomfort or pain as a result of using this toothbrush, discontinue use and contact your dentist.
- Do not clean this toothbrush in the dishwasher.
- The battery must be removed from this toothbrush before it is discarded. Ensure the battery power is fully drained before removing the battery. Dispose of the battery safely, according to local waste management regulations.
- Do not use outdoors or operate where aerosol spray products are being used, or where oxygen is being administered.

Questions? Contact support@aquasonic.com for assistance.

CONFIGURATION

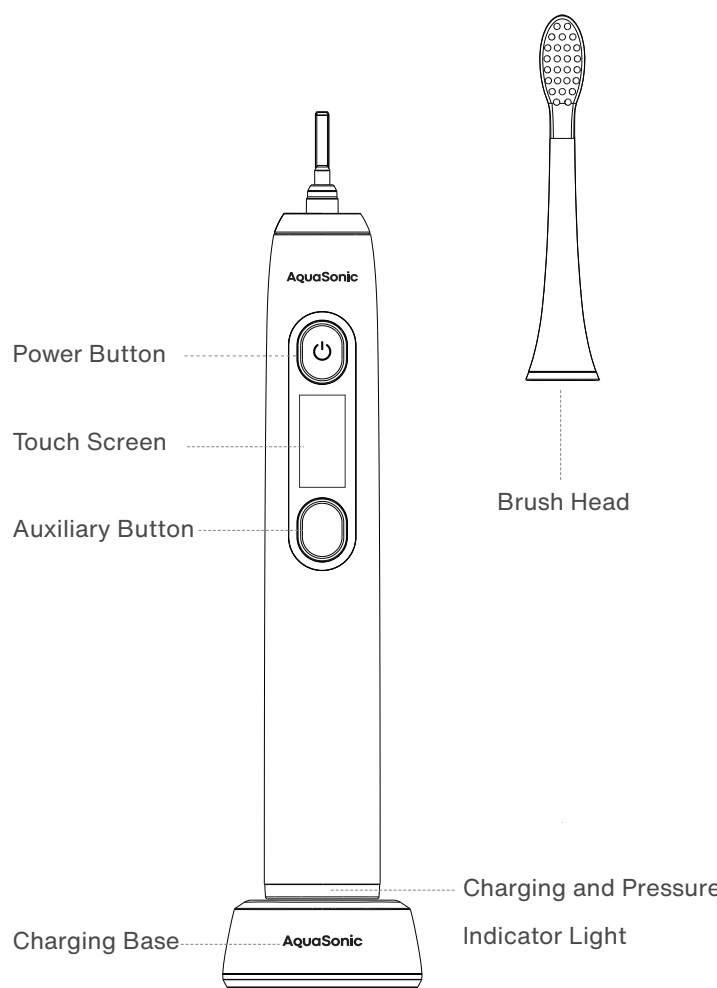


GETTING STARTED

Charging the Toothbrush

- Be sure the toothbrush is turned off before and during charging.
- Plug the power plug into an outlet, then place the toothbrush handle on the charger base.
- The charging indicator light will pulse while charging is in progress.

Please note:

- It takes at least 18 hours to fully charge the battery.
- The toothbrush must be charged for 18 hours:
 - prior to the first time using the toothbrush
 - if the toothbrush has not been used for 3+ months
- Charge your toothbrush overnight every few weeks, to keep the battery charged. Do not leave it charging at all times.
- While charging, the toothbrush may get warm, but this is normal.
- As the battery ages, its operating duration will shorten. If the battery only operates once or twice after being fully charged, it is nearing the end of its life cycle. Contact support@aquasonic.com for assistance.

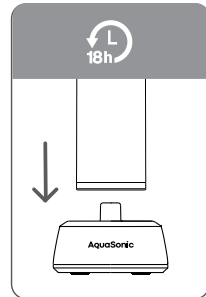

HOW TO USE THE APP

- Search "AquaSonic Clinical AQ" in the App store or scan the QR code to download and install the app on your phone.
- In your phone settings, make sure Bluetooth is on. Open the app and follow the prompts. Turn on your toothbrush to make it discoverable.
- When connecting your toothbrush and app for the first time, make sure they are within 4 inches of each other. When indicated, hold down the button for 15 seconds to complete the connection.
- When you exit the app, the toothbrush will automatically disconnect. When you open the app again, the app and toothbrush will automatically reconnect.



iOS



Android

HOW TO USE THE TOOTHBRUSH

- Slide the brush head onto the handle and push down firmly. There should be a small gap between the base of the brush head and the toothbrush handle (Fig 1).
- Wet the bristles and apply toothpaste.
- The first time you turn on the toothbrush, you'll be prompted to set the time. Swipe up and down until the time is correct, and push the auxiliary button to move to the next steps.
- To change modes, swipe up or down on the screen until the desired mode is displayed. There are 5 modes:
 - Clean** – standard brushing mode
 - Soft** – gentle mode for sensitive teeth
 - Whiten** – helps reduce the look of surface stains
 - Massage** – vibrating microbursts for gentle gum stimulation
 - Deep Clean** – targets stubborn discoloration and polishes teeth
- Press the power button to begin brushing. Brush evenly throughout the mouth, in 4 sections (Fig 2): 1-outside upper, 2-outside lower, 3-inside lower, 4-inside upper.
- Holding the bristles at a 45-degree angle to the gumline (Fig 3), guide the brush head slowly from tooth to tooth, pausing for a few seconds on each tooth. Brush the gums and the chewing surfaces, as well. Do not press too hard or scrub your teeth – simply let the toothbrush do all the work. This toothbrush is equipped with a pressure sensor, which will give an indication if you apply too much pressure.
- During brushing, the vibration will momentarily pause every 30 seconds, indicating to move to the next quadrant of the mouth. After 2 minutes, the toothbrush will shut off automatically. You may choose to set the brushing time to 2.5 or 3 minutes, if you habitually prefer a longer duration.
- After brushing, rinse the brush head and the toothbrush handle under running water (Fig 4), then wipe them dry.

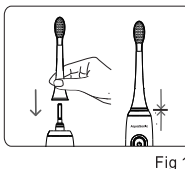


Fig 1

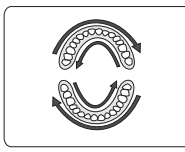


Fig 2

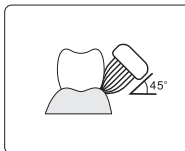


Fig 3

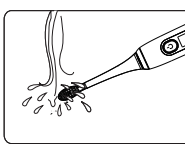


Fig 4



Fig 5

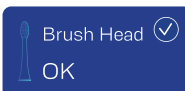


Fig 6


This toothbrush tracks the wear of your brush head. The app will display an indication when the brush head should be changed soon, and when a new brush head is overdue. Replace the brush head, then press and hold the auxiliary button for 8 seconds to reset the brush head status. (Fig 5, Fig 6).

CLEANING & MAINTENANCE

Be sure the toothbrush is powered off and the charger is unplugged before cleaning.


Toothbrush handle and brush head

- Remove the brush head.
- Wipe all surfaces with a soft, damp cloth. Do not submerge the handle in water.
- Rinse the brush head and bristles with water after each use.

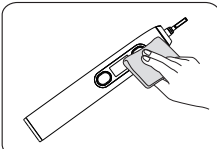

Charger

- Using a soft cloth, wipe down the charger base and power cord.
- Wipe the power plug periodically to ensure it remains in good working order. The accumulation of dust or dirt can cause an electrical shock, overheating, or fire.

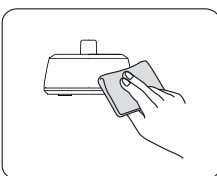

TROUBLESHOOTING




Problem	Possible cause	Solution
The toothbrush does not turn on	You have just purchased this product or it has not been used for more than 3 months	Charge the toothbrush handle for 18 hours
Weak or diminished vibration	The toothbrush is pressed too hard against your teeth	Place the bristles gently against your teeth
The toothbrush operates for only a short period of time, even immediately after being charged	The battery has reached the end of its life cycle	Contact AquaSonic for support
	Charging time was insufficient	Charge the toothbrush handle for 18 hours

DISPOSAL


AquaSonic Clinical AQ is powered by a lithium-ion battery, which can be recycled. At the end of this product's life, always remove the battery before discarding the handle, and recycle or dispose of the battery according to local waste management regulations. If you have any questions, please contact your local waste management service provider.

Removing and recycling the battery

- Before removing the battery, be sure the power is completely drained.
- Use a screwdriver to remove the screw cover and battery door.
- Holding the handle upside down, press down on the toothbrush shaft and remove the two hooks to release the internal components of the handle.
- Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
- Recycle or dispose of the battery properly.

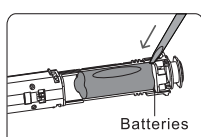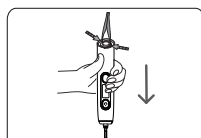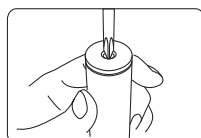

Batteries

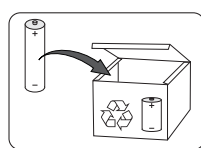

Be careful not to short-circuit the negative and positive poles of the battery. Do not reconnect the appliance to the mains after you have removed the battery.

DISCLAIMER

Users of this device do so at their own risk. Neither AquaSonic nor its retailers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting, directly or indirectly, from the use of this device. Further, AquaSonic reserves the right to revise this publication and to make changes from time to time in the contents thereof without obligation to notify any person of such revision or changes. Model may be changed for improvements without notice.

FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- This device may not cause harmful interference, and
- this device must accept any interference received, including interference that may cause undesired operation.

Please note

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

WARNING: Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment. The device must not be co-located or operating in conjunction with any other antenna or transmitter.



ISO9001 & ISO13485 Certified Manufacturer