
SEAGO 赛嘉



WiFi口腔镜
WiFi Oral Camera
SG-2002

/
Instruction Manual
使用说明书

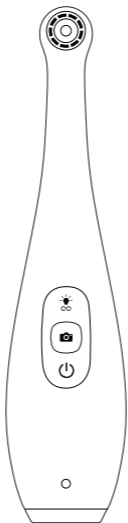
Thanks for choosing SEAGO.
Before use, please read the instructions carefully
and keep it at a safe place for future reference.

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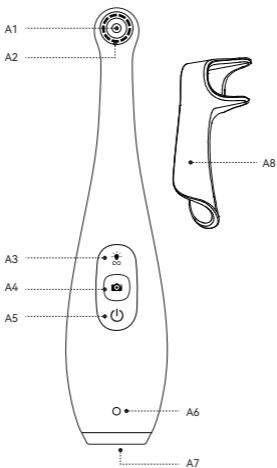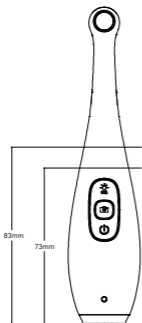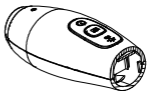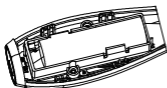
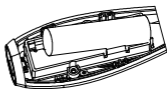
B-4

Dear customer,

Thank you very much for choosing SEAGO. We are committed to providing you with high-quality and cost-effective oral care products, assisting you in daily cleaning and tooth care, and protecting your oral health.

Safety Tips

- Do not disassemble, repair, or modify the product on your own.
- Do not hit and drop the product, and do not use it if damaged.
- Please check whether the local power supply voltage matches the product's specified voltage before connecting the power supply.
- Do not place the product in areas where it may easily fall.
- Please keep it out of the reach of infants and children. Do not use as a toy.
- This product is for indoor use only, and do not place in high-temperature or direct sunlight areas.

Product Information

Product Information: WIFI Oral Camera

Model: SG-2002

Input Voltage: DC5V

Input Power of Lamp: 1.5W

Product Materials: ABS / PC / PMMA

Charging Time: ≤4 hours

Product Size: 30×34.1×150.5mm

Use Conditions: 5-40°C

Lens Pixels: 300,000 pixels

Product Standard: GB4706.1-2005

Place of Origin: China

Component Names

A1	Camera	A5	Power Button
A2	LED Light	A6	Indicator Light
A3	LED Light Adjustment Button	A7	Charging Port
A4	Photo Button	A8	Bracket

Function Introduction

1. Screenshot and Photo Function: after connecting the mobile APP, align the head camera, press the photo button (camera icon), and take photos. The captured photos will be automatically stored in the APP album and the phone's gallery.
2. Video Recording Function: click on the video recording icon in the APP to start recording videos.
3. Camera Light Adjustment Function: press the LED light adjustment button on the product (light icon) to adjust the LED light brightness in three levels, suitable for use in different lighting environments.
4. WIFI Connection to Phone Function:
Scan the QR code below with your phone (Android system), or search "iVision" in the application market with your phone (IOS system), download and install the APP to connect to your mobile phone.

Android system

How to Use

1. After inserting the disposable film cover from the camera, use the transparent bracket for positioning.
2. Long press the power button for 3 seconds to turn on the product. The bottom blue indicator light starts flashing continuously, indicating that the WIFI function has been opened, and the phone can be connected.
3. Open the phone's WIFI settings, search for WIFI, and connect to the WIFI named "iVision-XXXX" (if unable to connect, please try ignoring the network or close the APP and reconnect). Once connected successfully, the product's blue indicator light will stay on, and the phone's APP interface will show the connection status.
4. After a successful connection, click on the connected product to view the image on the interface. If the phone prompts during use "Current WLAN is not available for internet access, do you want to switch to another network?", please select

"Not switch".

5. Before capturing images with the product, short press the LED light adjustment button to adjust the LED light brightness based on the surrounding environment. The gears are switched in the order of first gear → second gear → third gear → off → first gear.

6. While connected to the APP, short press the photo button on the product or click the photo icon on the APP interface to capture the current image, which will be automatically saved in the phone's album and the APP album.

Charging Instructions

1. TYPE-C charging: insert the TYPE-C charging cable into the bottom of the product. The red indicator light will illuminate, indicating that the product is charging. The red light goes out when the product is fully charged.

2. Overcharge protection: after the battery is fully charged, the charging will automatically shut off.

3. Over-discharge protection: when the battery level is critically low, the product will automatically power off.

4. Charge indication

Instructor	Battery Status
Charging indicator light not illuminated.	Less than 10%
Red charging indicator light illuminated.	10%~90%
Red charging indicator light off.	90%~100%

Recycling Instructions for Rechargeable Battery

This product uses a lithium battery. Before disposing of the product at the end of its life, be sure to remove the battery and recycle or dispose of it according to local waste management regulations. If in doubt, consult your local waste management authorities. To dismantle the battery, you will need the following tools: saw, wire cutters. Note: Only dismantle the battery when disposing of the product. Before dismantling, ensure that the battery is fully depleted. When using tools, follow basic safety rules and take appropriate precautions to avoid injury.

1) The picture shows cutting between distances 73-83mm (B-1).

2) Open the product by hand or with a tool (B-2).

3) Cut the two wires connected to the battery (B-3).

4) Remove the battery (B-4).

FCC Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

SAR Information Statement

Your Smart Visual Ear Cleaner is a radio transmitter and receiver. It is designed and manufactured not to exceed the emission limits for exposure to radio frequency (RF) energy set by the Federal Communications Commission of the U.S.

Government. These limits are part of comprehensive guidelines and establish permitted levels of RF energy for the general population. The guidelines are based on standards that were developed by independent scientific organizations through

periodic and thorough evaluation of scientific studies. The standards include a substantial safety margin designed to assure the safety of all persons, regardless of age and health. The exposure standard for Smart Visual Ear Cleaner employs a unit of measurement known as the Specific Absorption Rate, or SAR. The SAR limit set by the FCC is 4.0W/kg. * Tests for SAR are conducted with the phone transmitting at its highest certified power level in all tested frequency bands. Although the SAR is determined at the highest certified power level, the actual SAR level of the phone while operating can be well below the maximum value. This is because the phone is designed to operate at multiple power levels so as to use only the power required to reach the network. In general, the closer you are to a wireless base station antenna, the lower the power output. Before a phone model is available for sale to the public, it must be tested and certified to the FCC that it does not exceed the limit established by the government adopted requirement for safe exposure. The tests are performed in positions and locations (e.g., at the ear and worn on the body) as required by the FCC for each model. The highest SAR value for this model phone when tested for use at the ear is 0.448W/Kg and when worn on the body, as described in this user guide, is 0.165W/Kg (Body-worn measurements differ among phone models, depending upon available accessories and FCC requirements). While there may be differences between the SAR levels of various phones and at various positions, they all meet the government requirement for safe exposure. The FCC has granted an Equipment Authorization for this model phone with all reported SAR levels evaluated as in compliance with the FCC RFexposure guidelines. SAR information on this model phone is on file with the FCC and can be found under the Display Grant section of <http://www.fcc.gov/oet/fccid> after searching on FCC ID:2BC87SG-2002. Additional information on Specific Absorption Rates (SAR) can be found on the Cellular Telecommunications Industry Association (CTIA) web-site at <http://www.wow-com.com>. * In the United States and Canada, the SAR limit for mobile phones used by the public is 1.6 watts/kg (W/kg) averaged over one gram of tissue. The standard incorporates a substantial margin of safety to give additional protection for the public and to account for any variations in measurements.

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Dear customer, thank you very much for your care and support for SEAGO's oral care products! We appreciate your choice and purchase of our products! Our company promises that this product has passed our strict quality inspection, if there is any quality problem from the date of purchase, it will be returned within 7 days, exchanged within 15 days and warranted within 2 years. If you do not operate according to the manual or without the company's authorized units and individuals privately dismantled or repaired, does not belong to this scope.

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Fill out information below by customer

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Fill out information below by retailer

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