



BRUSHMO®
Brushmo BM01
Instruction Manual

Introduction


Congratulations on your purchase and welcome to Brushmo! To fully benefit from the support that Brushmo offers, register your product at www.brushmo.com

Intended use

Brushmo electric toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Brushmo electric toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

Contents

- 1 x Brushmo BM01 Handle
- 4 x Brush Heads
- 1 x Charger with USB Type-C
- 1 x Power Adapter

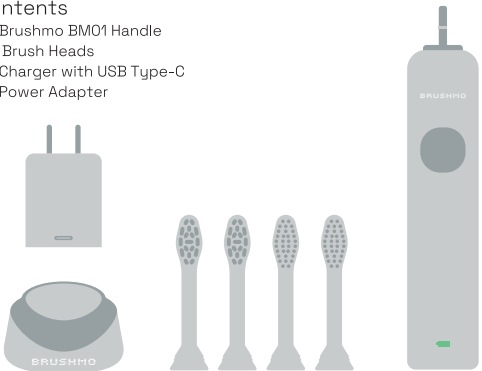


Getting Started

Attaching the brush head

- Align the brush head so the bristles point in the same direction as the front of the handle.
- Firmly press the brush head down onto the metal shaft until it stops.

Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

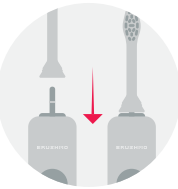

Charging your Brushmo toothbrush

- Plug the USB Type-C charging cord into a wall adapter that is connected to electrical outlet.
- Place the handle on the charger.

The light of the battery level indicator illuminates. This indicates that the toothbrush is charging.

Note: It can take up to 13 hours to fully charge the battery, but you can use the toothbrush before it is fully charged.

Using your Brushmo toothbrush

If you are using your toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience.

Brushing Instructions

- Wet the bristles and apply a small amount of toothpaste.
- Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.


Note: Keep the center of the brush in contact with the teeth at all times.

- Press the power on/off button to turn on the Brushmo Toothbrush.
- Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the 30 Seconds feature.

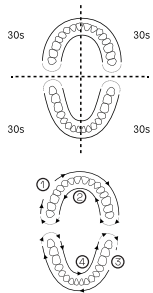

- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.





Your Brushmo toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Note: When the Brushmo toothbrush is used in clinical studies, the handle should be fully charged.

Battery Status

Battery status	Battery indicator
Full	Solid green (turn off after 30 seconds)
Charging	Flashing

- When handle is not placed on working charger.

Battery status	Battery indicator
75 - 100%	Solid green 
35 - 75%	2 green bars 
0 - 35%	1 green bar 
Empty	Flashing light 


Note: To save energy, the battery light will turn off when not in use.

BRUSHMO®

Since 2009

宽98.5MM
高91mm

Modes



5 Modes:

Clean:	For exceptional everyday clean
White:	To remove surface stains
Sensitive:	For sensitive teeth and gums
Gum Care:	Special attention to molars
Deep Clean:	For reaching deep between teeth

- Press the power button to turn on the toothbrush.
- Press a second time to switch modes.

Cleaning

You should clean your Brushmo toothbrush regularly to remove toothpaste and other residue. Failure to clean your toothbrush may result in an unhygienic product and damage to it may occur.

Warning: Do not clean product or accessories with sharp objects or in dishwasher, microwave, with chemicals or in boiling hot water.

Toothbrush handle

- Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.


Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

Brush head

- Rinse the brush head and bristles after each use.
- Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.



Charger

- Unplug the charger before you clean it.
- Wipe the surface of the charger with a damp cloth.

Storage

If you are not going to use the product for a long time, unplug it, clean it (see chapter "cleaning") and store it in a cool and dry place away from direct sunlight.

Replacement

Brush head

Replace Brushmo toothbrush heads every 3 months to achieve optimal results. Use only Brushmo replacement brush heads.

Disposal

- This product contains a rechargeable lithium-ion battery which must be disposed of properly.
- Contact your local town or city officials for battery disposal information. You can also call 1-800-822-8637 or visit www.cal2recycle.org for battery drop-off locations.



IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGERS

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not reach for a product that has fallen into water. Unplug immediately. After cleaning, make sure the toothbrush handle is completely dry before you connect it to the mains.

WARNINGS

If you've recently had surgery on your teeth or gums, consult your dentist before using the Brushmo toothbrush. If you experience excessive bleeding after using this device, please contact your dentist immediately. This toothbrush can be used by children and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not play with the appliance. The toothbrush complies with safety standards for electromagnetic devices. If you wear an implant or pacemaker and have questions about the device, please contact your doctor or the manufacturer of your implant or pacemaker.

SAFETY NOTE

Warning: For the purposes of re-charging the battery, only use the detachable supply unit provided with this appliance. This appliance contains batteries that are non-replaceable. Before using, please check that the voltage stated on the device complies with your local mains voltage. If the charger or connecting cable is damaged, please do not use or attempt to repair it yourself. The power cord cannot be replaced. If the cord is damaged, the entire device must be replaced. All repairs must be handled by a dealer or our customer service. Keep the charger away from water, and do not immerse in water or other liquids. Cleaning the charger: Disconnect from the mains before cleaning and reconnect only once the charger is completely dry. Never charge the device unattended.

ENVIRONMENTALLY FRIENDLY DISPOSAL

Do not dispose of this product with normal household waste at the end of its service life. Take it to a recycling center or official collection point. This helps protect the environment. The integrated rechargeable battery contains substances that may endanger the environment. So please take the used handset to an official recycling center or collection point.


BATTERY SAFETY INSTRUCTIONS

- Only use this toothbrush for its intended purpose and follow the general and battery safety instructions as described in this user manual. Any misuse can cause electric shock, burns, fire and other hazards or injuries.
- Charge, use and store the product at a temperature between 5 °C and 35 °C.
- To charge the battery, only use the USB cable provided with the product. Do not use an extension cord.

- Do not burn the toothbrush handle and its battery and do not expose them to direct sunlight or to high temperatures (e.g. in hot cars or near hot stoves). Batteries may explode if overheated.
- If the product becomes abnormally hot, gives off an abnormal smell, changes color or if charging takes much longer than usual, stop using and charging the product and contact Brushmo customer service.
- Do not place the toothbrush handle and its battery inside microwave ovens or on induction cookers.
- This product contains a rechargeable battery that is non-replaceable. Do not open the toothbrush handle to replace the rechargeable battery.
- When you handle the battery to remove them during final product disposal, make sure that your hands, the product and the battery are dry.
- To prevent the battery from heating up or releasing toxic or hazardous substances, do not modify, pierce or damage the toothbrush handle and battery, and do not disassemble, short circuit, overcharge or reverse charge the battery.
- To avoid accidental short-circuiting of the battery after removal, do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings). Do not wrap the battery in aluminum foil. Tape battery terminals or put the battery in a plastic bag before you discard it.
- If the battery is damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

BATTERY SAFETY INSTRUCTIONS

Within the EU, this symbol indicates that this product must not be disposed of with household waste. Old appliances contain valuable recyclable order materials. They should be recycled to prevent harm to the environment or human health through improper waste disposal. Therefore, please dispose of old devices through a proper collection system or send the device back to where you purchased it for proper disposal. They will then ensure that the device is recycled accordingly.



NOTE: This device complies with part 18 of the FCC Rules.

Information to the user.

1. The device has potential interference, but the interference is very small and meets the requirements of the FCC rule.
2. Equipment system maintenance is simple, please refer to the manual.
3. The user can take simple measures to correct the interference. Such as staying away from interference sources, turning off interference sources, etc.