
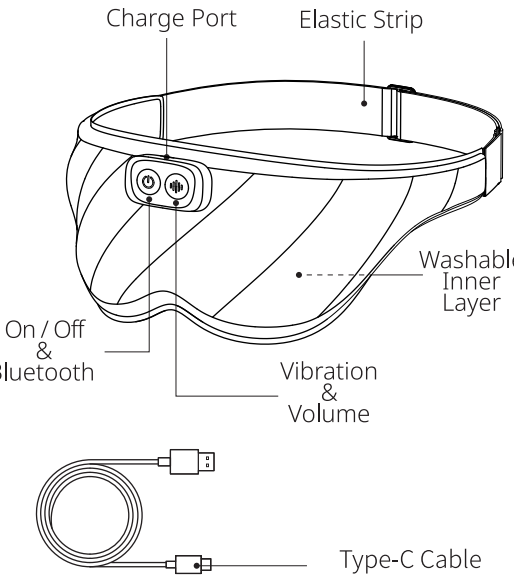

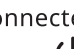
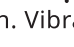














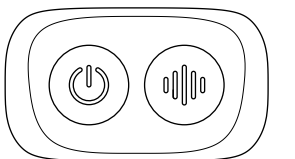



成品尺寸：80X110mm

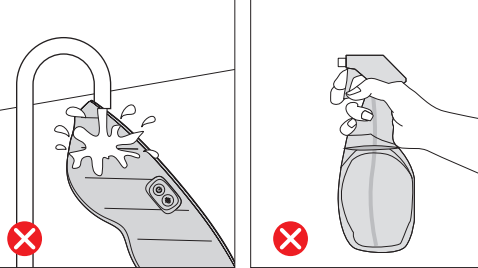 <p>User Manual</p> <p>Mynt Sleep Mask M13202A04</p>	<p>IMPORTANT SAFETY INSTRUCTIONS</p> <p>READ ALL INSTRUCTIONS BEFORE USING. Please follow the instructions and keep the user manual for future reference.</p> <p>USE THIS DEVICE ONLY AS DESCRIBED IN THIS MANUAL. Any misuse excludes any liability for damage, incorrect use may cause electric shocks, fires, burns, or other unexpected hazards.</p> <p>THE DEVICE CONFORMS TO RECOGNIZED TECHNICAL PRINCIPLES AND THE LATEST SAFETY REGULATIONS. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:</p>	<p>WARNING - TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, PLEASE</p> <ul style="list-style-type: none"> The massage device conforms to recognized technical principles and the latest safety regulations. DO NOT use near, on, by children and disabled persons without close supervision. DO NOT use it directly on open wounds, swollen skin, inflamed areas or skin eruptions. DO NOT use it on insensitive skin areas or in the presence of poor circulation. DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered. DO NOT setup or use it in, or around water/ other liquids, or places with high humidity. 	<ul style="list-style-type: none"> In case of malfunction, take off the device immediately. DO NOT use any other attachments not recommended by the seller. This appliance is designed for personal, non-professional use only. To reduce the risk of electric shock and warranty expiration, DO NOT attempt to disassemble, repair, or modify this product. Any possible repairs may only be carried out by authorized expert staff. For service, contact support@hellomynt.com. <p>CAUTION - READ ALL INSTRUCTIONS CAREFULLY BEFORE OPERATING</p> <ul style="list-style-type: none"> This product is not intended for medical use. It is intended only to provide a luxurious relaxing experience. DO NOT use this product as a substitute for medical attention. Consult your doctor prior to using this product, if <ul style="list-style-type: none"> You are pregnant You have a pacemaker You have any concerns regarding your health 	<p>At A Glance</p>  <p>Accessories: Type-C Cable×1, User Manual×1</p>	<ol style="list-style-type: none"> Long press  for 1.5 seconds to turn on/off the device and active bluetooth. Search for "Mynt_BT_Eye_Mask" and connect to the sleep mask through your smart phone. When connected to the sleep mask, short press  to play/pause the music. Double press to play next song. Long press  to turn on/off the vibration. Vibration doesn't require bluetooth connection. When connected to the sleep mask, short press  to adjust the volume. There're Mute-Low-Medium-High four settings. The default volume is at medium. <p>Notice: The vibration will automatically turn-off after 15 minutes of continuous use. All function will turn-off after 2 hours of continuous use.</p>	<p>3. Charge instructions</p> <ul style="list-style-type: none"> Fully charge the device before first use. To charge, connect the Type-C cable to the Power Adapter (not included). <p>Notice: Device cannot be switched on when charging.</p> <p>Care and Maintenance</p> <p>To Store</p> <ul style="list-style-type: none"> Place the product in its storage bag or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, do not wrap the power cord around the unit. Do not hang the unit by the charging cable. 	<p>Do not directly wash the functional layer with, or immerse in water.</p> <p>Never use detergents to clean this device as this may cause</p> <p>Specifications</p> <table border="1"> <tr><td>Product Model</td><td>M13202A04</td></tr> <tr><td>Charging Voltage</td><td>5V / 350mA</td></tr> <tr><td>Rated Power</td><td>0.3W</td></tr> <tr><td>Battery Type</td><td>Li-ion 300mAh 3.7V</td></tr> <tr><td>Battery Life</td><td>Approximately 1.5 Hours</td></tr> <tr><td>Battery Charging Time</td><td>Approximately 3 Hours</td></tr> </table>	Product Model	M13202A04	Charging Voltage	5V / 350mA	Rated Power	0.3W	Battery Type	Li-ion 300mAh 3.7V	Battery Life	Approximately 1.5 Hours	Battery Charging Time	Approximately 3 Hours	<p>FCC Statement</p> <p>This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.</p> <p>Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.</p> <p>Note: This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: (1) Reorient or relocate the receiving antenna. (2) Increase the separation between the equipment and receiver. (3) Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. (4) Consult the dealer or an experienced radio / TV technician for help.</p>	<p>Long Warranty Time for Quality-Related Issues We take care of all quality-related issues with a REPLACEMENT OR FULL REFUND – for 12 months, however the warranty can be extended(see below).</p> <p>FAQs</p> <ol style="list-style-type: none"> What isn't covered by the warranty? <ul style="list-style-type: none"> Purchases from unauthorized resellers Improperly operated devices Lost or stolen products Purchases made over limited warranty time (unless otherwise stated) Non quality-related issues (after 30 days of purchase) Free products When does the warranty begin? Warranty begins from the day you place your order. How do I claim the warranty? Refer to the specific FAQs, send an email to support@hellomynt.com with your receipt. What if I don't have any proof of purchase? If you made your purchase through an authorized reseller, contact the reseller for a copy of receipt. If the product was a gift, ask the giver for a copy of the receipt or claim the warranty on your behalf. Will the warranty be renewed if my product is replaced? No, it won't be renewed after a replacement has been provided. How do I extend the warranty to 24 months? Please scan the QR code below, then go to the Warranty tab for next steps. 	<p>FCC RF Radiation Exposure Statement:</p> <p>1.This Transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.</p> <p>2.This equipment complies with RF radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 5mm between the radiator and your body.</p> <p>FCC ID: 2BL7X-M13202A04</p>																
Product Model	M13202A04																																					
Charging Voltage	5V / 350mA																																					
Rated Power	0.3W																																					
Battery Type	Li-ion 300mAh 3.7V																																					
Battery Life	Approximately 1.5 Hours																																					
Battery Charging Time	Approximately 3 Hours																																					
 <p>For FAQs and more information, please visit: HelloMynt.com</p> <p>Model: M13202A04</p>	<p>DANGER - TO REDUCE THE RISK OF ELECTRIC SHOCK, PLEASE</p> <ul style="list-style-type: none"> Always unplug the product from the charger when not in use and before cleaning. DO NOT use a product that has fallen into water. NEVER use pins or other metallic fasteners with this product. Carefully examine the covering before each use. Discard the product if covering or power control cords show any sign of deterioration, such as blistering or cracking. 	<ul style="list-style-type: none"> DO NOT leave this device operating while unattended, especially if children are present. DO NOT insert objects into any openings. DO NOT place the device and charge cable near a heat source or open flame. DO NOT use pins or other metallic fasteners with this device. DO NOT ever touch the power plug, plug or unplug the charge cable with wet hands to avoid electrical shock. DO NOT use any damaged cables, plugs or loose sockets. If it is not working properly, please contact the Mynt support team support@hellomynt.com. 	<p>What's In The Box</p> <ul style="list-style-type: none"> Mynt Sleep Mask Type-C Cable User Manual Storage Bag 	<p>Before Using</p> <ul style="list-style-type: none"> Remove any eye accessories (e.g. eyeglasses, false eyelashes, etc.) before wearing. Position the device over eyes and fasten the elastic strip. <p>Operation</p>  <ol style="list-style-type: none"> Operation instructions  On/Off & bluetooth control button  Vibration & volume button 	<p>2. Indicators instructions</p> <p>There're several indicator lights on the two buttons that light up under certain situation:</p> <table border="1"> <thead> <tr> <th>Button</th> <th>LED Display</th> <th>Indicates</th> </tr> </thead> <tbody> <tr> <td rowspan="3"></td> <td>Orange light flashing</td> <td>Charging</td> </tr> <tr> <td>Green light on</td> <td>Fully charged</td> </tr> <tr> <td>White light on</td> <td>Power on Bluetooth connected</td> </tr> <tr> <td rowspan="3"></td> <td>White light flashing</td> <td>Next song</td> </tr> <tr> <td>Blue light on</td> <td>Vibration on</td> </tr> <tr> <td>Green light flashing</td> <td>Volume adjusted</td> </tr> </tbody> </table>	Button	LED Display	Indicates		Orange light flashing	Charging	Green light on	Fully charged	White light on	Power on Bluetooth connected		White light flashing	Next song	Blue light on	Vibration on	Green light flashing	Volume adjusted	<p>To Clean</p> <p>The mynt sleep mask can be divided into "Functional Layer" and "Washable Layer".</p> <ul style="list-style-type: none"> Unplug the charger from the device before cleaning. The washable layer can be washed in anyway convenient, including wash by hands or by washing machine. Clean the functional layer only with a soft, slightly damp sponge. Do not immerse the functional layer in any liquid to clean. Never use abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish, paint thinner to clean. 	<p>Trouble Shooting</p> <table border="1"> <thead> <tr> <th>Problems</th> <th>Possible Cause</th> <th>Solutions</th> </tr> </thead> <tbody> <tr> <td>Can not switch on</td> <td>Battery low</td> <td>Charge and try again</td> </tr> <tr> <td rowspan="2">Product suddenly stop working</td> <td>Battery low</td> <td>Charge and try again</td> </tr> <tr> <td>2-hour timer setting activated</td> <td>Long press  for 1.5 seconds to switch on again</td> </tr> </tbody> </table> <p>Notice</p> <p>This symbol means the device must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment. For more information on disposal and recycling of this device, contact your local municipality, disposal service, or the shop where you bought this product.</p>	Problems	Possible Cause	Solutions	Can not switch on	Battery low	Charge and try again	Product suddenly stop working	Battery low	Charge and try again	2-hour timer setting activated	Long press  for 1.5 seconds to switch on again	<p>Warranty</p> <p>For consumers, who are covered by consumer protection laws or regulations in their country of purchase or, if different, their country of residence, the benefits conferred by Mynt's Limited Warranty are in addition to all rights and remedies conveyed by such consumer protection laws and regulations, including but not limited to these additional rights.</p> <p>For purchases made through Mynt:</p> <p>30-Day Money-Back Guarantee for Any Reason You may contact marketplace to return your undamaged product and packaging within 30 days of purchase to receive a FULL REFUND for any reason. If the reason for return isn't quality-related, the customer must pay the return shipping costs.</p>	<p>Customer Service</p> <ul style="list-style-type: none"> 24-Month Warranty Lifetime Technical Support Support@HelloMynt.com Online support by scanning the QR code, warranty tab. 	
Button	LED Display	Indicates																																				
	Orange light flashing	Charging																																				
	Green light on	Fully charged																																				
	White light on	Power on Bluetooth connected																																				
	White light flashing	Next song																																				
	Blue light on	Vibration on																																				
	Green light flashing	Volume adjusted																																				
Problems	Possible Cause	Solutions																																				
Can not switch on	Battery low	Charge and try again																																				
Product suddenly stop working	Battery low	Charge and try again																																				
	2-hour timer setting activated	Long press  for 1.5 seconds to switch on again																																				