VIEWWORKS

VIVIX-S VW Service Manual

CE

© Vieworks. 2019 All rights reserved.

The specifications and related information in this manual may be changed without notice. Refer to Vieworks Download System (VDS) for the latest version of our manuals.

Under copyright laws, this manual should not be reproduced, in whole or in part, without the written permission of Vieworks.

Contents

1. Instruction	6
1.1 Document Guide	7
1.1.1 Caution	7
1.1.2 Target Audience.....	7
1.1.3 Symbols.....	7
1.1.4 Notations.....	7
1.1.5 Contact Us	7
1.2 Intended Use.....	8
1.2.1 Features.....	8
1.3 Product Use.....	9
1.3.1 Product Usage.....	9
1.3.2 Disclaimer.....	9
1.3.3 Product Disposal.....	10
1.3.4 Trademarks.....	10
2. Product.....	11
2.1 Product Components.....	12
2.1.1 Detector.....	12
2.1.2 Battery Pack.....	13
2.1.3 SCU (System Control Unit).....	14
2.1.4 Accessories.....	16
2.2 FXRD-3643VAW, FXRD-3643VAW PLUS.....	18
2.2.1 Specifications.....	18
2.2.2 Functions	19
2.2.3 Deco Sheet	21
2.2.4 Wireless Communication.....	22
2.2.5 Use Environment.....	22
2.2.6 Load Limit of Detector	23
2.3 FXRD-4343VAW, FXRD-4343VAW PLUS.....	24
2.3.1 Specifications.....	24
2.3.2 Functions	25
2.3.3 Deco Sheet	27
2.4 FXRD-2530VAW, FXRD-2530VAW PLUS.....	28
2.4.1 Specifications.....	28
2.4.2 Functions	29
2.4.3 Deco Sheet	31

2.5	Battery Pack (FXRB-04A)	32
2.5.1	Specifications	32
2.6	SCU Lite (FXRP-02A)	33
2.6.1	Specifications	33
2.6.2	Drawings	33
2.6.3	Functions	34
2.7	SCU mini (FXRS-04A)	35
2.7.1	Specifications	35
2.7.2	Drawings	35
2.7.3	Functions	36
2.8	Others	38
2.8.1	X-ray Generator (Recommended Exposure Condition)	38
2.8.2	Recommended Specifications of Workstation (PC)	38
2.8.3	Recommended Specifications of Grid	38
3.	System Configuration	39
3.1	Detector Connection Methods	40
3.1.1	Wireless Connection	40
3.1.2	Wired Connection	41
3.2	Diagram	43
3.2.1	Block Diagram	43
3.2.2	Wiring Diagram	43
3.3	System Configuration	44
3.3.1	SCU AP Mode	44
3.3.2	Tether Interface Mode	45
3.3.3	External AP Mode	46
3.3.4	AP Switching	47
3.3.5	Sync Wireless Settings	48
3.4	Generator Interface	49
3.4.1	AED (Auto Exposure Detection) Interface	49
3.4.2	DR Trigger Interface	51
3.5	Configuring DR Trigger Interface	53
3.5.1	Trigger Interface Way	53
3.5.2	Packet Trigger	53
3.5.3	Line Trigger	54
3.5.4	EXT_INF Port Pin Map	55
3.5.5	Input / Output Circuit	56
4.	Settings	57

4.1	Product Installation	58
4.1.1	Connecting SCU mini (FXRS-04A).....	58
4.1.2	Booting SCU mini (FXRS-04A)	59
4.1.3	Connecting SCU Lite (FXRP-02A)	60
4.1.4	Booting Up SCU Lite (FXRP-02A).....	61
4.1.5	Booting up the Detector	61
4.2	Device Setting	62
4.2.1	Software Installation	62
4.2.2	Setting Detector and SCU	62
4.3	Diagnosis of Devices	63
4.3.1	Diagnosing Images	63
4.3.2	Diagnosing Battery Pack	64
4.3.3	Diagnosing Communication Status.....	65
4.3.4	Diagnosing Communication Speed	68
4.3.5	Self-Diagnosis	69
5.	Maintenance.....	72
5.1	Product Initialization	73
5.1.1	SCU Initialization	73
5.1.2	Detector Initialization.....	74
5.1.3	Wireless Initialization of Detector	75
5.2	Detector Power Save Function	76
5.3	How to Replace the Battery Pack	78
5.3.1	Replacing the Battery of FXRD-3643VAW / FXRD-3643VAW PLUS / FXRD-4343VAW / FXRD-4343VAW PLUS.....	79
5.3.2	Replacing the Battery of FXRD-2530VAW / FXRD-2530VAW PLUS	81
6.	Troubleshooting	82
6.1	Troubleshooting.....	83
6.1.1	Troubleshooting Guide.....	83
6.1.2	Fail to Turn the Detector On	83
6.1.3	The Power Switch of SCU or Status LED is not worked normally	83
6.1.4	Communication Test is failed.....	84
6.1.5	There is no problem in communication connection, but no image is acquired.....	84
6.1.6	Getting Abnormal Images	85
6.1.7	Battery Pack or Installation Part of Battery is Getting Hot.....	85
7.	Information	86
7.1	Service Information	87

7.1.1 Product Lifetime87

7.1.2 Regular Inspection and Maintenance87

7.1.3 Repair87

7.1.4 Replacement Parts Support87

7.1.5 Consumables.....87

7.2 Warranty 88

7.3 Revision History 89

1. Instruction

This section gives basic information of this manual and products.

Document Guide

Intended Use

Product Use

1.1 Document Guide

This operation manual explains how to use the **VIVIX-S VW** detector models made by Vieworks, X-ray interface unit, and other peripheral equipment. Use this manual to set up the **VIVIX-S VW** detector as well as understand its various functions.

1.1.1 Caution

If the user is not fully acquainted with this manual, the product can be malfunctioned, or unsuspected problem can be happened due to carelessness. To prevent any medical accidents, the user should fully understand the instructions of this manual before operating this product.

1.1.2 Target Audience

This manual is intended for the users who set up and operate the **VIVIX-S VW** detector models.

1.1.3 Symbols

This product should be operated under the safety instructions with the warning or caution symbol in this manual. It is important for you to read and understand the contents to operate the products safely.

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1.1.4 Notations

Bold Types

Words in bold indicate products terms, or the sentences which are needed to transmit clear meaning to the customers.

1.1.5 Contact Us

For comments or inquiries regarding this document and relevant products, contact via email below.

Item	Contents
Department	Customer Support Team at Vieworks
E-mail	CustomerSupport@vieworks.com

1.2 Intended Use

VIVIX-S VW detectors are digital X-ray imaging solution. They acquire images by detecting X-rays that has been passed through the human body. When X-ray photons pass through the scintillator in the detector, the photons convert to visible ray, and the visible ray is converted to electronic signals through TFTs – thin film transistors (a-Si). Then the detector digitalizes X-ray images and transfers them to the PC (workstation) for diagnostic review using an image display monitor. Advanced digital image processing also allows efficient diagnosis, information management, and sharing of image information over the network.

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1.2.1 Features

- Since VIVIX-S VW detectors are compatible with a conventional film cassette, they enable to replace the analog radiographic diagnosis.
- VIVIX-S VW detectors provide high resolution digital images. The resolution varies depending on the models as follows.
 - **FXRD-3643VAW**, FXRD-3643VAW PLUS : 140µm pixel pitch sensor (Approx. 7.9 Mega Pixels)
 - **FXRD-4343VAW**, FXRD-4343VAW PLUS : 140µm pixel pitch sensor (Approx. 9.4 Mega Pixels)
 - **FXRD-2530VAW**, FXRD-2530VAW PLUS : 124µm pixel pitch sensor, (Approx. 5.2 Mega Pixels)
- The scintillator provides two types of CsI, and you can use one of the detectors.
- VIVIX-S VW detectors support IEEE 802.11n/ac to acquire images without a wired connection in anytime, anywhere.
- Checks the connecting status of the detector and battery remains through OLED Display.
- Makes direct wireless communication with the built-in wireless AP function. (Inside AP™)
- Makes quick application of various functions with two buttons on the side of the detector.
- Supports the stable and reliable AED (Auto Exposure Detection) function. (Anytime™)
- Designed with a thin and light structure, it is portable for easy shooting anytime, anywhere.
- Designed for easy carrying
- You can use the product with cable connection (tether interface cable) depending on the using environment.

1.3 Product Use

This chapter provides instructions about the use of the product, disposal and the liability limit of Vieworks.

1.3.1 Product Usage

- 1 Only a physician or a legally certified operator should use this product.
- 2 The equipment should be kept in a safe and operable condition by maintenance personnel.
- 3 Follow the guidelines in this manual when installing and using this product.
- 4 Use only computers and image display monitors recommended in this manual.
- 5 Use only the dedicated cables provided with this product.
- 6 For details about installation and use of the product, consult your sales representative or a distributor.

1.3.2 Disclaimer

- 1 In no event shall Vieworks be liable for damage or loss arising from a fire, earthquake, any action or accident by a third party, any intentional or negligent action by users.
- 2 In no event shall Vieworks be liable for damage or loss arising from any trial usage, or other usage under abnormal conditions.
- 3 In no event shall Vieworks be liable for personal physical harm or property damage that is sustained when the instructions of this manual are not followed.
- 4 In no event shall Vieworks be liable for any damage arising from moving, alteration, inspection or repair the product by a person other than an authorized service engineer by Vieworks.
- 5 In no event shall Vieworks be liable for loss of image data for any reason while using this product.
- 6 Roentgenography, image processing, image reading, and image data storage must be performed in accordance with the laws of the country or region in which the product is being used.
- 7 The user is responsible for maintaining the privacy of image data acquired from this product.
- 8 It is the responsibility of the attending physicians to provide medical care services. Vieworks will not be liable for faulty diagnoses.
- 9 Specifications, composition, and appearance of this product may change without prior notice.

1.3.3 Product Disposal

Disposal of this product in an unlawful manner may have a negative impact on health and on the environment. When disposing of this product, therefore, be absolutely sure to follow the procedure that conforms to the laws and regulations applicable in your area.

European Union (and EEA*) only

This symbol indicates that this product is not to be disposed with your household waste, according to the WEEE Directive (2012/19/EC) and your national law.

This product should be handed over to a designated collection point, e.g., on an authorized one-for-one basis when you buy a new similar product or to an authorized collection site for recycling electrical and electronic equipment (EEE). Improper handling of this type of waste

could have a negative impact on the environment and human health due to potentially hazardous substances that are generally associated with EEE. At the same time, your cooperation in the correct disposal of this product will contribute to the effective usage of natural resources. For more information on where you can drop off your waste equipment for recycling, please contact your local city office, waste authority, approved WEEE scheme, or your household waste disposal service.

*EEA: Norway, Iceland, and Liechtenstein

1.3.4 Trademarks

The name "Vieworks" and Vieworks logo are registered trademarks of Vieworks.

© Vieworks. 2019 All rights reserved.

2. Product

This section gives an instruction about the product components and their specifications.

Product Components

FXRD-3643VAW, FXRD-3643VAW PLUS

FXRD-4343VAW, FXRD-4343VAW PLUS

FXRD-2530VAW, FXRD-2530VAW PLUS

Battery Pack (FXRB-04A)

SCU Lite (FXRP-02A)

SCU mini (FXRS-04A)

Others

2.1 Product Components

2.1.1 Detector

Product Name	Image	Model/Description
		<ul style="list-style-type: none"> ※ Basic: Battery Pack (2ea) <ul style="list-style-type: none"> • FXRD-3643VAW (3.15kg) • FXRD-3643VAW PLUS (3.15kg) ※ Option: Battery Pack 1ea <ul style="list-style-type: none"> • FXRD-3643VAW (2.95kg) • FXRD-3643VAW PLUS (2.95kg)
<p>VIVIX-S VW Wireless Detector</p>		<ul style="list-style-type: none"> ※ Basic: Battery Pack (2ea) <ul style="list-style-type: none"> • FXRD-4343VAW (3.7kg) • FXRD-4343VAW PLUS (3.7kg) ※ Option: Battery Pack 1ea <ul style="list-style-type: none"> • FXRD-4343VAW (3.5kg) • FXRD-4343VAW PLUS (3.5kg)
		<ul style="list-style-type: none"> ※ Basic: Battery Pack (1ea) <ul style="list-style-type: none"> • FXRD-2530VAW (1.9kg)

- The design of the decor sheet attached to the detector can be different for each customer.
- The battery pack is built into the detector. The battery pack can be replaced at the end of its life.
- FXRD-3643VAW, FXRD-3643VAW PLUS, FXRD-4343VAW, and FXRD-4343VAW PLUS are equipped two battery packs, and only one can be installed as an option upon request.
- FXRD-2530VAW and FXRD-2530VAW PLUS can only be equipped with one battery pack.

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2.1.2 Battery Pack

Name	Image	Model/Description
Battery Pack		<ul style="list-style-type: none"> • FXRB-04A (3,400mAh, 185g)

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- AC Power Cable (2m)
- 110V (1EA) / 220V (1EA)

- SCU is not a standard item provided to customers. If necessary, the user can optionally purchase the desired SCU according to the intended use.
- SCU Lite does not support wireless connections. Only use it as a wired connection.
- The SCU mini supports both wired and wireless connections. If you only need a wireless connection, you can use a separate external AP or detector AP.
- SCU n'est pas un article standard fourni aux clients. Si nécessaire, l'utilisateur peut éventuellement acheter le SCU souhaité en fonction de l'utilisation prévue.
- SCU Lite ne prend pas en charge les connexions sans fil. Utilisez-le uniquement comme connexion filaire.
- Le SCU mini prend en charge les connexions filaires et sans fil. Si vous n'avez besoin que d'une connexion sans fil, vous pouvez utiliser un point d'accès externe ou un point de détection séparé.

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Velcro tape
(25mm x 50mm / 6ea)

- Attached to the corners of bucky and the friction surface of a detector.

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2.2 FXRD-3643VAW, FXRD-3643VAW PLUS

2.2.1 Specifications


Item	Specifications
Model	<ul style="list-style-type: none"> FXRD-3643VAW FXRD-3643VAW PLUS
Image Sensor	<ul style="list-style-type: none"> TFT: a-Si (Amorphous Silicon)
X-ray Scintillator Type	<ul style="list-style-type: none"> FXRD-3643VAW : CsI type A FXRD-3643VAW PLUS : CsI type B
Pixel Pitch	<ul style="list-style-type: none"> 0.14mm (140μm)
Field of View	<ul style="list-style-type: none"> 36cm x 43cm (14" x 17")
Active Area (H x V)	<ul style="list-style-type: none"> 358.4mm x 430.08mm
Active Array	<ul style="list-style-type: none"> 2560 x 3072 pixels
Effective Area	<ul style="list-style-type: none"> 355.04mm x 426.72mm
Effective Array	<ul style="list-style-type: none"> 2536 x 3048 pixels
Grayscale	<ul style="list-style-type: none"> 16 bit
Spatial Resolution	<ul style="list-style-type: none"> Min. 3.5 lp/mm
Image Acquisition Time (Wired)	<ul style="list-style-type: none"> Up to 3 sec. (Set exposure time to 500ms / except exposure time)
Image Acquisition Time (Wireless)	<ul style="list-style-type: none"> Up to 4.5 sec. (IEEE802.11n, MIMO 2x2, 5GHz, 40MHz) Up to 3 sec. (IEEE802.11ac, MiMO 3x3, 5GHz, 80MHz) (Exposure time is set to 500ms, except exposure time)
Cycle Time	<ul style="list-style-type: none"> Min. 4 sec. (with optimal wired / wireless environment, exposure time is set to 500ms, excluding software processing time)
X-ray Synchronous Control	<ul style="list-style-type: none"> AED (Auto Exposure Detection) DR Trigger (External Line Trigger)
Power Supply	<ul style="list-style-type: none"> Powered from SCU via tether interface cable: DC 24V, Max. 1.0A Powered via AC-DC adapter: DC 18V, Max. 4.44A Powered via 2 battery packs: DC 9 ~13.2V, Max. 78.54Wh
Power Consumption	<ul style="list-style-type: none"> Max. 24W (without battery charged) Max. 80W (when charging battery)
Operating Time (Early life of battery)	<ul style="list-style-type: none"> One battery pack <ul style="list-style-type: none"> ▫ 7 hours (image acquisition every 100 seconds), 8 hours (standby) Two battery packs <ul style="list-style-type: none"> ▫ 14 hours (image acquisition every 100 seconds), 16 hours (standby)
Dimensions (H x W x D)	<ul style="list-style-type: none"> 384.0mm x 460.0mm x 15.0mm
Weight (including battery pack)	<ul style="list-style-type: none"> With one battery pack : 2.95kg With two battery packs : 3.15kg
Image Transmission	<ul style="list-style-type: none"> Wired: Gigabit Ethernet (1000BASE-T) via Power over Ethernet (PoE) Wireless: IEEE 802.11n / ac (2.4GHz / 5GHz), 3 antennas
Data Transmission Speed (Wired)	<ul style="list-style-type: none"> Max. 1Gbps
Data Transmission Speed (Wireless)	<ul style="list-style-type: none"> Max. 300Mbps (IEEE802.11n, MIMO 2x2, 5GHz, 40MHz) Max. 1300Mbps (IEEE802.11ac, MIMO 3x3, 5GHz, 80MHz)

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
2.2.2 Functions

Side

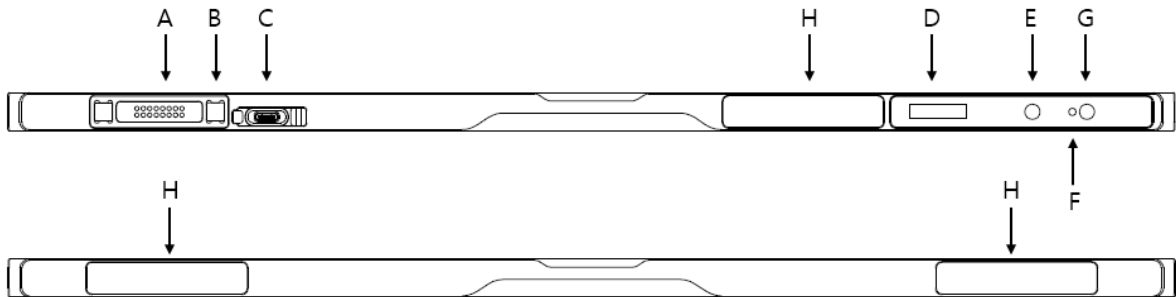
Name	Description
A Tether Interface Connector	Tether interface cable connector. <ul style="list-style-type: none"> • Used for wired connection between a detector and SCU.
B Magnet for fixing the tether interface	Used for fixing a tether interface cable
C AC-DC Adapter Connector	<ul style="list-style-type: none"> • Connector for fastening the AC-DC adapter • Used for fast battery charging
D OLED Display	<ul style="list-style-type: none"> • Displays battery status • Displays wired / wireless connection status • Displays sleep mode status
E AP Button	<ul style="list-style-type: none"> • Button for changing AP setting (Detector AP mode ON / OFF, AP switching when detector is in STATION mode) • Changes OLED screen
F Power Indicator LED	<ul style="list-style-type: none"> • Displays system power status • Displays system boot status
G Power Button	<ul style="list-style-type: none"> • System power on/off • Changes OLED screen
H Antenna for Wireless LAN	Antennas for wireless communication (3ea)

Rear

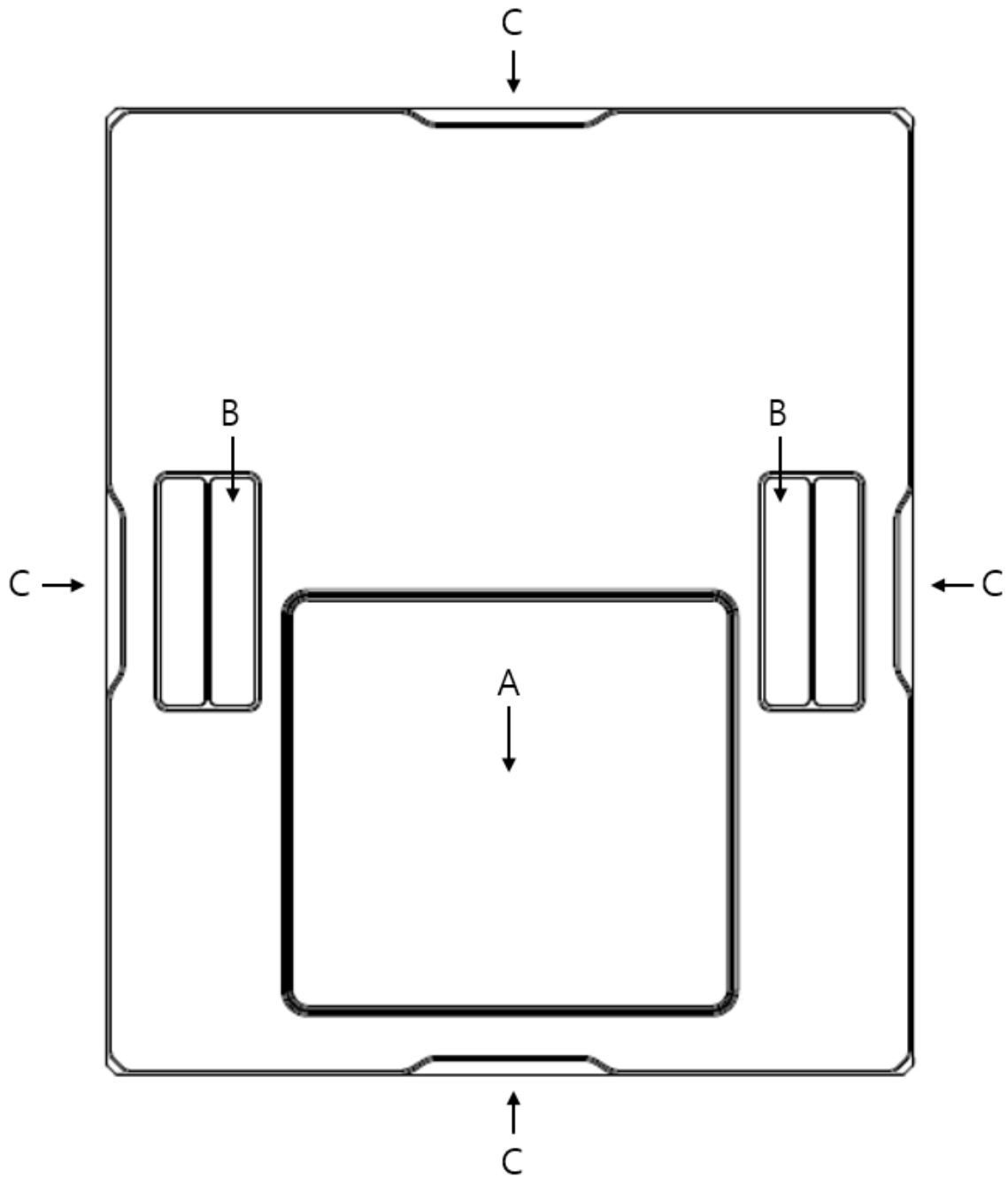	Name	Description
A	Battery Pack Cover	The cover needs to be opened and closed when replacing the battery pack.
B	Handle	A handle for carrying a detector
C	Lift Structure	Used when the detector is placed on flat surface

2.2.3 Deco Sheet

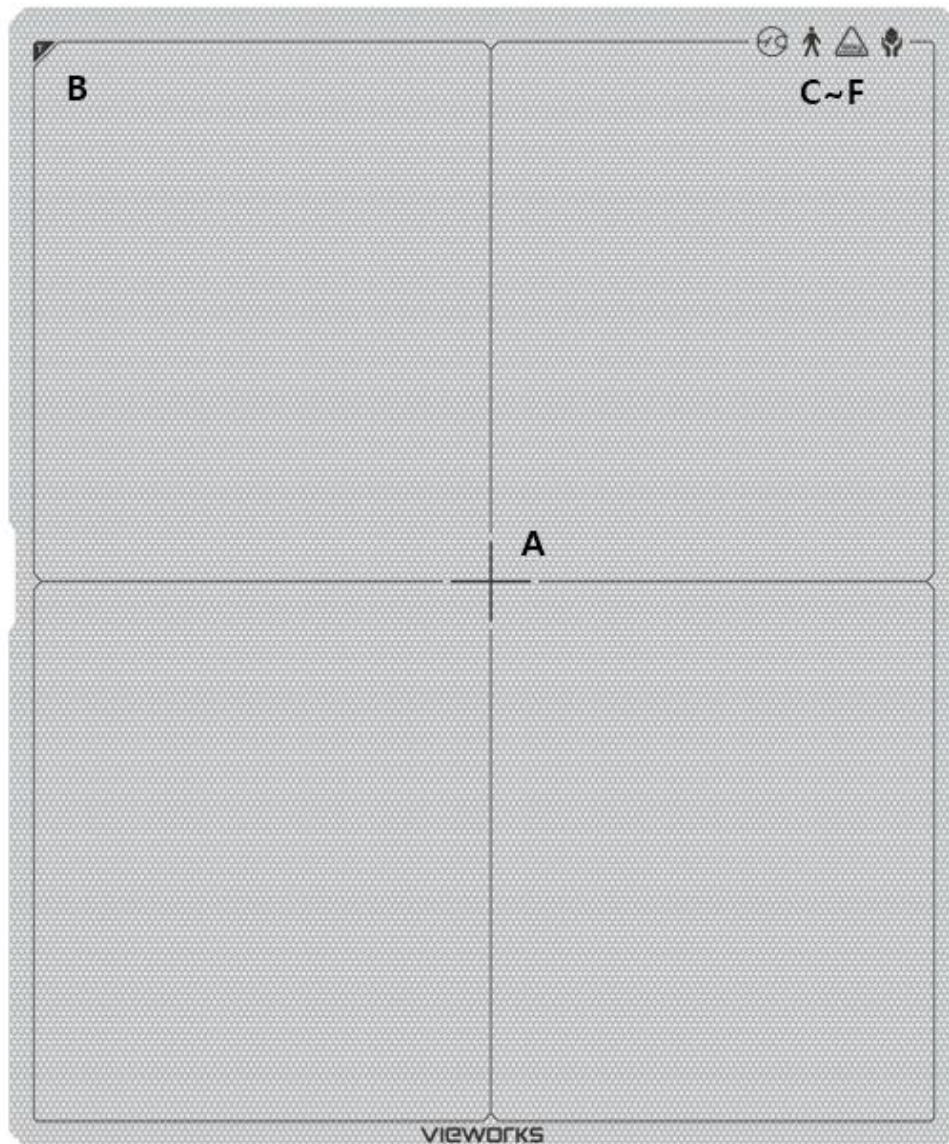	Indication Info.	Description
A	Center of the detector	Indicates the central position of detector.
B	Image starting point	Indicates the starting point of an original image.
C~F	Certification logo	Indicates the certification logos relating to a medical device.

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2.2.5 Use Environment

Item	Operation	Storage & Transportation
Temperature	0 ~ +40°C	-15 ~ +55°C
Humidity	5 ~ 90% (Non-condensing)	5 ~ 90% (Non-condensing)
Atmospheric pressure	700 ~ 1060 hPa	500 ~ 1060 hPa
Shock	20G	30G
Vibration	2G	5G



- The operating environment specifications of **VIVIX-S VW** detectors (FXRD-3643VAW, FXRD-3643VAW PLUS, FXRD-4343VAW, FXRD-4343VAW PLUS, FXRD-2530VAW, FXRD-2530VAW PLUS) are the same.
- Les spécifications d'environnement d'exploitation des détecteurs VW VIVIX-S (FXRD-3643VAW, FXRD-3643VAW PLUS, FXRD-4343VAW, FXRD-4343VAW PLUS, FXRD-2530VAW, FXRD-2530VAW PLUS) sont identiques.

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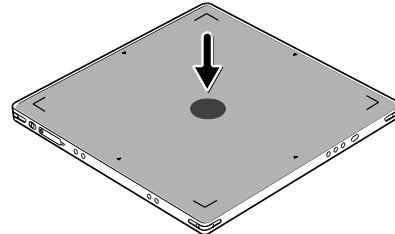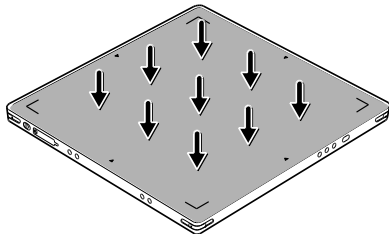2.3 FXRD-4343VAW, FXRD-4343VAW PLUS

2.3.1 Specifications


Item	Specifications
Model	<ul style="list-style-type: none"> FXRD-4343VAW FXRD-4343VAW PLUS
Image Sensor	<ul style="list-style-type: none"> TFT: a-Si (Amorphous Silicon)
X-ray Scintillator Type	<ul style="list-style-type: none"> FXRD-4343VAW: CsI A type FXRD-4343VAW PLUS: CsI B type
Pixel Pitch	<ul style="list-style-type: none"> 0.14mm (140μm)
Field of View	<ul style="list-style-type: none"> 43cm x 43cm (17" x 17")
Active Area (H x V)	<ul style="list-style-type: none"> 430.08mm x 430.08mm
Active Array	<ul style="list-style-type: none"> 3072 x 3072 pixels
Effective Area	<ul style="list-style-type: none"> 426.72mm x 426.72mm
Effective Array	<ul style="list-style-type: none"> 3048 x 3048 pixels
Grayscale	<ul style="list-style-type: none"> 16 bit
Spatial Resolution	<ul style="list-style-type: none"> Min. 3.5 lp/mm
Image Acquisition Time (Wired)	<ul style="list-style-type: none"> Up to 3 sec. (Set exposure time to 500ms / except exposure time)
Image Acquisition Time (Wireless)	<ul style="list-style-type: none"> Up to 4.5 sec. (IEEE802.11n, MIMO 2x2, 5GHz, 40MHz) Up to 3 sec. (IEEE802.11ac, MiMO 3x3, 5GHz, 80MHz) (Exposure time is set to 500ms, except exposure time)
Cycle Time	<ul style="list-style-type: none"> Min. 4 sec. (with optimal wired / wireless environment, exposure time is set to 500ms, excluding software processing time)
X-ray Synchronous Control	<ul style="list-style-type: none"> AED (Auto Exposure Detection) DR Trigger (External Line Trigger)
Power Supply	<ul style="list-style-type: none"> Powered from SCU via tether interface cable: DC 24V, Max. 1.0A Powered via AC-DC adapter: DC 18V, Max. 4.44A Powered via 2 battery packs: DC 9 ~13.2V, Max. 78.54Wh
Power Consumption	<ul style="list-style-type: none"> Max. 24W (without battery charged) Max. 80W (when charging battery)
Operating Time (Early life of battery)	<ul style="list-style-type: none"> One battery pack <ul style="list-style-type: none"> ▫ 6.5 hours (image acquisition every 100 seconds), 8 hours (standby) Two battery packs <ul style="list-style-type: none"> ▫ 13 hours (image acquisition every 100 seconds), 16 hours (standby)
Dimensions (H x W x D)	<ul style="list-style-type: none"> 460.0mm x 460.0mm x 15.0mm
Weight (including battery pack)	<ul style="list-style-type: none"> With one battery pack : 3.5kg With two battery packs : 3.7kg
Image Transmission	<ul style="list-style-type: none"> Wired: Gigabit Ethernet (1000BASE-T) via Power over Ethernet (PoE) Wireless: IEEE 802.11n / ac (2.4GHz / 5GHz), 3 antennas
Data Transmission Speed (Wired)	<ul style="list-style-type: none"> Max. 1Gbps
Data Transmission Speed (Wireless)	<ul style="list-style-type: none"> Max. 450Mbps (IEEE802.11n, MIMO 2x2, 5GHz, 40MHz) Max. 1300Mbps (IEEE802.11ac, MIMO 3x3, 5GHz, 80MHz)

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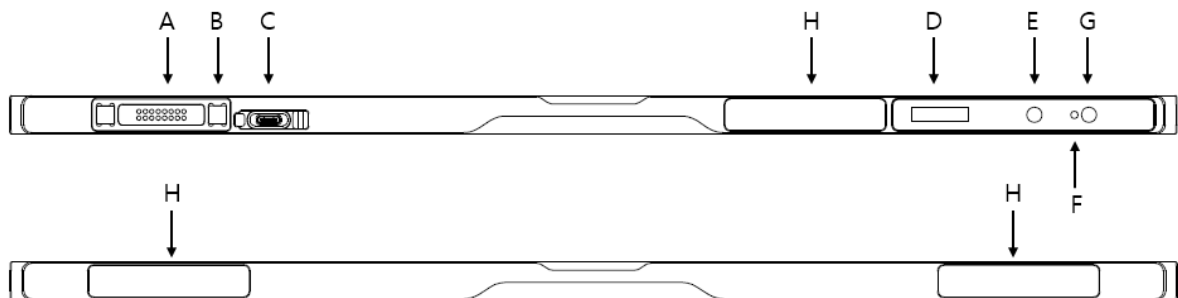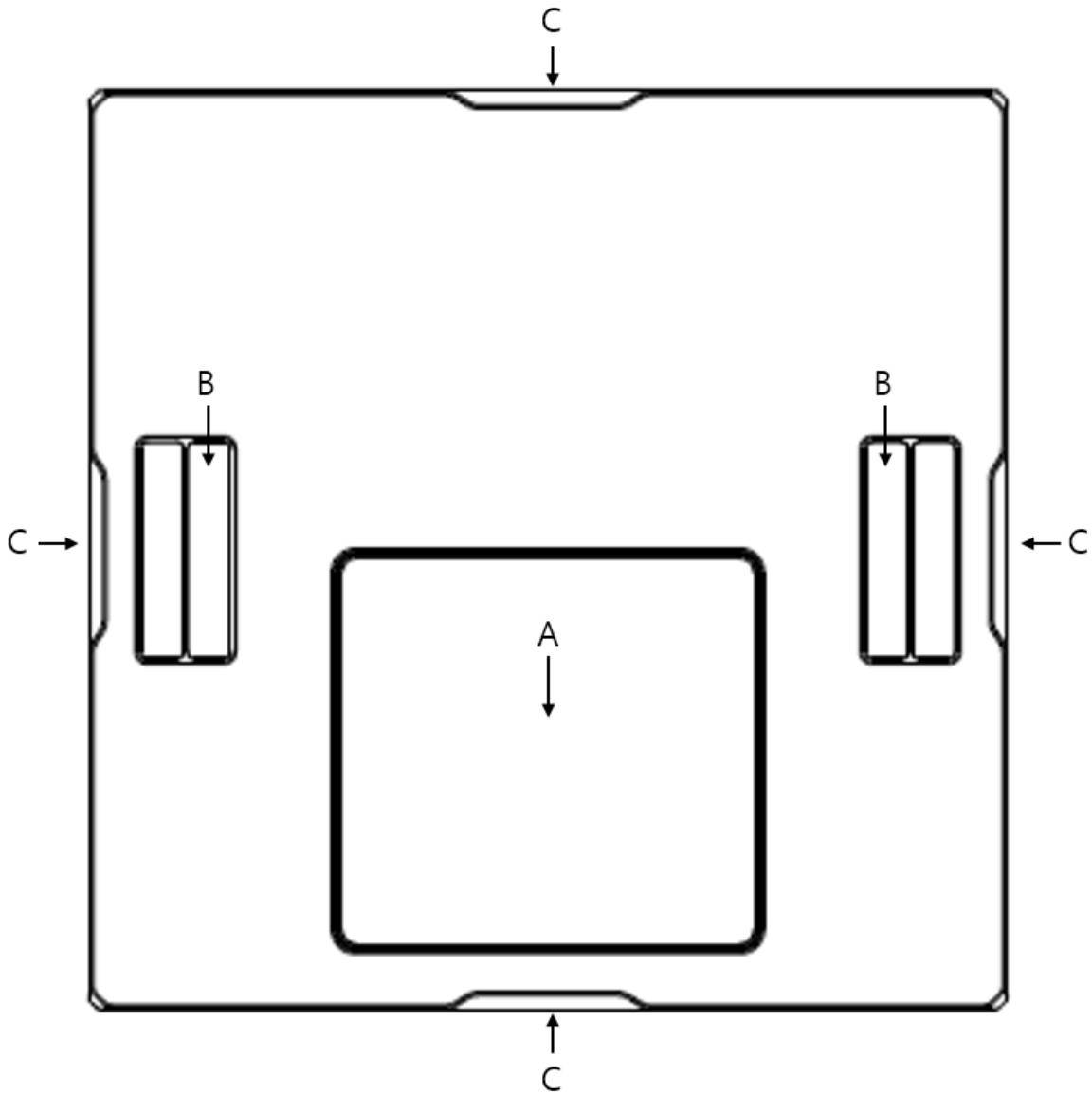
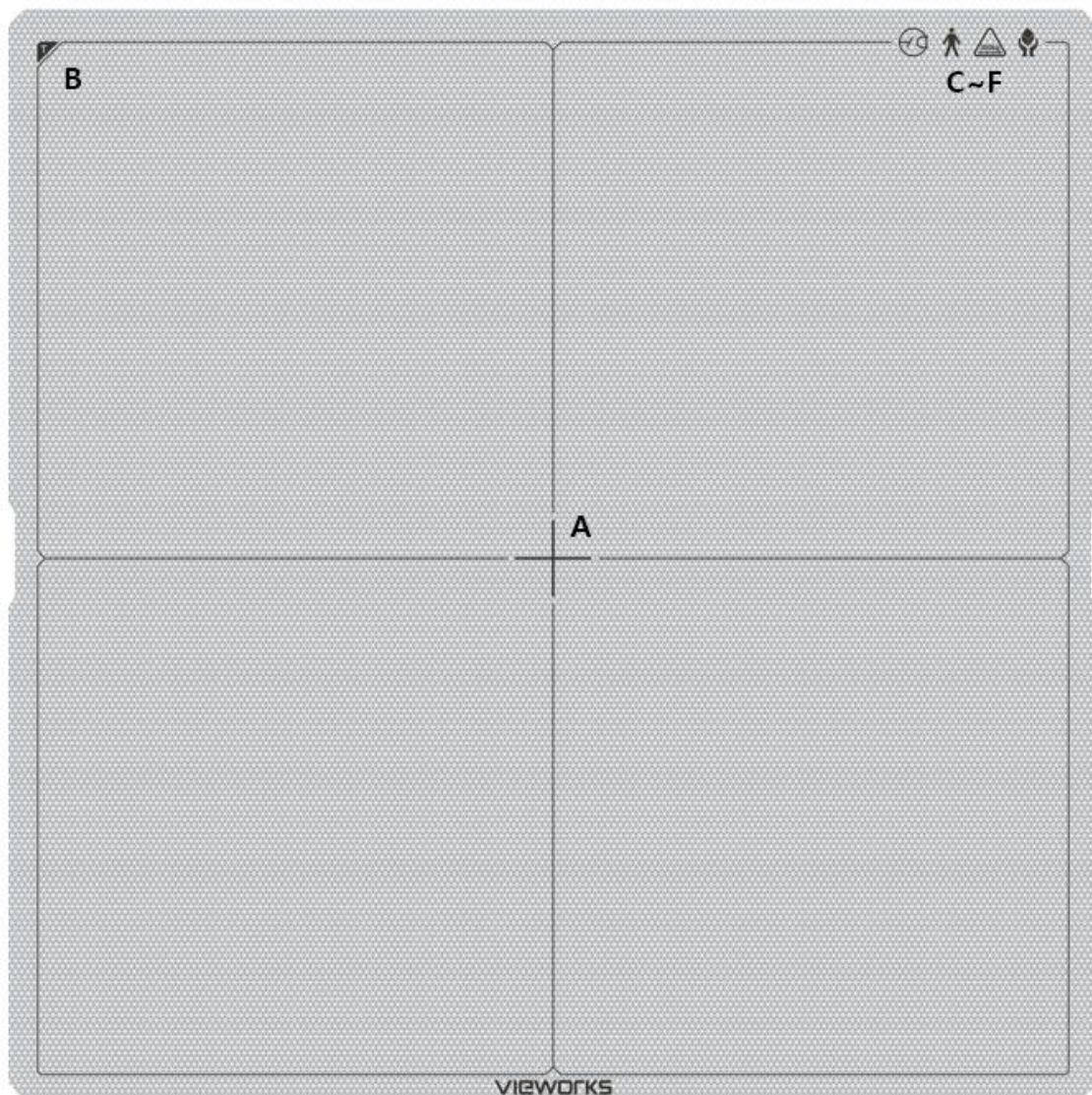
2.4 FXRD-2530VAW, FXRD-2530VAW PLUS

2.4.1 Specifications


Item	Specifications
Model	<ul style="list-style-type: none"> FXRD-2530VAW / FXRD-2530VAW PLUS
Image Sensor	<ul style="list-style-type: none"> TFT: a-Si (Amorphous Silicon)
X-ray Scintillator Type	<ul style="list-style-type: none"> FXRD-2530VAW: CsI A type FXRD-2530VAW PLUS: CsI B type
Pixel Pitch	<ul style="list-style-type: none"> 0.124mm (124μm)
Field of View	<ul style="list-style-type: none"> 25cm x 32cm (10" x 12")
Active Area (H x V)	<ul style="list-style-type: none"> 253.95mm x 317.44mm
Active Array	<ul style="list-style-type: none"> 2048 x 2560 pixels
Effective Area	<ul style="list-style-type: none"> 251.0mm x 314.5mm
Effective Array	<ul style="list-style-type: none"> 2024 x 2536 pixels
Grayscale	<ul style="list-style-type: none"> 16 bit
Spatial Resolution	<ul style="list-style-type: none"> Min. 4.0 lp/mm
Image Acquisition Time (Wired)	<ul style="list-style-type: none"> Up to 3 sec. (Set exposure time to 500ms / except exposure time)
Image Acquisition Time (Wireless)	<ul style="list-style-type: none"> Up to 4.5 sec. (IEEE802.11n, MIMO 3x3, 5GHz, 40MHz) Up to 3 sec. (IEEE802.11ac, MiMO 3x3, 5GHz, 80MHz) (Exposure time is set to 500ms, except exposure time)
Cycle Time	<ul style="list-style-type: none"> Min. 4 sec. (with optimal wired / wireless environment, exposure time is set to 500ms, excluding software processing time)
X-ray Synchronous Control	<ul style="list-style-type: none"> AED (Auto Exposure Detection) DR Trigger (External Line Trigger)
Power Supply	<ul style="list-style-type: none"> Powered from SCU via tether interface cable: DC 24V, Max. 0.625A Powered via AC-DC adapter: DC 18V, Max. 2.78A Powered via 2 battery packs: DC 9 ~13.2V, Max. 39.27Wh
Power Consumption	<ul style="list-style-type: none"> Max. 15W (without battery charged) Max. 50W (when charging battery)
Operating Time (Early life of battery)	<ul style="list-style-type: none"> One battery pack <ul style="list-style-type: none"> 7 hours (image acquisition every 100 seconds), 8 hours (standby)
Dimensions (H x W x D)	<ul style="list-style-type: none"> 287.0mm x 350.0mm x 15.0mm
Weight (including battery pack)	<ul style="list-style-type: none"> 1.9kg
Image Transmission	<ul style="list-style-type: none"> Wired: Gigabit Ethernet (1000BASE-T) via Power over Ethernet (PoE) Wireless: IEEE 802.11n / ac (2.4GHz / 5GHz), 3 antennas
Data Transmission Speed (Wired)	<ul style="list-style-type: none"> Max. 1Gbps
Data Transmission Speed (Wireless)	<ul style="list-style-type: none"> Max. 450Mbps (IEEE802.11n, MIMO 3x3, 5GHz, 40MHz) Max. 1300Mbps (IEEE802.11ac, MIMO 3x3, 5GHz, 80MHz)

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
2.4.2 Functions

Side

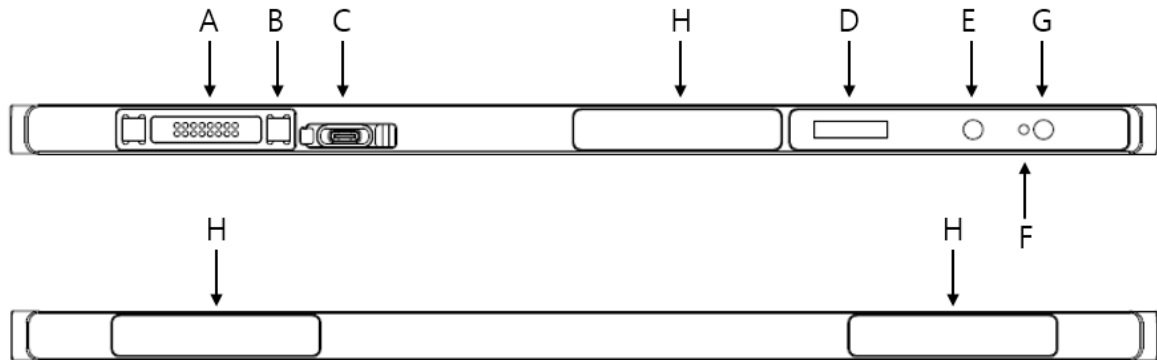
Name	Description
A Tether Interface Connector	Tether interface cable connector. <ul style="list-style-type: none"> Used for wired connection between a detector and SCU.
B Magnet for fixing the tether interface	Used for fixing a tether interface cable
C AC-DC Adapter Connector	<ul style="list-style-type: none"> Connector for fastening the AC-DC adapter Used for fast battery charging
D OLED Display	<ul style="list-style-type: none"> Displays battery status Displays wired / wireless connection status Displays sleep mode status
E AP Button	<ul style="list-style-type: none"> Button for changing AP setting (Detector AP mode ON / OFF, AP switching when detector is in STATION mode) Changes OLED screen
F Power Indicator LED	<ul style="list-style-type: none"> Displays system power status Displays system boot status
G Power Button	<ul style="list-style-type: none"> System power on/off Changes OLED screen
H Antenna for Wireless LAN	Antennas for wireless communication (3ea)

Rear

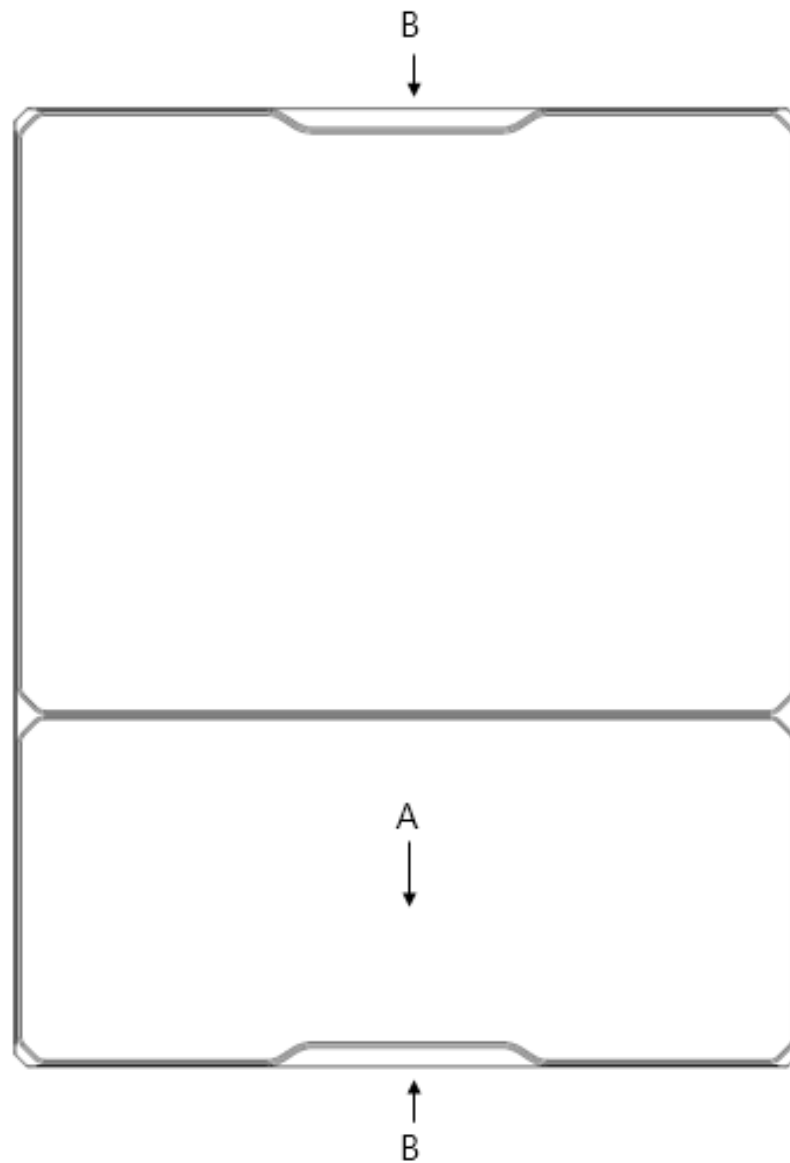	Name	Description
A	Battery Pack Cover	The cover needs to be opened and closed when replacing the battery pack.
C	Lift Structure	Used when the detector is placed on flat surface

2.4.3 Deco Sheet

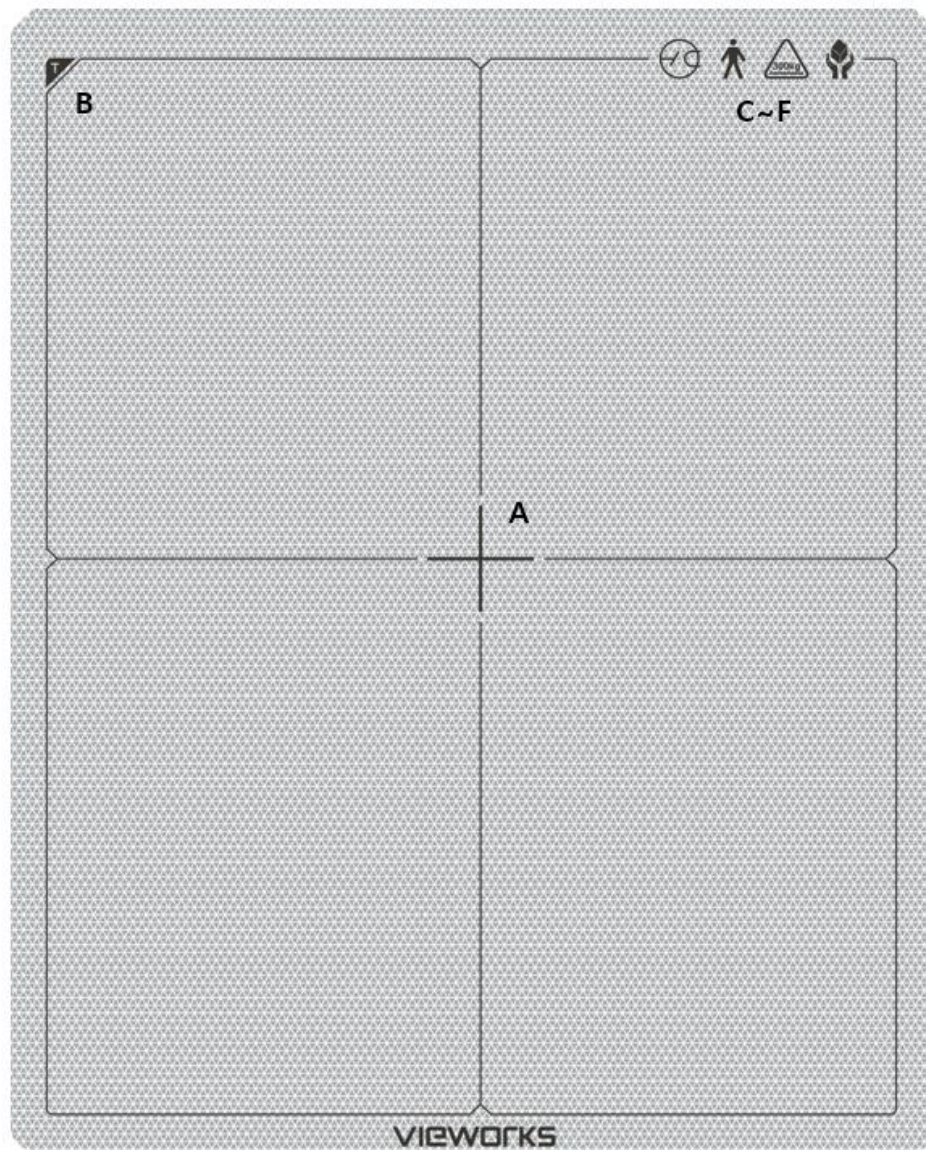	Indication Info.	Description
A	Center of the detector	Indicates the central position of detector.
B	Image starting point	Indicates the starting point of an original image.
C~F	Certification logo	Indicates the certification logos relating to a medical device.

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2.5 Battery Pack (FXRB-04A)

2.5.1 Specifications

Item	Specifications
Model	• FXRB-04A
Type	• Lithium Ion Polymer
Nominal Voltage	• DC +11.55V
Nominal Capacity	• 3,400mAh
Number of Cells	• 3S1P (3 Series 1 Parallel)
Battery Life	• Approx. 800 times (One time standard: Fully charged and then fully discharged)
Dimensions (H × W × D)	• 189.0mm × 89.0mm × 6.65mm (max.)
Weight	• 185g (max.)

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
2.6 SCU Lite (FXRP-02A)

SCU Lite is a device that transmits images by wired connection between the detector and workstation.

2.6.1 Specifications


Item	Specifications
Model	<ul style="list-style-type: none"> FXRP-02A
Power Supply	<ul style="list-style-type: none"> Input: DC +24V 1.0A Max.
Cable Connection Port	<ul style="list-style-type: none"> Gigabit Ethernet port (1ea) PoE (Power over Ethernet) Port (1ea)
Dimensions (H × W × D)	<ul style="list-style-type: none"> 108mm × 109mm × 29.5mm
Weight	<ul style="list-style-type: none"> 330g

2.6.2 Drawings

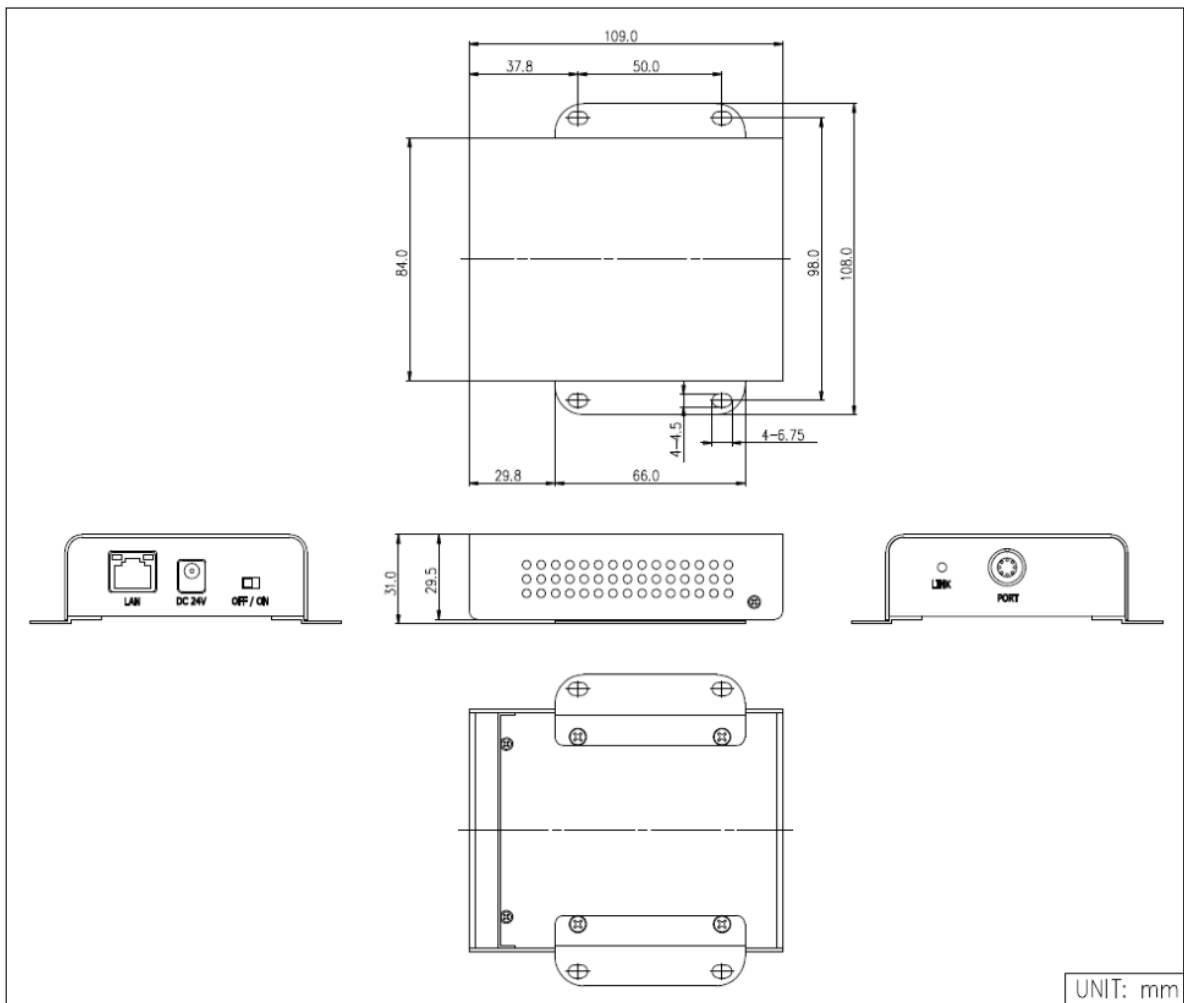
2.6.3 Functions

Front Side

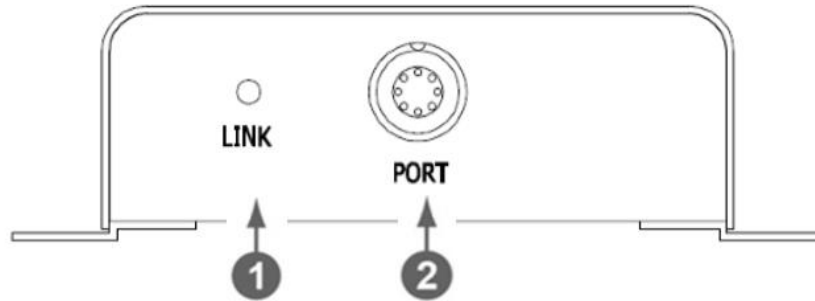No.	Name	Description
1	LINK LED	Indicates the status of PoE port. <ul style="list-style-type: none"> • Green : 1 Gbps • Yellow : 100 Mbps
2	POE Port	PoE (Power over Ethernet) port (1000BASE-T) <ul style="list-style-type: none"> • Communication between the detector and SCU Lite. • Supplies power to the detector.

Rear Side

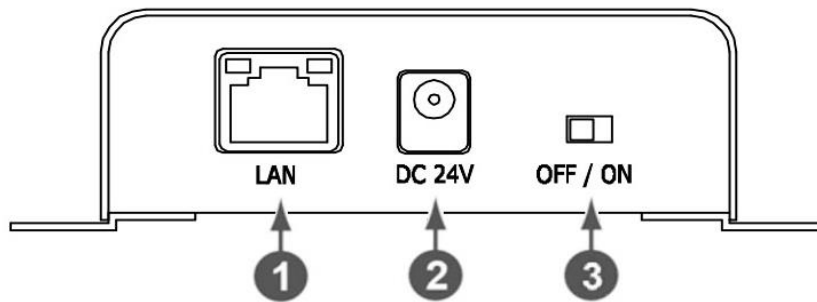No.	Name	Description
1	LAN port	Gigabit Ethernet port (1000BASE-T) <ul style="list-style-type: none"> • Communication between the workstation and SCU Lite.
2	DC power input port	DC +24V <ul style="list-style-type: none"> • Supplies power to SCU Lite.
3	Power Switch	Turns on/off the power of SCU Lite.

2.7 SCU mini (FXRS-04A)


The SCU mini is located between the detector, the workstation and the X-ray generator to synchronize images and signals.

The SCU mini has a wired or wireless connection to the detector and a wired connection to the workstation. It can also be used by wired connection with the X-ray generator (DR Trigger).

2.7.1 Specifications


Item	Specifications
Model	• FXRP-04A
Power Supply	• Input: DC +24V, 2A Max.
Cable Connection Port	• Gigabit Ethernet port (3ea) • PoE (Power over Ethernet) Port (1ea)
Wireless Communication	• IEEE 802.11n (2.4 GHz / 5 GHz)
Dimension (H × W × D)	• 210.0mm × 170.0mm × 45.0mm
Antenna	• 140mm (2ea, dual band)
Weight	• 1.2kg

2.7.2 Drawings

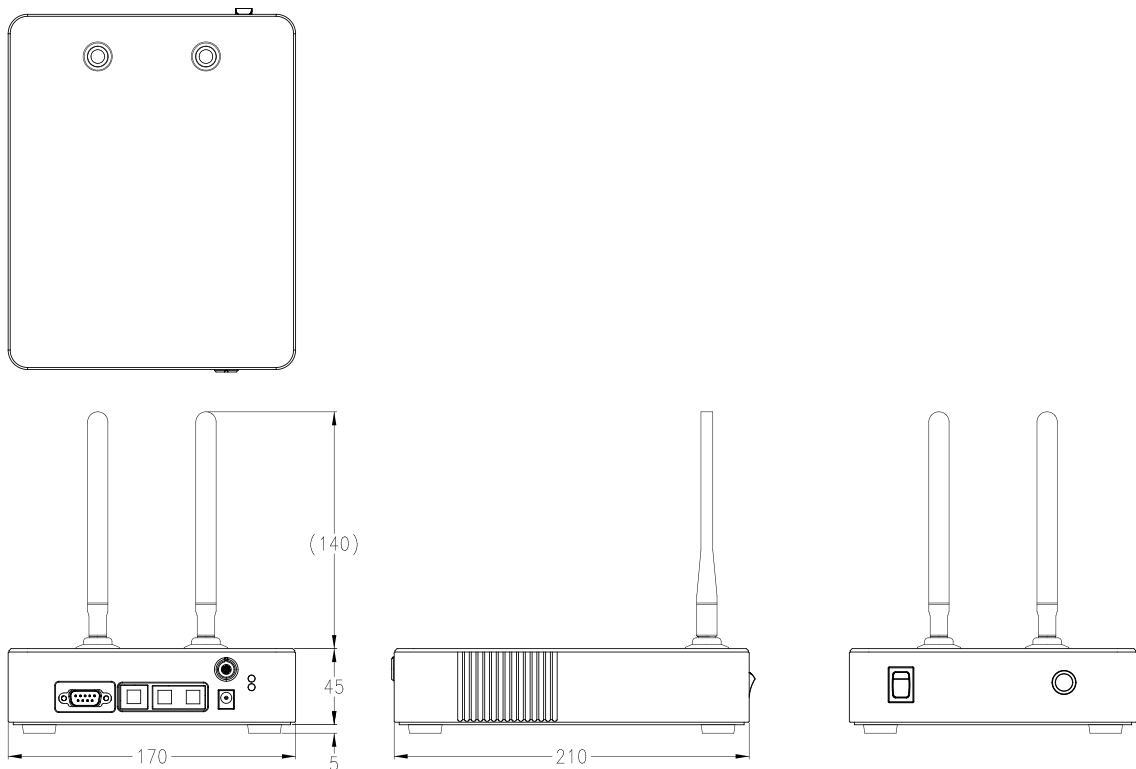
2.7.3 Functions

Front Side

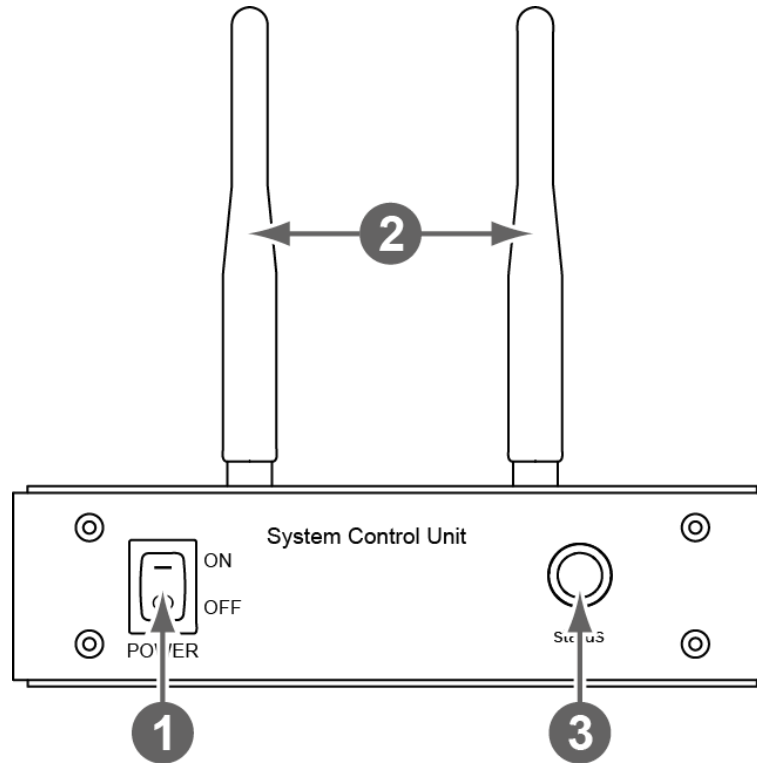No.	Name	Description
1	Power switch	Power switch of SCU mini. (ON / OFF)
2	Antenna	Assists communications between the detector and SCU mini. Indicates status of SCU mini operation and connection.
3	Status LED	<ul style="list-style-type: none"> • Blinking green: Booting • Green: Completed to boot up • Blue: The detector is connected and ready to communicate.

Rear Side

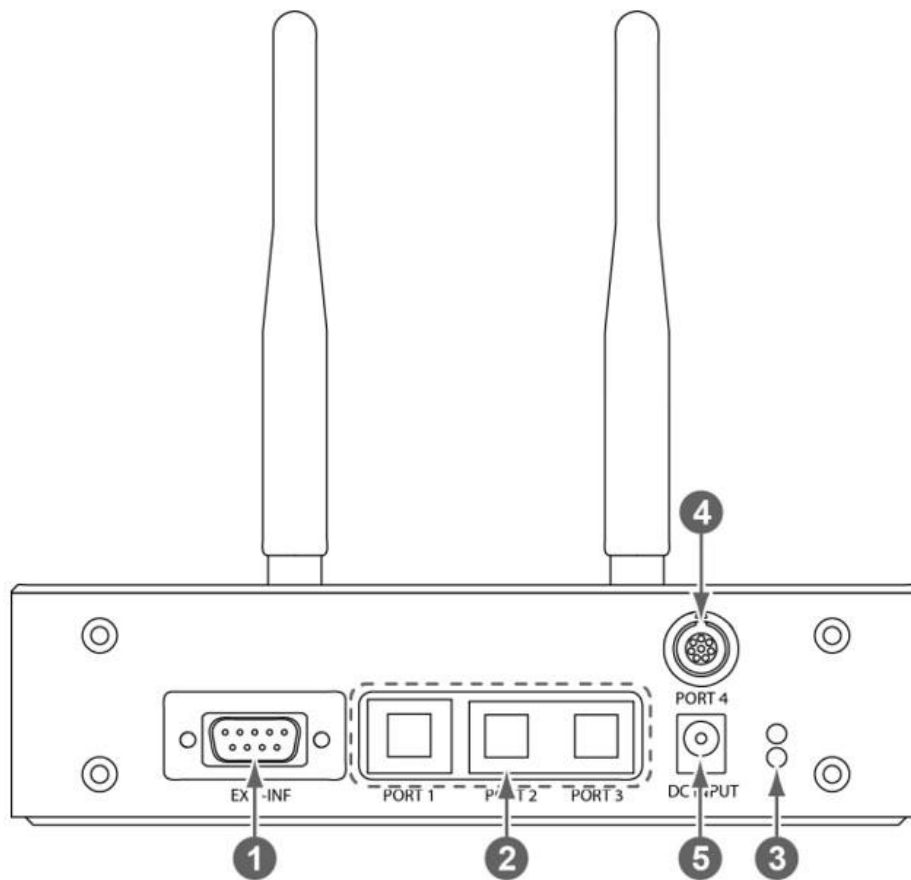No.	Name	Description
1	EXT_INF	X-ray generator interface connector (D-SUB 15pin, Female)
2	LAN port (Port 1, 2, 3)	Gigabit Ethernet port (1000BASE-T) <ul style="list-style-type: none"> • Port 1: Communication between the workstation and SCU mini. • Port 2, 3: For configuring multiple detectors.
3	PoE status lamp	Indicats the status of PoE port. <ul style="list-style-type: none"> • Green: 1 Gbps • Orange: 100 Mbps
4	PoE port	PoE (Power over Ethernet) port (1000BASE-T) <ul style="list-style-type: none"> • Communication between the detector and SCU mini. • Supplies power to the detector.
5	DC power input port	DC +24V <ul style="list-style-type: none"> • Supplies power to SCU mini.

2.8 Others

2.8.1 X-ray Generator (Recommended Exposure Condition)

Item	Recommended condition
X-ray energy range	<ul style="list-style-type: none"> • 40kVp ~ 150kVp
Reliability (Lifetime Dose)	<ul style="list-style-type: none"> • 74Gy or more (35 uGy x 365days x 24hrs. x 60min. x 60sec. / 15sec.)

2.8.2 Recommended Specifications of Workstation (PC)

Item	Recommended specification
Operating System	<ul style="list-style-type: none"> • VIVIX Setup <ul style="list-style-type: none"> ▫ Microsoft Windows 7 SP1 (32bit / 64bit) ▫ Microsoft Windows 8 / 8.1 / 8.2 SP1 (32bit / 64bit) ▫ Microsoft Windows 10 (32bit / 64bit)
CPU	Intel® Core™ i7-8700 CPU (or compatible CPU)
Memory	16GB or higher
Hard Disk	1TB or higher
LAN Card (only for detector communication)	<ul style="list-style-type: none"> • 1000 BASE-T Gigabit Ethernet card • Gigabit (for detector only) • Intel® PRO 1000 Series (Gigabit LAN Card for network interface) • Min.requiremetns <ul style="list-style-type: none"> ▫ Speed: 1Gbps or more ▫ Jumbo Frames: 9K ▫ Receive Descriptors: 2K (1024 or higher)
Monitor	<ul style="list-style-type: none"> • Minimum: 1280 × 768 or higher • Recommended: 1680 x 1050 or higher
CD-ROM	CD or DVD R/W

2.8.3 Recommended Specifications of Grid

Item	Recommended specification
SID	100cm / 130cm / 150cm / 180cm
Ratio	8.1 / 10:1 / 12:1
Frequency	215 line/inch

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3. System Configuration

This section gives information about the various connection / configuration ways among the detector, SCU, workstation and X-ray generator.

You can figure out the most suitable way of system configuration through this instruction.

Detector Connection Methods

Diagram

System Configuration


Generator Interface

Configuring DR Trigger Interface

3.1 Detector Connection Methods

You can organize wired or wireless connection between the **VIVIX-S VW** detector and **SCU**. You can also choose other connection ways suitable for the use environment in case of need.

3.1.1 Wireless Connection

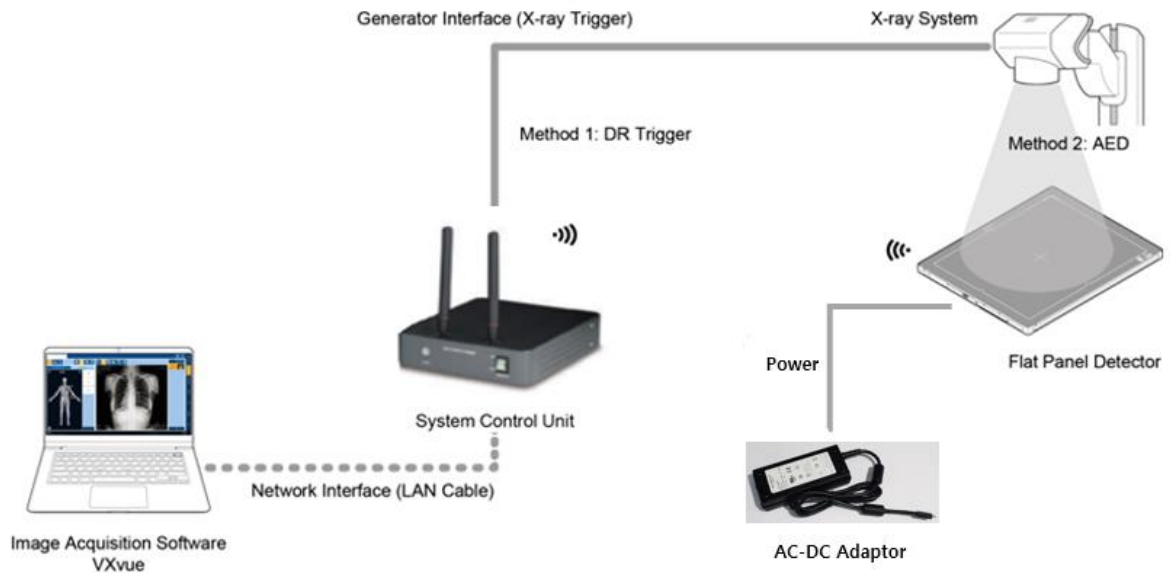- The VIVIX-S VW Detector and PC (Workstation) connect wirelessly to transfer video and data to a PC.
- Using the SCU mini (FXRS-04A), you can use the DR Trigger as well as AED (Auto Exposure Detection) with the generator interface.
- Instead of using the SCU mini (FXRS-04A) as an AP, you can use an external AP. In this case, DR Trigger cannot be used.
- If your PC (Workstation) supports Wi-Fi, you can connect directly to your PC using the detector as an AP without a separate AP (SCU, external AP). In this case, DR Trigger cannot be used.
- In a wireless connection, the detector is powered by a battery basically. However, if you connect an AC-DC adaptor, the detector receives power from the adaptor and the battery charges quickly.


- We recommend you use the wireless connection method for the following cases.
 - When using the unfixed detector.
 - When the wired communication is in trouble by the interruption of cables.
- The wireless connection is more freely to operate devices than the wired connection.
- Use the wireless connection with a laptop computer to enhance mobility.
- Nous vous recommandons d'utiliser la méthode de connexion sans fil dans les cas suivants.
 - Lors de l'utilisation du détecteur non fixé.
 - Lorsque la communication filaire est en difficulté par l'interruption des câbles.
- La connexion sans fil permet plus facilement d'utiliser des appareils que la connexion filaire.
- Utilisez la connexion sans fil avec un ordinateur portable pour améliorer la mobilité

- Be sure to operate the wireless communication between the detector and SCU within a maximum of **8m**.
- Use of multiple WLAN devices within the same frequency band may interfere with each wireless communication and cause a decline in transmission speed.
- Do not cover or block the wireless module of the detector. Also, do not put any shielding materials between the detector and SCU. Otherwise, the transmission speed or operable distance may be reduced.
- Veuillez à utiliser la communication sans fil entre le détecteur et le SCU dans un délai maximal de 8 m.
- L'utilisation de plusieurs périphériques WLAN dans la même bande de fréquence peut interférer avec chaque communication sans fil et provoquer une diminution de la vitesse de transmission.
- Ne couvrez pas et ne bloquez pas le module sans fil du détecteur. De même, ne placez aucun matériau de protection entre le détecteur et le SCU. Sinon, la vitesse de transmission ou la distance de fonctionnement peut être réduite.

3.1.2 Wired Connection

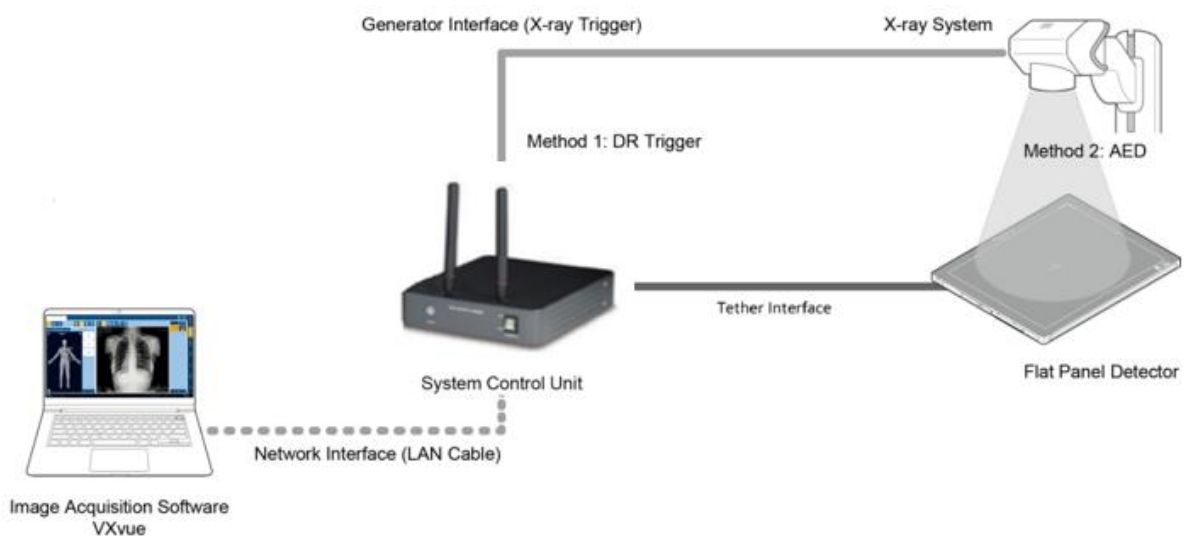
- The VIVIX-S VW detector is wired to the PC (Workstation) via the SCU and transmits image and data to the PC.
- Using the SCU mini (FXRS-04A), you can use the DR Trigger as well as AED (Auto Exposure Detection) with the generator interface.
- You can use SCU Lite (FXRP-02A) instead of the SCU mini (FXRS-04A). In this case, DR Trigger cannot be used.
- When the detector is wired, the battery charges while the detector is powered by the tether interface cable.

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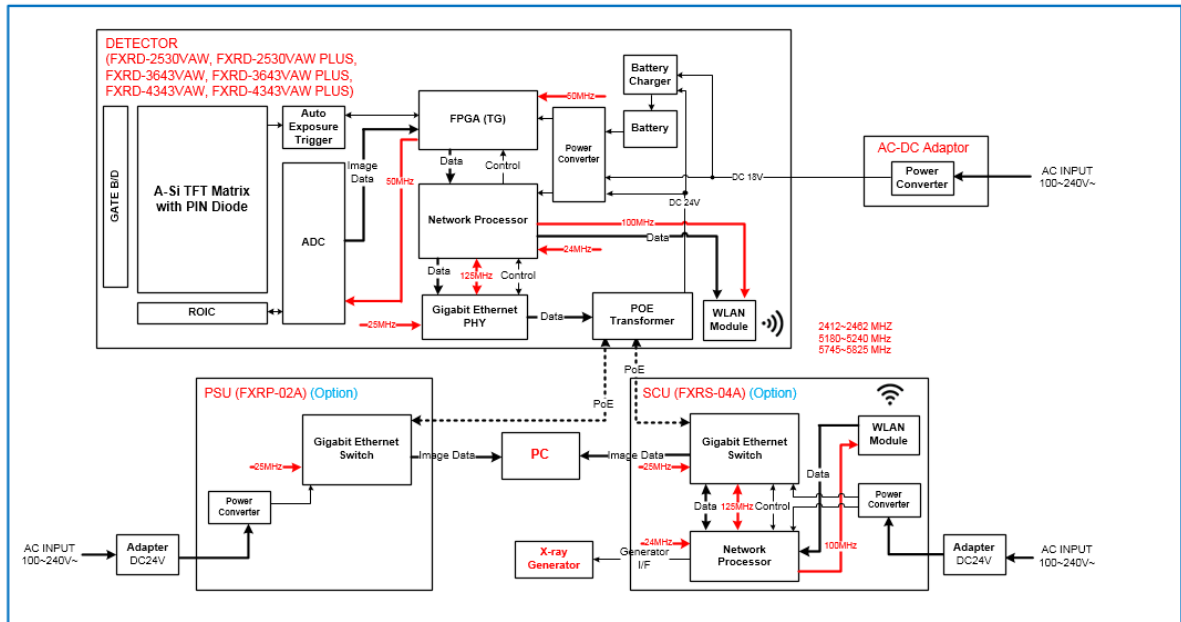
3.2.2 Wiring Diagram

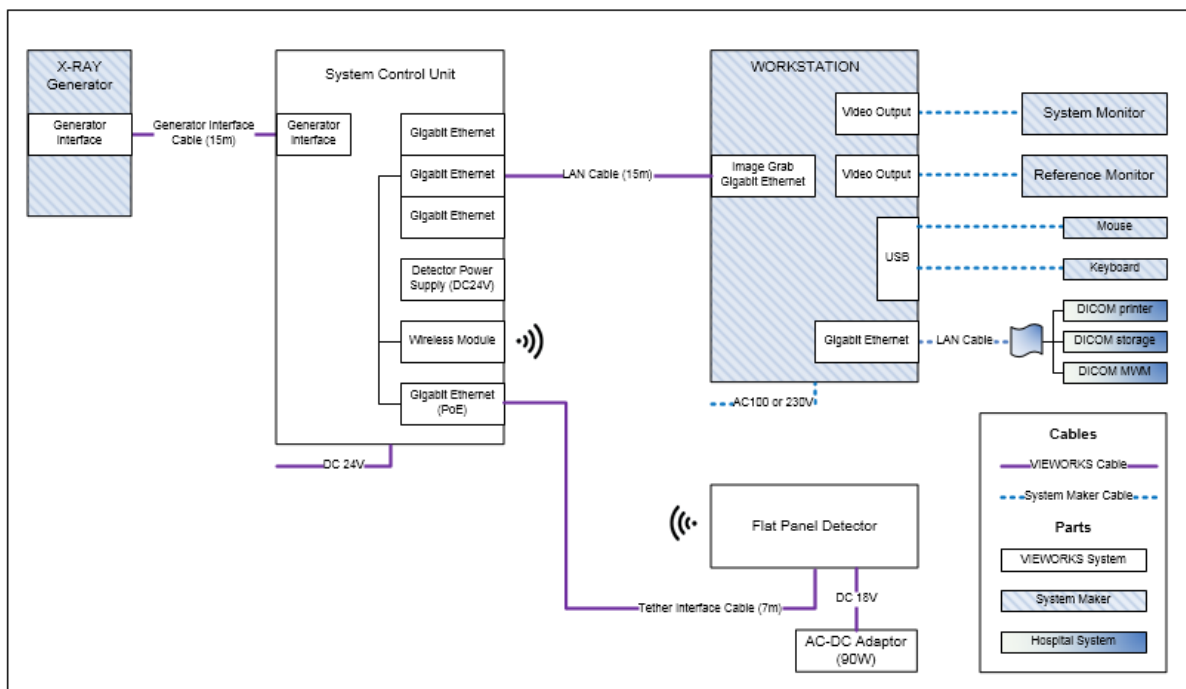
3.3 System Configuration

This section gives information about the configuration modes of **VIVIX-S VW** detectors, SCU, workstation and the generator.

3.3.1 SCU AP Mode

In SCU AP mode, the detector and the SCU mini (FXRS-04A) are connected wirelessly, and the SCU mini and PC are connected by wired.

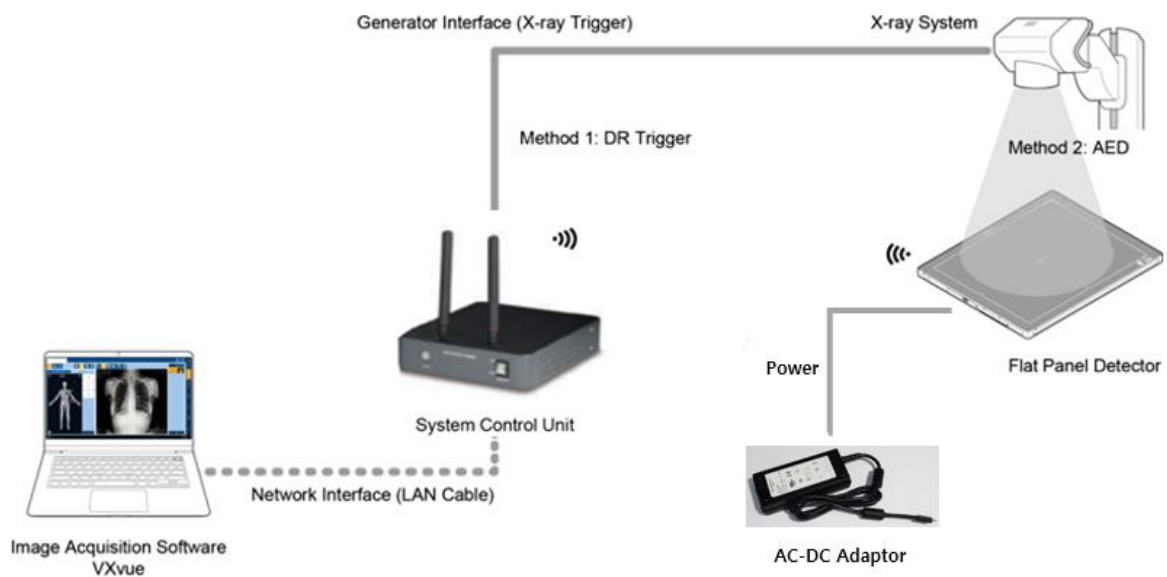SCU mini & Detector

- SCU and the detector are connected wirelessly. In this case, SCU operates as AP (Access Point) and detector as STATION.
- Basically, the detector is powered by a battery. However, when you connect the AC-DC adapter to the detector, it receives power from the adapter and the battery charges quickly.

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SCU mini & PC (Workstation)

- SCU and PC (Workstation) are connected with the LAN cable.

SCU mini & X-ray Generator


- SCU and the X-ray generator can be connected with a generator interface cable. (DR trigger mode)

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3.3.2 Tether Interface Mode

Tether interface mode is the wired connection method that connecting the detector with SCU as the center.

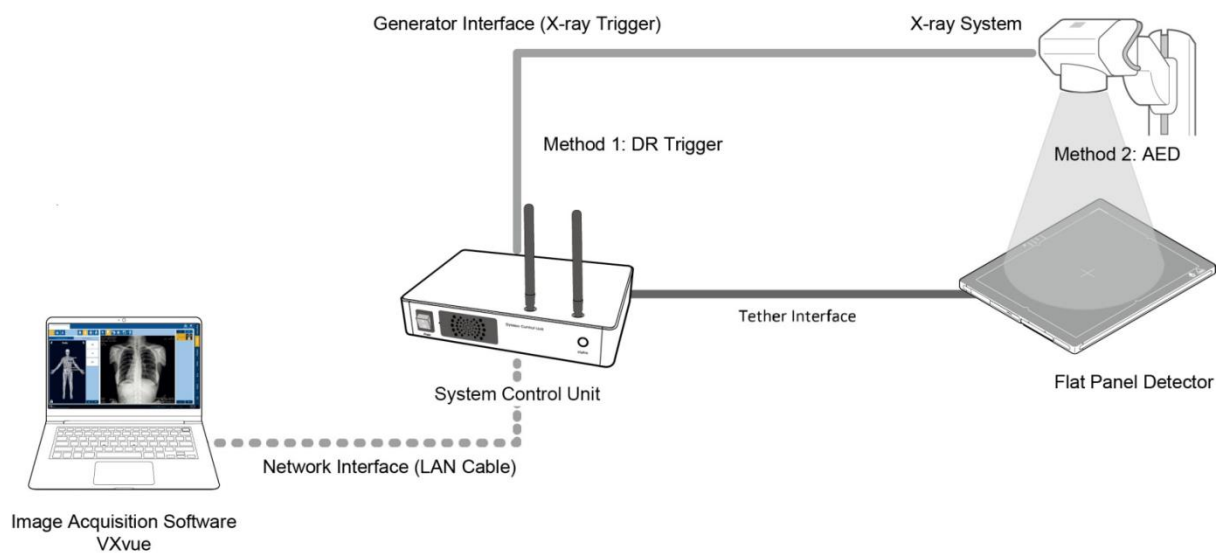SCU & Detector


- SCU and the **VIVIX-S VW** detectors are connected with the tether interface cable.

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3.3.3 External AP Mode

External AP mode is the wireless connection method that adding the external AP device instead of SCU.

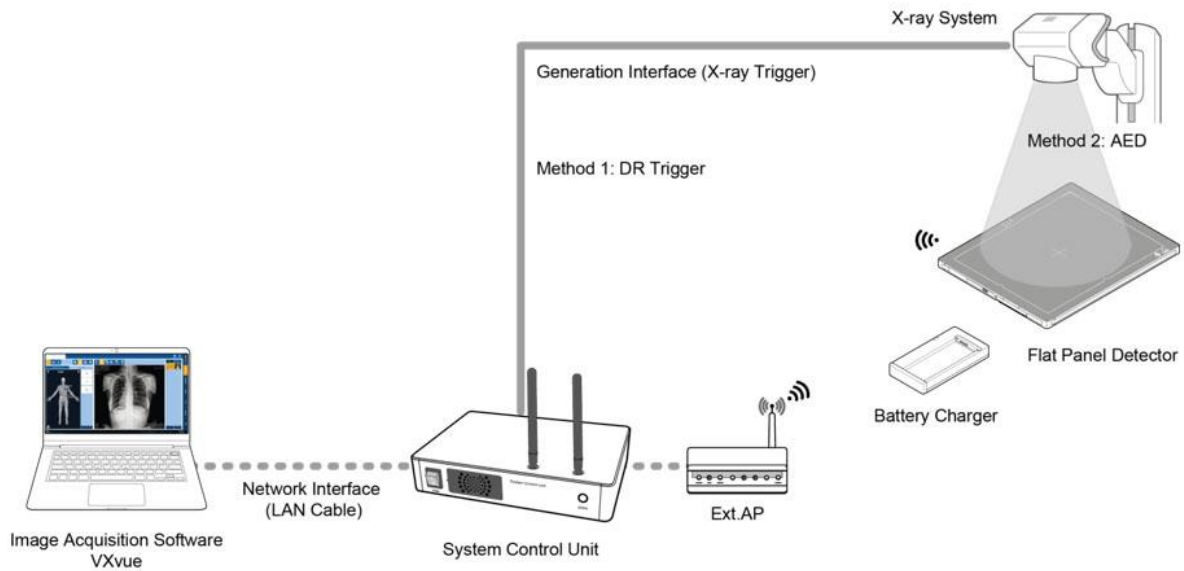External AP & SCU

- The external AP and SCU are connected with a LAN cable.

External AP & Detector

- The external AP and **VIVIX-S VW** detectors are connected wirelessly. In this case, the external AP and the detector operate as **AP** (Access Point) and **STATION** respectively.

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
3.3.4 AP Switching

Using this function, you can conveniently change the external AP when the detector is STATION.

- Ex.) If you want to use one detector in both ROOM #1 and ROOM #2.

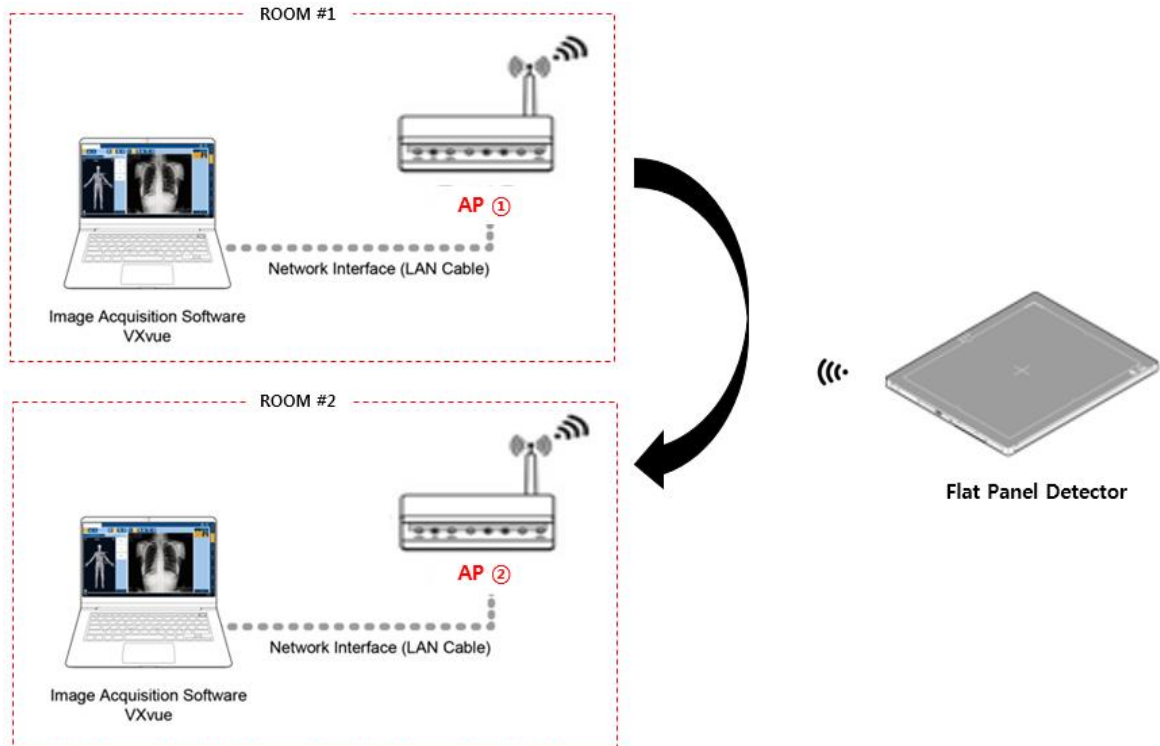- 1 When the detector is in STATION, you can change the external AP by pressing the AP button for about 3 seconds.
 - That is, the detector changes to the same SSID and key as those of the external AP to be changed.
 - AP switching disconnects the existing AP and automatically changes to the AP with the largest signal strength among the list of registered APs. Therefore, the SSID and KEY of the AP to be used must be registered before switching the AP.
- 2 During AP switching, the OLED display shows Changing Settings, and you can check the changed settings after the switching.

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3.3.5 Sync Wireless Settings

You can easily switch from tether interface mode to SCU AP mode.

- 1 In the tether interface mode, press the AP button for about 3 seconds to change the detector to the same SSID and KEY as the wired SCU.
- 2 During synchronization, the OLED display shows 'Changing Settings'. After switching, you can check the changed settings.

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

3.4 Generator Interface

VIVIX-S VW detectors provide **AED** interface and **DR Trigger** interface as a generator interface method to acquire images by detecting X-ray.

Mode	Description
AED	The detector detects X-ray exposure from the generator automatically and then performs image acquisition without any cable connection.
DR Trigger	<p>The detector and generator receive and send their signal to each other for image acquisition.</p> <ul style="list-style-type: none"> • SCU mini (FXRS-04A) and X-ray generator should be connected with a generator interface cable. • DR Trigger interface is not supported in case you connect SCU Lite (FXRP-02A) with the X-ray generator.

3.4.1 AED (Auto Exposure Detection) Interface

If VIVIX-S VW detectors are used as the **AED** interface, you can acquire images without connecting the generator to the detector with a generator interface cable.

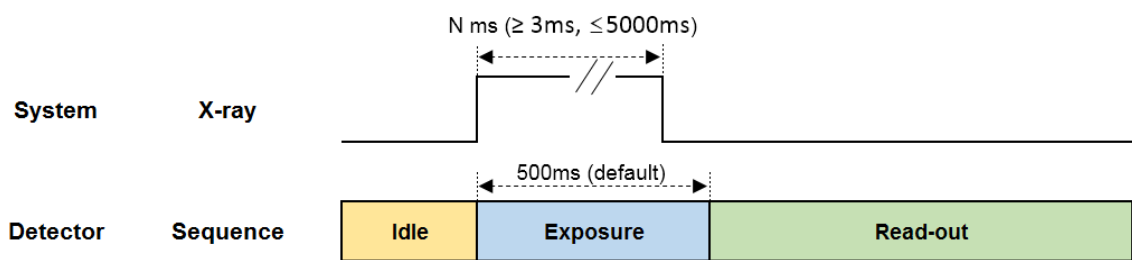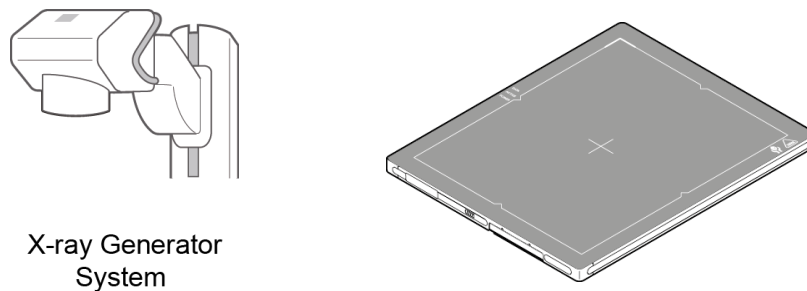Signal processing steps

- 1 X-ray generator exposes X-rays.
- 2 The detector acquires the image and sends the image data.

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

3.4.2 DR Trigger Interface

Connect the detector and X-ray generator with a generator interface cable, and then acquire images by sending and receiving their signal.

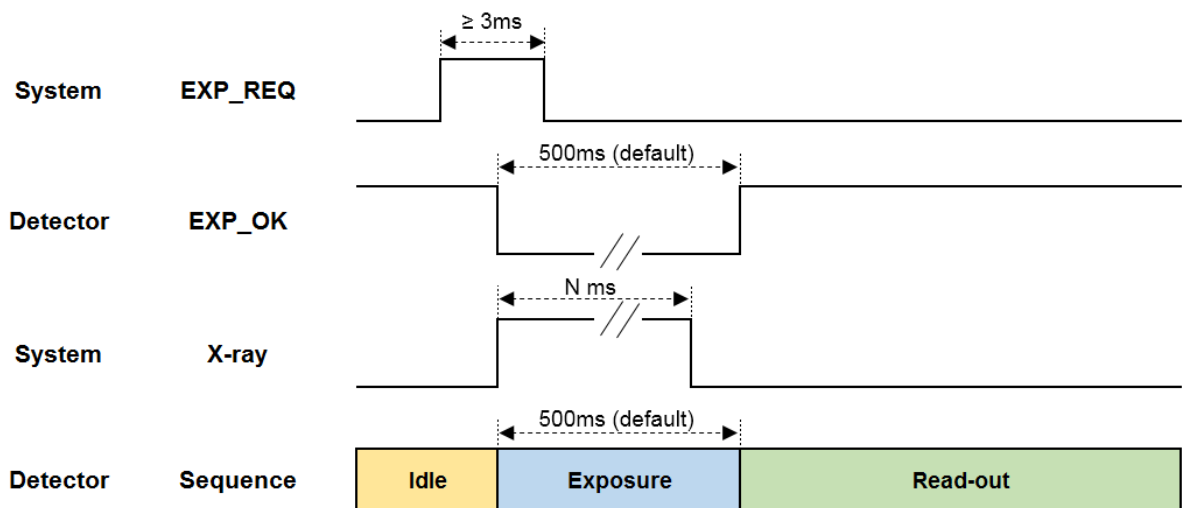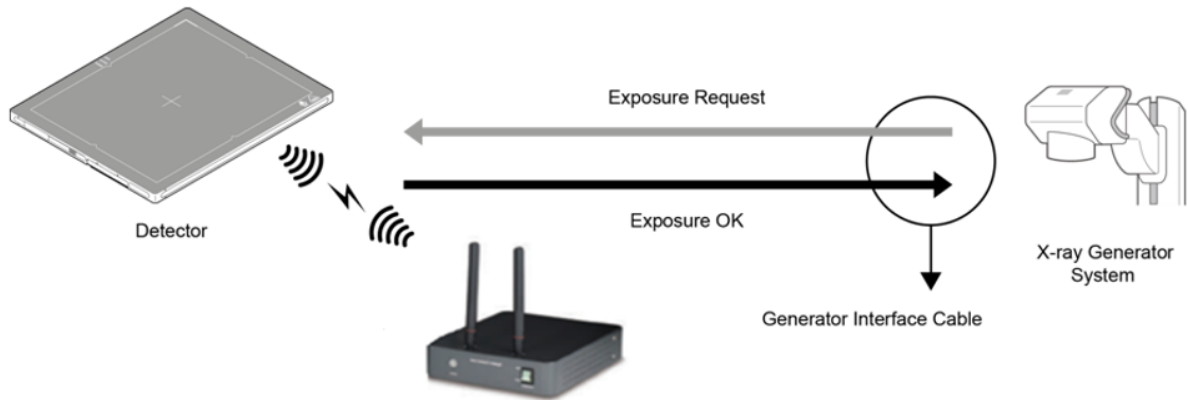Signal Processing Steps

- 1 The detector receives **EXP_REQ** signal from the X-ray generator.
- 2 After the detector completes to prepare image acquisition, it sends **EXP_OK** signal to the X-ray generator.
- 3 The X-ray generator confirms **EXP_OK** signal and generates X-rays.
- 4 The detector acquires images and then transmits the image data.

- EXP_REQ is a shooting request signal sent by the X-ray generator to the detector.
- EXP_OK is a shooting ready signal sent by the detector to the X-ray generator.
- The default X-ray exposure time set in the detector is 500ms and can be changed by the user's settings.
- EXP_REQ est un signal de demande de prise de vue envoyé par le générateur de rayons X au détecteur.

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3.5 Configuring DR Trigger Interface

To configure the **DR Trigger** interface, SCU mini (FXRS-04A) and the generator should be connected with a generator interface cable. Connect the one end of generator interface cable to **EXT_INF** port of SCU, and then connect the other end to the generator.

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
3.5.1 Trigger Interface Way

The generator interface cable is connected differently depending on the trigger interface ways.

Trigger Interface	Description
Packet Trigger	Handles exposure signal by exchanging packet.
Line Trigger	Handles exposure signal by exchanging electrical signal.

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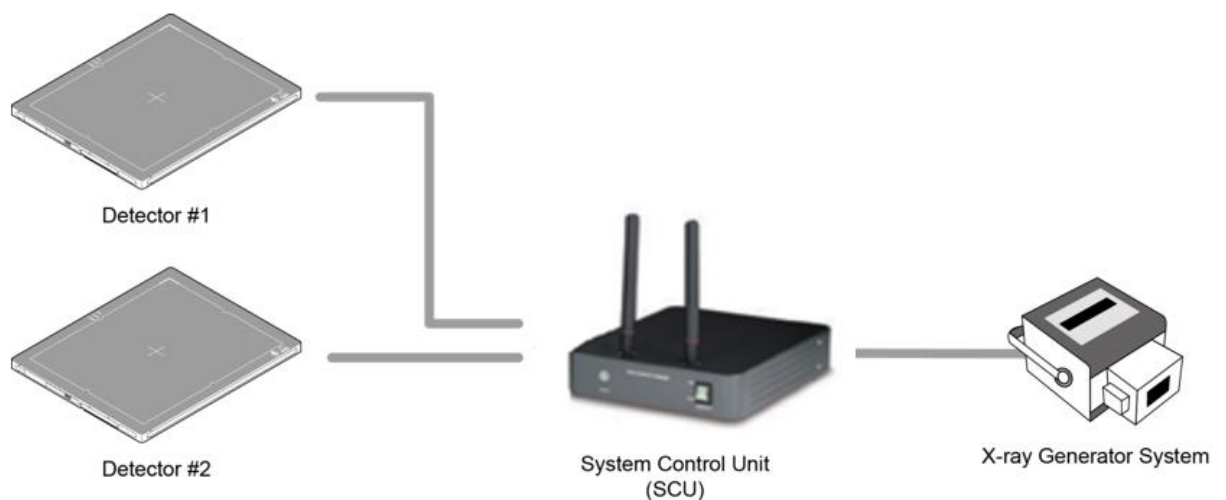
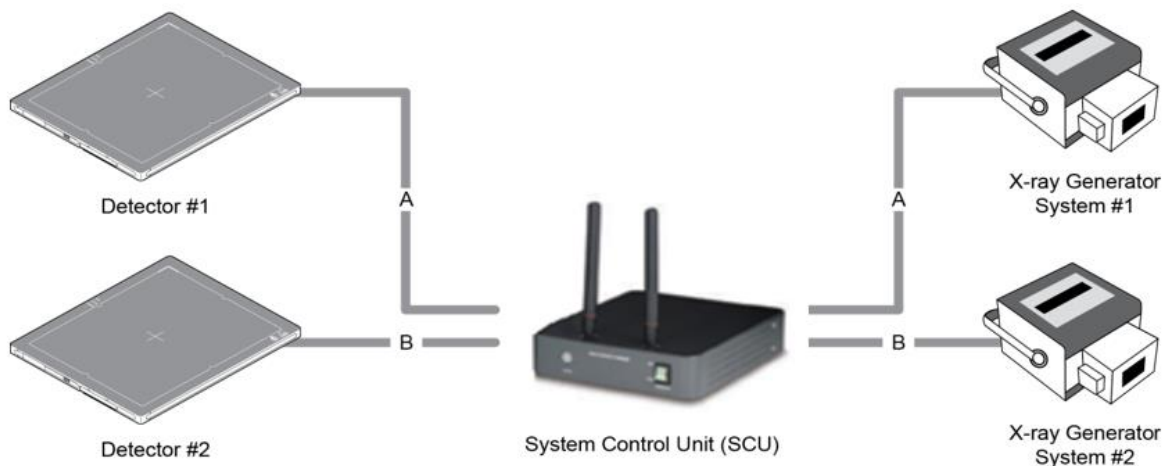
3.5.4 EXT_INF Port Pin Map

EXT_INF port is D-SUB 15 pin female connectors of SCU mini (FXRS-04A), and used for the generator interface.

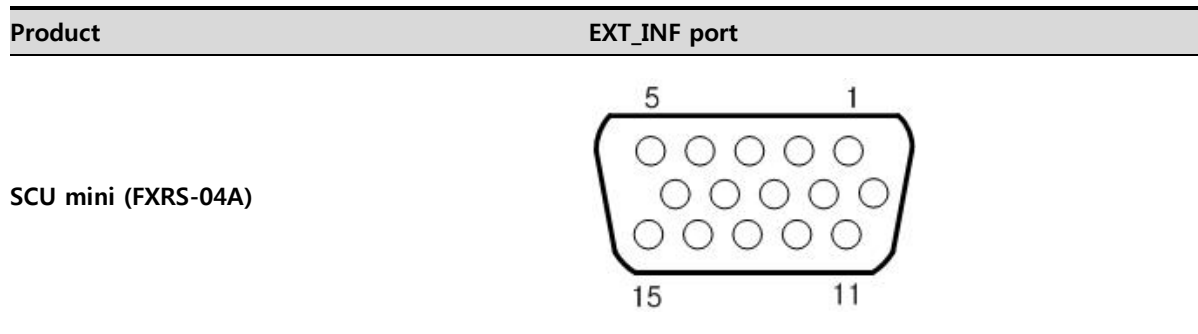EXT_INF1 port pin map definition (1 ~ 15)

No.	Signal name	I/O	Type	Color	Ch.	Description
1	EXP_REQ+_A	Input	Contact	Red	A	Receives EXP_REQ
2	EXP_REQ-_A	Input	Contact	Black	A	Returns signal from EXP_REQ+_A
3	EXP_REQ_TTL_A	Input	TTL	Orange	A	Receives EXP_REQ
4	EXP_REQ_GND_A	Input	TTL	Gray	A	Returns signal from EXP_REQ_TTL_A
5	EXP_OK_POWER_A/B	Input	-	Yellow	A/B	Power of TTL signal
6	EXP_OK+_A	Output	-	Green	A	Sends EXP_OK
7	EXP_OK-_A	Output	-	Brown	A	Returns signal from EXP_OK+_A
8	EXP_OK+_B	Output	-	Blue	B	Sends EXP_OK
9	EXP_OK-_B	Output	-	Pink	B	Returns signal from EXP_OK+_B
10	Reserved	-	-	-	-	Reserved for test only.
11	EXP_REQ+_B	Input	Contact	White	B	Receives EXP_REQ
12	EXP_REQ-_B	Input	Contact	Purple	B	Returns signal from EXP_REQ+_B
13	EXP_REQ_TTL_B	Input	TTL	White/Red	B	Receives EXP_REQ
14	EXP_REQ_GND_B	Input	TTL	White/Black	B	Returns signal from EXP_REQ_TTL_B
15	Reserved	-	-	-	-	Reserved pin for test only.

- There are **Contact** type (OPEN/ CLOSE) and **TTL** type (ON/ OFF) on the signal in/output pin.
- TTL type information
 - ON : VCC
 - OFF : GND
 - Current: 5 mA ~ 10 mA
 - Voltage: 12 V ~ 24 V
- Il existe un type de contact (OPEN / CLOSE) et un type TTL (ON / OFF) sur la broche signal entrée / sortie.


- Informations de type TTL
 - ON: VCC
 - OFF: GND
 - Courant: 5 mA ~ 10
 - Tension: 12 V ~ 24 V

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
3.5.5 Input / Output Circuit

Exposure Request Input Circuit (EXP_REQ)

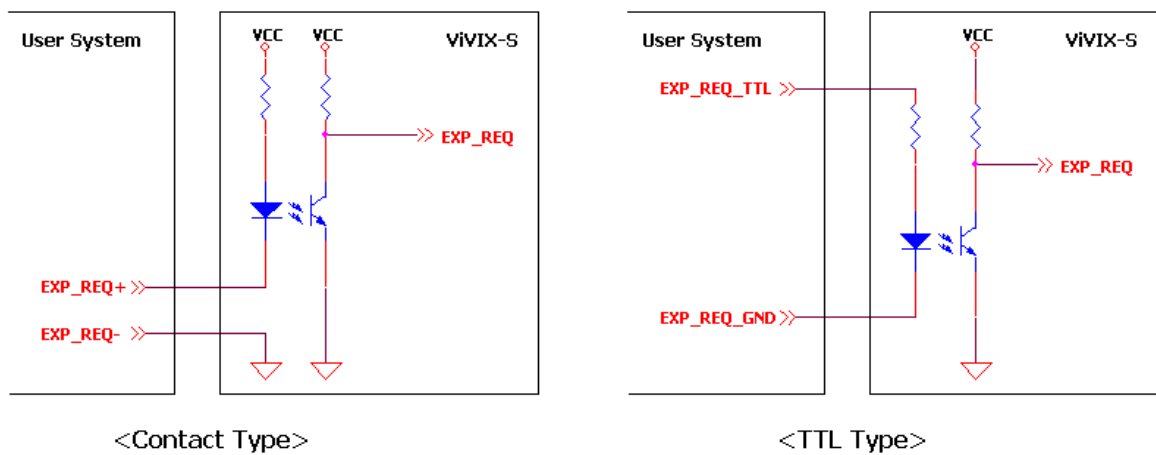Exposure Respond Output Circuit (EXP_OK)

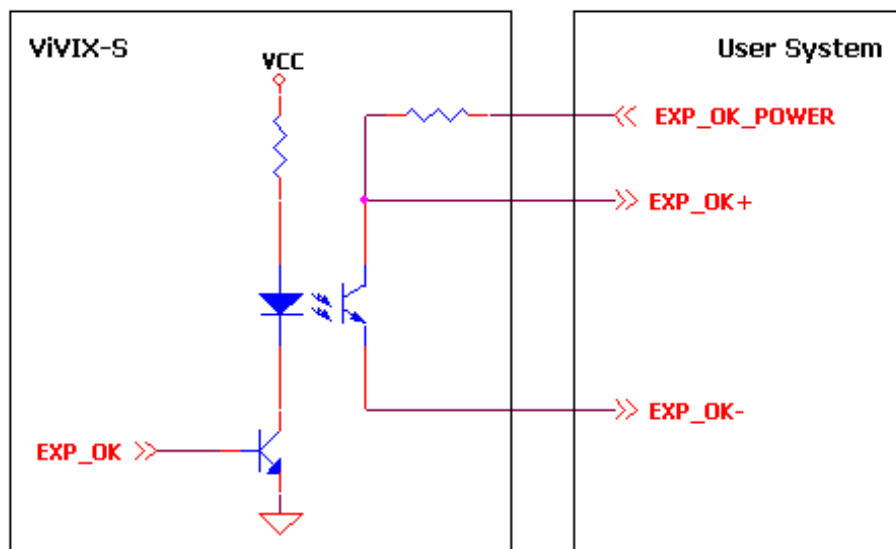4. Settings

This section gives information about the way of installation, setting, and diagnosis to use a detector with SCU.

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

4.1.1 Connecting SCU mini (FXRS-04A)

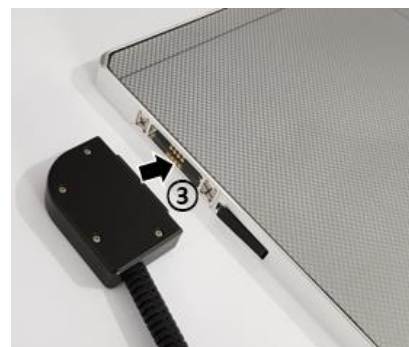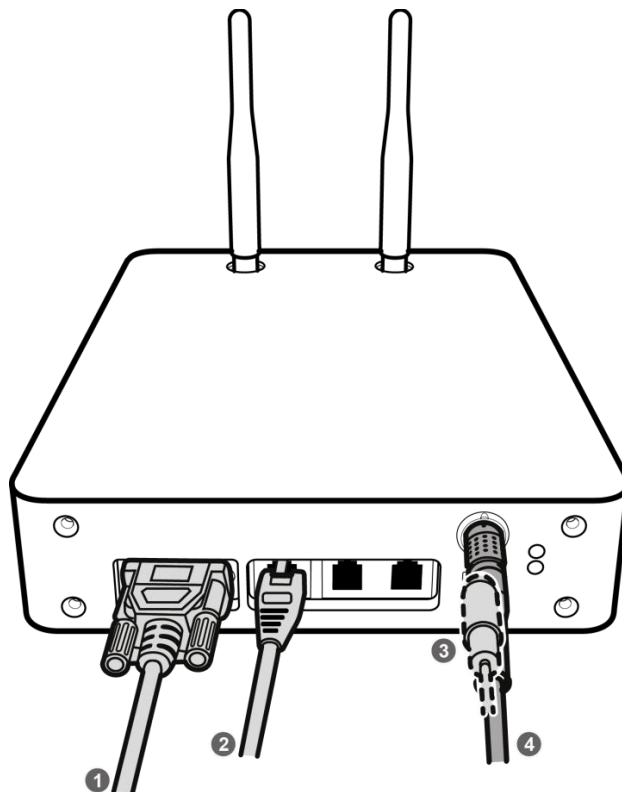- 1 If you use the **DR Trigger** interface, connect the one end of generator interface cable to the **EXT_INF** port of SCU, and the other to the X-ray generator.

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2 Connect one end of the LAN cable to one of the LAN ports of SCU mini, and the other to the LAN Card Connector of workstation assigned for data transfer. Stand the antenna of SCU mini upright.


3 To transmit image data using a tether interface, connect the one end of the tether interface cable to the port of SCU mini and the other to the detector.

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4 To supply power, connect the DC power cable to the DC power input port of SCU mini.

4.1.2 Booting SCU mini (FXRS-04A)



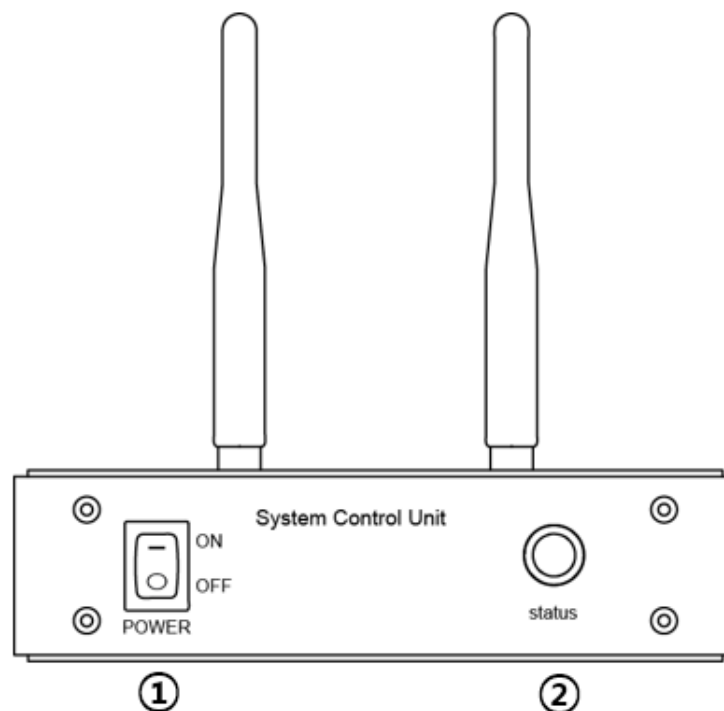
- 1 Turn on the power switch on the front of the SCU mini. (ON).
- 2 Make sure the status LED is lit green.

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4.1.3 Connecting SCU Lite (FXRP-02A)

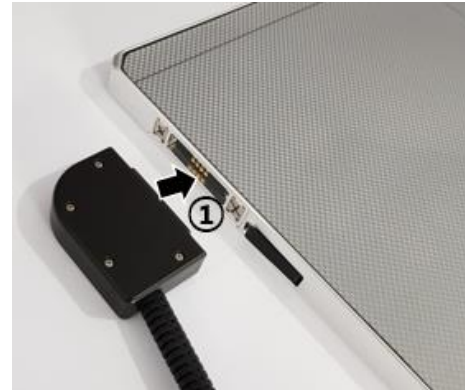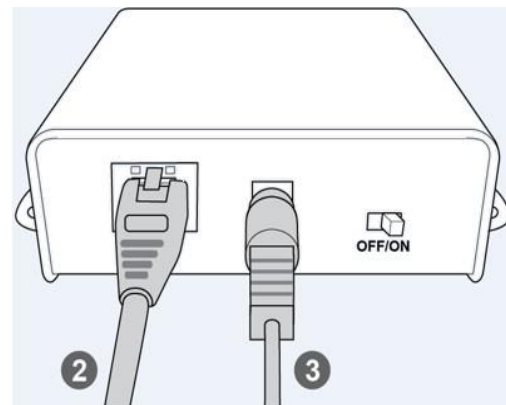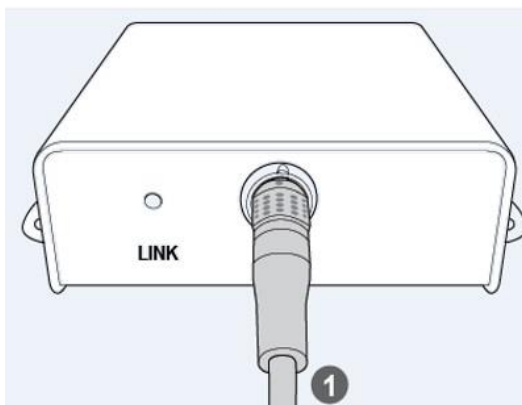
- 1 Connect one end of the tether interface cable to the detector and connect the other end to **PoE port** of SCU Lite.
 - The tether interface cable can be connected to the detector in any direction.
- 2 Connect one end of the LAN cable to the LAN port on the SCU Lite and the other end to the LAN port for data transfer on the workstation.
- 3 Connect the AC-DC power adaptor to the power input port of SCU Lite.

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4.1.4 Booting Up SCU Lite (FXRP-02A)

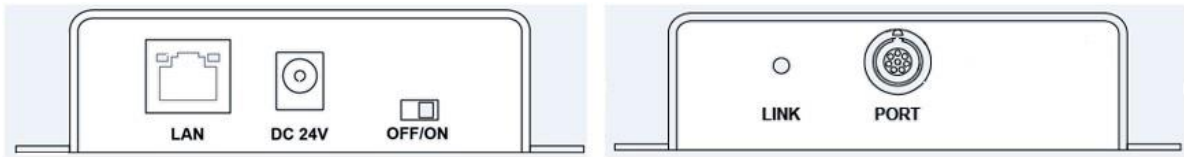
- 1 Turn on the power switch at rear side of SCU Lite (ON).
- 2 If you connect a detector to SCU Lite by wire, check if the LINK LED of SCU Lite lights up green.

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4.1.5 Booting up the Detector

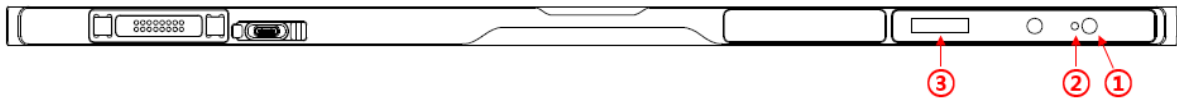- 1 Press the power button on the detector for about 1 second.
- 2 Check that the detector's power LED is lit green.
- 3 Check the battery status (remaining, charged) and communication connection status (wired / wireless, detector AP / STATION, wireless connection strength, etc.) on the detector's OLED display.

- If the power LED flashes green, power is applied, and the detector is booting up.
- When the power indicator LED is lit green, the detector has finished booting.
- If you press and hold the power button for about 3 seconds, the power indicator LED and the OLED display turn off, then the system is powered off. However, the battery pack can be charged when the SCU is powered via a tether interface cable or an AC-DC adapter is attached.
- Si le voyant d'alimentation clignote en vert, l'alimentation est appliquée et le détecteur démarre.
- Lorsque le voyant d'alimentation est allumé en vert, le détecteur a fini de démarrer.
- Si vous maintenez le bouton d'alimentation enfoncé pendant environ 3 secondes, le voyant d'alimentation et l'écran OLED s'éteignent, puis le système est mis hors tension. Toutefois, la batterie peut être chargée lorsque le SCU est alimenté via un câble d'interface d'attache ou lorsqu'un adaptateur CA-CC est connecté.

4.2 Device Setting

4.2.1 Software Installation

- 1 After connecting all devices, prepare the following softwares to set, calibrate and operate the detector / SCU.

Software	Description
VIVIX Device Driver (VDD)	Image filter driver for acquiring images from a detector.
VIVIX Setup	A program for setting and managing the detector / SCU.

- 2 Install **VIVIX Device Driver** and **VIVIX Setup** in sequence.

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- 3 Configure environment for the workstation.

4.2.2 Setting Detector and SCU

- 1 After executing **VIVIX Setup**, access to the detector and SCU to set each device properly.
- 2 Perform detector calibration to acquire images suitable for the installation environment.
- 3 Take radiographic images to check if the shooting is conducted normally.

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4.3 Diagnosis of Devices

Execute the **VIVIX Setup** program to check if there is any problem to operate the detector / SCU after installing and setting devices.

Diagnosis Items


Items	Description
Image	Diagnoses the acquired images.
Battery Pack	Diagnoses the condition of a battery pack.
Wired/Wireless Communication Status	Diagnoses the status of wired/wireless communication.
Wired/Wireless Communication Speed	Diagnoses the speed of wired/wireless communication.
Self-diagnosis	Diagnoses defects of a detector by self-diagnosis.

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4.3.1 Diagnosing Images

- 1 Execute **VIVIX Setup** and move to the **Image** dialog.
- 2 Take an image and check if it has any problem.
- 3 Take a dark image and check if it has any problem.
- 4 Check the effective area and whole area of the image.

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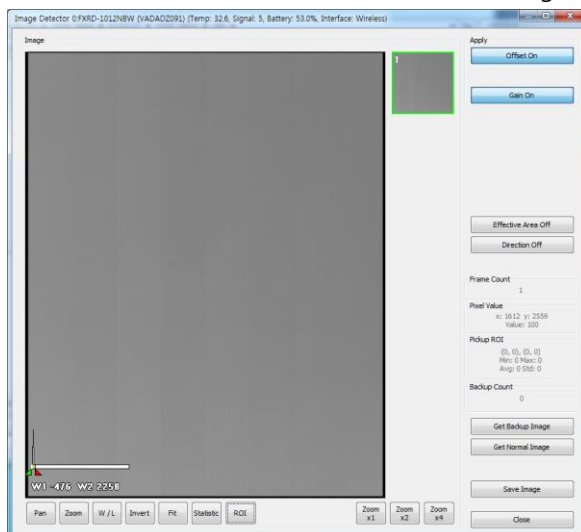4.3.2 Diagnosing Battery Pack

Check from OLED display of the detector

- 1 Check the mounting status / charging status of a battery, and battery levels on the OLED display.

Icon	Description
	No battery pack
	Charging (Capacity changes depending on remaining charge.)
	Charging complete or discharged. (Capacity changes depending on remaining charge.)

Check from VIVIX Setup

- 1 Execute **VIVIX Setup** and go to the **Information** tab of the **Diagnosis** dialog.
- 2 Check the mounting status, charging status and remains of a battery.

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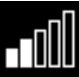
4.3.3 Diagnosing Communication Status

Checking on the OLED display of the detector

- 1 Check the detector communication status on the initial screen of OLED display.

Category	Item	Icon	Description
Wireless	AP		Operates in AP Mode
	STATION		Identifier of external AP
			Level 5 (Link Quality: 66~70, Very Good)
			Level 4 (Link Quality: 56~65, Good)
			Level 3 (Link Quality: 41~55, Normal)
			Level 2 (Link Quality: 31~40, Bad)
			Level 1 (Link Quality: 1~30, Very Bad)
			Level 0 (Link Quality: 0, Unknown)
Wired	Tether Interface		Link speed: 1Gbps


			Link speed: 100Mbps
Operating Status	Changing Settings		Changing communication settings

2 On the initial display of the OLED display, pressing the power button or AP button for about 1 second will switch the OLED display to the second screen. On the second screen, check the SSID and KEY important for wireless connection.

Category	Item	Icon	Description
Communication Status	Identifier of wireless connection		Displays SSID of the detector's current state (AP or STATION)
			Displays KEY of the detector's current state (AP or STATION)


Wired Communication Diagnosis

- 1 Connect the detector with a cable.
- 2 Execute **VIVIX Setup** and go to the **Maintenance** tab of the **Maintenance** window.
- 3 Check the status of wireless communication from the **Network** item.

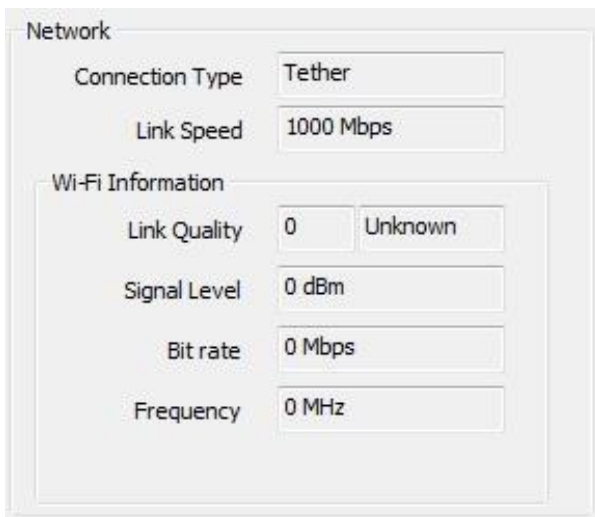Wireless Communication Diagnosis

- 1 Connect the detector wirelessly.
- 2 Execute **VIVIX Setup** and go to the **Maintenance** tab of the **Maintenance** window.
- 3 Check the status of wireless communication from the **Network** item.

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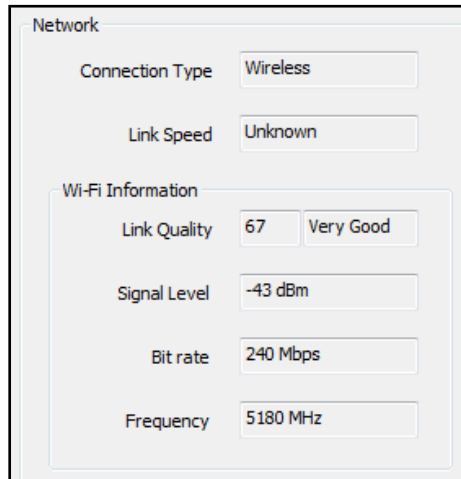
- In case of using the detector under wireless communication, be sure to check the communication status before starting to use the detector. If the status is bad, the speed of acquiring images will be very slow or you may fail to acquire images.
- Be sure to check the surrounding wireless communication to prevent communication interference. If wireless communication module in the detector has any problems, contact the service engineer in Vieworks.
- Si vous utilisez le détecteur en communication sans fil, assurez-vous de vérifier l'état de la communication avant de commencer à utiliser le détecteur. Si l'état est mauvais, la vitesse d'acquisition des images sera très lente ou vous risquez de ne pas acquérir les images.

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4.3.4 Diagnosing Communication Speed

- 1 Execute **VIVIX Setup** and go to the **Maintenance** tab of the **Maintenance** dialog.
- 2 Click on the **Download** button in **Image Transmission Time** and check image transmission speed of the detector.
- 3 Click on the **Start** button in **Throughput Measurement** to check the data traffic per setting time.

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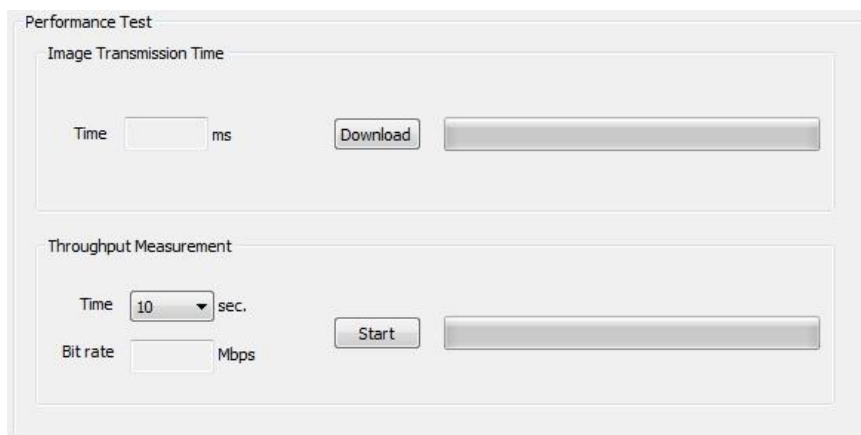
4.3.5 Self-Diagnosis

- 1 Execute **VIVIX Setup** and go to the **Self Diagnosis** tab in the **Diagnosis** dialog.
- 2 Check the desired items to diagnose from the **Test Case** list.
- 3 Click on the **Execute Test** button to perform self-diagnosis.
- 4 Check the status and result of diagnosis for each item in the **Result** dialog.

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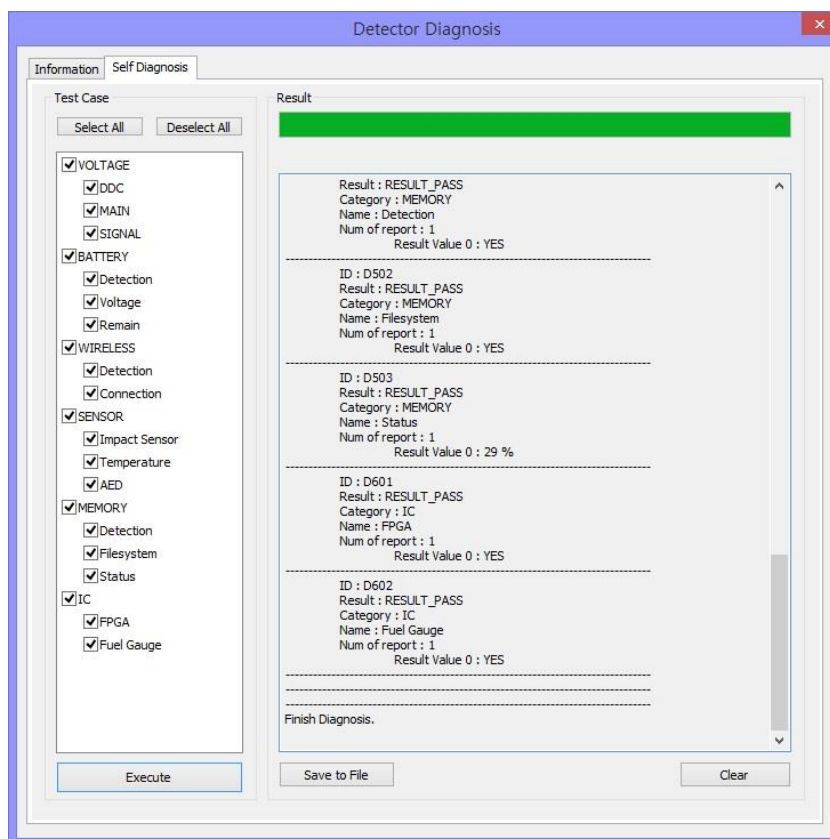Self-diagnosis items of detector and measures

Voltage

Item	Form	Expected problem	Measures
SOURCE	Decision	Problem with input power	Replace the battery pack.
			Replace the tether interface cable.
			Replace AC-DC adapter.
MAIN	Decision	Problem with system power	Contact a service engineer.

Battery

Item	Form	Expected problem	Measures
Detection	Decision	The battery is not attached.	Check if a battery is attached or not.
		A defective circuit is connected to a battery pack.	Contact a service engineer.
Voltage	Information	N/A	N/A
Remain	Information	N/A	N/A

Wired

Item	Form	Expected problem	Measures
Detection	Decision	Defective ethernet PHY	Contact a service engineer.
Connection	Decision	The cable is connected incorrectly.	Reconnect the tether interface cable.
		Defective cable	Replace the tether interface cable.
		Defective ethernet PHY	Contact a service engineer.

Wireless

Item	Form	Expected problem	Measures
Detection	Decision	Defective wireless module	Contact a service engineer.
Connection	Decision	Inconsistent environment of the wireless communication.	Check wireless communication environment.
		Defective wireless module	Contact a service engineer.

Sensor

Item	Form	Expected problem	Measures
Impact Sensor	Decision	Defective shock sensor	Contact a service engineer.
Temperature	Decision	Defective temperature sensor	Contact a service engineer.
AED	Decision	Defective AED sensor	Contact a service engineer.

Memory

Item	Form	Expected problem	Measures
Detection	Decision	Impossible to save backup images.	Contact a service engineer.
		Impossible to save logs.	Contact a service engineer.
		The calibration data is inapplicable.	Contact a service engineer.
File system	Decision	Impossible to save backup images.	Contact a service engineer.
		Impossible to save logs.	Contact a service engineer.
Status	Information	N/A	N/A

IC

Item	Form	Expected problem	Measures
FPGA	Decision	Impossible to take images from a detector.	Contact a service engineer.
Fuel Gauge	Decision	Impossible to check the remaining of a battery pack.	Contact a service engineer.

Self-diagnosis items of SCU and measures**Wired**

항목	Form	Expected problem	Measures
Detection	Decision	Defective ethernet PHY	Contact a service engineer.

Wireless

Item	Form	Expected problem	Measures
Detection	Decision	Defective wireless module	Contact a service engineer.

Memory

Item	Form	Expected problem	Measures
Detection	Decision	Impossible to save logs.	Contact a service engineer.
File system	Decision	Impossible to save logs.	Contact a service engineer.
Status	Information	N/A	N/A

IC

Item	Form	Expected problem	Measures
Switching IC	Decision	Impossible to connect the detector and PC.	Contact a service engineer.
Current Controller	Decision	Impossible to block overcurrent when using the wired mode.	Contact a service engineer.

5. Maintenance

This section gives information about maintenance of the product.

Product Initialization

Detector Power Save Function

How to Replace the Battery Pack

5.1 Product Initialization

If connection status of the system is not stable or setting value is not correct, you can initialize the products.

5.1.1 SCU Initialization

- 1 Click on the **Configuration** button of SCU after running the **VIVIX Setup** program.
- 2 Click on the **Factory Reset** button in the **Configuration** dialogue.
- 3 Wait for SCU to be initialized and rebooted automatically.
- 4 Make sure the SCU is initialized.

Default value of SCU initialization

Item	Default Value
Network	
IP Address	169.254.2.100
Subnet Mask	255.255.0.0
Gateway	169.254.2.100
AP	
AP On/Off	ON
Frequency	5GHz
Country	KR
Band	40MHz
Channel	+36
SSID	vivix
Key	1234567890
Trigger	
Method	Packet
Polarity	Auto

5.1.2 Detector Initialization

- 1 Click the **Maintenance** button after executing the **VIVIX Setup** program.
- 2 Go to **Maintenance** window > **Maintenance** tab and click the **Factory Reset** button.
- 3 Wait for the detector until it is rebooted automatically.
- 4 Click the **Configuration** button and check if detector initialization is complete.

Default value of detector initialization

Item	Default Value
Network	
IP Address	169.254.1.10
Subnet Mask	255.255.0.0
Gateway	169.254.2.100
Wireless Network	
SSID	vivix
Key	1234567890
AP Button	AP switching
Wireless Only	OFF
Power Mode	
Sleep	OFF
Shut down	OFF
Power Control	By Detector
AP	
Enable	OFF
Frequency	5GHz
Country	KR
Band	40MHz
Channel	+36
SSID	vivix_ap
Key	1234567890
Image Acquisition	
Exposure Mode	AED
Auto Offset Refresh Setting	
Use offset refresh	Uncheck
Time Interval	30Min
Temperature Interval	3°C
Number of shot	5
Image Transmission	
Timeout	20sec
Option	None

5.1.3 Wireless Initialization of Detector

- 1 Connect the detector wirelessly without connecting the tether interface cable.
- 2 Press and hold both the power button and AP button over 3 seconds.
- 3 Initialization will be conducted when the green LED is flickering. After that, the detector will reboot automatically.
- 4 When booting is complete, make sure the detector is initialized.

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5.2 Detector Power Save Function

You can extend the battery life by using the power saving function of the detector.

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Types of Power Save Mode

Item	Description
Normal	The detector can be operated and can take images at any time.
Sleep	The detector cannot be operated. User can take an image by disabling the sleep mode.
Shut Down	The detector has been turned off. User can take an image after the detector is rebooted.

Setting Items of Power Save Function

Item	Description
	Sets whether you use the sleep mode of the detector or not.
Sleep	<ul style="list-style-type: none"> • If the detector is not used for the set time, it switches to Sleep mode if the Sleep function is set. (Off / After 10min / After 15min / After 20min / After 25min / After 30min)
	Sets whether you use the shut down function of the detector or not.
Shut Down	<ul style="list-style-type: none"> • If the detector is not used during the set time, the detector is turned off. (Off / After 30min / After 60min / After 90min / After 120min)

- Refer to **VIVIX Setup Operation Manual** for the detailed information about the power save mode.

Entry Condition of Power Save Mode

Item	Description
Normal	-
Sleep	The detector turns to sleep mode if not used for the setting time.
Shut Down	The detector is turned off if not used for the setting time under the sleep mode. If the Sleep function is not used, the detector will shut down if it is not used for a set time.

Checking Indication in Power Save Mode

Item	Description
Normal	The power indicator LED is on and the OLED display shows the initial screen.
Sleep	<ul style="list-style-type: none"> The power indicator LED is on and the OLED display shows a sleep mode icon. You can check the sleep status in VIVIX Setup or VXvue. VIVIX SDK notices sleep state.
Shut Down	Power indicator LED and OLED display are off.

Shows Sleep Mode on OLED Display

Item	Icon	Description
Sleep		Enter sleep mode
Changing Settings		Canceling sleep mode

Disabling Power Save Function

Item	Description
Normal	-
Sleep	<ol style="list-style-type: none"> You can turn off the sleep state in VIVIX Setup or VXvue. You can call the function for turning off the sleep state from VIVIX SDK. You can cancel the sleep state by pressing the detector power button or the AP button.
Shut Down	Reboot the detector by pressing a power button on the detector.

5.3 How to Replace the Battery Pack

The battery pack is built into the detector. FXRD-3643VAW / FXRD-3643VAW PLUS / FXRD-4343VAW / FXRD-4343VAW PLUS can install up to two battery packs, and FXRD-2530VAW / FXRD-2530VAW PLUS can install only one battery pack.

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
- Battery packs lose capacity over time.
 - The battery pack can be replaced at the end of its life.
- The life of a battery pack is the number of charge and discharge cycles until it reaches 80% or less of its initial capacity (nominal capacity).
- The battery pack used for VIVIX-S VW detector has a service life of approximately 800 times. (One time standard: Fully charged and then fully discharged)
- Les batteries perdent de leur capacité avec le temps.
 - La batterie peut être remplacée en fin de vie.
- La durée de vie d'un bloc batterie est le nombre de cycles de charge et de décharge jusqu'à ce qu'il atteigne 80% ou moins de sa capacité initiale (capacité nominale).
- La batterie utilisée pour le détecteur VIVIX-S VW a une durée de vie d'environ 800 fois. (Norme ponctuelle: complètement chargée puis complètement déchargée.)


5.3.1 Replacing the Battery of FXRD-3643VAW / FXRD-3643VAW PLUS / FXRD-4343VAW / FXRD-4343VAW PLUS

How to remove the battery pack

- 1 Remove the deco sheet attached to the battery pack cover.
- 2 Loosen the bolts that secure the battery pack cover using a hexagon wrench, and then open the battery pack cover.

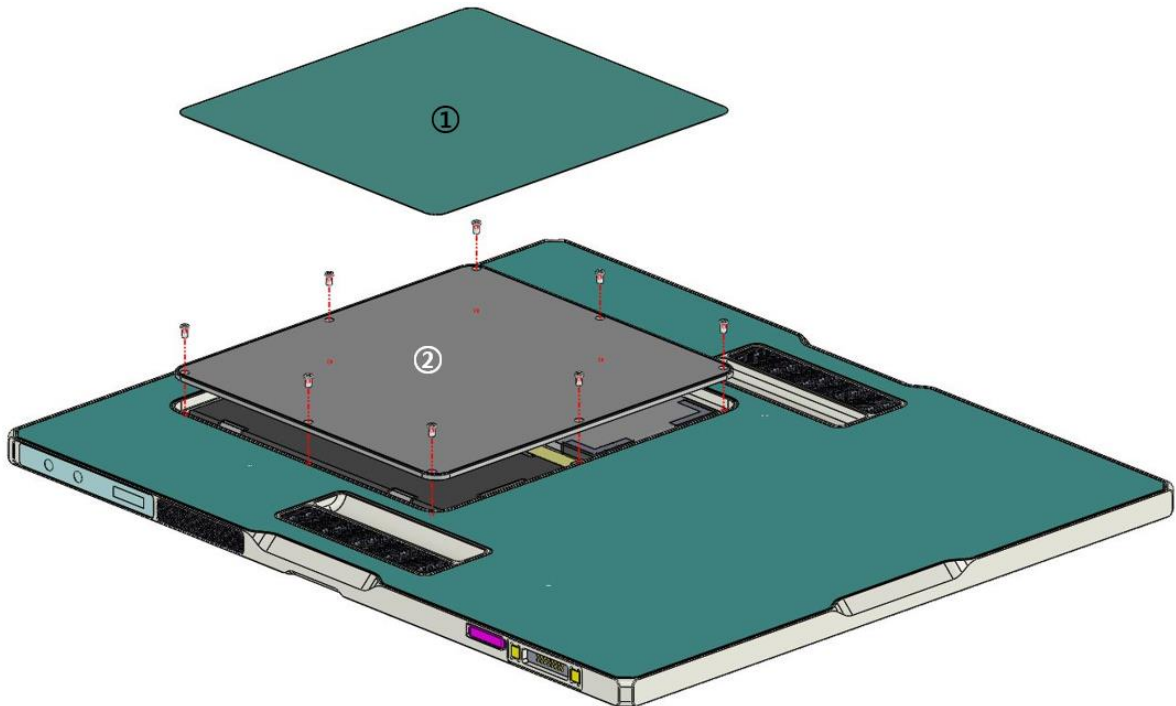- 3 Take the battery pack out of the detector.
- 4 Disconnect the cable connected to the battery pack.

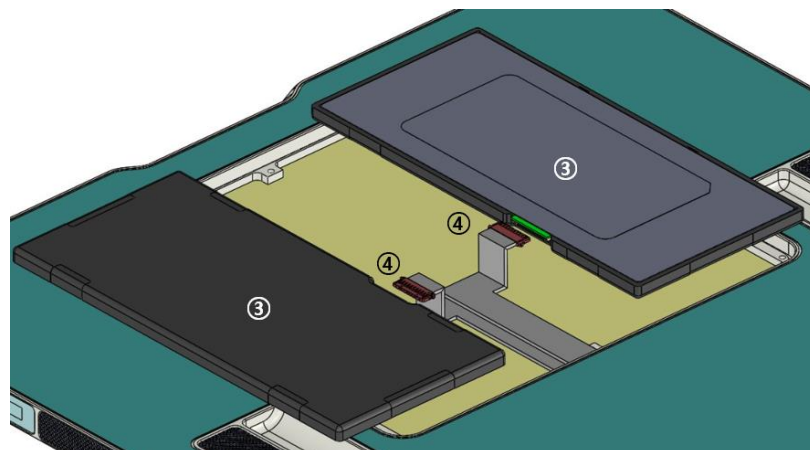- Before removing the battery pack from the detector, turn off the detector. If you press and hold the power button for about 3 seconds, the power indicator LED and the OLED display turn off, then the system is powered off.
- Disconnect wired power before removing the battery pack from the detector. If a tether interface cable or AC-DC adapter is attached, remove it from the detector.
- Avant de retirer la batterie du détecteur, éteignez le détecteur. Si vous maintenez le bouton d'alimentation enfoncé pendant environ 3 secondes, le voyant d'alimentation et l'écran OLED s'éteignent, puis le système est mis hors tension.
- Déconnectez l'alimentation câblée avant de retirer la batterie du détecteur. Si un câble d'interface d'attache ou un adaptateur CA-CC est connecté, retirez-le du détecteur.

How to attach the battery pack


- 1 Connect the battery pack to the cable.
- 2 Put the battery pack into the detector.
- 3 Close the battery pack cover and use a hex wrench to tighten the bolts that secure the battery pack.
- 4 Attach the new deco sheet to the battery pack cover.

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
5.3.2 Replacing the Battery of FXRD-2530VAW / FXRD-2530VAW PLUS

- 1 Remove the deco sheet attached to the battery pack cover.
- 2 Loosen the bolts that secure the battery pack cover using a hexagon wrench, and then open the battery pack cover.

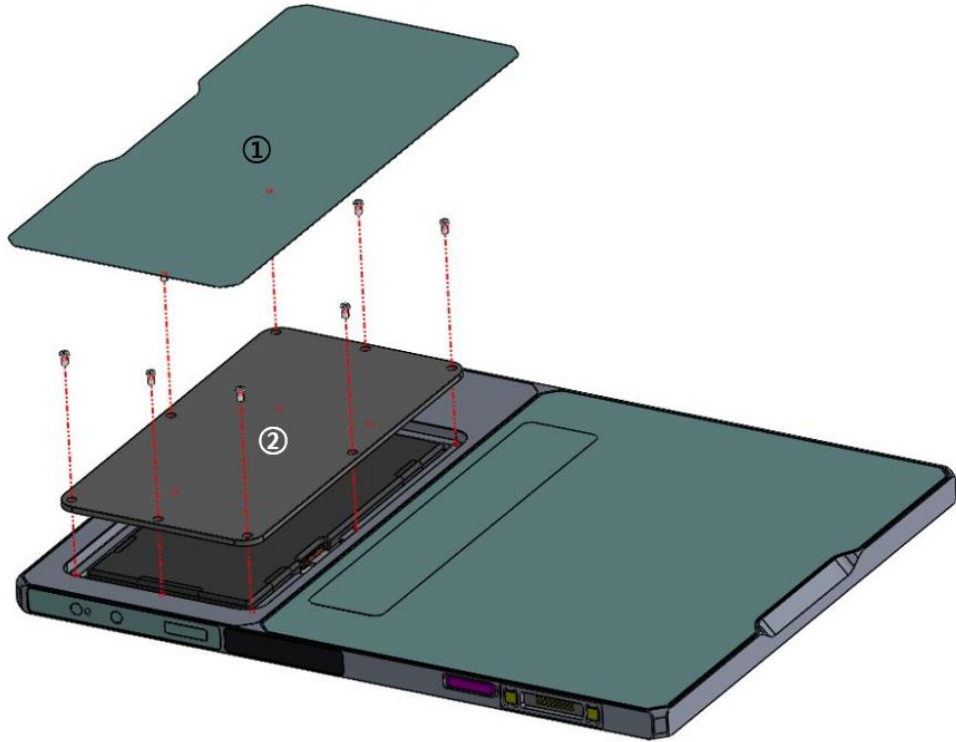- 3 Take the battery pack out of the detector.
- 4 Disconnect the cable connected to the battery pack.

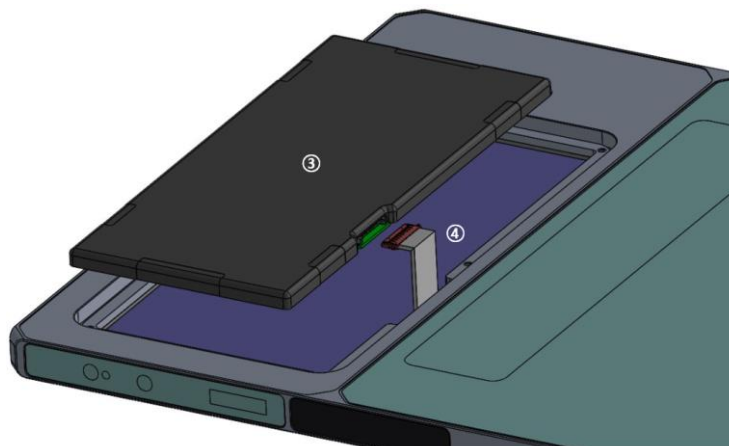How to attach the battery pack

- 1 Connect the battery pack to the cable.
- 2 Put the battery pack into the detector.
- 3 Close the battery pack cover and use a hex wrench to tighten the bolts that secure the battery pack.
- 4 Attach the new deco sheet to the battery pack cover.

6. Troubleshooting

6.1 Troubleshooting

6.1.1 Troubleshooting Guide

When you encounter problems while using the equipment, search for the table below for the problem or error messages and try the solutions. If the problem persists, turn off the detector and consult your sales representative or a distributor. Please refer to the details of the following symptoms or error messages.

- Troubleshooting must be performed by service engineer who is authorized by Vieworks. If an unqualified person performs troubleshooting on the system resulting in damaging the detector, software or hardware, then the Vieworks or its representative is not responsible for the detector repair regardless of remain warranty. For more detailed information, refer to <7.1 Service Information> and <7.2 Warranty>.
- Le dépannage doit être effectué par un technicien de maintenance agréé par Vieworks. Si une personne non qualifiée effectue un dépannage sur le système entraînant une détérioration du détecteur, du logiciel ou du matériel, Vieworks ou son représentant n'est pas responsable de la réparation du détecteur, quelle que soit la garantie restante. Pour plus d'informations, reportez-vous aux sections <7.1 Informations sur la maintenance> et <7.2 Garantie>.

6.1.2 Fail to Turn the Detector On

Category	Description
Symptom	<ul style="list-style-type: none"> • Failed to turn the power of the detector.
Expected Causes	<ul style="list-style-type: none"> • The battery pack is discharged. • The battery pack is attached wrongly. • The battery pack is broken. • The detector is broken.
Solutions	<ol style="list-style-type: none"> 1 Charge a battery pack 2 Attach a battery pack again. 3 Replace with another battery pack. 4 Replace with another detector.

6.1.3 The Power Switch of SCU or Status LED is not worked normally

Category	Description
Symptom	<ul style="list-style-type: none"> • The power switch of SCU is not working. • The status LED of SCU is not responding.
Expected Causes	<ul style="list-style-type: none"> • Problems with the connection of power cable. • Power cable is broken. • AC-DC adapter is broken. • LAN cable is broken. • SCU is broken.

Solutions	1 Check the cable connections (AC-DC adapter, power cable, LAN cable) of the SCU.
	2 Turn on and off SCU.
	3 Replace the power cable.
	4 Replace the AC-DC adapter.
	5 Replace the LAN cable.
	6 Replace SCU.

6.1.4 Communication Test is failed

Category	Description
Symptom	<ul style="list-style-type: none"> • The communication connection fails or is lost. • A transmission error occurs and abnormal images or data are sent.
Expected Causes	<ul style="list-style-type: none"> • Network connection problem • Network setting problem • PC environment setting problem • Wireless environment environment problem • SCU (or external AP) failure • Detector failure
Solutions	<ol style="list-style-type: none"> 1 Check the connection of network cable between workstation and SCU (or external AP). 2 Check if the accurate network cable is used or not. (CAT 5E or 6) 3 Set the network information of workstation, SCU (or external AP) and detector again. 4 Set the workstation environment again such as firewall setting. 5 Check surrounding wireless communication environment. If there is an environment that prevents wireless communication, remove the obstacle or change the wireless communication settings to prevent the interference. 6 Reboot the detector and SCU. 7 Replace SCU (or external AP). 8 Replace the detector.

6.1.5 There is no problem in communication connection, but no image is acquired

Category	Description
Symptom	<ul style="list-style-type: none"> • The image cannot be sent to PC even though the X-ray is exposed.
Expected Causes	<ul style="list-style-type: none"> • Problems with the generator interface • Problems with the generator interface connection • Generator interface cable failure • SCU failure • Detector failure
Solutions	<ol style="list-style-type: none"> 1 Check the generator interface setting of the detector. (AED mode, DR Trigger mode)

-
- 2 In DR trigger mode, check the generator interface connection.
 - 3 Replace the generator interface cable.
 - 4 Replace the SCU.
 - 5 Replace the detector.
-

6.1.6 Getting Abnormal Images

Category	Description
Symptom	<ul style="list-style-type: none"> • The image sent to the PC is abnormal.
Expected Causes	<ul style="list-style-type: none"> • Problems with calibration data • Problems with setting exposure time • Wrong operating environment • Detector failure.
Solutions	<ol style="list-style-type: none"> 1 Perform the calibration again. 2 Make sure that the detector is set to be applied calibration data. 3 Set the exposure time to be longer than the X-ray exposure time. 4 Check the operating environment (temperature, humidity, shock, vibration, X-ray energy range). 5 Replace the detector.

6.1.7 Battery Pack or Installation Part of Battery is Getting Hot

Category	Description
Symptom	<ul style="list-style-type: none"> • The fully charged battery pack runs down quickly.
Expected Causes	<ul style="list-style-type: none"> • Reduced performance (capacity) with prolonged use of the battery pack • Using the battery pack in low or high temperature environments
Solutions	<ol style="list-style-type: none"> 1 If the battery pack is to be used for a long time, replace it with a new battery pack. (The battery pack is a consumable.) 2 Use the battery pack at room temperature. In cold or hot environments, the battery pack's capacity decreases.

7. Information

This section gives overview information for service and warranty of the product.

Service Information

Warranty

Revision History

7.1 Service Information

7.1.1 Product Lifetime

The estimated product lifetime may be up to seven (7) years under the appropriate regular inspection and maintenance.

7.1.2 Regular Inspection and Maintenance

In order to ensure the safety of patients, operating personnel and third parties, and to maintain the performance and reliability of the equipment, be sure to perform regular inspection at least once a year. If necessary, clean up the equipment, make adjustments, or replace consumables.

There may be cases where overhaul is recommended depending on the conditions. Contact your sales representative or distributor for regular inspections or maintenance.

7.1.3 Repair

If a problem cannot be solved even after taking the measures indicated in Troubleshooting and contact your sales representative or a distributor for repairs. Please refer to the name label and provide the following information.

- Model name
 - **FXRD-3643VAW / FXRD-3643VAW PLUS**
 - **FXRD-4343VAW / FXRD-4343VAW PLUS**
 - **FXRD-2530VAW / FXRD-2530VAW PLUS**
- Serial number
 - 9 digit-number on the product label
- Explanation of problem
 - Describe as detailed as possible.

7.1.4 Replacement Parts Support

Performance parts (parts required to maintain the functioning of the product) of this product will be stocked for seven years after discontinuance of production, to allow for repair.

7.1.5 Consumables

The following consumable can deteriorate because of its characteristics and structure. For purchase of consumables, contact your sales representative or distributor.

- Battery pack: **FXRB-04A**

7.2 Warranty

Vieworks Co., Ltd. follows the contract with each customer regarding the product warranty period. If any such product proves defective during this warranty period, Vieworks at its option, either will repair the defective product without charge for parts and labor, or will provide a replacement in exchange for the defective product. In order to obtain service under this warranty, Customer must notify Vieworks of the defect before the expiration of the warranty period and make suitable arrangements for the performance of service. Customer shall be responsible for packaging and shipping the defective product to the service center designated by Vieworks with shipping charges prepaid.

Vieworks shall pay for the return of the product to customer if the shipment is to a location within the country in which Vieworks designated service center is located. Customer shall be responsible for paying all shipping charges, duties, taxes, and any other charges for products returned to any other locations.

This warranty shall not apply to any defect, failure, or damage caused by improper or inadequate maintenance and care. Vieworks shall not be obligated to furnish service under this warranty to repair damage resulting from attempts by personnel other than Vieworks or its representatives to install, repair, or service this product, to repair damage resulting from improper use or connection to incompatible equipment or power source; or to service a product that has been modified or integrated with other products when the effect of such modification or integration increases the time or difficulty of servicing the product.

THIS WARRANTY IS GIVEN BY VIEWORKS WITH RESPECT TO THIS PRODUCT IN LIEU OF ANY OTHER WARRANTIES, EXPRESS OR IMPLIED. VIEWORKS AND ITS VENDOR DISCLAIM ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. VIEWORKS RESPONSIBILITY TO REPAIR OR REPLACE DEFECTIVE PRODUCTS IS THE SOLE REMEDY PROVIDED TO THE CUSTOMER FOR BREACH OF THIS WARRANTY. VIEWORKS AND ITS VENDORS WILL NOT BE LIABLE FOR ANY INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES IRRESPECTIVE OF WHETHER VIEWORKS OR THE VENDOR HAS ADVANCE NOTICE OF THE POSSIBILITY OF SUCH DAMAGES.

There are no warranties which extend beyond the description mentioned in this document

FCC Compliance

- This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of **FCC Rules**. These limits are designed to provide reasonable protection against harmful interference in a residential installation.
- Operation is subject to the following conditions.
 - This device may not cause harmful interference.
 - This device must accept any interference received, including interference that may cause undesired operation.
- This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measure.
 - Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from where the receiver is connected.
 - Consult the distributor or an experienced radio/TV technician for help.

For product available in the USA/Canada market, only channel 1~11 can be operated. Selection of other channels is not possible.

RF Exposure Statement :

The equipment complies with FCC RF radiation exposure limits set forth for an uncontrolled environment. This device and its antenna must not be co-located or operation in conjunction with any other antenna or transmitter.

This Mobile device has been evaluated for and shown compliant with the FCC RF Exposure limits. The unit of measurement for RF exposure is Specific Absorption Rate (SAR). The FCC SAR limits for mobile is 1.6W/Kg per 1g of tissue. For the SAR values specific to your device, see the Quick Start Guide or User Manual.

7.3 Revision History

Version	Date	Descriptions
1.0	2019-08-22	<ul style="list-style-type: none">Initial Release

VIEWWORKS

Vieworks Co., Ltd.

(Gwanyang-dong) 41-3, Burim-ro 170beon-gil, Dongan-gu, Anyang-si,

Gyeonggi-do, 14055 Republic of Korea

Telephone: +82-70-7011-6161

Fax: +82-31-386-8631

Homepage: <http://www.vieworks.com>

**European representative: Obelis s.a**

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1030 Brussels, BELGIUM

Tel: +(32) 2.732.59.54

Fax: +(32) 2.732.60.03

E-mail: mail@obelis.net