



iDream 5

HEAD MASSAGER

USER MANUAL

Thank you for purchasing Breo

Please read this manual carefully before using and keep it properly for future reference



Instruction

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capability or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- details regarding safe disposal of used batteries;
- the battery must be removed from the appliance before it is scrapped;
- the appliance must be disconnect from the supply mains when removing the battery;
- the battery is to be disposal of safely.

5124004141

Instruction



Correct Disposal of this product. This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device. Please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.


Contents

Product function	1
Inside the package	2
Product Deployment	3
Operation	4-5
Precautions	6
Troubleshooting	7


Product function



Point Massage

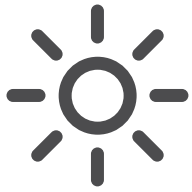

Heating

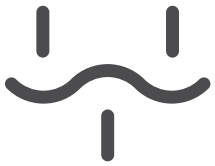



Air Pressure



Music

Inside the package

Your massager includes main unit, power adapter, earphone, MICRO SD CARD, user manual

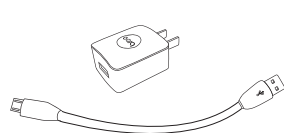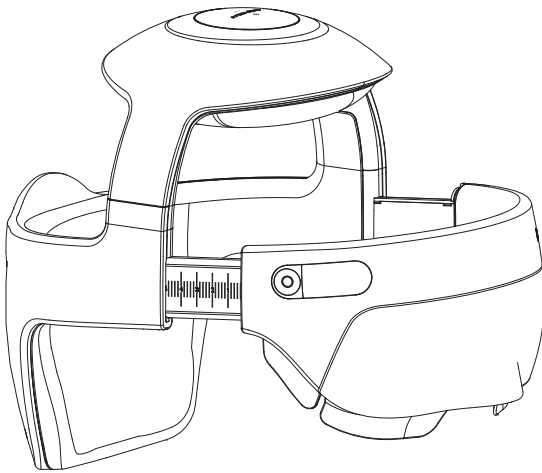

Power adapter




Earphone



User manual

Product Deployment

Main unit

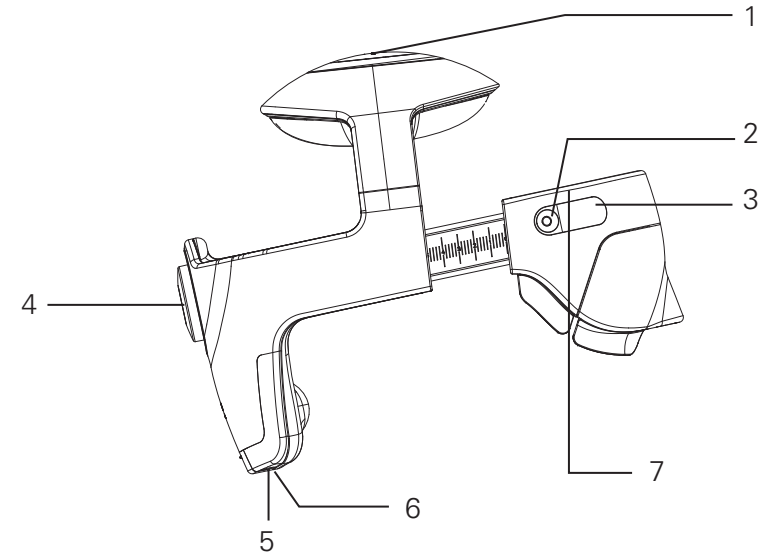

1. Massage head height adjust bottom
2. On/off /mode switch
3. OLED display
4. Size adjust Knob
5. Earphone jack
6. DC power interface
7. Disassemble button for eye part

Operation

Preparation before use

1. Check if there are missing parts of the package.
2. Make sure the main unit with enough power.
3. Make sure there is no damage of the adaptor.

Charge the unit

1. Open the cover of ⑥DC power interface , insert the adaptor.
2. When charging, the battery status on ③OLED display will flashing,when fully charged,it will stop flashing.

Attention

1. Make sure the main unit is power off when charging; Charging time is about 4 hours.
2. Do not use the product when charging.The main unit doesn't work while charging
3. Download APP for customizing massage.

How to wear

1. Press the button in the middle of ④ Size adjust Knob and wear the main unit.
2. Turn ④ Size adjust Knob to adjust the size to fit your head and eyes.

Operation

3. The top side size could be adjusted after machine on.
Short press ① Massage head height adjust botton for shorter the size(maximum 3 times), and long press for higher size.
4. The eye part could be disassembled before you wear the main unit by ⑦ Disassemble button for eye part.

Operation& Mode Selection

Short press ②On/off /mode switch, MED mode will on. Short press ② On/off /mode switch to switch the massage mode.

Operation& Mode Selection

MED	HARD
Air pressure	Air pressure
Heating (39-42 ℃)	Heating (39-42 ℃)
Music	Music
15 min	Kneading
	15min

Shutdown

Long press ② On/off /mode switch.

Precautions

1. Do not use the head massager when driving.
2. Do not use the massager in places too hot and wet.
3. Check before using the massager : No hairpin, hair accessory on your head and no wear contact lenses.
4. Wear the massager before switch on to avoid overpressure damage to the air bag.
5. The massager is suitable for adults of 56cm~64cm head circumference.
6. For those who suffer from heavy inflammation, blood stasis or brain surgery, do not use the massager.
7. For those who suffer from heavy slow action, poor physical and mental function or disturbance of intelligence, use under the supervision of companion, and stop using immediately if there is any discomfort.
8. For those who suffer from apoplexy, high blood pressure and low blood pressure, please use with caution and consult professional physician, and stop using immediately if there is any discomfort.
9. Use the stabilized DC5V adaptor, and use safely.
10. The product is not suitable for diagnosis, treatment and other commercial purposes.
11. Use sterile cotton with neutral detergent for cleaning, rather than water.
12. The massager has heating function, those who are insensitive to temperature shall pay attention.
13. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Troubleshooting

If there is any problem when you use the device, please follow the instructions below to check. Do not disassemble or repair the device yourself.

Faults	Possible reasons	Solutions
Unable to boot	Low battery	If yes, please charge the main unit
		If not, check the connection of charging cable
Air pressure or vibration is over low	Low battery	Please charge the machine

Technical parameters

DC	5V---
Power	5W
Ideal environment	Temperature: 5°C~40 °C; Relative Humidity: ≤ 80% Atmospheric Pressure: 860hpa~1060hpa
Dimension	L254mmXW210 mmXH242mm
Net weight	1075g

Storage

1. If the storage period exceeds 3 months, please check and recharge it.
2. Do not store close to floor, wall or roof.
3. Keep in good ventilation condition and away from direct sunshine and corrosive gases.

FCC Statement

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

IC Statement

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes : (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radio électrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.