


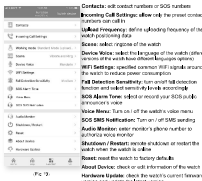

 <p>Smart A Fall Detection GPS Watch for Elderly BPC010</p>	<p>Introduction</p> <p>ProVista Care Smart A is a fall detection alert and position tracking device to provide seniors, especially those with Alzheimer's and cardiovascular diseases, an extra safety to lead their lives.</p> <p>Contact Us:</p> <p>If you have any questions about using ProVista Care of if you'd like any further details, please feel free to contact us or visit our website before for an up-to dated information, downloads and frequently asked questions.</p> <p>Email your questions to: info@provista.com</p> <p>Call us at: +852 65435320</p> <p>Website: www.provista.com</p> <p>Store nearest you can visit at: www.provista.com</p>	<p>In The Box</p> <p>Your Smart A2 pack contains the following items:</p> <ol style="list-style-type: none"> 1 x Smart A2 Watch 1 x 2-MHz rechargeable USB cable 1 x Manual 	<p>Watch Installation & Operation</p> <p>1. Charging</p> <p>Before you charge the watch with the 2x magnetic USB cable, when the watch is connected to the power, the screen will show light up and show a battery icon, (pic.1)</p> <p>2. SIM Card</p> <p>First open the SIM card cover. Slide the stem to left and lift. Insert a Nano-SIM card into the tray on the right direction. Then cover the stem and slide to the right. And close up the cover. (pic.2)</p> <p>(SIM card size service check is available. Data consumption about 30-100MB/month)</p> <p>3. Turn on / off</p> <p>Press the start button before turn on the watch. (pic.3)</p> <p>Turn on: press the power button for 3 seconds, and click the "Start/Stop" icon on the watch screen.</p> <p>Turn off: press the power button for 3 seconds, and click the "Start/Stop" icon on the watch screen.</p> <p>Force on/turn on: press the power button for 8 seconds.</p>	<p>4. Light up the watch screen</p> <p>Press the "Start/Stop" button to light up the watch screen, (pic.4)</p> <p>5. Phone Call / SOS Call</p> <p>Phone call: click the screen of the watch to select contact and tap and hold the call button to talk. (pic.5)</p> <p>Accord call: press the SOS button to accept call.</p> <p>Click the call: press the ring to button to decline call.</p> <p>SOS call: press the SOS button for 3 seconds to call emergency numbers. SOS call number and SOS number should be preset by user on the APP. (pic.10)</p> <p>6. Chat</p> <p>Send voice message: click the screen to [Chat] and click [add] icon, then press and hold the screen to talk, and release to send the message.</p>	<p>7. Heart Rate Measurement</p> <p>Place the watch on your wrist correctly. And slide the watch screen to [Heart Rate] and click to start measurement. The watch will display on the watch and connect to the APP.</p> <p>8. Step Counter</p> <p>Click the screen to [Steps] to check daily steps. The data will synchronize to the APP. The function should be preset on the APP. (pic.11) (pic.12)</p> <p>9. Weather</p> <p>Click the screen to [Weather] to check the weather and conditions in the city of the watch.</p>	<p>APP Setup</p> <p>1. Download APP</p> <p>Before use, please download a smart phone for remote monitoring (Android 4.4 and above, iOS 8 and above)</p> <p>Download from App Store or Google Play by searching "ProVista Care Smart A" or scan the APP QR code on the gift box or watch screen to download.</p> <p>2. Register</p> <p>Enter the APP and click "Register" in the mobile phone number. This verification option will be required. Set a password (20 characters), QIN (8 characters), (pic.13)</p>	<p>3. Add Device</p> <ol style="list-style-type: none"> 1. Click "Add Device" on the APP Home page. (pic.14) 2. Scan the QR code on the watch screen. "QR Code - Watch" (pic.15) 3. Enter your nickname for the watch and click "The go". (pic.16) 4. Enter the phone number of the watch, and the nickname of the watch. 5. Select time zone of the watch and click "Next". 	<p>(pic.16) After successfully set up the watch, the watch status will show on the homepage. (pic.17)</p> 
---	---	---	--	--	---	---	--	---

<p>Map (pic.18)</p> 	<p>Function (pic.19)</p> 	<p>Me (pic.21)</p> <p>Check on APP APP account information, and log out.</p> 	<p>FAQS</p> <p>Q. What kind of SIM card should be used for the device? And how many data consumption is needed?</p> <p>A. Nano-SIM card. Data service should be enabled. Data consumption is about 30-100MB/month.</p> <p>Q. Can I use the Device and APP without the watch?</p> <p>A. Yes. Check server in higher resolution. Please see the Read the User Guide and make sure the network support GSM.</p> <p>Q. Can I add multiple phones with the APP?</p> <p>A. Yes. Click the watch screen in Contacts page, then select and click the "number" you want to call. The number should be preset on the APP.</p> <p>Q. How long can I use the watch? Is it all through battery working to maintain the APP?</p> <p>A. The power consumption of the watch varies from the setting, such as GPS frequency. We also recommend charging the watch daily.</p> <p>Q. Why would the fall alarm triggered?</p> <p>A. When Device fall more than 30m and stay 3 seconds, will trigger the alarm.</p>	<p>Technical Information</p> <ul style="list-style-type: none"> ● Charging: 4.5V, DC 0.1-0.2V/1.000mAh (Magnetic USB Port) ● Battery: 3.6V non-rechargeable polymer battery on battery ● Built-in MC: 2540MHz, 32KHz, 32.768KHz, 32.768KHz ● Frequency: GPS (GPS), WLF, LBS ● Network Frequency: GSM/GPRS (850, 900, 1800, 1900), GPRS (Class 12) ● Network Standard: GPRS/GSM 1/2G ● Data Service: GPRS ● Call Time: About 4 hours ● Standby Time: 2 days ● Sleep Battery Indicator ● Waterproof: IP64 	<p>(DISCLAIMER AND IMPORTANT STATEMENT)</p> <p>We warrant the accuracy of the operating manual, but do not guarantee that there is no mistake in it.</p> <p>We do not warrant the accuracy of the software and hardware of the product without your authorization.</p> <p>The manual is part in whole, shall not be copied, reprinted, reproduced or stored in the retrieval system in any form, or stored into any database without the written authorization from ProVista Care (BPC010) or ProVista Group.</p> <p>You have no liability to any loss or damage of the product (Watch, Accessories, Software, Cable) and the third person. (BPC010) caused by fault or careless operation.</p> <p>We have no responsibility for any claims brought about by improper or careless installation or handling of the product.</p> <p>We have no liability for any loss or damage caused by the wireless or cable, network, server, phone or information download or upload.</p> <p>All the information provided by our website is for reference purpose, not for professional medical purpose. And we are not responsible for any errors, omissions or losses caused by the products.</p> <p>All the information provided here is the original one on the official website. It is for reference only. We are not liable to any mistake or correction.</p>	<p>FCC STATEMENT</p> <p>The equipment does not contain any FCC-certified components. The equipment is not intended for use in the United States or in any other country where FCC certification is required. The equipment is not intended for use in the United States or in any other country where FCC certification is required. The equipment is not intended for use in the United States or in any other country where FCC certification is required.</p> <p>Specific Absorption Rate (SAR) Declaration Statement</p> <p>The SAR value of this device is 0.0018 W/kg. The SAR value of this device is 0.0018 W/kg. The SAR value of this device is 0.0018 W/kg. The SAR value of this device is 0.0018 W/kg.</p>	<p>These FCC SAR requirements are achieved from the measurements of the equipment in accordance with the FCC Part 15.207 and 15.209. The SAR value of this device is 0.0018 W/kg. The SAR value of this device is 0.0018 W/kg. The SAR value of this device is 0.0018 W/kg. The SAR value of this device is 0.0018 W/kg.</p> <p>ProVista</p> <p>ProVista (Far East) Ltd.</p> <p>Website: www.provista.com</p> <p>Reg. No. 602 8322 000</p> <p>Email: info@provista.com</p> <p>Head Office: Unit 2, 2/F, COS Center, 35 Sun Yu Street, Kowloon, Hong Kong</p> <p>China: ProVista, Building Park, Xiangye Village, Shaoxing Town, Dongguan, Guangdong, China</p>	<p>ProVista</p> <p>ProVista (Far East) Ltd.</p> <p>Website: www.provista.com</p> <p>Reg. No. 602 8322 000</p> <p>Email: info@provista.com</p> <p>Head Office: Unit 2, 2/F, COS Center, 35 Sun Yu Street, Kowloon, Hong Kong</p> <p>China: ProVista, Building Park, Xiangye Village, Shaoxing Town, Dongguan, Guangdong, China</p>
--	--	--	---	---	---	---	---	--