
Quick Reference Booklet And User Guide

Version 1.4

SafeKare™ Model CRU101A

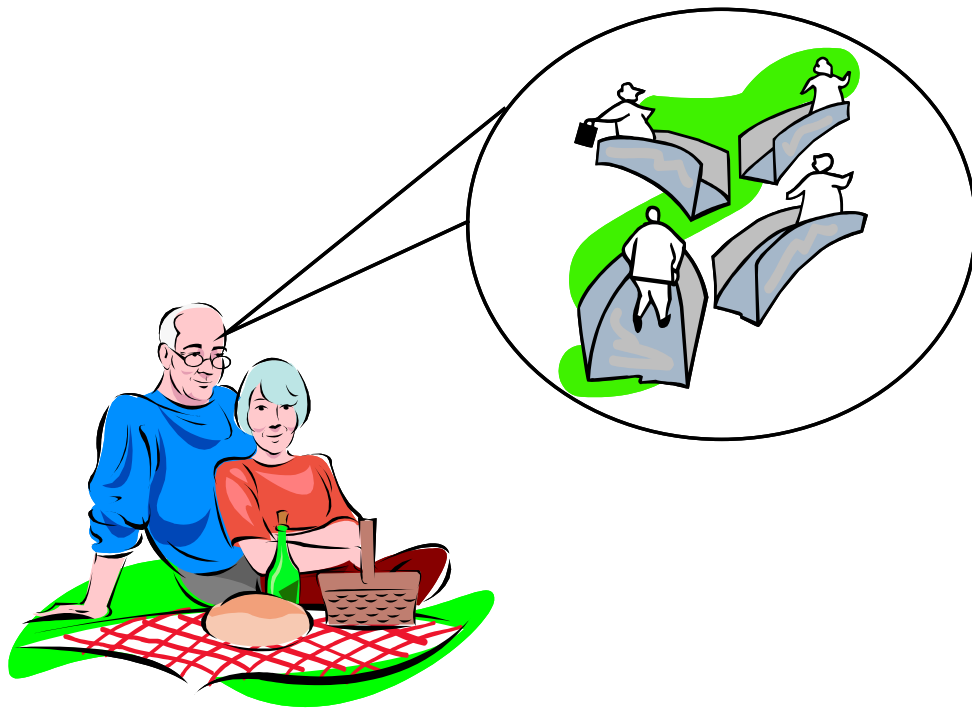SafeKare™

Alzheimer's disease or different forms of Dementia can affect all who have to deal with the care of individual family members stricken with this unfortunate condition. A major problem in providing home care for patients, is knowing that the patient is within the caregiving control area. All too often, for whatever reason, the patient will leave the control area without the caregiver's knowledge. This can be a dangerous situation for the patient, as well as a frustrating problem for the caregiver.

DT Communications Inc. has developed a simple method of alerting a busy caregiver whenever the patient leaves a defined control area.

Our company has designed a unit that can be used both inside and outside the home or caregiving facility. The unit will alert the caregiver if the patient leaves or moves outside the range of the SafeKare™ system.

Congratulations on the purchase of your SafeKare™ Caregiver Assistant. This information sheet will guide you through the Quick Start setup and operation of your SafeKare™ system. Refer to your User's Manual for details on setting up and operating the caregiver system.

“DO NOT” power on the Caregiver Receiver CRU101A until the battery has been charged. See Section 5 of the user’s manual for help if necessary.

Table of Contents

1 UNPACKING THE SYSTEM COMPONENTS 4

1.1 CHARGING THE CAREGIVER RECEIVER..... 5

1.2 UNPACKING THE BADGE/FOB 5

1.3 FOB OPTION: ACTIVATING THE PATIENT FOB..... 6

1.4 BADGE OPTION: ACTIVATION OF THE PATIENT BADGE 7

1.5 SYSTEM LEARN MODE..... 7

2 USING THE SAFEKARE CAREGIVER SYSTEM..... 8

2.1 CHECKING THE SYSTEM RANGE 8

2.2 PLACING THE BADGE OR FOB UNIT ON THE PATIENT..... 9

2.3 ENSURING CORRECT USAGE OF THE SAFEKARE™ SYSTEM..... 10

1 Packing List

Before unpacking this caregiver kit, take out the packing list. As you unpack the kit, check off the items listed on the packing list

Example packing list


	<i>Sent</i>	<i>Received</i>
Caregiver Receiver, CRU101A	<u> X </u>	<u> </u>
Wall Power Module, WPU101A	<u> X </u>	<u> </u>
Hanging Patient Fob, FOB101A	<u> X </u>	<u> </u>
Coin style Battery #CR2032	<u> X </u>	<u> </u>
Hanging Patient Badge, PB101A (Optional)	<u> </u>	<u> </u>
Base Stand for Caregiver Unit, BS101A (Optional)	<u> </u>	<u> </u>
Belt Clip For Receiver Unit BC101 (Optional)	<u> </u>	<u> </u>
Charging Plug For Car/Vehicle CP101 (Optional)	<u> </u>	<u> </u>
READY TO USE Option (Optional)	<u> </u>	<u> </u>

1.1 Unpacking the System Components

Unpack the box containing your SafeKare™ system. Ensure that all items are present and accounted for. Check your packing list to see what items were included in your system. The packing list for a basic system is shown below.

If a part of the system you ordered is missing, consult the User's Manual for instructions on contacting the company.

Figure 1. Location Of Receiver Functions

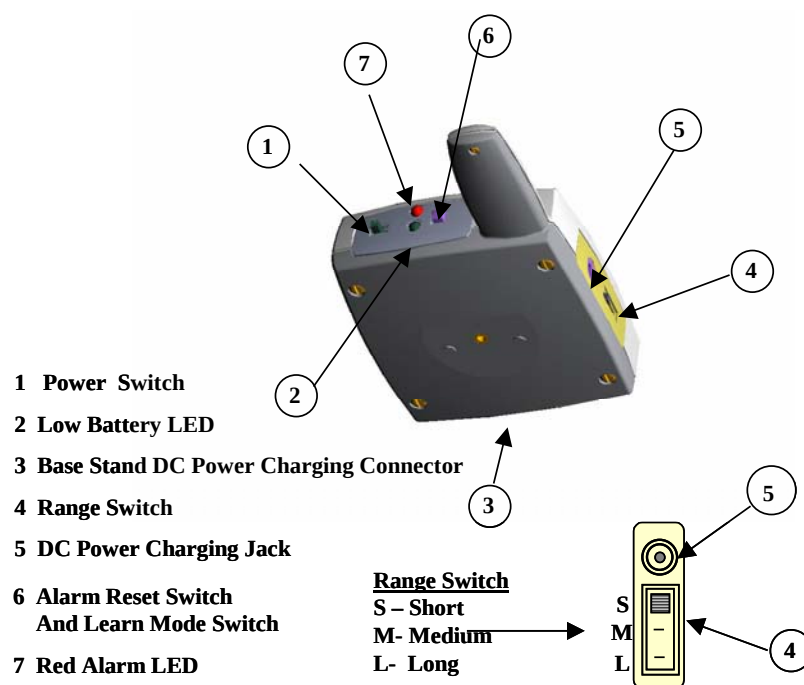Before your SafeKare system is used, the battery in the caregiver receiver must be fully charged.

The locations of the switches, indicators, and power jack, for the CRU101A receiver, are shown in Figure 1. Become familiar with each function location.

1.2 Charging the Caregiver Receiver

Plug the AC adapter into the wall, and then plug the other end into the caregiver receiver. See Figure 2 for reference.

A full charge will be reached in approximately 6-8 hours. The unit will not be damaged if left charging for longer than 8 hours.

After the caregiver receiver has been charged, the receiver is ready for use.

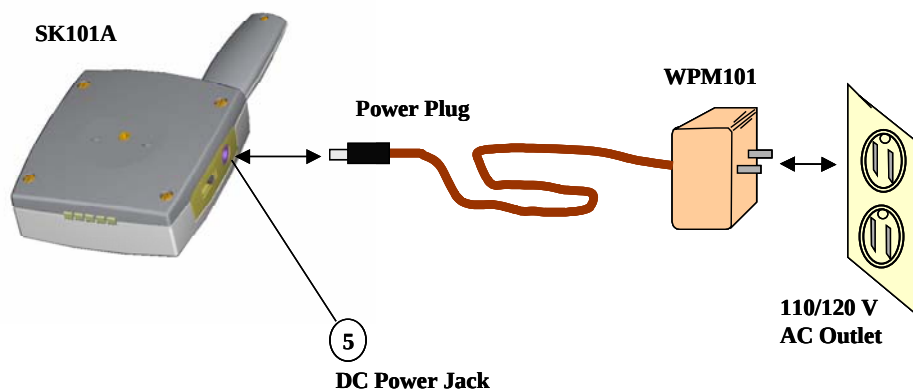
Figure 2. Charging the Battery

1.3 Unpacking the Badge/Fob

Remove the badge or fob from the metal case.

Do not dispose of this metal box. It may be needed to store or ship the fob or badge.

1.4 Fob Option: Activating the Patient Fob

Place the patient fob FOB101A face down on a flat table like surface. Remove the screw located in the center of the fob and remove the cover. See Figure 3.

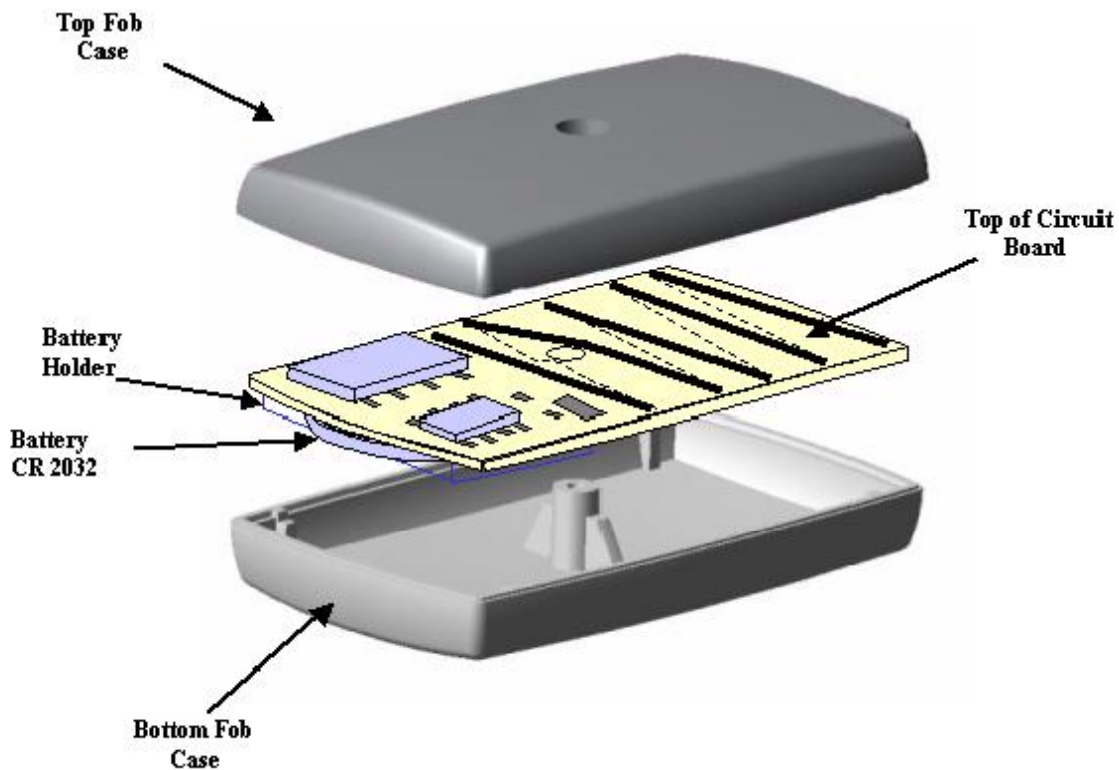
Figure 3 Removing The Back of The Fob

The fob battery sits on the backside of the circuit board. Remove the battery from the protective plastic pouch that came in your kit. Insert the battery as shown in Figure 4. Note, the battery must be placed into the holder with the minus (-) side of the battery against the circuit board. The plus (+) side of the battery faces away from the circuit board. Replace the circuit board after the battery is installed. Replace the cover and screw the back together. The fob is now ready to use.

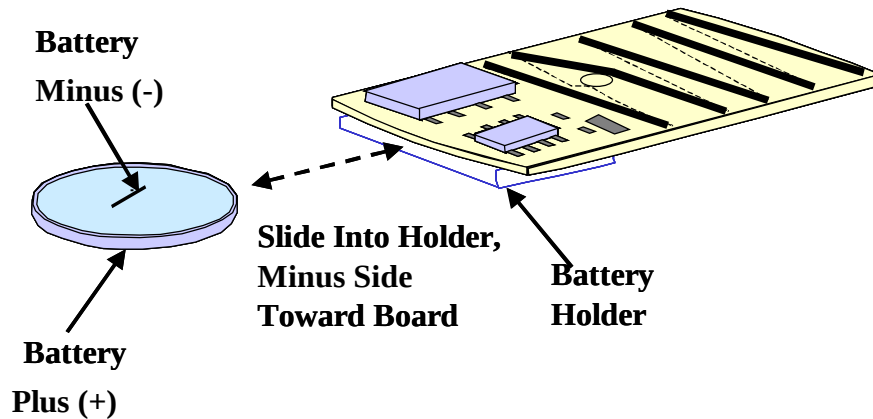
Figure 4. Inserting the Battery

1.5 Badge Option: Activation of the Patient Badge, OPTIONAL Equipment

The patient badge comes with a preinstalled battery. The battery has not been factory activated in order to preserve battery life. On the backside of the badge, locate the circle about ½ inches in diameter with the inscription START located near the middle of the badge. Hold the badge firmly with one hand and with the other hand use a finger to press the circle firmly for 6 to 7 seconds. The badge should now be active and ready for use

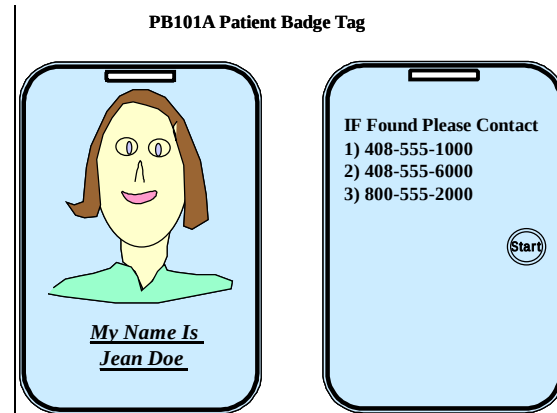
Figure 5, PB101A Hanging Badge

1.6 System Learn Mode

Before the caregiver receiver can operate with the patient badge or fob, it must first learn the badge or fob ID code. Place the badge or fob 2 – 3 feet away from the caregiver receiver. Refer to figure 6. Turn on the caregiver receiver (refer to figure 1). Wait for receiver to initialize, indicated when the two LED's stop flashing. Locate the “Reset-Learn” button on the caregiver receiver. Hold caregiver receiver in the upright position and press down on the ‘Reset/Learn” button switch. Hold the button down until an audible beep is heard from the caregiver receiver. Within 3 to 15 seconds a second audible sound (a long beep) should be heard. This indicates the receiver has learned the badge/fob code. If the receiver starts to alarm, repeat the learn function. If the alarm condition continues, consult the Trouble shooting section of the User’s Manual

Press RS/Learn Switch, Hold Down Until “Beep”

Figure 6. Setting Up Receiver Learn Mode

2 Using the SafeKare Caregiver System

A detailed description of using the CRU101A system is described in the User's Manual. Refer to the User's Manual if you have any problems using the "Quick Start User Guide." A short system check of the receiver operating at different range settings is described in the next section. For a detailed description see Section 6 of the User's Manual. Note the locations where the receiver starts to alarm as you move away from the badge or fob. This will be an indication of how the unit will operate in your home environment.

2.1 Checking the System Range

Refer to Figure 1 for the following instructions. Turn the caregiver receiver on, and place the badge or fob unit on a table or similar object. Set the switch located on the side of the receiver to the "S" position (short range setting).

Walk very slowly away from the badge or fob and record the distance at which the caregiver receiver begins to alarm (a short beep, followed by a continuous beep 12 seconds after the short beep).

Switch the receiver to the "M" position (medium range setting). The beeping should stop within 12 to 15 seconds. Continue moving *slowly* away from the badge or fob. Again, record the distance at which the caregiver receiver begins to alarm (a short beep).

Switch the receiver range to the "L" position (long range setting). The beeping should stop in about 12 to 15 seconds. Continue walking *slowly* away from the fob. Again, record the distance at which the receiver alarms, a short "beep" followed by a "continuous beeping" after 12 seconds. The receiver is now in full alarm mode.

Walk back towards the badge/fob unit. The receiver should stop alarming when you are approximately 50 feet away from the badge or fob.

! WALK VERY SLOWLY during the system range check

It is important to move slowly during the test. The unit responds within a 12 second time period. If you move too quickly, you will not be able to get an accurate measurement of the operating distance in your house or facility. The 12 second time delay is a requirement placed on the system by the FCC, but in a practical sense, this delay cuts down on the number of receiver false alarms. This time delay will not degrade the performance of the SafeKare™ System.

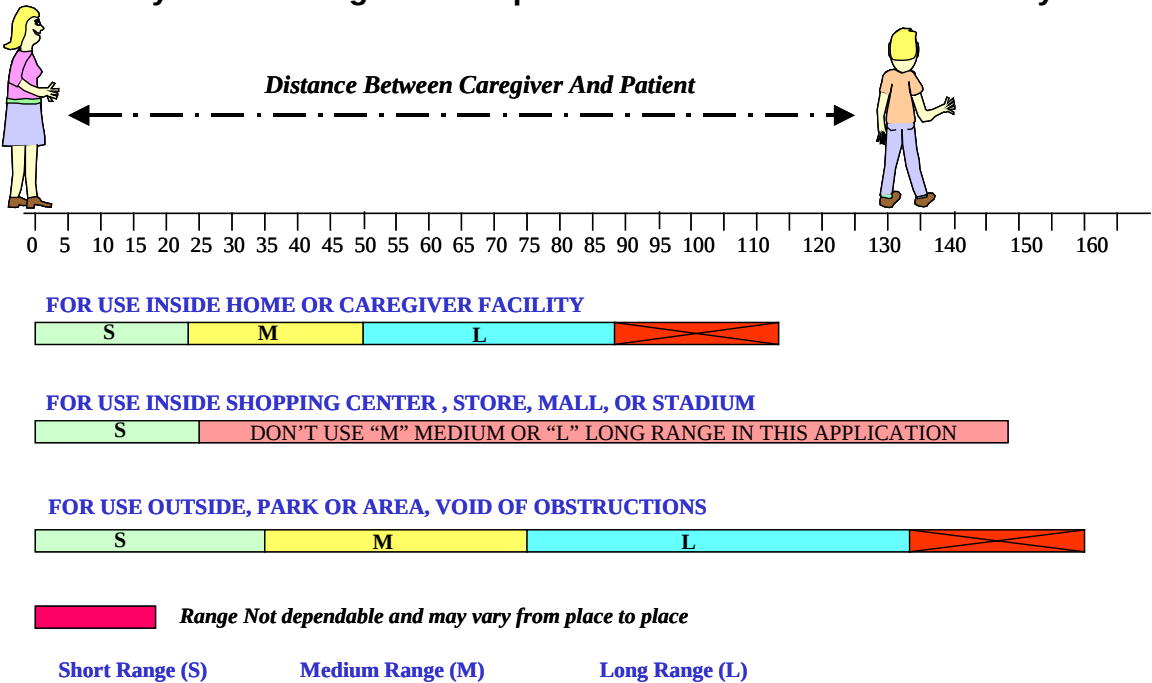
Figure 7. Distance Range Check

2.2 Placing the Badge or Fob Unit on the Patient

The range of the patient badge/fob is affected by how the patient fob or badge is worn. An example of the various ways the badge or fob can be worn is shown in Figure 8. Note that the badge or fob unit is always shown in the hanging position. Consult the User's Manual for more details.

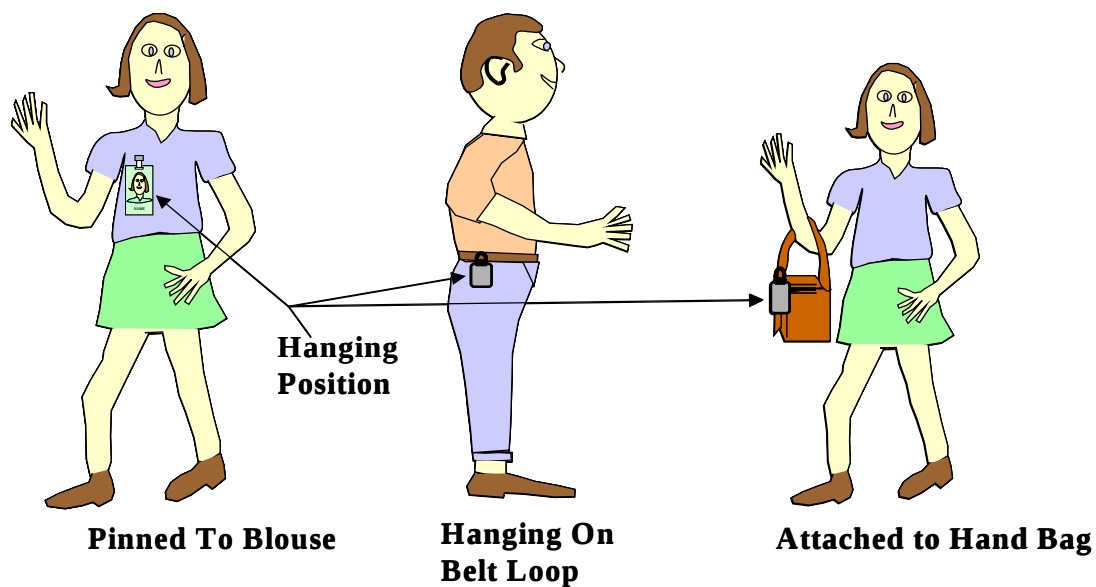
Figure 8. Placing the Badge or Fob on the patient.

After completing the range check, the caregiver unit is ready to use. To begin, it is suggested that the unit be switched to the “M” (medium range setting). After you have become familiar with using the caregiver system, you can make adjustments to suit your environment

2.3 Ensuring Correct Usage of the SafeKare™ System

Some tips on using the system to ensure the most effective usage.

- Never turn off the caregiver receiver when using the system.
- Charge the caregiver receiver at night, but keep it turned on if still monitoring a patient at night.
- Never let the battery run completely dead. This could damage the battery.
- Check on the patient periodically to ensure that the badge or fob is still attached to the patient.
- When using the system outside of the home, use the short range “S” mode.
- Consult the User’s Manual if you have problems.

Notice, FCC Compliance

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING

The user of these device is cautioned that: any unauthorized modifications or changes to these electrical devices, not expressly approved by the party responsible for FCC compliance of these devices, is prohibited and could void the user's authority to operate these device.

WARNING

This Device, the SafeKare 101A System, is an “AIDE” in the supervision of a patient and should not, for any reason, be used to replace active caregiver physical monitoring of a patient.

WARNING

The user or operator of the CRU101A caregiver receiver devices is warned that all warranties for the receiver are VOID if the receiver is opened for any reason.

NOTICE FOR USER:

The USER’S MANUAL should be read for a better understanding of the operation of your SafeKare™ system. Refer to the User’s Manual for any trouble, and for contacting the company, if the need arises.

DT Communications, Inc.
245-M Mt. Hermon Rd.
PMB #305
Scotts Valley, CA. 95066

Email: info@alzsafe.com

Web Site: www.alzsafe.com

NOTES

SafeKare™ is a registered trademark of DT Communications

NOTES:

SafeKare[™] is a registered trademark of DT Communications

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